

# ***Spartan'10***

*for Windows, Macintosh and Linux*

## Tutorial and User's Guide



**Wavefunction, Inc.**

18401 Von Karman Avenue, Suite 370  
Irvine, CA 92612 U.S.A.  
[www.wavefun.com](http://www.wavefun.com)

**Wavefunction, Inc., Japan Branch Office**

3-5-2, Kouji-machi, Suite 608  
Chiyoda-ku, Tokyo, Japan 102-0083  
[japan@wavefun.com](mailto:japan@wavefun.com) • [www.wavefun.com/japan](http://www.wavefun.com/japan)

Copyright © 2010 by Wavefunction, Inc.

All rights reserved in all countries. No part of this book may be reproduced in any form or by any electronic or mechanical means including information storage and retrieval systems without permission in writing from the publisher, except by a reviewer who may quote brief passages in a review.

ISBN978-1-890661-41-4

Printed in the United States of America

# Acknowledgements

---

The *Spartan'10* Tutorial and User's Guide was prepared by Warren Hehre and Sean Ohlinger, who take both credit for its utility and blame for its limitations and inaccuracies. Several members of Wavefunction contributed significantly, with Philip Klunzinger, Bernard Deppmeier, Andy Driessen and Jeffrey Johnson warranting special mention. As with past manuals, sincere thanks goes to Pamela Ohsan for turning a sloppy manuscript into a finished book.



# Scope of this Guide

---

This guide provides a general reference for *Spartan*'10 for Windows, Macintosh and Linux. It is divided into 22 chapters grouped into three sections, along with several appendices.

**Section I (Introduction, Chapter 1)** introduces *Spartan* as a tool for exploring organic, bioorganic, inorganic and organometallic chemistry by way of molecular mechanics and quantum chemical calculations, together with an array of graphical models for conveying the results of these calculations. The section ends with a brief discussion of new capabilities in *Spartan*'10 as well as changes and improvements over previous versions.

**Section II (Getting Started, Chapters 2 to 11)** describes the overall operating environment of *Spartan*, and then provides an extensive set of hands-on tutorials. This section is the place to start for new users of the program and should also be examined by users of previous versions of *Spartan*.

**Section III (Features and Functions, Chapters 12 to 22)** describes in detail the functions available from the menus and dialogs incorporated into the graphical user interface for *Spartan*. The focus is on graphical input and manipulation of structure, input of other required information and text, spectral and graphical output resulting from molecular mechanics and quantum chemical calculations, and on use of databases of previously-calculated structures, energies, properties and spectra accessible from *Spartan*. This section is intended as a general reference to *Spartan*'10.

What this guide *does not do* is document the performance and cost (in computation time) of the different molecular mechanics and quantum chemical models available in *Spartan*'10, or recommend specific models or combinations of models for use on chemical problems. Nor does it show the utility of graphical models in presenting and interpreting the results of the calculations. These topics are touched on in *Topics* available from the **Activities** menu and are covered in

depth in *A Guide to Molecular Mechanics and Quantum Chemical Calculations* included as a PDF under the **Help** menu and available from Wavefunction as a hardbound volume. The guide also provides a collection of illustrative examples.

Appendices provide an overview of the program's overall architecture as well as its present capabilities and limitations (**A**), a directory of functions under its menus (**B**), a listing of commonly-used options (**C**), a listing of units (**D**), the proper citation for the program (**E**), instructions for installing the Cambridge Structural Database (**F**), directions for making databases from *Spartan* calculations (**G**), examples of pharmacophore input (**H**), input of experimental infrared, UV/visible and NMR spectra (**I**), specification of steps and empirical parameters involved for calculation of T1 heats of formation (**J**), specification of empirical parameters for correction of  $^{13}\text{C}$  chemical shifts obtained from B3LYP/6-31G\* and EDF2/6-31G\* calculations (**K**), directions for installing a network HASP (**L**) and a listing of NMR chemical shift standards (**M**). Additional materials relating to several of these appendices may be found as PDFs under the **Help** menu.

An up-to-date version of this *Tutorial and User's Guide* is available on Wavefunction's website (directly accessible from the **Help** menu).

# Table of Contents

<b>Section I</b>	<b>INTRODUCTION .....</b>	<b>1</b>
<b>Chapter 1</b>	<b><i>Spartan</i>'10 .....</b>	<b>9</b>
<b>Section II</b>	<b>GETTING STARTED .....</b>	<b>17</b>
<b>Chapter 2</b>	<b>Operating Spartan .....</b>	<b>19</b>
	Opening and Exiting <i>Spartan</i> .....	19
	Menus .....	19
	Keystroke Equivalents .....	21
	Icons .....	21
	Tabs .....	22
	Mouse and Keyboard Operations .....	23
	Selection .....	23
	Manipulation .....	24
	Selecting Molecules, etc .....	25
	Dialogs .....	26
	Database Access .....	27
	Stereo Displays .....	27
	Changing Colors and Setting Preferences .....	27
	Monitoring and Terminating Jobs .....	27
<b>Chapter 3</b>	<b>Basic Operations .....</b>	<b>29</b>
<b>Chapter 4</b>	<b>Organic Molecules .....</b>	<b>47</b>
	Acrylonitrile .....	49
	Cyclohexanone .....	57
	Camphor .....	60
	Ethinamate .....	63
	3-Cyano-4-methylcyclohexenyl Radical .....	65
	Androsterone .....	68
<b>Chapter 5</b>	<b>Groups of Organic Molecules .....</b>	<b>71</b>
	Dienophiles in Diels-Alder Cycloadditions .....	73
	Addition vs. Substitution .....	79
	Allyl Vinyl Ether .....	81
	Internal Rotation in Dimethylperoxide .....	84

	Hydration of Carbonyl Compounds.....	88
	Acidities of Carboxylic Acids.....	92
	Positional Selectivities of Substituted Naphthalenes.....	95
	Tautomers of Nucleotide Bases .....	98
<b>Chapter 6</b>	<b>Spectra of Organic Molecules.....</b>	<b>101</b>
	Infrared Spectrum of Methyl Formate.....	103
	Searching Spartan's Infrared Spectral Database..	106
	Proton NMR Spectrum of 2-Norbornene.....	110
	<sup>13</sup> C NMR Spectrum of Caulophylline.....	115
	<sup>13</sup> C NMR Spectrum of <i>cis</i> -1,2-Dimethylcyclohexane .....	117
	Stereochemical Assignments from <sup>13</sup> C Spectra....	119
	<sup>13</sup> C Chemical Shifts Depend on Conformation....	121
<b>Chapter 7</b>	<b>Organic Reactions.....</b>	<b>123</b>
	Ene Reaction of 1-Pentene.....	125
	S <sub>N</sub> 2 Reaction of Bromide and Methyl Chloride...	129
	Carbene Additions to Alkenes.....	133
	Stereospecific Diels-Alder Reactions .....	137
	Thermodynamic vs. Kinetic Control .....	140
	Activation Energies of Diels-Alder Reactions....	143
<b>Chapter 8</b>	<b>Medicinal Chemistry .....</b>	<b>147</b>
	Anticipating Blood-Brain Transport.....	149
	Terfenadine. A Potassium Channel Blocker?.....	152
	Morphine. Structure vs. Pharmacophore .....	155
<b>Chapter 9</b>	<b>Polypeptides to Proteins .....</b>	<b>163</b>
	Polyglycine .....	165
	Gleevec. Protein-Bound vs. Free Conformer.....	169
	Gleevec. Making a Pharmacophore from PDB ...	173
<b>Chapter 10</b>	<b>Inorganic and Organometallic Molecules.....</b>	<b>177</b>
	Sulfur Tetrafluoride.....	179
	Reactivity of Silicon-Carbon Double Bonds .....	182
	Benzene Chromium Tricarbonyl.....	185
	Ziegler-Natta Polymerization of Ethylene.....	187



<b>Chapter 11</b>	<b>“Dry Labs”: Using the Spartan Spectra and Properties Database .....</b>	<b>189</b>
	Isomeric C <sub>5</sub> H <sub>8</sub> Dienes .....	191
	Using Infrared Spectroscopy to Identify an Unknown Ethyl Benzoate Derivative .....	193
	Infrared Spectra of Short-Lived Molecules .....	196
	Using <sup>13</sup> C NMR to Distinguish Structural Isomers.....	198
	Using <sup>13</sup> C NMR Spectra to Distinguish Stereoisomers.....	201
	Finding Stable Enols.....	203
	Electrophilic Reactivity of Polycyclic Aromatics.....	205
<b>Section III.</b>	<b>FEATURES AND FUNCTIONS .....</b>	<b>207</b>
<b>Chapter 12</b>	<b>The File Menu .....</b>	<b>209</b>
	New.....	209
	Open.....	209
	Close .....	210
	Save.....	211
	Save As .....	211
	New Molecule.....	211
	Delete Molecule.....	212
	Append Molecule(s) .....	212
	Access PDB Online .....	212
	Print/Print Output/Print Spreadsheet .....	213
	Embedded Data.....	214
	Start/Stop QuickTime Recording.....	215
	Exit.....	215
<b>Chapter 13</b>	<b>The Edit Menu .....</b>	<b>217</b>
	Undo .....	217
	Cut.....	217
	Copy.....	217
	Paste.....	217
	Find.....	218
	Find Next .....	218
	Center.....	218

	Clear.....	218
<b>Chapter 14</b>	<b>The Model Menu.....</b>	<b>219</b>
	Wire.....	220
	Ball and Wire .....	220
	Tube .....	222
	Ball and Spoke .....	222
	Space Filling.....	223
	Line .....	223
	Hide.....	224
	Global Model .....	224
	Coupled.....	224
	Hydrogens.....	225
	Labels.....	225
	Ribbons .....	225
	Ramachandran Plot.....	225
	Hydrogen Bonds .....	226
	CFD's .....	226
	Configure .....	228
	Configure Labels .....	228
	Configure Objects.....	229
	Configure Ribbons .....	229
	Configure CFD's .....	230
<b>Chapter 15</b>	<b>The Geometry Menu.....</b>	<b>233</b>
	Measure Distance.....	233
	Measure Angle .....	233
	Measure Dihedral.....	233
	Freeze Center .....	234
	Set Torsions.....	236
	Set Similarity Centers .....	238
	Similarity Based on Structure .....	238
	Similarity Based on CFD's.....	239
	Similarity Based on a Pharmacophore .....	240
	Constrain Distance.....	240
	Constrain Angle .....	240
	Constrain Dihedral.....	240

Define Point .....	243
Define Plane .....	243
Define CFD .....	244
Align .....	244
Align Based on Structure .....	245
Align Based on CFD's .....	246
Align to a Pharmacophore .....	247
<b>Chapter 16 The Build Menu .....</b>	<b>249</b>
Model Kits .....	249
Organic Model Kit .....	250
Groups .....	251
Rings .....	252
More .....	252
Clipboard .....	253
Inorganic Model Kit .....	254
Ligands .....	256
Peptide Model Kit .....	257
Specification of Amino Acid Sequence .....	258
Specification of Macroscopic Structure .....	258
Termination .....	259
Nucleotide Model Kit .....	260
Specification of Base Sequence .....	261
Specification of Helical Structure .....	262
Substituent Model Kit .....	262
Library Preparation .....	264
Attachment .....	266
List Generation .....	266
Accessing ChemDraw .....	267
General Molecule Building Functionality .....	268
Multiple Fragments .....	268
Bond Rotation/Bond Stretching .....	268
Atom/Fragment Replacement .....	269
Chirality Inversion .....	269
Building/Editing Menu Functions .....	270
View .....	270
Add Fragment .....	270

Delete .....	271
Make Bond .....	271
Break Bond.....	272
Minimize .....	272
<b>Chapter 17 The Setup Menu .....</b>	<b>273</b>
Calculations .....	274
Calculate.....	276
Task.....	277
State .....	280
Type of Calculation.....	281
Basis Sets.....	287
Pseudopotentials .....	288
Solvent .....	288
Start From.....	289
Subject to.....	290
Compute .....	291
IR .....	291
Raman.....	292
NMR.....	292
UV/vis.....	293
QSAR.....	293
Total Charge .....	294
Multiplicity.....	294
Print .....	294
Orbitals and Energies.....	294
Thermodynamics .....	294
Vibrational Modes .....	295
Atomic Charges .....	295
Options .....	295
Converge .....	295
Similarity .....	295
Global Calculations.....	299
Surfaces.....	299
Common Surfaces and Property Maps.....	306
Surface.....	307
Property .....	308

Resolution .....	309
IsoValue.....	309
Global Surfaces .....	310
Submit.....	310
<b>Chapter 18 The Display Menu.....</b>	<b>313</b>
Output .....	313
Properties .....	314
Molecule Properties, QSAR Descriptors and Thermodynamics.....	315
Molecule .....	316
QSAR.....	317
Thermodynamics .....	318
Molecule Utilities .....	318
Atom Properties .....	319
Bond Properties.....	320
Constraint Properties.....	320
Point and Plane Properties .....	322
Surface Properties .....	322
Curve and Plot Properties.....	323
Substituent Properties.....	327
CFD Properties.....	328
Regression Properties.....	329
“Selected” Style .....	329
Orbital Energies .....	330
Surfaces.....	330
Spectra .....	330
IR.....	331
Raman .....	331
NMR.....	335
UV/vis .....	339
Formulas .....	341
Spreadsheet.....	346
From the Add Dialog.....	347
From Post (P) Buttons.....	350
From the Clipboard .....	351
User-Defined Expressions .....	352

Numerical Data .....	352
Plots .....	355
XY Plots .....	356
XYZ Plots.....	357
Similarities .....	359
Reactions.....	360
<b>Chapter 19 The Search Menu .....</b>	<b>363</b>
Define Search Query.....	363
Define Reaction Query.....	364
Open/Close Databases .....	367
Cambridge Structural Database (CSD) .....	367
Spartan Molecular Database (SMD) and Spartan Spectra and Properties Database (SSPD).....	372
Accessing and Processing Data from SMD and SSPD.....	375
Molecule Replacement.....	375
Substructure and Other Modes of Searching SMD and SSPD.....	377
Plotting Data from an SMD or SSPD Search.....	381
Regression Analysis of Data Resulting from an SMD or SSPD Search .....	383
Relationships Involving Chemical Reactions.....	384
Retrieving Data from an SMD or SSPD Search.	387
Spartan Reaction Database (SRD) .....	387
Substructure Searching.....	388
Retrieving Data from an SRD Search.....	390
Spartan Infrared Database (SIRD) .....	391
Spectra Searching .....	391
Filters .....	393
Retrieving Data from an SIRD Search .....	396
Guess Transition State .....	397
Identify Tautomers .....	398
Extract Ligands .....	400

<b>Chapter 20</b>	<b>The Options Menu .....</b>	<b>403</b>
	Preferences.....	403
	Settings.....	404
	Molecule.....	407
	VDW Radii.....	408
	Paths.....	409
	Hosts.....	410
	Jobs.....	411
	Parallel.....	412
	X-Platform.....	413
	Toolbars .....	414
	Colors.....	414
	Fonts/Graphics Fonts .....	414
	Monitor .....	415
	Calculator.....	416
	Model Kit.....	416
	File Toolbar.....	416
	Geometry Toolbar .....	416
	Build Toolbar .....	416
	Search Toolbar .....	416
	Cascade .....	416
	Tile .....	416
<b>Chapter 21</b>	<b>The Activities Menu .....</b>	<b>417</b>
	Tutorials.....	417
	Topics.....	417
	Lookup in Wikipedia .....	418
<b>Chapter 22</b>	<b>The Help Menu .....</b>	<b>419</b>
	Help.....	419
	About .....	420
<b>Appendix A</b>	<b>Capabilities and Limitations.....</b>	<b>421</b>
	Molecular Mechanics Models.....	421
	Semi-Empirical Models .....	421
	Hartree-Fock Models .....	421
	Density Functional Models.....	422
	Møller-Plesset Models .....	422

	Advanced Correlated Models .....	424
	Thermochemical Recipes.....	424
	Excited-State Models .....	425
	Basis Sets and Pseudopotentials .....	425
	Solvent Models .....	426
	Properties and Spectra .....	426
	Graphical Models.....	427
	Similarity Analysis.....	428
<b>Appendix B</b>	<b>Menus .....</b>	<b>429</b>
	<i>Spartan</i> Screen .....	429
	Contextual.....	435
<b>Appendix C</b>	<b>Commonly-Used Program Options.....</b>	<b>437</b>
	Conformational Search .....	437
	Similarity Analysis.....	438
	Geometry Optimization .....	438
	Molecular Mechanics Calculations.....	439
	Quantum Chemical Calculations .....	439
	Semi-Empirical Calculations .....	441
	Property and Spectra Calculations.....	441
	Miscellaneous and Printing.....	443
<b>Appendix D</b>	<b>Units .....</b>	<b>445</b>
<b>Appendix E</b>	<b>Citation .....</b>	<b>447</b>
<b>Appendix F</b>	<b>Installing the Cambridge Structural Database.....</b>	<b>449</b>
<b>Appendix G</b>	<b>Constructing Custom Databases .....</b>	<b>451</b>
<b>Appendix H</b>	<b>Pharmacophore Input .....</b>	<b>453</b>
<b>Appendix I</b>	<b>Input of Experimental Spectra .....</b>	<b>455</b>
<b>Appendix J</b>	<b>T1 Thermochemical Recipe .....</b>	<b>457</b>
<b>Appendix K</b>	<b><sup>13</sup>C Chemical Shift Corrections.....</b>	<b>459</b>
<b>Appendix L</b>	<b>Installing a Network HASP.....</b>	<b>463</b>
<b>Appendix M</b>	<b>NMR Chemical Shift Standards.....</b>	<b>465</b>



**Index..... 467**

# Section I

## Introduction

---

Molecular mechanics calculations and quantum chemical calculations play an ever-increasing role in modern chemistry. Primary functions are to supply information about structures, relative stabilities and other properties of isolated molecules. Because of their inherent simplicity, molecular mechanics calculations on complex molecules are widespread throughout the chemical community. Quantum chemical calculations, including Hartree-Fock molecular orbital calculations, but especially calculations that take account of electron correlation, are much more time demanding. Only recently, have fast enough computers become widely available to make their application routine among mainstream chemists.

Quantum chemical calculations may also be called upon to furnish information about the mechanisms and product distributions of chemical reactions, either directly by calculations on transition states, or indirectly based on the Hammond Postulate, by modeling the steric and electronic demands of the reactants. Quantitative calculations, leading directly to information about the geometries of transition states, and about reaction mechanisms in general, are becoming more and more common, while qualitative models are still needed for systems that are too large or too complicated to be subjected to the more rigorous treatments. Finally, quantum chemical calculations can supply information to complement existing experimental data or to replace it altogether, for example, atomic charges for QSAR analyses, and intermolecular potentials for molecular mechanics and molecular dynamics calculations.

*Spartan'10* (“*Spartan*”) has been designed to address the ever increasing role that calculations play in chemistry and related fields. It represents a continuing collaboration between Wavefunction, Inc., and Q-Chem, Inc. Q-Chem codes supplement and extend the traditional strengths of *Spartan* as an easy to learn and use tool for

molecular mechanics, semi-empirical and Hartree-Fock molecular orbital calculations, as well as a wide variety of graphical models, with a full range of density functional models and a broad selection of wavefunction-based post-Hartree-Fock models. All models have been implemented using what we believe are the most robust algorithms currently available, and have been tuned for high performance on Intel and AMD processors including multi-core processors.

*Spartan is intended to be utilized by chemists*, not only computational chemists who are already familiar with the capabilities of molecular mechanics and quantum chemical methods, but also experimental chemists who may have little or no prior experience, but who want to use calculations much in the same way as experimental techniques such as NMR spectroscopy. This ambitious goal is directly reflected in the program's overall design criteria: "convenient access both to a full range of modern molecular mechanics and quantum chemical models and tools for interpretation of the results of calculations", and clearly distinguishes *Spartan* from other molecular modeling packages.

*Spartan* is supported under Windows (XP, XP 64, Vista and 7, 32 and 64 bit), Macintosh (OS X 10.5 and 10.6), Linux (Linux Kernel 2.6 or later). The Linux version (only) may be used in a server mode with Windows, Macintosh or Linux acting as a front end, and can also perform calculations locally. The appearance and functionality of the different versions of *Spartan* are nearly identical.

New features in *Spartan'10* include the following:

### **Infrastructure**

**Tabs.** Each open molecule is now given a tab, any or all of which may either be turned "on" or "off". All tabs "on" mimics *Spartan's* previous display paradigm (all open molecules are displayed). Turning tabs "off" simplifies screen clutter.

**64-Bit Implementation (Windows only).** The graphical user interface now supports 64 bit chip architectures (64 bit versions of back-end codes were released with *Spartan'08*). This allows *Spartan* documents in excess of 2 GB, meaning that very large molecule lists

may now be handled. Choice between 32 and 64 bit architectures is automatic and made at the time of installation.

**File Compression.** Files generated from *Spartan* (.spartan files) have been significantly reduced in size due to elimination of the part of the wavefunction that is not used, to file compression and reduction of precision of data used for construction of graphical displays. This reduction ranges from a factor of five to between one and two orders of magnitude.

**Monitor.** The *Spartan* job monitor has been extended to allow real-time visualization of molecular structure. This means that progress of a geometry/transition-state geometry optimization can be followed.

**Color Palette.** A new color palette has been provided.

**Customizable Icons.** *Spartan's* icons may now be individually selected for display.

## Graphical User Interface

**Simplified Graphics Requests.** The most commonly requested graphics (HOMO, LUMO, electron density, and electrostatic potential, local ionization potential and LUMO maps) may now be requested by the above names instead of as combinations of surfaces and properties.

**Multi-Dimensional Energy Profiles.** Multi-dimensional energy profiles may now be constructed with arbitrary limits on the number of steps in each dimension. Plots from two-dimensional profiles may be made.

**Ramachandran Plot.** A Ramachandran plot may be made for a protein structure retrieved from PDB.

**Formula Editor.** A dialog has been provided for construction of quantities (properties or combination of properties) of individual molecules or of molecules related by chemical reaction. These formulas are used in spreadsheet calculations or as queries for database searches (mining) and subsequent statistical analysis and plotting.

**QSAR Descriptors.** Descriptors for QSAR type analyses have been collected in one place, and may now be conveniently transferred for data analysis. A series of new descriptors based on the electron density surface and the electrostatic potential and local ionization potential maps have been added.

## Computation

**RI-CIS(D) Gradients.** Gradients from RI-CIS(D) models are available allowing efficient calculation of the geometries of excited states.

**Multi-Core Parallel.** Vibrational frequency calculation has been parallelized for Hartree-Fock and density functional models, resulting in a speed improvement of approximately 1.8 for a two-core machine and 3.2 for a four-core machine. RI-MP2 has been parallelized. The underlying Hartree-Fock part (which completely dominates the calculation for less than 500 basis functions) gives speed improvements of approximately 1.9 for a two-core machine and 3.5 for a four-core machine. Speed improvements for the MP2 part are not as large. An important practical consequence is that the speed of T1 thermochemical calculations now improve significantly on multi-core machines.

**PM6.** The PM6 semi-empirical model has been implemented. This offers some improvement over PM3. PM6 becomes the second semi-empirical method available in *Spartan* that may be applied to transition metals.

**T1 Extensions.** The T1 thermochemical recipe has been extended to include molecules with silicon and phosphorus. Extension to boron is underway. A procedure for reliably estimating the best (lowest) T1 energy for a flexible molecule has been specified and implemented inside the graphical user interface.

## Graphical Models

**Hydride Potential.** The potential felt by hydride anion may be calculated and either displayed as a surface or as a property mapped

on an electron density surface. The latter is intended as an indicator of attraction of a molecule for a nucleophile.

**Orbital Energy Diagram.** An orbital energy diagram (spanning up to 20 occupied valence molecular orbitals and 2 unoccupied molecular orbitals) may now be displayed for the selected molecule. *Clicking* on a line (energy level) in this diagram generates on-the-fly and displays the orbital (assuming that the wavefunction is available). Orbitals generated in this way are not kept when the molecule is closed.

**Slicing Planes.** Graphical surfaces may now be sliced allowing the inside or graphical surfaces (in particular, property maps) to be seen.

**Areas of Banded Property Maps.** This allows the area of a band on a property map to be determined. It is useful in quantifying the results of graphical models, for example, which face of a molecule is more likely to undergo electrophilic or nucleophilic attack.

## Spectroscopy

**NMR.** The NMR spectra dialog has been completely reorganized allowing (independent or simultaneous) display of six different spectra:  $^1\text{H}$ ,  $^{13}\text{C}$ , DEPT, COSY, HMBC and HSQC. Empirically-estimated H-H coupling constants may or may not be included in  $^1\text{H}$ , COSY, HMBC and HSQC plots and all plots may be averaged using Boltzmann weights (if available). Experimental  $^1\text{H}$ , and  $^{13}\text{C}$  and DEPT spectra (available from a public on-line database or accessed from a local file) may be displayed on top of calculated spectra.

A new scheme for correcting  $^{13}\text{C}$  chemical shifts has been implemented for spectra obtained from the EDF2/6-31G\* model (only). The supplements the scheme provided in *Spartan'08* (based on topology and applicable only to spectra obtained from the B3LYP/6-31G\* model). The new scheme (based on topology and on bond orders to atoms connected to carbon) provides shifts with an (RMS) error of approximately 1.7 ppm (down from approximately 2.5 ppm for the first-generation correction scheme and approximately 6 ppm for uncorrected shifts).

## Infrared

**Raman.** Raman spectra may be calculated from Hartree-Fock and density functional models. Calculated Raman spectra may be fit to experimental spectra the using an overall frequency scale parameter and a width parameter (analogous to what is done for infrared spectra).

## Databases

**Spartan Molecular Database (SMD).** The full Spartan Molecular Database (provided with maintenance) continues to expand. Approximately 3,000 transition-metal inorganic and organometallic compounds have been added to SMD obtained from calculations using the B3LYP/6-31G\* model (or the analogous B3LYP/LACVP\* model for second and third-row metals). This model has been shown to provide a good account of equilibrium geometries for organometallic compounds (based primarily on comparisons with X-ray crystal structures). Experimental thermochemical data for organometallic compounds are almost non-existent, and detailed assessment of the B3LYP for reaction energies is not possible. However, the B3LYP/6-31G\* model appears to properly account for observed products in cases where multiple products are possible. All entries correspond to optimized geometries, and infrared spectra have been provided for most entries. (A sample of approximately 5,000 molecules is provided with all copies of *Spartan*.)

**Spartan Spectra and Properties Database (SSPD).** A new database of IR and NMR spectra and of molecular and atomic properties and QSAR descriptors has been developed and is available as an option to *Spartan'10*. (A sample of approximately 5,000 molecules is provided with all copies of *Spartan*.) At first release, it contains ~60,000 molecules and is expected to double in size by the end of 2011. SSPD entries are based on equilibrium geometries from the EDF2/6-31G\* density functional model using the best conformation obtained from the T1 thermochemical recipe. In addition to spectra, properties and QSAR descriptors, SSPD entries include the T1 heat of formation, the vibrational modes (allowing animation of the motions associated with individual frequencies) and the wavefunction

(allowing on-the-fly generation and display of graphical surfaces) including the molecular orbitals and property maps such as the electrostatic potential map.

**Data Mining.** Data mining allows access to all information relating to molecular structure, properties and spectra contained in SMD and SSPD (presently >600K entries in total). The primary objective is to allow relationships among different quantities or combinations of quantities to be examined across entire series of molecules. Information retrieved from mining may be presented in the form of histograms or XY plots and may be used as the basis for linear regression. Data mining opens up a wealth of information to exploration. Data mining is accomplished in a series of steps. A **structural query** is devised and all entries containing this query are retrieved. For example, were the structural query the molecule phenol with an open valence at the *para* position, all *para*-substituted phenols available in SMD and/or SSPD would be retrieved. Next, one or more **property queries** are formulated. These might simply refer to individual quantities such as the dipole moment or the charge on a particular atom. More generally, they may reference a combination of properties, for example, the charge on one atom times the charge on another atom. By associating reaction arrows with a structural query, property queries can be extended to refer to one of the products of a chemical reaction (fully defined given the structural query and the reaction arrows) or to a property of the reaction itself, most conspicuously, the energy of the reaction. Finally, one or more properties queries and brought together with the retrieved data from SMD and/or SSPD either for the purpose of both the results of the search and are combined.

Results from mining SMD or SSPD may be displayed in terms of histograms or scatter plots. In addition, multiple regression analysis may be carried out to uncover relationships among properties.

**Spartan Infrared Database (SIRD).** The interface to the Spartan Infrared Database has been completely reworked. Matching to an unknown (experimental) spectrum may now be accomplished with or without functional group, formula and substructure filters.



## Documentation

**Tutorials, Topics and Problems.** The full set of *Spartan* tutorials and a selection of discussion topics and problems may be opened inside of *Spartan*. Tutorials are indexed by chapter (in this manual).

**Wikipedia Access.** A Wikipedia page may be brought up from *Spartan*. Access is provided both from the **Activities** menu where the query is entered explicitly and from the **Database Preview** dialog where the query is the name of the molecule.

Changes from *Spartan'10* and *Spartan'08* include the following:

**Attachment Points/Reaction Arrows.** Specification of attachment points (used in database searches and reaction arrows (used both in database searches and transition-state guessing) have been separated from these functions.

**Database Dialogs.** A *Spartan* document (including a blank document in the case of a name search or spectra search) is required for entry into the databases. This is a consequence of being able to save search and data mining results.

**T1 Database.** The T1 database is no longer available. Rather, T1 heats of formation are included with all entries in the SMD and SSPD databases.

**On-Line PDB Access.** PDB access has been moved from the **Search** menu to the **File** menu.

Since the release of *Spartan'08*, *Spartan* users have diligently reported program bugs via [support@wavefun.com](mailto:support@wavefun.com). All reported bugs have been addressed. Wavefunction would like to thank our customers for their help in making *Spartan* a better program. User-feedback is greatly appreciated!

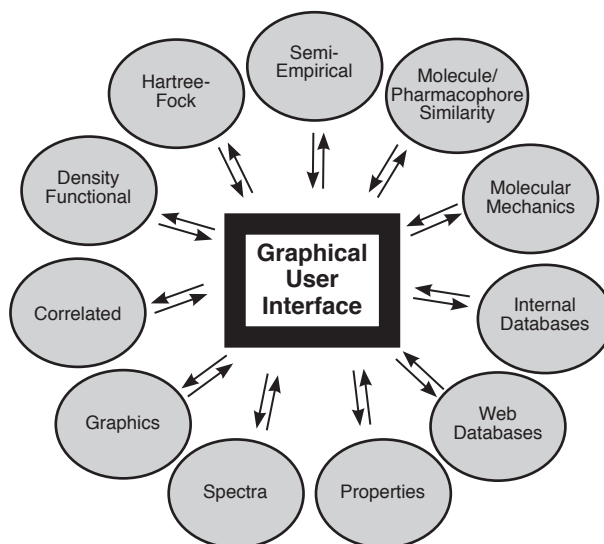
# Chapter 1

## *Spartan'10*

*This chapter describes the architecture of **Spartan'10**, focusing on the connectivity of graphics, database and computational components. Available molecular mechanics and quantum chemical methods are enumerated, and their overall role suggested.*

**Spartan'10** (“**Spartan**”) is a cross-platform application. With very few exceptions, **Spartan's** operations and features are identical between Windows, Mac and Linux versions. The following nomenclature will be used throughout this manual to indicate the few cases where features are not available: **NAL**; Not Available Linux, **NAM**; Not Available Macintosh and **NAW**; Not Available Windows.

**Spartan** comprises a series of independent modules tightly connected via a graphical user interface that is highly functional, yet elementary and uncluttered. It has been designed not only to greatly reduce the drudgery and possibility for human error associated with the preparation of input, but also to guide the interpretation of output. The interface is perhaps best viewed as an interactive and intuitive window into a full range of modern computational techniques.



Included in the interface are builders for organic, inorganic and organometallic molecules, polypeptides and polynucleotides, as well as for creating substituted molecules, and a procedure for guessing transition states. Access to ChemDraw<sup>1</sup> is provided without having to exit *Spartan's* interface. The interface provides access by substructure to the Spartan Molecular Database (SMD)<sup>2</sup>, a collection of over 150,000 calculated structures and associated properties of molecules, each obtained from as many as nine different theoretical models (Hartree-Fock models with 3-21G, 6-31G\* and 6-311+G\*\* basis sets, B3LYP and MP2 models with the 6-31G\* and 6-311+G\*\* basis sets, the EDF1 model with the 6-31G\* basis set and the G3(MP2) model. New in *Spartan'10* is a collection of approximately 3,000 transition-metal inorganic and organometallic molecules obtained from the B3LYP/6-31G\* model. The Spartan Reaction Database (SRD) has been significantly enhanced and now contains transition states for approximately 1,500 organic and organometallic reactions each obtained at one or more theoretical models. It may be accessed either by substructure searching (as SMD) or automatically, based on the structure of the reactant and a set of “curly arrows” defining the reaction.

The Spartan Spectra and Properties Database (SSPD)<sup>2</sup> appears for the first time with *Spartan'10*. It comprises IR and NMR spectra and diverse atomic and molecular properties together with QSAR descriptors for ~60,000 molecules (set to double by the end of 2011). Database entries include the vibrational modes, allowing animation of atomic motions associated with lines in the IR spectrum, and the wavefunction, allowing on-the-fly generation and display of graphical models. Entries are based on equilibrium structures obtained from the EDF2/6-31G\* density functional model and (where T1 has been defined) assuming the lowest-energy conformer assigned by the T1 thermochemical recipe.

The Spartan Infrared Database (SIRD)<sup>2</sup> comprises almost all the molecules with molecular weights < 400 amu in SSPD. It differs from SSPD in that it is searched by matching to an unknown (experimental) infrared spectrum.

New to *Spartan*'10 is capability to “mine” SMD and SSPD, not only for information about individual molecules, but also about relationships among molecules in balanced chemical reactions. Basic statistical analysis and plotting capabilities have been provided, including a new capability for singular value decomposition (finding the most important of a starting set of variables).

Also provided is an interface to the Cambridge Structural Database, comprising nearly 500,000 experimental X-ray crystal structures<sup>3</sup>, and automatic retrieval from the PDB database<sup>4</sup> of nearly 70,000 protein and related crystal structures. Finally, *Spartan* allows for import of a variety of data files, most important among them being PDB (protein structures) and SDF (structure drawings). Thus, it is not only possible to construct molecules from scratch, but also to access existing structures.

*Spartan*'s interface provides the gateway to a range of modern computational methods, including molecular mechanics models, semi-empirical and Hartree-Fock molecular orbital models, and a variety of correlated models including density functional and Møller-Plesset models. None of these models is likely to be ideal for every application<sup>5</sup>. While the most sophisticated quantum chemical models may yield excellent results, they will likely be too time consuming for routine application, and it will usually be necessary to contend with lesser treatments. *Spartan*'s interface facilitates mixing and matching different molecular mechanics and quantum-chemical models. Results from one model may easily be passed on for further analysis with another (more rigorous) model.

The simplest computational methods in *Spartan* are molecular mechanics models using the SYBYL and MMFF force fields, with aqueous energy corrections available for the latter (MMFFaq). These are available to determine equilibrium geometries and conformations of molecules comprising upwards of several thousand atoms. These are the only computational techniques that are applicable to biopolymers.

In addition, databases of diverse conformers obtained from MMFF molecular mechanics and intended to be used for similarity analysis are available<sup>6</sup>.

Quantum chemical models are required to account for the geometries of transition states as well as for reaction and activation energies. The simplest of these, semi-empirical molecular orbital models, can be routinely applied to systems with 100-200 atoms. Supported in *Spartan* are the MNDO model (with extensions for second-row and heavier main-group elements), the AM1 model, the PM3 model (with parameters for transition metals) and the RM1 model (a reparameterization of AM1). The PM6 semi-empirical molecule has been made available with the introduction of *Spartan'10*. While semi-empirical models have proven to be quite successful for determining equilibrium geometries including the geometries of transition-metal compounds, they have not proven to be reliable for the calculation of the reaction energies.

Hartree-Fock molecular orbital models are a mainstay of quantum chemical techniques, in particular, for equilibrium and transition-state structure determination, and for thermochemical comparisons. The models available in *Spartan* may be routinely applied to molecules with up to one hundred atoms. Hartree-Fock models are, however, *not adequate* for thermochemical comparisons where bonds are broken or formed, making them suspect for the calculation of activation energies. Also, Hartree-Fock models provide a poor account of the geometries of molecules incorporating transition metals. So-called correlated models are required. Several classes of correlated models are available in *Spartan*: density functional models, Møller-Plesset models, coupled cluster models and quadratic configuration interaction (QCI) models. Density functional models may be routinely applied to molecules of the same size as with Hartree-Fock models. On the other hand, the simplest Møller-Plesset model (MP2) is limited to molecules comprising thirty or fewer atoms. The recently introduced RI-MP2 model is nearly an order of magnitude faster than MP2 for energy calculations and 2-3 times faster for equilibrium and transition-state geometry calculations, but provides nearly identical results.

Higher order (MP3 and MP4) Møller-Plesset, coupled-cluster and QCI models are available for energy calculations only. They are limited in practice to molecules with less than 10-20 heavy (non-hydrogen) atoms.

Hartree-Fock and correlated models may be used with a variety of all-electron Gaussian basis sets. Hartree-Fock models are often quite successful with split-valence basis sets, while correlated models require polarization basis sets, and in some cases basis sets incorporating diffuse functions. Dual-basis set calculations, which approximate changes in energy and wavefunction resulting from extension from a small basis set to a large basis set, are also available. They can lead to an order of magnitude performance enhancement over calculations that directly use large basis sets. Pseudopotentials are provided for use with both Hartree-Fock and correlated models for calculations on molecules incorporating elements heavier than Kr (for which all-electron bases sets are not available).

A series of recipes (combinations of different methods) are available to estimate heats of formation. G2, G3 and G3 (MP2) recipes, while applicable only to very small systems, have been shown to yield heats of formation in excellent accord with experiment, and are benchmarks against which other calculations may be judged. The T1 recipe closely reproduces heats of formation from the G3 (MP2) recipe but requires 2-3 orders of magnitude less computation time. T1 is used as a basis for choosing the lowest-energy conformer for almost all entries in the Spartan Spectra and Properties Database, and the T1 heat of formation is included as a property.

*Spartan* provides singles configuration interaction (CIS) techniques for energies, equilibrium and transition-state geometry and conformation calculations on molecules in electronically-excited states. In addition, energy calculations (only) on excited states with the CIS (D) and RI-CIS(D) models as well as a range of density functional models may be performed.

*Spartan* provides access to several common spectral quantities, in particular infrared spectra (molecular mechanics, semi-empirical, Hartree-Fock, density functional and MP2 models), Raman spectra (Hartree-Fock and density functional models), NMR spectra (Hartree-Fock and density functional models) and UV/visible spectra (CIS, CIS(D) and density functional models). These are available both as numerical data (vibrational frequencies, chemical shifts, etc.) as well

as spectral plots. Spectra may be directly input to *Spartan* and internet connections to public IR, NMR and UV/visible databases are available allowing direct comparison of calculated and experimental spectra. The infrared spectra fitting procedure first provided in *Spartan'08* has been significantly improved in *Spartan'10* and a new interface has been implemented. This accesses the Spartan Infrared Database, a collection of 60,000 calculated infrared spectra.

The empirical scheme introduced in *Spartan'08* to correct  $^{13}\text{C}$  chemical shifts obtained from the B3LYP/6-31G\* model has been supplemented by a new scheme. This applies only to the EDF2/6-31G\* density functional model and reduces the rms error in  $^{13}\text{C}$  shifts to 1.2-1.8 ppm, down from 2.0-2.6 ppm for the original correction scheme and from >6 ppm for uncorrected shifts. We believe that  $^{13}\text{C}$  NMR spectra obtained from the corrected EDF2/6-31G\* model are now of sufficient quality to be useful in assigning experimentally determined NMR spectra. A collection of ~60,000 calculated proton and  $^{13}\text{C}$  NMR spectra is available in the Spartan Spectra and Properties Database.

*Spartan* provides tools to quantify both the similarity of molecular structures and of chemical environments characteristic of these structures. Also available is the ability to identify molecules that fit into a specific chemical environment (a so-called pharmacophore). Finally, pharmacophores may be extracted from small molecules bound to proteins and nucleotides in the PDB.

*Spartan* provides a variety of graphical tools to assist in interpreting the results of calculations as well as similarity analyses. These include structure models, and also molecular orbitals, electron and spin densities, local ionization potentials and electrostatic potentials that can be displayed as surfaces, slices and property maps. New to *Spartan'10* is an orbital energy diagram with subsequent *clicking* on an individual energy level in the diagram leading to display of the associated molecular orbital (assuming that a wavefunction is available). Also new to *Spartan'10* is the hydride potential map, providing the energy of interaction of hydride anion as a function of its position on an electron density surface. *Spartan* provides the ability

to distinguish accessible and inaccessible regions on a density surface and on property maps based on this surface. New to *Spartan*'10 is the ability to measure the area of a banded region on a property map and slice into a graphical surface or property map allowing the inside to be seen.

Graphical models may be manipulated in real time, and as many as desired can be displayed simultaneously. Animation files may be made, and can be used to depict conformational changes or chemical reactions.

- 
1. ChemDraw is not included with *Spartan*, but may be obtained from CambridgeSoft ([www.cambridgesoft.com](http://www.cambridgesoft.com)). **NAL**, **NAM**. While all versions of *Spartan* are able to read ChemDraw files, seamless access is presently available only in the Windows version.
  2. Subsets of SMD, SIRD and SSPD are provided with *Spartan*. The full SMD database is provided under maintenance contract. The full SSPD and SIRD databases (SIRD uses the entries in SSPD) may be licensed separately from Wavefunction.
  3. The Cambridge Structural Database is not included with *Spartan*, but is available by subscription from the Cambridge Crystallographic Data Centre ([www.ccdc.cam.ac.uk](http://www.ccdc.cam.ac.uk)).
  4. PDB web reference: <http://www.rcsb.org>.
  5. Full discussion and assessment of the specific molecular mechanics and quantum chemical models available in *Spartan* is provided in: W.J. Hehre, *A Guide to Molecular Mechanics and Quantum Chemical Calculations*, Wavefunction, Irvine, 2003. This is available as a PDF under the **Help** menu. See also: W.J. Hehre, L. Radom, P.v.R. Schleyer and J.A. Pople, *Ab Initio Molecular Orbital Theory*, Wiley, New York, 1986; Y. Shao, L.F. Molnar, Y. Jung, J. Kussmann, C. Ochsenfeld, S.T. Brown, A.T.B. Gilbert, L.V. Slipchenko, S.V. Levchenko, D.P. O'Neill, R.A. DiStasio Jr., R.C. Lochan, T. Wang, G.J.O. Beran, N.A. Besley, J.M. Herbert, C.Y. Lin, T. Van Voorhis, S.H. Chien, A. Sodt, R.P. Steele, V.A. Rassolov, P.E. Maslen, P.P. Korambath, R.D. Adamson, B. Austin, J. Baker, E.F.C. Byrd, H. Dachsel, R.J. Doerksen, A. Dreuw, B.D., Dunietz, A.D. Dutoi, T.R. Furlani, S.R. Gwaltney, A. Heyden, S. Hirata, C-P. Hsu, G. Kedziora, R.Z. Khallilulin, P. Klunzinger, A.M. Lee, M.S. Lee, W.Z. Liang, I. Lotan, N. Nair, B. Peters, E.I. Proynov, P.A. Pieniazek, Y.M. Rhee, J. Ritchie, E. Rosta, C.D. Sherrill, A.C. Simmonett, J.E. Subotnik, H.L. Woodcock III, W. Zhang, A.T. Bell, A.K. Chakraborty, D.M. Chipman, F.J. Keil, A. Warshel, W.J. Hehre, H.F. Schaefer, J. Kong, A.I. Krylov, P.M.W. Gill and M. Head-Gordon, *Phys. Chem. Chem. Phys.*, **8**, 3172 (2006).
  6. Databases comprising sets of diverse conformers pharmaceutical candidates as well as Maybridge and Life Chemicals collections may be licensed from Wavefunction.





# Section II

## Getting Started

---

The chapters in this section provide a brief introduction to the graphical user interface for *Spartan*'10, and following this, a series of tutorials. The objective of the latter is to provide hands on experience, and in doing so, illustrate the way in which molecular mechanics, quantum chemical calculations and similarity analyses may be set up, performed and interpreted. This is the place to start for new users. Users of previous versions of *Spartan* are also advised to review this section as there are a variety of new capabilities.

No attempt has been made to illustrate the full range of *Spartan*'s capabilities. For this, the reader should consult **Section III** of this guide. Focus is on the use of *Spartan* to calculate equilibrium and transition-state geometries, to search conformation space, to align molecules and evaluate the similarities of molecules, to evaluate reaction thermochemistry and activation energetics and to obtain NMR, infrared and UV/visible spectra with emphasis on NMR. Likewise, no attempt has been made to assess the performance of the various molecular mechanics and quantum chemical models.

Graphical models are illustrated in several of the tutorials that follow, and connections between specific models and chemical observables are pointed out, for example, between electrostatic potential maps and electrophilic character.

Discussion of both the behavior of different molecular mechanics and quantum chemical models and of graphical modeling techniques is provided in *Topics* available under the **Activities** menu. A more complete account is provided in *A Guide to Molecular Mechanics and Quantum Chemical Calculations*, included as a PDF under the **Help** menu and available from Wavefunction as a hardbound volume.

*Operating Spartan (Chapter 2)* outlines the basic workings of *Spartan*'s graphical user interface. It should be reviewed before

starting the tutorials and referred to as needed. The tutorials that make up the bulk of this section are divided across nine chapters: *Basic Operations*, *Organic Molecules*, *Groups of Organic Molecules*, *Spectra of Organic Molecules*, *Organic Reactions*, *Medicinal Chemistry*, *Polypeptides to Proteins*, *Inorganic/Organometallic Molecules* and “*Dry Labs*”: *Using the Spartan Spectra and Properties Database*. The initial three cover a number of important general features and should be completed first, and in order. The next five chapters cover more specialized topics and may be completed in any order. The final chapter points to an entirely different mode of using *Spartan* by drawing on the large amount of information in the Spartan Spectra and Properties Database.

Estimates for required computer time (in minutes, assuming a single processor on an Intel I7/920 chip) are indicated at the top of the first page of each tutorial. Users with dual and quad-core processors and who access *Spartan*'s parallel Hartree-Fock, density functional and RI-MP2 codes can expect significant reductions in these times.

5 mins

Where two times are provided, the second corresponds to the total time including the optional parts of the tutorial.

# Chapter 2

## Operating *Spartan*

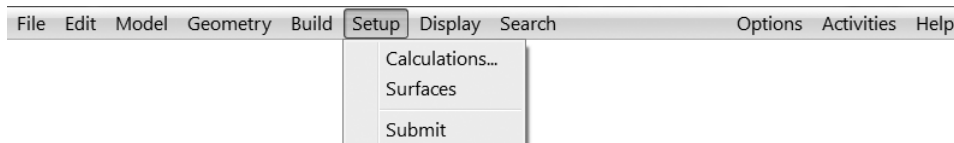
*This chapter describes the general operating features of Spartan'10. It should be read prior to starting the tutorials.*

### Opening and Exiting *Spartan*

To open under Windows, *click* on the **Start** button, then *click* on **All Programs**, and finally *click* on **Spartan'10** (or *double click* on the **Spartan** icon on your desktop). To open under Linux, bring up a terminal and *type* **spartan**. To open under Macintosh, *double click* on the **Spartan'10** icon in the Applications Folder. To exit, select **Exit** from the **File** menu (Select **Quit Spartan 10** from the **Spartan 10** menu on Mac), or *click* the **Close** button (⌵) at the top right (top left for Mac) of the **Spartan** interface.

### Menus

Program functions may be accessed using pull-down menus under the headings in the menu bar, for example, the **Setup** menu.



**File** Allows you to create a new molecule or read in a molecule that you have previously saved, or a file created by another program (for example, an SD file created by ChemDraw) to print text and graphics, embed external files (for example, Word and Excel files) into **Spartan** documents, to access PDB files online and to make QuickTime movies.

**Edit** Allows you to transfer information to and from the clipboard, to undo the last operation, to find text strings and molecule fragments, to center molecules on screen, and to clear the active molecule by deleting it.


<b>Model</b>	Allows you to control the style of your model, to display hydrogen bonds and chemical function descriptors, and to couple or decouple molecules in a multi-molecule file. Allows you to display of a Ramachandran plot for a protein structure brought in from PDB.
<b>Geometry</b>	Allows you to measure and constrain bond lengths, angles and dihedrals, define points and planes, setup frozen atoms, select conformational degrees of freedom, select centers for similarity analysis and align molecules.
<b>Build</b>	Allows you to build and edit molecules and to create lists of substituted molecules. Provides seamless access to ChemDraw™. (NAL, NAM)
<b>Setup</b>	Allows you to specify the task to be performed and the theoretical model to be employed for this task, to specify graphical surfaces and property maps and to submit jobs for calculation.
<b>Display</b>	Allows you to display text output, molecular and atomic properties, QSAR descriptors, orbital energy diagrams, surfaces and property maps and infrared, Raman, NMR and UV/visible spectra, as well as to access experimental IR, NMR and UV/visible spectra over the internet. Allows you to present data in a spreadsheet and make plots from and perform regression analysis on these data, to access the results of a similarity analysis and to compute reaction energies based either on user data or from entries in the Spartan Molecular Database (SMD) or Spartan Spectra and Properties Database (SSPD).
<b>Search</b>	Allows you to specify queries for the Spartan Molecular Database (SMD), the Spartan Spectra and Properties Database (SSPD) and the Spartan Reaction Database (SRD). Allows you to search and mine SMD and SSPD for calculated structures, properties and spectra, and to search SRD for calculated transition-state structures. Allows you to search the Cambridge Structural Database (CSD) of experimental X-ray crystal structures, and to extract ligands from the Protein Data Bank (PDB). Allows you to match “unknown” (experimental) infrared spectra to calculated spectra in the Spartan Infrared Database (SIRD). Allows you to guess a transition-state geometry

	based on a library of reactions and to identify tautomers.
<b>Options</b>	Allows you to set display standards, specify databases, assign compute servers, monitor executing jobs and specify queues and customize icons.
<b>Activities</b>	Allows you to display tutorials and topics inside of <i>Spartan</i> and open up a Wikipedia search page.
<b>Help</b>	Provides access to information on <i>Spartan's</i> general operation, the <i>Spartan'10 Tutorial and User's Guide</i> , <i>A Guide to Molecular Mechanics and Quantum Chemical Calculations</i> , and a number of computational FAQ's.

A complete listing of menu functions is provided in **Appendix B**.

## Keystroke Equivalents (NAM)

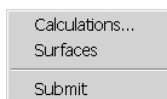
Keystroke equivalents for menus are designated by the underlined letter: **File**, **Edit**, **Model**, **Geometry**, **Build**, **Setup**, **Display**, **Search**, **Options**, **Activities** and **Help**. To access a menu, *press* the **Alt** key (**Option** key on Mac) and the appropriate letter key (case insensitive), for example, **Alt S** to access the **Setup** menu.



Each menu entry is designated by an underlined letter. To access the entry, *press* the appropriate letter key, for example, *press* the **C** key to access the **Calculations** dialog.

## Icons

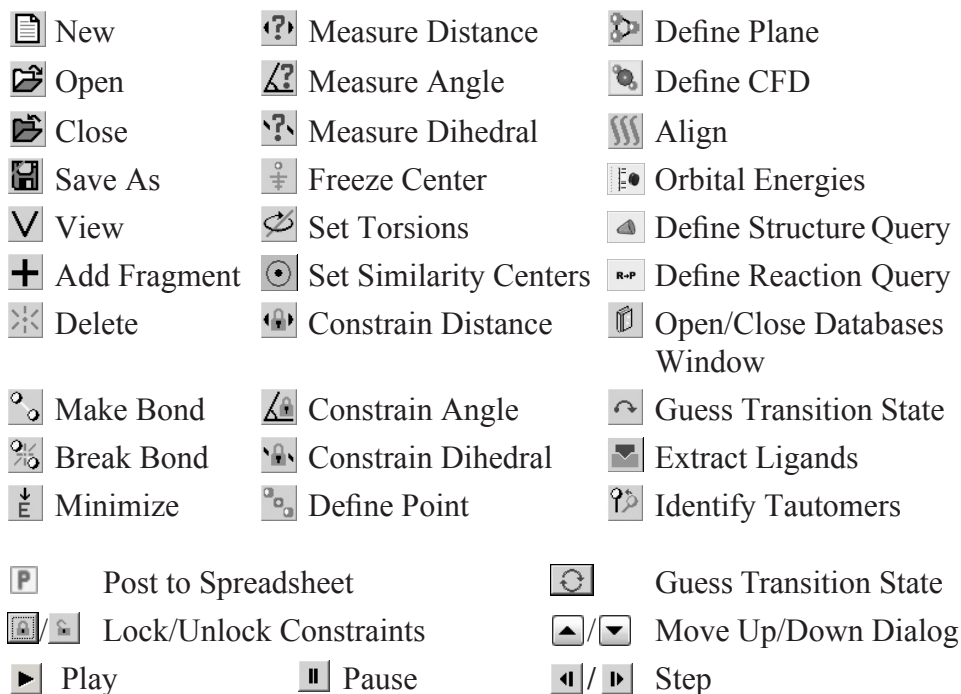Icons provide convenient access to selected functions under the **File**, **Build**, **Geometry** and **Search** menus, as well as other specialized functions. Icons may be selected or deselected for display from the **Toolbars Preferences** dialog (**Preferences** under the **Options** menu; **Chapter 20**).



## Tabs

*Spartan'10* assigns a tab to each open document. These appear in a single row along the bottom of the screen in the order that the documents were created or read. Backward and forward step keys ( and ) at the far right provide access to tabs outside of those displayed. Documents are selected for display by *clicking* on its tab. To the left of each tab is a check box, which if checked indicates that the document will be displayed on screen even if it is not the selected document. The behavior of previous versions of *Spartan* (all tabs checked) can be preserved by selecting **Pin New Documents** from the **Settings Preferences** dialog (**Options** menu; **Chapter 20**). Alternatively, one or only a few tabs can be checked leading to less clutter on the screen.

## Mouse and Keyboard Operations

The following functions are associated with the two-button mouse.

keyboard	button	
	left	right
–	selection, X/Y rotate, atom/fragment exchange <sup>a</sup> , fragment insertion <sup>a</sup>	X/Y translate
<b>Shift</b>	range selection, Z rotate	scaling <sup>b</sup>
<b>Ctrl</b>	global X/Y rotate <sup>c</sup>	global X/Y translate
<b>Ctrl + Shift</b>	multiple selection, global Z rotate <sup>c</sup>	scaling <sup>b</sup>
<b>Ctrl</b> (add fragment mode)	fragment X/Y rotate chirality invert <sup>a,d</sup>	fragment X/Y translate
<b>Ctrl + Shift</b> (add fragment mode)	fragment Z rotate absolute configuration invert <sup>a,d</sup>	scaling <sup>b</sup>
<b>Ctrl + Shift</b> <sup>d</sup> (Mac & Linux)	bond rotation	bond stretching
<b>Alt</b> (Windows), <b>option</b> (Mac)	group selection <sup>e</sup> , bond rotation <sup>f</sup>	bond stretching

a) Requires *double clicking*.  
 b) Scaling is always applied to all open molecules and all fragments. The center mouse wheel will also scale molecules.  
 c) Global rotations can be either molecule or screen centered. This is controlled by **Global Rotate** in the **Miscellaneous Preferences** dialog (**Preferences...** under **Options** menu).  
 d) For Macintosh, **command** key replaces **Ctrl** key for chirality inversion.  
 e) In **View** mode only.  
 f) In **Add Fragment** mode with bond selected (red arrow visible).

Mouse/keyboard operations may be broadly separated into two categories: selection (picking) and manipulation (translation/rotation).

**Selection.** The left button is used for selection of objects on screen and/or of menu items. Left and right buttons together are used to define a selection box for copying/cutting to the clipboard, as well as for multiple model selection. Together with the **Shift** key, the left button allows for selection over a range. Together with the **Ctrl** (**Control**) key, the left button allows for multiple selection. Both range and multiple selection apply not only to text items in lists, but



to atoms and bonds in molecules as well. Together with the **Alt** key (**option** key for Mac), the left button allows for selection of an entire group (detached molecular fragment). In add fragment mode, *double clicking* the left button leads to atom or group exchange. *Double clicking* on an atom with the **Ctrl** key depressed leads to inversion in chirality of the atom and *double clicking* on an atom with both **Ctrl** and **Shift** keys depressed inverts the absolute configuration of the molecule. Once an initial fragment, group or ring has been drawn, *double clicking* on the background will insert it alongside (but not bonded to) whatever fragments currently exist on screen.

**Manipulation.** The left button is used for rotation and the right button is used for translation and scaling of objects on screen. With no keys depressed, the left mouse button gives rise to rotation about the X and Y (screen) axes, while the right mouse button gives rise to translation in the X and Y (screen) directions. Together with the **Shift** key, the left mouse button gives rise to rotation about the Z direction and the right mouse button gives rise to scaling. The center (scroll) wheel on the mouse may also be used for scaling.

The **Ctrl** key in conjunction with the left or right mouse buttons and (optionally) the **Shift** key, signifies a change in focus away from the default for the purpose of rotations and translations. Outside of add fragment mode, the default is focus on a single molecule (the selected molecule). Use of the **Ctrl** key changes focus to the entire set of molecules on screen, meaning that rotations and translations are carried out globally. In add fragment mode, the default is focus on the full set of fragments that make up the molecule being constructed, and rotations and translations refer to this set of fragments as a whole. Use of the **Ctrl** key changes focus to a single fragment (the selected fragment), and rotations and translations now refer only to this fragment.

In add fragment mode, use of the **Alt** key (**option** key for Mac, **Ctrl** and **Shift** keys for Linux) with the left mouse button allows for rotation about a selected bond and, with the right mouse button, for changing the length of the selected bond.

Additional keys control various *Spartan* functions.






<b>3</b>	<b>3</b> selects red-cyan stereo display. Pressing again returns to non-stereo display.
<b>Page Up, Page Down, Home, End</b>	Moves up ( <b>Page Up</b> ), down ( <b>Page Down</b> ), to the top ( <b>Home</b> ) and to the bottom ( <b>End</b> ) of the set of open molecules. Also, moves up and down pages in the <b>Output</b> dialog.
<b>Insert</b> ( <b>option</b> for Mac)	In add fragment mode only, inserts a new fragment on screen. This is accomplished by selecting the fragment from the model kit, holding down the <b>Insert</b> key and <i>clicking</i> on screen. Insertion may also be accomplished by <i>double clicking</i> on the background following selection of a fragment.
<b>Delete</b>	Deletes a fragment, free valence, CFD, reaction arrow, the contents of a selection box, spectrum, curve, or plot. This is accomplished by holding down the <b>Delete</b> key and <i>clicking</i> on the fragment, etc.
<b>Enter</b> ( <b>return</b> for Mac)	Required following text or data entry into spreadsheet or dialogs.

### Selecting Molecules, etc.

Two or more molecules may be simultaneously displayed in *Spartan's* window (see **Tabs** earlier in this chapter). However, only one molecule may be selected. The selected molecule has access to all capabilities (molecule building, job setup and submission and text and graphical display and manipulation), while non-selected molecules may only be displayed as static images. The exceptions involve scaling and the use of the **Ctrl** key.

Selection of a molecule occurs by *clicking* on its structure model or on any of its associated graphical surfaces. This will also result in deselection of any previously selected molecule. Molecular properties for the selected molecule are available in the **Properties** dialog (**Display** menu; **Chapter 18**). Atom, bond, substituent and surface display properties, pharmacophores and group descriptors, as well as information about geometrical constraints may be

accessed by subsequent *clicking* on an atom, bond, etc., associated with the selected molecule, that are then highlighted (colored gold). *Clicking* on the selected atom, bond etc. resets the display to molecular properties. *Clicking* on another molecule results in display of molecular properties for that molecule. Information about plots (including spectra) and changes to plot style and range are also available from the **Plot Properties** dialog.

Where the molecule belongs to a document with more than a single member, selection from among the different members in the document may be made using either the  and  buttons or the scroll bar at the bottom left of the screen. Alternatively, if the spreadsheet for the document is open on screen (**Spreadsheet** under the **Display** menu; **Chapter 18**), selection can be made by *clicking* on the molecule label at the left of the spreadsheet. *Clicking* on  at the bottom left of the screen animates the display of molecules in the document, that is, steps through them sequentially. This is useful for displaying a progression of graphical surfaces along a reaction coordinate. *Clicking* on  (that replaces ) stops the animation.

Two or more molecules from the same document may be displayed at once (although only one may be selected). Molecules are marked for display by *checking* the box immediately to the left of the molecule label in the spreadsheet.

## Dialogs

Dialogs are either modal (as indicated by ... following their reference in a menu) or non-modal. Modal dialogs must be dismissed before program operation can continue, while non-modal dialogs may be kept open on screen. Only one modal dialog may be open on screen, while several different non-modal dialogs and/or several copies of the same non-modal dialog (each copy referring to a different molecule in a document or different documents) may be open on screen.

## Database Access

*Spartan* comes with sample copies of both the Spartan Molecular Database (SMD) and the Spartan Spectra and Properties Database (SSPD), each covering ~5,000 molecules. Full versions (~150,000 molecules for SMD and ~60,000 molecules for SSPD) are available from Wavefunction. SMD accesses up to nine different quantum chemical models for each molecule, whereas SSPD is limited to a single model (the EDF2/6-31G\* density functional model). SSPD is based on a more reliable model for selecting best conformer and includes infrared and NMR spectra as well as the wavefunction. Either (or both) SMD and SSPD may be searched by substructure (**Open Databases** under the **Search** menu; **Chapter 19**), but explicit molecular replacement requires that either SMD or SSPD be selected. This is done from a menu at the bottom of *Spartan*'s screen. Most of the tutorials that follow make use of SMD, and the menu should be set to SMD. The tutorials in **Chapter 11** make use of SSPD and here the menu should be set to SSPD.

## Stereo Displays

*Spartan* supports red-cyan stereo. Red/blue glasses must be worn. To enter stereo mode, *press* the **3** key. *Press* again to return to non-stereo mode.

## Changing Colors and Setting Preferences

**Colors and Preferences...** under the **Options** menu (**Chapter 20**) allows for changing default background and graphical object colors, and for setting (and resetting) program defaults, respectively. Further control of color is available in the various **Style** dialogs accessible from **Properties** under the **Display** menu; **Chapter 18**.

## Monitoring and Terminating Jobs

**Monitor** under the **Options** menu (**Chapter 20**) allows for monitoring of executing jobs as well as for terminating jobs.




# Chapter 3

## Basic Operations

---

*This tutorial introduces a number of basic operations in **Spartan** required for molecule manipulation, property query and spectra and graphics display. It is a must for new **Spartan** users, and should at least be skimmed by experienced users to expose them to new features of **Spartan**'10. Specifically it shows how to: i) open molecules, ii) view different models and manipulate molecules on screen, iii) measure bond distances, angles and dihedral angles, iv) display energies, dipole moments, atomic charges and infrared, NMR and UV/visible spectra, and v) display graphical surfaces and property maps. Spreadsheet operations are not illustrated, no molecules are built and no quantum chemical calculations are performed.*

1. Start **Spartan**. Click (left mouse button) on **File** from the menu bar that appears at the top of **Spartan**'s main window. Click on **Open...** from the **File** menu that appears. Alternatively, click on the  icon at the top of the screen. A file browser appears.

Move to the *tutorials* directory\*, click on *basic operations* and click on **Open** (or double click on *basic operations*). A single file containing *ethane*, *acetic acid dimer*, *propene*, *ammonia*, *hydrogen peroxide*, *acetic acid*, *water*, *cyclohexanone*, *camphor*, *3-aminobenzophenone*, *ethylene*, *benzene*, *aniline* and *cyclohexenone* will be opened. A ball-and-spoke model for the first molecule (*ethane*) will be displayed, and its name appears at the bottom right of the screen. The appearance of the name means that ethane is in the Spartan Molecular Database.\*\*

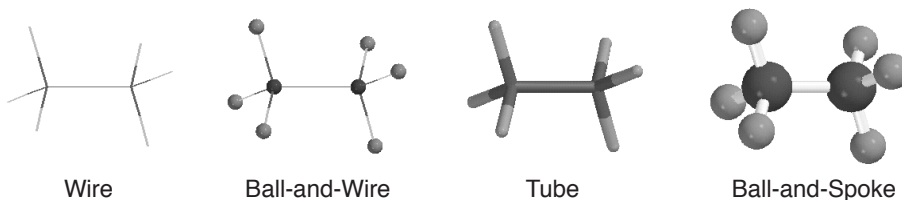
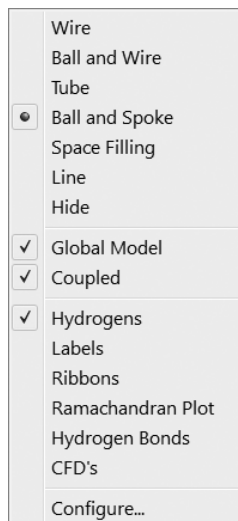
2. Practice rotating (*move* the mouse while holding down the left

---

\* For Windows, this is found in *Program Files/Wavefunction/Spartan10* or *Program Files(x86)/Wavefunction/Spartan10*. For Linux, this is found in the directory where **Spartan** was installed. For Macintosh, this is located in the Spartan.10.folder on the installation CDROM.

\*\* Make certain that **SMD** and not **SSPD** is selected from the menu at the bottom of the screen.


button) and translating (*move* the mouse while holding down the right button). *Click* on **Model** from the menu bar.



One after another, select **Wire**, **Ball and Wire**, **Tube** and finally **Ball and Spoke** from the **Model** menu. All four models for *ethane* show roughly the same information. The wire model looks the most like a conventional line formula. It uses color to distinguish different atoms, and one, two and three lines between atoms to indicate single, double and triple bonds, respectively.



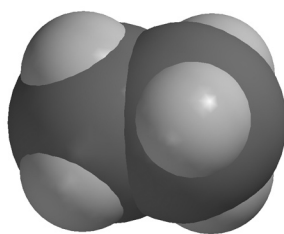
The ball-and-wire model is identical to the wire model, except that atom positions are represented by small spheres, making it easy to identify atom locations. The tube model is identical to the wire model, except that bonds are represented by solid cylinders. The tube model is better than the wire model in conveying three-dimensional shape. The ball-and-spoke model is a variation on the tube model; atom positions are represented by colored spheres, making it easy to see atom locations.


Select **Space Filling** from the **Model** menu.



Space-Filling

The space-filling model is different from the other models in that bonds are not shown. Rather, each atom is displayed as a colored sphere that represents its approximate size. Thus, the space-filling model for a molecule provides a measure of its size. While lines between atoms are not drawn, the existence (or absence) of bonds can be inferred from the extent to which spheres on neighboring atoms overlap. If two spheres substantially overlap, then the atoms are almost certainly bonded, and conversely, if two spheres barely overlap, then the atoms are not bonded. Intermediate overlaps suggest weak bonding, for example, hydrogen bonding.

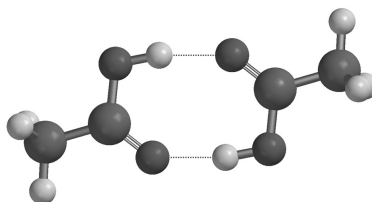3. *Click* once on the right arrow key  at the bottom left of the screen. This will move to the next molecule in the document, ***acetic acid dimer***. Its name will appear at the bottom of the screen. If you make a mistake, use the backward  or forward  step keys to get to ***acetic acid dimer*** in the document. Switch to a space-filling model and look for overlap between the (OH) hydrogen on one acetic acid molecule and the (carbonyl) oxygen on the other. Return to a ball-and-spoke model and select **Hydrogen Bonds** from the **Model** menu.




Ball-and-Spoke model for acetic acid dimer  
with hydrogen bonds displayed




The two hydrogen bonds, that are responsible for holding the acetic acid molecules together, will be drawn.




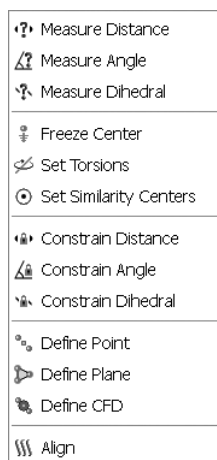
Use the **3** key to toggle between stereo 3D and regular display. To view in 3D you will need to wear red/blue glasses.\*

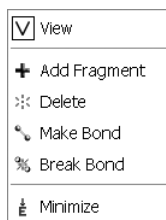4. Distances, angles, and dihedral angles can easily be measured with *Spartan* using **Measure Distance**, **Measure Angle**, and **Measure Dihedral**, respectively, from the **Geometry** menu.





- a) **Measure Distance:** This measures the distance between two atoms. *Click* once on  to move to the next molecule, *propene*, and then select **Measure Distance** from the **Geometry** menu (or *click* on the  icon at the top of the screen). *Click* on a bond or on two atoms (the atoms do not need to be bonded). The distance (in Ångstroms) will be displayed at the bottom of the screen. Repeat the process for several atoms. When you are finished, select **View** from the **Build** menu (or *click* on the  icon at the top of the screen).






- b) **Measure Angle:** This measures the angle around a central atom. *Click* once on  to move to the next molecule,

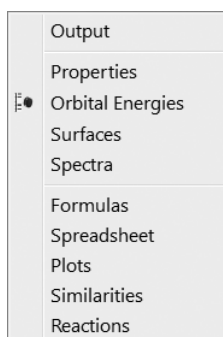
---


\* Active stereo will be implemented in *Spartan* as soon as a standard television protocol is established.


*ammonia*, and then select **Measure Angle** from the **Geometry** menu (or *click* on the  icon at the top of the screen). *Click* first on H, then on N, then on another H. Alternatively, *click* on two NH bonds. The HNH angle (in degrees) will be displayed at the bottom of the screen. *Click* on  when you are finished.

c) **Measure Dihedral**: This measures the angle formed by two intersecting planes, one containing the first three atoms selected and the other containing the last three atoms selected. *Click* once on  to move to the next molecule, *hydrogen peroxide*, then select **Measure Dihedral** from the **Geometry** menu (or *click* on the  icon at the top of the screen) and then *click* in turn on the four atoms (HOOH) that make up hydrogen peroxide. The HOOH dihedral angle will be displayed at the bottom of the screen. *Click* on  when you are finished.

5. Energies, dipole moments and atomic charges among other calculated properties, are available from **Properties** under the **Display** menu.



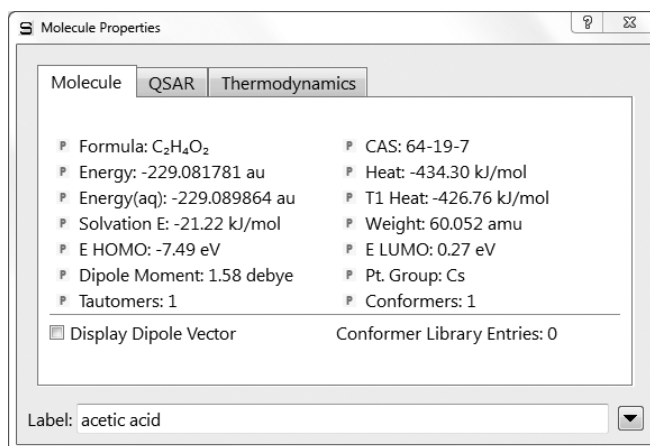a) **Energy**: *Click* once on  to move to the next molecule, *acetic acid*, and then select **Properties** from the **Display** menu. The **Molecule Properties** dialog appears. It is divided into two parts designated by tabs. **Molecule** provides a number of basic quantities while **QSAR** provides quantities that may be used as QSAR descriptors. Make certain that the **Molecule** tab is selected.





This provides the energy for acetic acid in atomic units (**Energy** in au). Also provided is an estimate of the energy in water (**Energy(aq)** in au) and the difference between the two (**Solvation E** kJ/mol).

- b) **Dipole Moment:** The magnitude of the dipole moment (**Dipole Moment** in debyes) is also provided in the **Molecule Properties** dialog. A large dipole moment indicates large separation of charge. You can attach the dipole moment vector,  $\text{+} \rightarrow$  where the + side refers to the positive end of the dipole, to the model on the screen, by *checking* the box to the left of **Display Dipole Vector** near the bottom of the dialog.

The vector will not be displayed if the dipole moment is zero. The dipole moment will not be reported if the molecule is charged because in this case it depends on the location and orientation of the molecule in space.

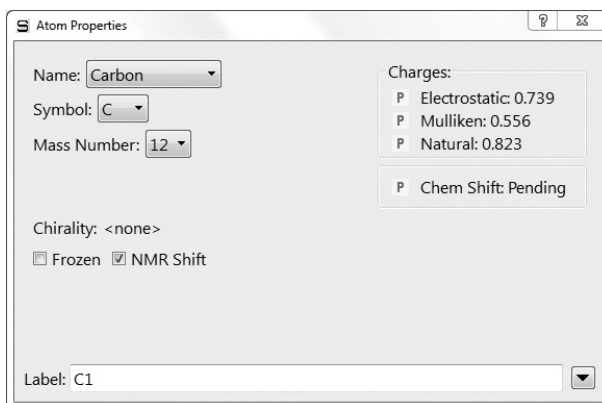
- c) **Atomic Charges:** To display the charge on an atom, *click* on it with the **Molecule Properties** dialog on the screen. The **Atom Properties** dialog replaces the **Molecule Properties** dialog.




Three different atomic charges, **Electrostatic**, **Mulliken** and **Natural**, are given in units of electrons. A positive charge indicates a deficiency of electrons on an atom and a negative charge, an excess of electrons. Repeat for other atoms. Confirm that the positively-charged atom(s) lie at the positive end of the dipole moment vector. When you are finished, close the dialog by *clicking* on  at the top of the dialog.


The three sets of atomic charges for acetic acid are different because the charge on an atom in a molecule cannot be uniquely defined, let alone measured. The problem is not the nuclear charge which is equal to the atomic number, but rather with saying how many electrons “belong” to a particular nucleus. The different calculated charges correspond to different ways of counting the number of electrons associated with an atom.


- d) **Infrared Spectra:** Molecules vibrate (stretch, bend, twist) even if they are cooled to absolute zero. This is the basis of infrared spectroscopy, where absorption of energy occurs when the frequency of a particular molecular motion matches the frequency of the light. Infrared spectroscopy is important for identifying molecules as different functional groups vibrate at noticeably different and characteristic frequencies.

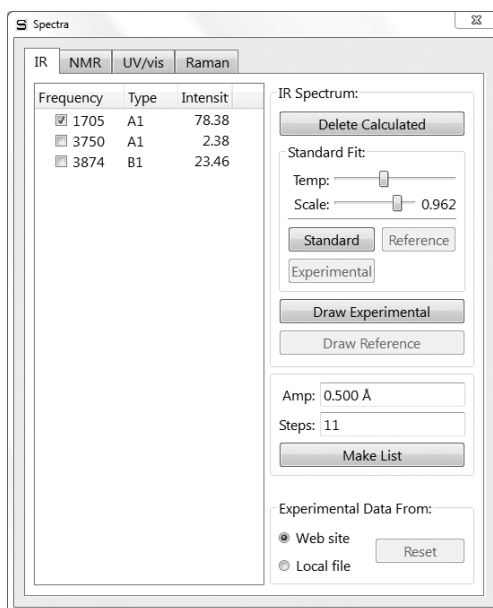*Click* once on  to move to the next molecule in the document, *water*. To animate a vibration, select **Spectra** from the **Display** menu and *click* on the **IR** tab. This leads to the **IR Spectra** dialog.



This displays the three vibrational frequencies for the water molecule, corresponding to bending and symmetric and antisymmetric stretching motions. One after the other, *click* in the box to the left of each frequency and examine the motion. Turn off the animation by *clicking* in the box of the selected frequency when you are finished.

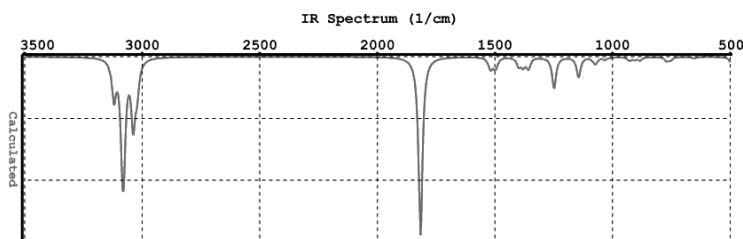
*Click* once on  to move to *cyclohexanone*, the next molecule in the list. The **Spectra** dialog now lists its 45 frequencies corresponding to the 45 vibrational motions of this molecule with 17 atoms. Examine each in turn (*click* on the entry in the dialog) until you locate the frequency corresponding to the CO (carbonyl) stretch. Next, *click* on **Draw Calculated** at the top of the dialog. The infrared spectrum of cyclohexanone appears.



You can move the spectrum around the screen by first *clicking* on it to select it (it will turn yellow) and then *moving* the mouse while holding down the right button. You can size it with the scroll wheel or by *moving* the mouse up and down while holding down both the **Shift** key and the right button.

Identify the line in the spectrum associated with the C=O stretch (a small red circle moves from line to line as you step through the frequencies in the **Spectra** dialog). Note that this line is isolated and that it is very intense, making it easy to find.

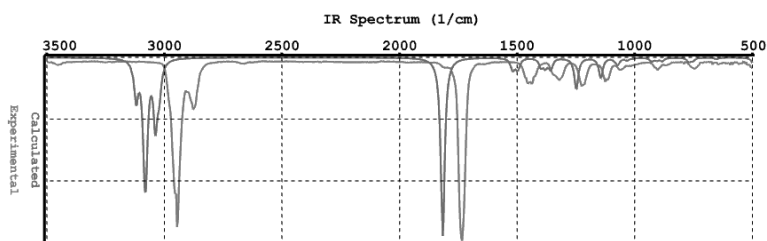
If your computer is connected to the internet, you can draw the experimental IR spectrum for cyclohexanone on top of the calculated spectrum. Select **Web Site** under **Experimental Data From:** at the bottom of the dialog and *click* on **Draw Experimental** in the middle of the dialog. If you are not online or have trouble downloading the experimental spectrum, select **Local File** under **Experimental Data From:**, find *cyclohexanoneIR.dx* in the *tutorials* directory\* and *click* on **Draw Experimental**.



Note that the lines in the calculated spectrum are consistently of higher frequency. To scale the calculated spectrum, bringing it into closer agreement with the experimental spectrum, *click* on it and move the slider bar to the right of **Scale** near the top of the dialog. Calculated and experimental spectra are automatically fit by selecting **Experimental** under **Calculated Fit** at the center of the **IR Spectra** dialog.


---

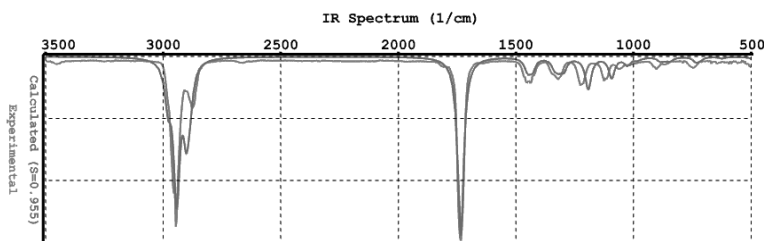
\* **All Files** needs to be selected in the **Open** dialog.




You can remove the plot by *clicking* on both **Delete Calculated** and **Delete Experimental** in the **Spectra** dialog. (These buttons have replaced **Draw Calculated** and **Draw Experimental**, respectively.)

- e) **NMR Spectra:** Along with mass spectrometry, NMR spectroscopy is the most powerful tool available with which to assign molecular structure. Many nuclei exhibit NMR spectra, but proton and  $^{13}\text{C}$  are by far the most important.

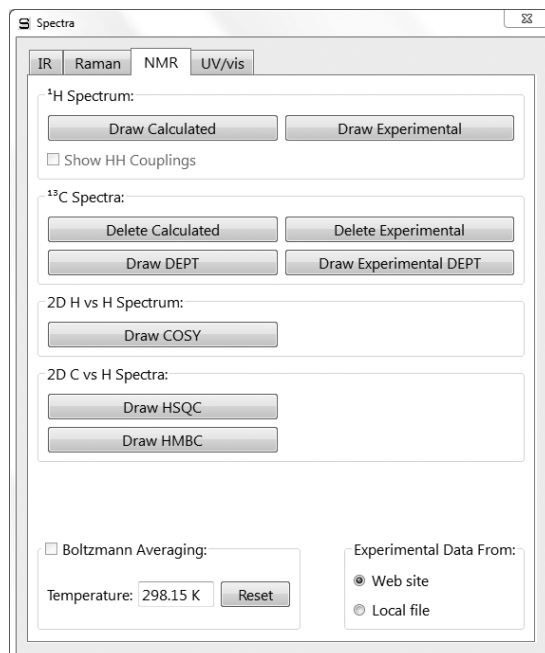*Click* once on  to move to the next molecule in the document, *camphor*. With the **Spectra** dialog on screen, *click* on the **NMR** tab to bring up the **NMR Spectra** dialog.




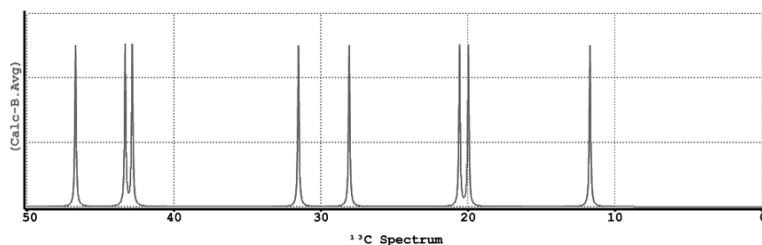Make certain tht **Boltzmann Averaging** is not checked and *click* on **Draw Calculated** under  $^{13}\text{C}$  **Spectrum** to show the calculated  $^{13}\text{C}$  spectrum.



This comprises nine lines, in the range of 50 to 0 ppm (there is a tenth line corresponding to the carbonyl carbon at 218 ppm). Note that the scale used here is different than the default range (150 to 0 ppm). Scale adjustment will be discussed in a later tutorial.

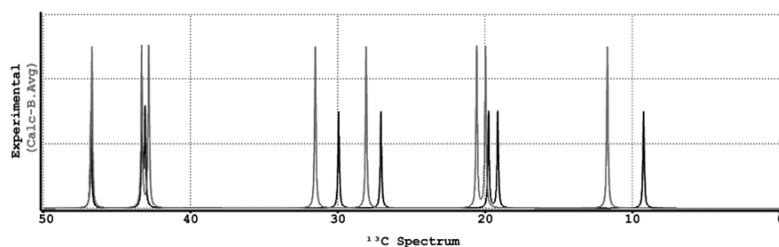If you are connected to the internet, select **Web Site** under **Experimental Data From:** at the bottom right of the dialog and *click* on **Draw Experimental** (under  $^{13}\text{C}$  Spectrum) to superimpose the experimental spectrum onto the calculated one. If you are not online or have trouble downloading the experimental spectrum, select **Local File** under **Experimental Data From:**, find *camphorNMR.dx* in the *tutorials* directory\* and *click* on **Draw Experimental**.



Remove the spectra by *clicking* on **Delete Calculated** and **Delete Experimental** in the NMR Spectra dialog. (These buttons have replaced **Draw Calculated** and **Draw Experimental**, respectively.)

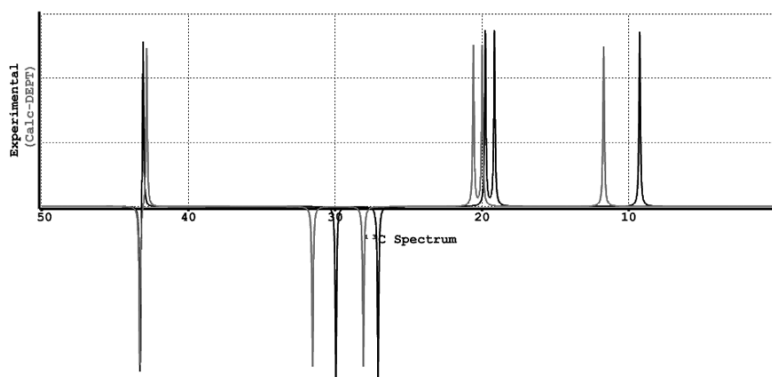
You can obtain numerical values for the individual chemical shifts by bringing up a properties dialog (**Properties** from the **Display** menu) and *clicking* on the atom of interest. Note that proton shifts in addition to  $^{13}\text{C}$  shifts are available.

---

\* **All Files** needs to be selected in the **Open** dialog.



Finally, draw both calculated and measured DEPT spectra. There is no experimental file, so before you *click* on **Draw Experimental DEPT**, select **Local File** under **Experimental Data From:**, find *camphorDEPT.dx* in the *tutorials* directory\* and *click* on **Draw Experimental**.




Note that both spectra now contain only seven lines. The lines corresponding to the three quaternary carbons (including the carbonyl carbon that is off scale) are missing. The three lines pointing down correspond to the methylene carbons while four lines pointing up correspond to the methine carbon and the three methyl groups.


- f) **UV/visible Spectra:** Absorption of light in the visible or ultraviolet range of the electromagnetic spectrum leads to electronic excitation from the ground-state to excited-states and (in the case of absorption in the visible), is responsible for a molecule's color. UV/visible spectroscopy not only offers a valuable fingerprint but is also an important screen to identify molecules that may be damaged by exposure to light. *Click* once on  to move to the next molecule, **3-aminobenzophenone**. Bring up the **Spectra** dialog and *click* on the **UV/vis** tab to give the **UV/vis Spectra** dialog.

---

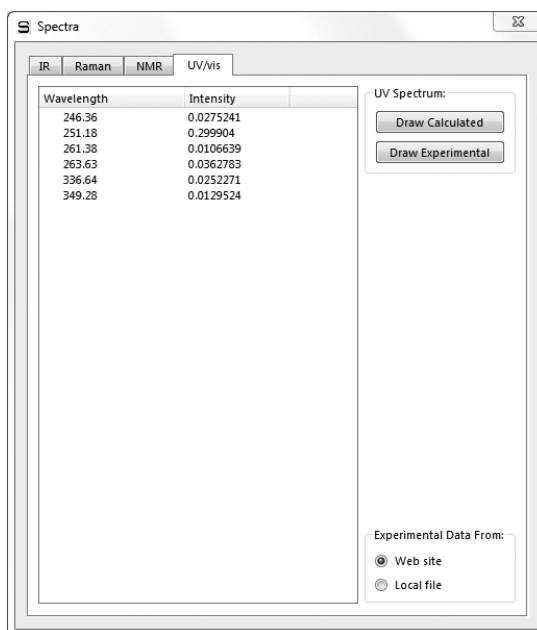\* **All Files** needs to be selected in the **Open** dialog.



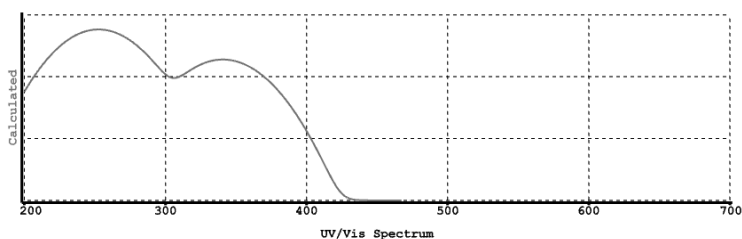
Click on **Draw Calculated** to show the calculated spectrum.



If you are on the internet, select **Web Site** under **Experimental Data From:** at the bottom of the dialog and *click* on **Draw Experimental**. If you are not online or have trouble downloading the experiment spectrum, select **Local File** under **Experimental Data From:**, find **3-aminobenzophenone.dx** in the *tutorials* directory\* and *click* on **Draw Experimental**. In either case, the experimental spectrum will be superimposed onto the calculated one.

---

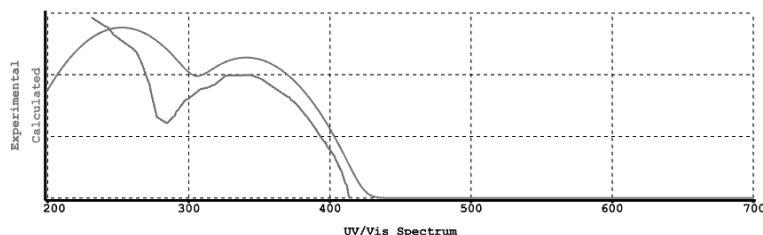
\* **All Files** needs to be selected in the **Open** dialog.




The two spectra are visually similar at least in a qualitative sense. Note, however, that any vibrational structure revealed in the experimental spectrum will not be mirrored in the calculated spectrum. Remove the spectra by *clicking* on **Delete Calculated** and **Delete Experimental** in the **UV/vis Spectra** dialog. (These buttons have replaced **Draw Calculated** and **Draw Experimental**, respectively.)

Even though the spectra look similar, they are sufficiently different that the theory would not likely be successful in accounting for the “color” of the molecule. Where it is successful is in anticipating changes in color resulting from changes in structure.


6. *Spartan* permits display, manipulation and query of a number of important graphical quantities resulting from quantum chemical calculations. Most important are the *electron density* (that reveals how much space a molecule actually takes up), the *bond density* (that reveals chemical bonds), and key *molecular orbitals* (that provide insight into both bonding and chemical reactivity). In addition, the *electrostatic potential map*, an overlaying of the electrostatic potential (the attraction or repulsion of a positive charge for a molecule) onto the electron density, is valuable for describing overall molecular charge distribution as well as anticipating sites of electrophilic addition. Another indicator of electrophilic addition is provided by the *local ionization potential map*, an overlaying of the energy of electron removal (ionization) onto the electron density. Finally, an indicator of nucleophilic addition is provided by the *|LUMO| map*, an overlaying of the absolute value of the lowest-unoccupied molecular orbital (the LUMO) onto the electron density.

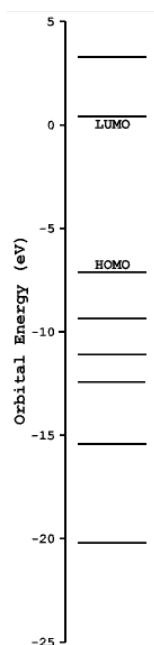Click once on  to move to the next molecule in the list, *ethylene*. Select **Orbital Energies** from the **Display** menu (or *click* on icon). An orbital energy diagram for ethylene will appear at the left of the screen. This will contain all six occupied valence molecular orbitals and two unoccupied molecular orbitals.





Click on the energy level in the diagram labelled *HOMO*. In a second, the familiar  $\pi$  bond in ethylene will appear. Note that the graphic has “blue and red” regions. These correspond to positive and negative values of the orbital (the absolute sign is arbitrary). Examine the other occupied orbitals (by *clicking* on their respective energy levels in the diagram) as well as the lowest-unoccupied molecular orbital (the *LUMO*).

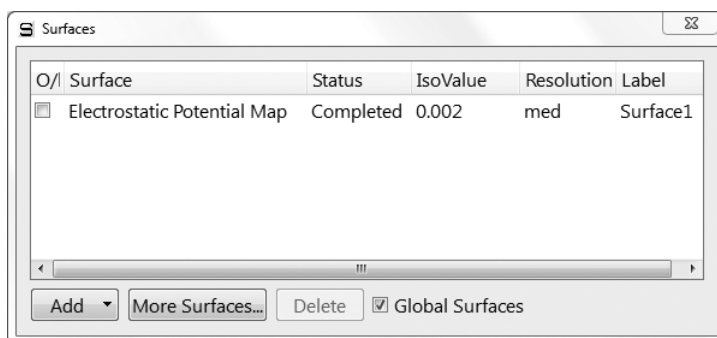
Click once on  to move to the next molecule in the list, *benzene*. Select **Surfaces** from the **Display** menu. The **Surfaces** dialog appears.




Select ***Electrostatic Potential Map*** inside the **Surfaces** dialog (*click* inside the box to the left of the name). An electrostatic potential map for benzene will appear. *Click* on the map. The **Style** menu will appear at the bottom right of the screen. Select **Transparent** from this menu. Making the map transparent allows you to see the molecular skeleton underneath. Go back to a **Solid** display (**Style** menu) in order to clearly see color differences. The surface is colored red in the  $\pi$  system (by convention, indicating negative potential and the fact that this region is attracted to a positive charge), and blue in the  $\sigma$  system (by convention, indicating positive potential and the fact that this region is repelled by a positive charge). Bring up the **Properties** dialog (**Display** menu) and *click* on the surface. Remove the checkmark from the box to the left of **Bands** in the **Surface Properties** dialog to replace the series of color bands (discrete display) by a continuous display.

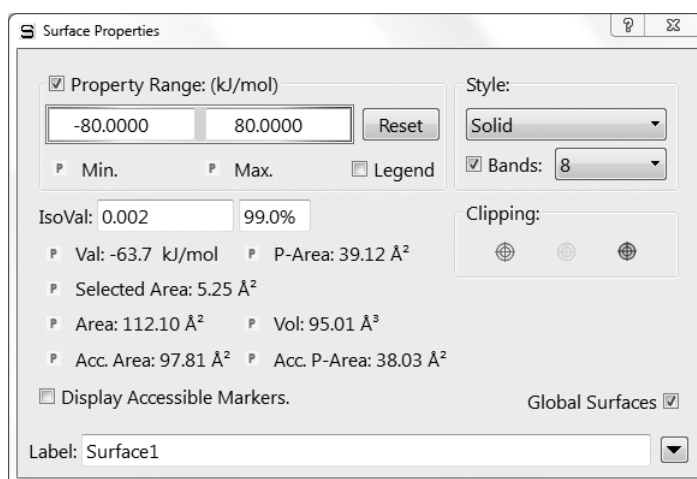


Click once on  to move to the next molecule in the list, *aniline*, and select **Local Ionization Potential Map** inside the **Surfaces** dialog. The graphic that appears, a local ionization potential map. By convention, red regions on the density surface indicate areas from which electron removal (ionization) is relatively easy, meaning that they are subject to electrophilic attack. These are easily distinguished from regions where ionization is relatively difficult (by convention, colored blue). Note that the *ortho* and *para* ring carbons are more red than the *meta* carbons, consistent with the known directing ability of the amino substituent.

Click once on  to move to the next molecule in the list, *cyclohexenone*, and select **LUMO** inside the **Surfaces** dialog. The resulting graphic portrays the lowest-energy empty molecular orbital (the LUMO) of cyclohexenone. This orbital is delocalized onto several atoms and it is difficult to tell where exactly a pair of electrons (a nucleophile) will attack the molecule.

A clearer portrayal is provided by a LUMO map, that displays the (absolute) value of the LUMO on the electron density surface. By convention, the color blue is used to represent maximum value of the LUMO and the color red, minimum value. First, remove the LUMO from your structure (select **LUMO** in the **Surfaces** dialog) and then turn on the LUMO map (select **|LUMO|Map** in the dialog). Note that there are two blue regions, one directly over the carbonyl carbon and the other over the  $\beta$  carbon. This is entirely consistent with known chemistry. Enones may either undergo carbonyl addition or conjugate (Michael) addition.

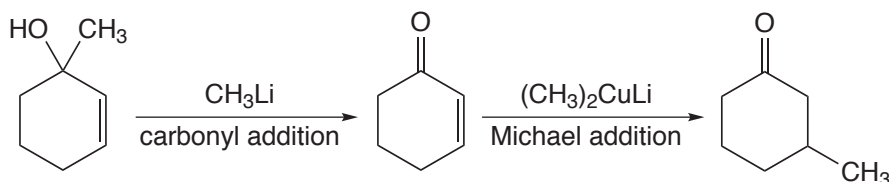

- When you are finished, close the document by selecting **Close** from the **File** menu or alternatively by *clicking* on the  icon at the top of the screen.




# Chapter 4
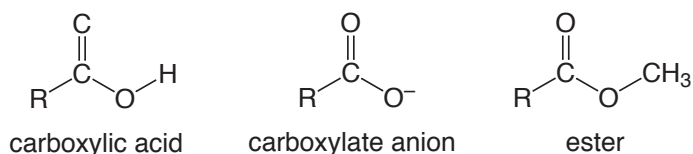
## Organic Molecules

---

*This chapter shows how to: i) construct organic molecules from atomic fragments, functional groups and rings, ii) carry out molecular mechanics and quantum chemical calculations and iii) calculate and interpret graphical models. It also illustrates use of the Spartan Molecular Database\* of calculated structures, energies and molecular properties as well as **Spartan's** interface to the Cambridge Structural Database. Discussion of the calculation and display of infrared and NMR spectra is deferred until **Chapter 6**.*

The simplest building blocks incorporated into **Spartan's** organic model kit are atomic fragments. These constitute specification of atom type, for example, carbon, and local environment, for example, tetrahedral. However, much of organic chemistry is organized around functional groups, collections of atoms the structure and properties of which are roughly the same in every molecule. The organic model kit also incorporates a small library of functional groups that can easily be extended or modified. For example, the carboxylic acid group may be modified to build a carboxylate anion (by deleting a free valence from oxygen), or an ester (by adding tetrahedral carbon to the free valence at oxygen).




---

\* **Spartan'10** includes subsets of approximately 5,000 molecules from the Spartan Molecular Database (SMD) which comprises more than 150,000 molecules along with their structures, energies, spectra and properties calculated with up to 9 quantum chemical models, and from the Spartan Spectra and Properties Database (SSPD), which comprises more than 60,000 molecules, includes the wavefunction allowing on-the-fly generation of graphical surfaces and property maps, and most of which provide both infrared and NMR spectra.



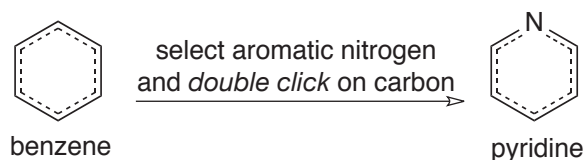
Rings are also common components in organic molecules, and the organic model kit incorporates a library of commonly-encountered hydrocarbon rings that can easily be modified by atom replacement. For example, pyridine may be built starting from benzene by selecting aromatic nitrogen from the list of atomic fragments, and then *double clicking* on one of the carbons. This results in the substitution of aromatic nitrogen for aromatic carbon.




Functional groups may also be modified in this manner.


In addition to illustrating the construction of organic molecules, the tutorials in this chapter also illustrate the way in which molecular mechanics and quantum chemical calculations are set up and their results examined and interpreted. A variety of graphical models are introduced and illustrated. Finally, examples of the use of the Spartan Molecular Database of calculated molecular structures, energies and properties and of the Cambridge Structural Database of experimental X-ray crystal structures are provided.


## Acrylonitrile

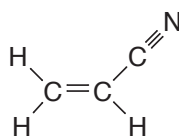

Acrylonitrile provides a good opportunity to illustrate the basics of molecule building in *Spartan*, as well as the steps involved in carrying out and analyzing a quantum chemical calculation.

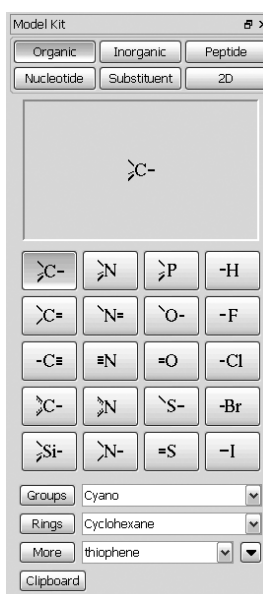1. Click with the left mouse button on **File** from the menu bar. Then click on **New** from the menu that appears (or click on the  icon at the top of the screen). The organic model kit appears.



At the center of the kit is a library of atomic fragments. Click on trigonal planar  $sp^2$  hybridized carbon  from the fragment library. A model of the fragment appears at the top of the model kit. Bring the cursor anywhere on screen and click (left button). Rotate the carbon fragment (*drag* the mouse while holding down the left button) so that you can clearly see both the double free valence (=) and the two single free valences (-).

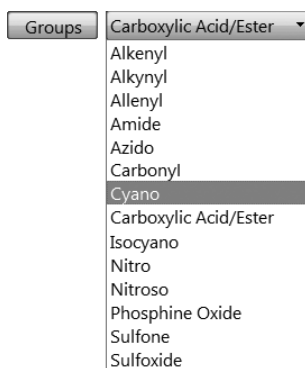
*Spartan's* model kits connect atomic fragments (as well as groups, rings and ligands) through free valences. Any free valences that

remain upon exiting a model kit are automatically converted to hydrogen atoms; it is not necessary to explicitly add hydrogens to open valences.

2.  $sp^2$  carbon is still selected. *Click* on the double free valence. The two fragments are connected by a double bond, leaving you with ethylene. The name “ethylene” will appear at the bottom right of the screen. If you make a mistake and *click* instead on the single free valence, select **Undo** from the **Edit** menu. You can also start over by selecting **Clear** from the **Edit** menu.



*Spartan's* organic model kit allow only the same type of free valences to be connected, for example, single to single, double to double, etc.



3. *Click* on **Groups** in the model kit, and select **Cyano** from the functional groups available from the menu.




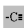*Click* on any of the four single free valences on ethylene (they are equivalent). This bonds the cyano group to ethylene, leaving you with acrylonitrile.\* Its name will now appear at the bottom right of the screen.

---

\* You could also have built acrylonitrile without using the **Groups** menu. Starting from scratch (**Clear** from the **Edit** menu), first build ethylene as above, then select  $sp$  hybridized carbon  from the model kit and then *click* on one of the free valences on ethylene. Next, select  $sp$  hybridized nitrogen  from the model kit and *click* on the triple free valence on the  $sp$  carbon. Alternatively, you could have built the molecule entirely from groups. Starting from scratch, *click* on **Groups**, select **Alkenyl** from the menu and *click* anywhere on screen. Then select **Cyano** from the menu of functional groups and *click* on one of the free valences on ethylene. In general, molecules can be constructed in more than one way.

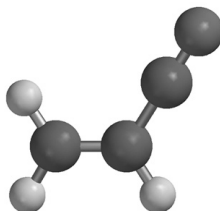
- Select **Minimize** from the **Build** menu (or *click* on the  icon at the top of the screen). The final molecular mechanics energy (36.2 kJ/mol) and symmetry point group ( $C_s$ ) are provided at the bottom right of the screen.
- Select **View** from the **Build** menu (or *click* on the  icon at the top of the screen). The model kit disappears, leaving only a ball-and-spoke model of acrylonitrile on screen.



ball-and-spoke model

This model can be rotated, translated and zoomed by using the mouse in conjunction with keyboard functions. To rotate the model, *drag* the mouse while holding down the left button; to rotate in the plane of the screen also hold down the **Shift** key. To translate the model, *drag* the mouse with the right button depressed. To zoom the model (translation perpendicular to the screen), use the center mouse wheel (scroll wheel) if available, or hold down the **Shift** key in addition to the right button while *dragging* the mouse up (zoom in) or down (zoom out).

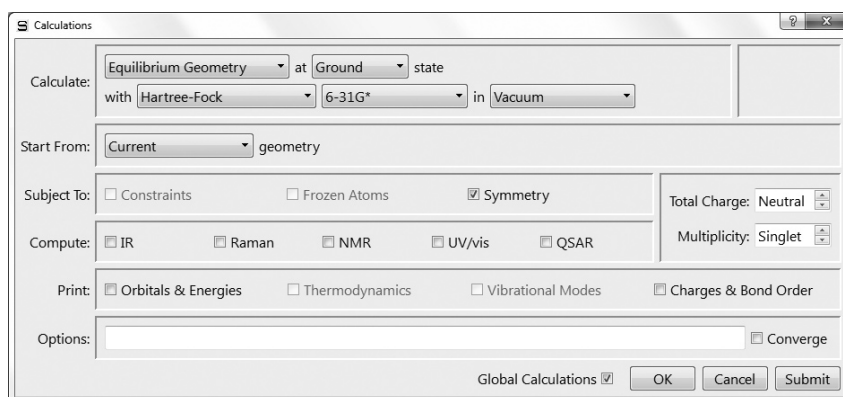
- Select **Calculations...** from the **Setup** menu, and perform the following operations in the **Calculations** dialog which appears.




- a. Select **Equilibrium Geometry** from the top menu to the right of **Calculate**. This specifies optimization of equilibrium geometry. **Ground** and not **First Excited** should appear in the menu to the left of **State**.
- b. Select **Hartree-Fock** and then **3-21G** from the two bottom menus to the right of **Calculate**. This specifies a Hartree-Fock calculation using the 3-21G split-valence basis set. Keep the default **in Vacuum** setting.

When you finish, *click* on **OK** to remove the dialog.

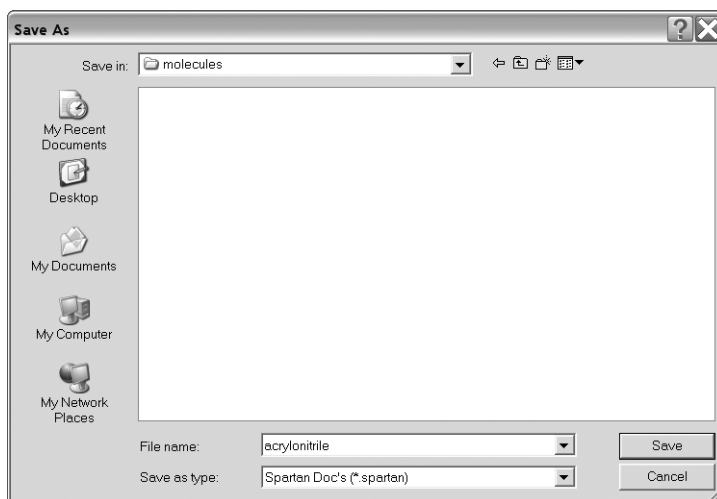7. Select **Submit** from the **Setup** menu.\* A file browser appears.



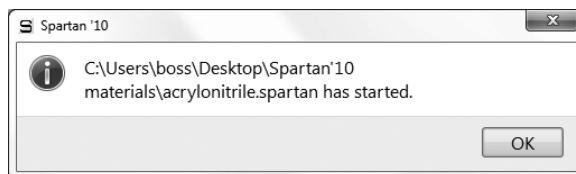The name *acrylonitrile* will be presented to you in the box to the right of **File name**. Either use it or type in whatever name you like and then *click* on **Save**.\*\* You will be notified that the calculation has been submitted.




- 
- \* You could also have *clicked* on **Submit** inside the **Calculations** dialog.
  - \*\* Proper names will automatically be provided for you to accept, modify or replace whenever the molecule exists in the Spartan Molecular Database or Spartan Spectra and Properties Database (whichever is selected) and where the document being submitted contains only one molecule. Otherwise the names “spartan1”, “spartan2” will be provided.


Click on **OK** to remove the message from the screen.

After a molecule has been submitted, and until the calculation has completed, you are not permitted to modify any dialogs or other information associated with it.

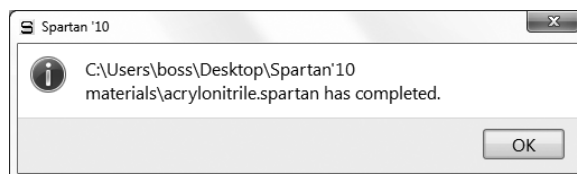8. You will be notified when the calculation has completed.



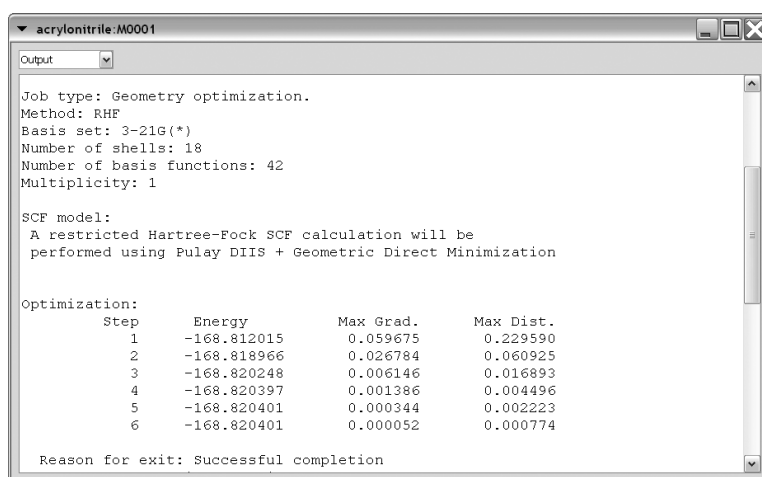
Click on **OK** to remove the message from the screen. Select **Output** from the **Display** menu. A window containing text output for the job appears.




You can scan the output from the Hartree-Fock calculation by using the scroll bar at the right of the window or by *clicking* (left button) on or inside the output window and using the scroll wheel on your mouse. The information at the top of the dialog includes the task, basis set, number of electrons, charge and multiplicity, as well as further details of the calculation. Below this is the symmetry point group of the molecule that was maintained during the optimization.

Eventually, a series of lines appear, under the heading “Optimization”. These tell the history of the optimization process. Each line (or “Step”) provides results for a particular geometry. Ideally, the energy will monotonically approach a minimum value for an optimized geometry. If the geometry was not optimized satisfactorily an error message, such as: “Optimization has exceeded N steps – Stop”, will be displayed following the last optimization cycle. If this were the case, you would have been notified that the job had failed, rather than seeing the “completed” message dialog.

Near the end of the output is the final energy\* (-168.82040 atomic units for acrylonitrile with the 3-21G basis set), and the computation time. *Click* on  at the top of the output dialog to close it.

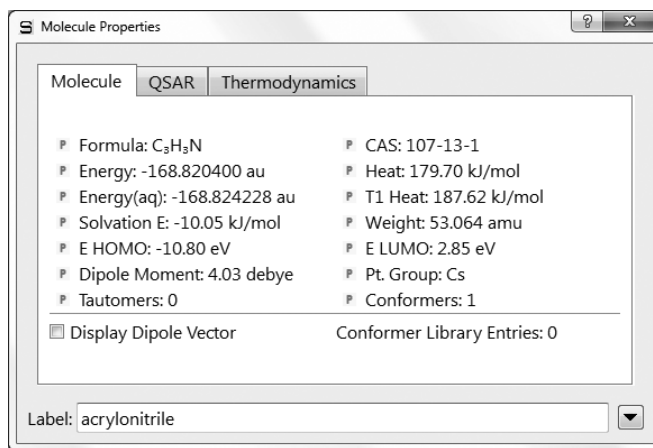You may examine the total energy and dipole moment among other calculated properties without having to go through the output. Select **Properties** from the **Display** menu to bring up the **Molecule Properties** dialog (make certain that the **Molecule** tab and not the **QSAR** tab is selected).




To see the dipole moment vector (indicating the sign and direction of the dipole moment), *check* the box to the right of **Display Dipole Vector**. Wire, ball-and-wire or tube models are best for this display.


\* See **Calculations...** (**Setup** menu; **Chapter 17**) for a discussion of how energy relates to the more familiar quantities, heat of formation and strain energy.



*Uncheck* the box to remove the dipole moment vector.

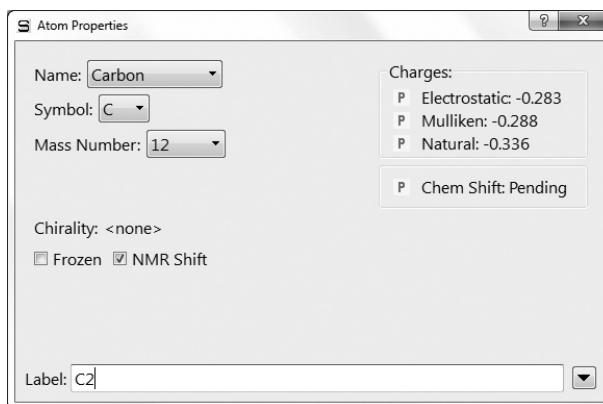*Click* on an atom. The (**Molecule Properties**) dialog will be replaced by the **Atom Properties** dialog.



Among other things, this provides three different sets of atomic charges: **Electrostatic**, **Mulliken** and **Natural**. To obtain the charge on another atom, simply *click* on it. Inspect all the atomic charges on acrylonitrile (by *clicking* on the appropriate atoms). When you are finished, *click* on  at the top of the **Atom Properties** dialog to close it.


9. Select **Surfaces** from either the **Setup** or **Display** menu. *Click* on **Add** (at the bottom of the **Surfaces** dialog that results) and select **Electrostatic Potential Map** from the menu.\* This requests an

---

\* Surfaces selected from this menu will be obtained at the resolution set in the last access to the **Add Surfaces** dialog. This is normally set to medium. Both lower and higher resolution images are available. The former provides rough images that require significantly less computer time to generate. The latter provides production quality images that require significantly more computer time and disk space. To specify resolution, *click* on **More Surfaces** instead of **Add** in the **Surfaces** dialog (or select **More Surfaces...** from the menu that follows from *clicking* on **Add** in the **Surfaces** dialog). This gives rise to the **Add Surfaces** dialog from which you can select the desired resolution from the **Resolution** menu. Computer time and storage can also be reduced by *checking* **Fixed** to the right of **Isovalue** in the **Add Surface** dialog. This results in constant isovalue surfaces that cannot be changed and cannot be used to generate “Slice” type displays.



electrostatic potential map (an electron density surface onto which the value of the electrostatic potential will be mapped). *Click* on **OK**. A line *Electrostatic Potential Map* appears at the top of the dialog. If you make a mistake, *click* on this line (*Electrostatic Potential Map*) and then *click* on **Delete** at the bottom of the dialog.

10. The graphics calculation will run automatically (without needing to resubmit the job) following your request.\* When it has completed (a few seconds), select *Electrostatic Potential Map* by *clicking* in the selection box in the **Surfaces** dialog. The surface itself corresponds to the electron density and provides a measure of the overall size and shape of acrylonitrile. The colors indicate values of the electrostatic potential on this surface; by convention, colors toward red correspond to negative potential (stabilizing interaction between the molecule and a positive charge), while colors toward blue correspond to positive potential. The nitrogen (the most electronegative atom) is red and the hydrogens (the most electropositive atoms) are blue.\*\*
11. Select **Close** from the **File** menu (or *click* on ) to remove acrylonitrile from the screen.\*\*\* Also, close any dialogs that may be open.


---

\* If it does not, you need to select **Auto-Gen Graphics** in the **Settings Preferences** dialog (**Preferences** under the **Options** menu; **Chapter 20**). Note that automatic generation of graphics is suppressed irrespective of the preferences setting when the request applies to a list with more than 10 molecules.


\*\* The default display is a set of discrete color bands. You can alter the number of bands by bringing up (**Molecule**) **Properties** dialog (**Properties** from the **Display** menu), *clicking* on the map to bring up the **Surface Properties** dialog and altering the number to the right of **Bands**. You can change to a continuous display by removing the checkmark to the left of **Bands**.


\*\*\* While *Spartan* permits as many molecules as desired on screen at a given time, it will be less confusing for first-time users to keep only a single molecule open at a time.


## Cyclohexanone

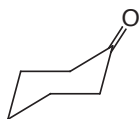

Cyclohexanone provides the opportunity to illustrate additional building features as well as graphical models for investigation of the stereoselectivity of an important class of organic reactions.

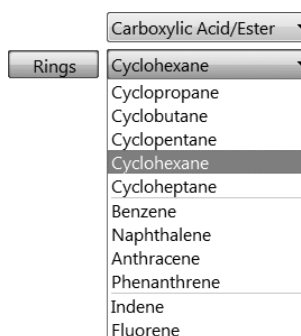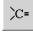
1. Click on  to bring up the organic model kit. Click on **Rings** and select **Cyclohexane** from the menu of rings.



Click anywhere on screen and cyclohexane will appear. Select  $sp^2$  carbon  from the model kit. Double click on any carbon (*not free valence*). The  $sp^3$  hybridized center will be replaced by an  $sp^2$  hybridized carbon.

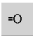*Spartan's* organic model kit allows replacement of fragments, subject both to the usual valence rules and to the availability of free valences. For example, replacement of an  $sp^3$  carbon by an  $sp^2$  carbon requires that at least two free valences be available. This is possible for cyclohexane (above), but would not have been allowed for the substituted ring carbon in methylcyclohexane.



Select  $sp^2$  oxygen  from the model kit. Click on the double free valence on the  $sp^2$  carbon. You have made cyclohexanone


(its name will appear at the bottom right of the screen). *Click* on **E** to produce a final structure with  $C_s$  symmetry, and then on **V** to remove the model kit.

2. Select **Calculations...** from the **Setup** menu. Specify **Equilibrium Geometry** from the top menu to the right of **Calculate**, and **Hartree Fock, 3-21G** and **Vacuum** from the three bottom menus to the right of **Calculate**. *Click* on **OK**.
3. Cyclohexanone undergoes nucleophilic attack at the carbonyl carbon, and it is reasonable to expect that the molecule's lowest-unoccupied molecular orbital (the LUMO) will be localized here. To visualize the LUMO, enter the **Surfaces** dialog (**Setup** or **Display** menu). *Click* on **Add** and select **LUMO** from the menu. Also request an electron density surface onto which the (absolute) value of the LUMO has been mapped in color (a so-called LUMO map). *Click* on **Add** and select **|LUMO| Map** from the menu. Submit the job (**Submit** from the **Setup** menu) with the name *cyclohexanone* (the selection that is offered).

Graphical displays may either be requested at the time of geometry calculation (as in the present example), or at anytime after (as in the previous example). In the latter case, it is not necessary to submit the job. *Spartan* keeps track of what has been done and will not perform unnecessary calculations.

4. After the calculations have completed, reenter the **Surfaces** dialog and *click* on **LUMO**. As expected, this is a  $\pi^*$  orbital primarily localized on the carbonyl group.
5. Inside the **Surfaces** dialog, turn "off" the display of the LUMO (*click* to deselect **LUMO**). Next, *click* on **|LUMO| Map** to display the electron density surface onto which the (absolute) value of the LUMO has been mapped. By convention, colors toward red indicate small (absolute) values of the LUMO (near zero), while colors toward blue indicate large (absolute) values of the LUMO. Note the blue spot directly over the carbonyl carbon. This corresponds to the maximum value of the LUMO and is where nucleophilic attack will occur.

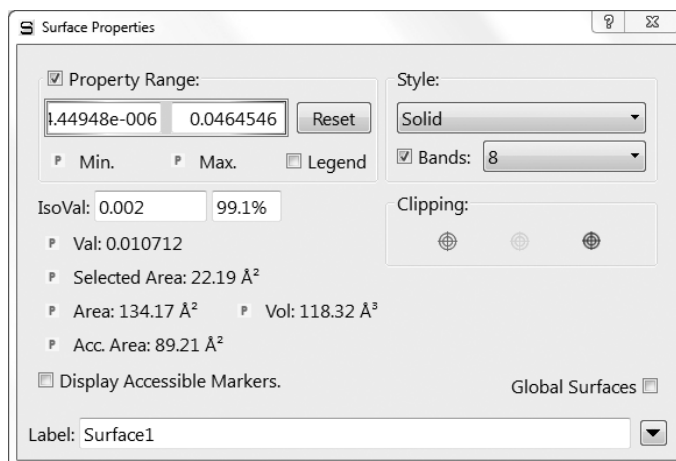6. You will see that the blue spot over the *axial* face of the carbonyl carbon is bigger (and “more” blue) than that over the *equatorial* face. This indicates preferential attack by nucleophiles onto the *axial* face. Quantify the difference by measuring the (absolute) value of the LUMO on these two faces. Bring up the (**Molecule**) **Properties** dialog (**Display** menu) and *click* anywhere on the LUMO map. The **Surface Properties** dialog appears.




*Check* the box next to **Legend** to display the property range on screen. (Remove the checkmark to remove the legend.) To translate the legend, *click* on the legend to select, then hold down the right mouse button and move the mouse. The legend is useful when making qualitative comparisons of property values. Turn the map such that you can clearly see the *axial* face of the carbonyl carbon, and *click* on the area of maximum blue. The (absolute) value of the LUMO at the surface point you have selected is provided in the dialog to the right of **Val**. Note the value, and then turn the map over such that you can now see the *equatorial* face of the carbonyl carbon, and *click* on the region of maximum blue on this face\*. Do these values support your qualitative conclusions from viewing the image?

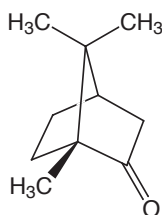7. Turn off the LUMO map (*click* to deselect). Keep cyclohexanone on screen, it will be used in the next tutorial.




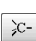



---

\* *Click* on the background to return to the **Molecule Properties** dialog.



## Camphor



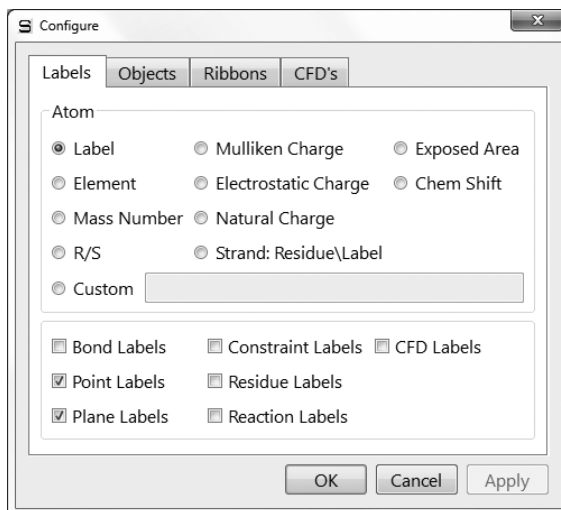
Camphor illustrates how one molecule may be used as the starting point for another. It also illustrates a search on the Cambridge Structural Database (CSD).<sup>\*</sup> No calculations are involved in this example.

1. If you do not already have cyclohexanone on screen (see previous tutorial), select **Open...** from the **File** menu (or *click* on ) and *double click* on *cyclohexanone* in the dialog which results. Select **Save As...** from the **File** menu (or *click* on the  icon at the top of the screen), supply the name *camphor* and *click* on **Save**.
2. Select **Add Fragment** from the **Build** menu (or *click* on the  icon at the top of the screen) to bring up the organic model kit. *Click* on  $sp^3$  carbon  from the model kit and then *click* on the *axial* free valence two positions removed from the carbonyl carbon. You have made *axial* 3-methylcyclohexanone. Its name will appear at the bottom right of the screen.
3. Select **Make Bond** from the **Build** menu (or *click* on the  icon at the top of the screen). *Click* on one of the free valences of the  $sp^3$  carbon which you have just added and then on the *equatorial* free valence on the opposite side of the six-membered ring (adjacent to the carbonyl group). A bond will be drawn. *Click* on  to produce a refined (intermediate) structure. It is 2-norbornone, the name of which will appear at the bottom right of the screen. Finish making camphor by adding three methyl groups (*click* on  from the model kit and, one after the other, *click* on the three appropriate free valences; refer to image above


<sup>\*</sup> CSD will need to be licensed and installed in order to complete this tutorial. For information refer to the CCDC website ([www.ccdc.cam.ac.uk](http://www.ccdc.cam.ac.uk)). For installation instructions, see **Appendix F**. If you don't have access to CSD, a trial license may be obtained from the Cambridge Crystallographic Data Centre ([www.ccdc.cam.ac.uk](http://www.ccdc.cam.ac.uk)).


for CH<sub>3</sub> positions). Finally, *click* on  to give a minimized camphor geometry with C<sub>1</sub> symmetry. Its name will appear at the bottom right of the screen. *Click* on  to remove the model kit.

4. Select **Configure...** from the **Model** menu to bring up the **Configure** dialog.




*Check R/S* under **Atom** in the dialog that results, and *click* on **OK**. R/S labels now appear (only) on the two bridgehead carbons. Both should be R (or both S). You can change from one to the other by *double clicking* on any atom while holding down both the **Ctrl** and **Shift** keys. You can turn the labels off by selecting **Labels** from the **Model** menu.

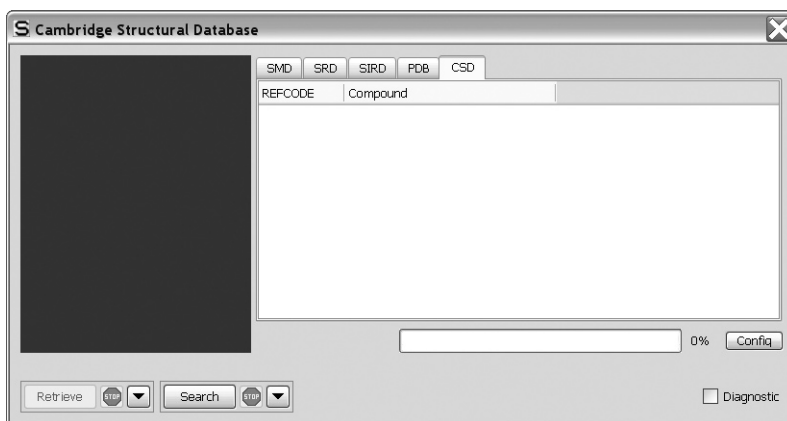
5. Move the molecule to the far right of the screen. Then, select **Open Databases** from the **Search** menu (or *click* on the  icon at the top of the screen). *Click* on the **CSD** tab at the top of the **Databases** dialog that results. This brings up the **CSD** dialog.



Click on one of the free valences on the methyl group at the bridgehead carbon adjacent to the carbonyl group. The free valence is replaced by an orange cone. This restricts the search to molecules that derive from substitution of camphor at this carbon.

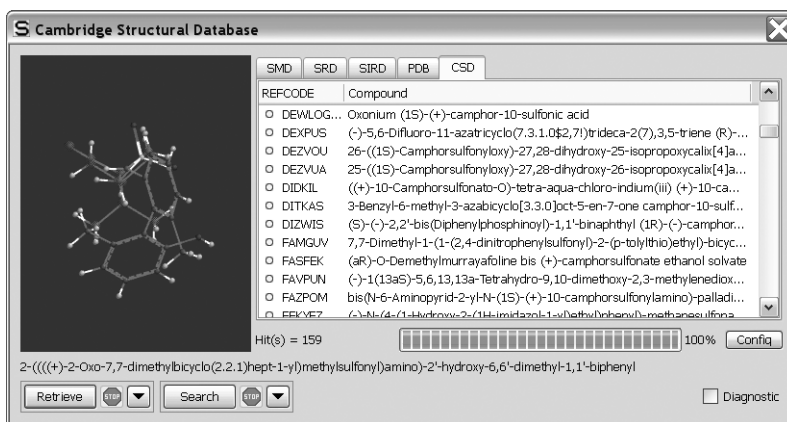- Click on the **Search** button at the bottom of the dialog. Hits will appear in the text box at the right of the **CSD** dialog.




- Wait for the search to complete. Click on any line marked by a filled yellow circle.\* Click on any one of these. Its name will be highlighted and a structure model will appear in the window at the top of the dialog. This model can be manipulated in the usual way.\*\*
- Close camphor and remove the **CSD** dialog (click on ).

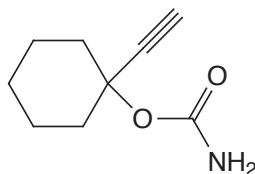\* Unfilled circles correspond to entries for which structural data are unavailable. The literature reference is available.


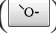
\*\* One (or more) structures may be transferred to *Spartan's* file system using the **Retrieve** button. Full discussion is provided in **Chapter 19**.


## Ethinamate





Acrylonitrile, cyclohexanone and camphor are examples of molecules that adopt only a single conformation. Much more common are molecules that can adopt more than one conformation, usually due to rotation about single bonds or puckering of rings. The sedative ethinamate offers both rotation about the ring oxygen bond and the possibility for the carbamate functionality to occupy either the *equatorial* or *axial* position on the ring.

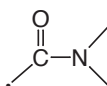
This tutorial and the one immediately following, illustrate the simplest approach to dealing with conformationally-flexible molecules, that is, finding the lowest-energy conformer (the global minimum). Tutorials in later chapters address the more difficult problem of identifying all reasonable low-energy conformers.


1. Bring up the organic model kit () . Click on **Rings**, select **Cyclohexane** from the menu and *click* anywhere on screen. Click on  $sp^3$  oxygen () in the model kit and add to either an *equatorial* or *axial* free valence on cyclohexane. Click on **Groups** and select **Amide** from the menu. Make certain that free valence on the amide group to be used for bonding is associated with carbon and not with nitrogen. This is indicated by a • on the icon at the top of the model kit.



You can move the bonding position by repeatedly *clicking* on the icon. When you are satisfied, add the amide to oxygen to make the carbamate functionality. Finally, select **Alkynyl** from the **Groups** menu and add to the same carbon on the cyclohexane ring as occupied by the carbamate group (using the remaining *axial* or *equatorial* free valence). Click on  and then on  to remove the model kit.



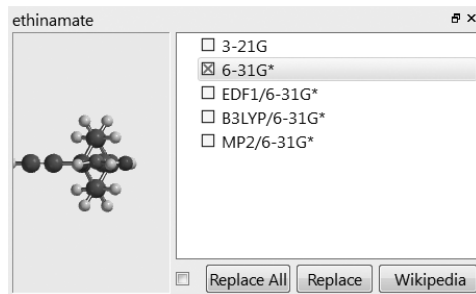- Bring up the **Calculations** dialog (**Setup** menu). Select **Equilibrium Conformer** from the top menu to the right of **Calculate** and **Molecular Mechanics** and **MMFF** from the two bottom menus. You have requested that *Spartan* examine all possible conformers\*, but then select only the most stable conformer (according to the MMFF model). *Click* on **Submit** at the bottom of the dialog. The name *ethinamate* will be suggested. Accept it or provide an alternative of your choosing.
- When completed, the preferred conformer will replace the initial conformer (just like a calculated equilibrium geometry replaces a starting geometry). Which group, alkynyl or carbamate, occupies the *equatorial* position?
- The name “ethinamate” appears at the bottom right of the screen. This indicates that the molecule exists in the Spartan Molecular Database (SMD) of structures and properties. *Click* on  immediately to the left of the name. This brings up the **Preview dialog\*\***.



Listed at the right are theoretical models for which ethanimatate is available, while a ball-and-spoke model of the selected SMD entry is displayed at the left. Select **6-31G\*** and *click* on **Replace**. The entry in SMD will replace your model. Is the conformation of ethanimatate in the database entries the same as you have found?

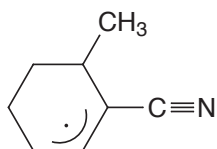
- Close *ethinamate* as well as any open dialogs.

---


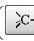



\* This system is small, and the possible conformers will be systematically examined. This is not practical for molecules with a large number of degrees of conformational freedom. Here, *Spartan* employs Monte-Carlo methods.

\*\* Make certain that SMD and not SSPD is selected from the menu at the bottom of *Spartan*'s screen.

### 3-Cyano-4-methylcyclohexenyl Radical



This tutorial illustrates how different models may be combined to achieve a desired result with minimum computation cost. MMFF molecular mechanics will be employed to establish equilibrium conformation, the semi-empirical AM1 model to establish equilibrium geometry and the Hartree-Fock 3-21G model to provide a basis for graphical calculations.

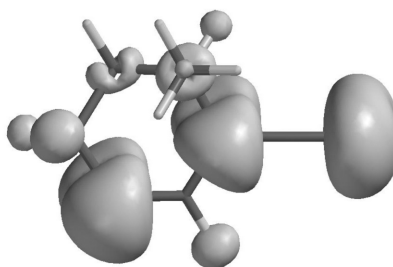
1. Build 3-cyano-4-methylcyclohexene. Start with **Cyclohexane** from the **Rings** menu, *click* on , and then *click* on *axial* free valences on any two adjacent carbons. Add **Cyano** from the **Groups** menu to one of the *equatorial* free valences on a carbon adjacent to the double bond. Finally, add a methyl group () to the appropriate ring position using either the *equatorial* or *axial* free valence.
2. Select **Delete** from the **Build** menu (or *click* on the  icon at the top of the screen), and *click* on the (*axial*) free valence on the carbon to which the cyano group is attached. *Click* on  to produce a refined geometry with C<sub>1</sub> symmetry. *Click* on  to remove the model kit.
3. Select **Calculations...** (**Setup** menu). Specify **Energy** from the top menu to the right of **Calculate**, **Hartree-Fock** and **3-21G** from the two bottom menus. Next, specify **MMFF Conformer** from the **Start from** menu and **AM1** from the menu which then appears to the right. You have requested that a Hartree-Fock energy is to be obtained using an AM1 geometry, which in turn is to be based on the best conformation according to MMFF.

This molecule has one unpaired electron. It is necessary to change **Multiplicity** from **Singlet** to **Doublet**. *Click* on **OK**.

4. Select **Surfaces** (**Setup** or **Display** menu). *Click* on **Add** and select **Spin** from the menu. *Click* on **Add** and select **Spin**

**Density Map** from the menu. You have requested two different representations of spin distribution. The first presents spin density as a surface of constant value, while the second color maps the value of the spin density onto an electron density surface. Finally, request the singly-occupied molecular orbital. *Click* on **Add** one more time and select **aHOMO** (the highest occupied molecular orbital of  $\alpha$  spin) from the menu.

5. Submit the job. Name it *3-cyano-4-methylcyclohexenyl radical*\* When completed, examine the conformation. Is the methyl group *equatorial* or *axial*? Next, enter the **Surfaces** dialog. Select *spin* to display the spin density surface. Note that the spin density is delocalized over two of the ring carbons and onto the cyano group.




6. Again enter the **Surfaces** dialog. First, remove the spin density surface. Then, display the spin density map. Select *Spin Density Map*. Note that the areas of maximum spin (colored blue) exactly match those where the surface is large in the previous image. Select **Properties (Display menu)** and *click* anywhere on the spin density map. This brings up the **Surface Properties** dialog. **Range** provides the range of spin density values. Note, in particular, that negative values of the spin density are possible.
7. Enter the **Surfaces** dialog, turn off the display of the spin density map and then select *aHOMO* to turn on the display of the  $\alpha$  HOMO. Note that the singly-occupied molecular orbital is nearly identical to the previously-displayed image of the spin.

---

\* SMD contains very few radicals and/or charged species so a chemical name will only rarely be provided.

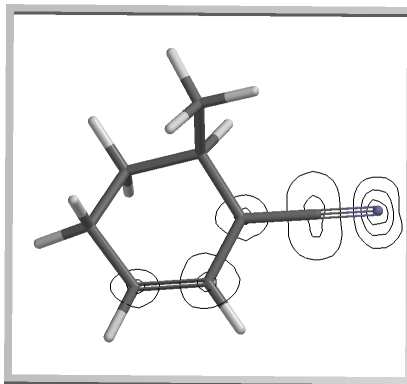8. Turn off the display of the  $\alpha$  HOMO. Then, *click* on **More Surfaces...** at the bottom of the (**Surfaces**) dialog, and select **Slice** from the **Surfaces** menu and **spin** from the **Properties** menu. *Click* on **OK**. A new line **Slice, spin** appears in the window at the top of the dialog.\* Select it.
9. A plane (a slice of spin density) surrounded by a frame appears in the middle of the model on screen. *Click* inside the frame to select. The frame will turn gold. Position the cursor outside the frame, then *press* both the **Shift** key and right button and move the mouse up and down (or use the scroll wheel). This will zoom the plane. You can also translate and rotate the plane independently of the molecule using the usual mouse operations. Alternatively, you can move the molecule and plane together by first *clicking* on the molecule (the frame will now turn white) and then using the mouse. For all operations, be certain to keep the cursor positioned outside of the frame. Size and orient the slice as you wish.\*\*




10. Remove **3-cyano-4-methylcyclohexenyl radical** and any remaining dialogs from the screen.

---

\* No further calculations are involved as this is just another way of representing the spin density.


\*\* You can change the display style from **Contours** to **Solid** or **Transparent** using the **Style** menu at the bottom right of the screen. This will appear only when the slice is selected.

## Androsterone

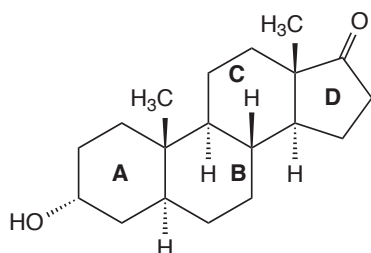


The steroid androsterone illustrates the use of ChemDraw directly from *Spartan*. Also, the molecule provides a good opportunity to compare space-filling models and electron density surfaces as metrics of molecular size and shape.

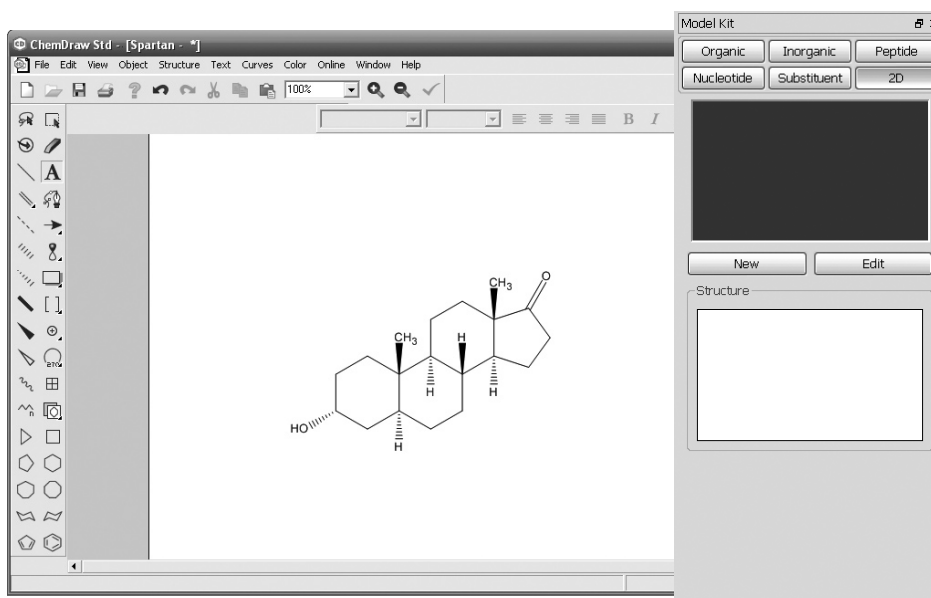1. Click on and then *click* on the **2D** tab. Click on the **New** button to open ChemDraw.\* Draw androsterone (pay close attention to stereochemistry) and close ChemDraw. You will now see the ChemDraw image of androsterone at the bottom of the 2D model kit and a 3D model at the top. Click in the workspace. The name *androsterone* will appear in the lower right hand corner. Click on (V).




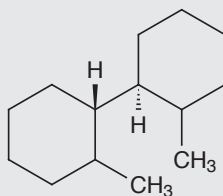
\* The 2D model kit in *Spartan* will automatically open ChemDraw if it is available on your machine (NAM, NAL). If you do not have access to ChemDraw, contact CambridgeSoft ([www.cambridgesoft.com](http://www.cambridgesoft.com)). While *Spartan* for Macintosh and Linux cannot automatically access ChemDraw, both programs can open 2D molecule files in .CDX or .SKC formats.


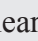
If you are not running *Spartan* for Windows, or if you do not have ChemDraw installed on your machine, follow the directions below to build androsterone.



Click on . Select **Cyclohexane** from the **Rings** menu, and *click* anywhere on screen to create the A ring. **Cyclohexane** is still selected. The icon at the top of the model kit should indicate **eq** (*equatorial*) rather than **ax** (*axial*); attachment needs to be made using an *equatorial* free valence. If not, *click* on the icon to switch to **eq**. *Click* on an *equatorial* free valence of the A ring to join it to the C ring. You have made bicyclohexane. Its name will appear at the bottom right of the screen.


Rotate around the bond connecting the A and C rings to position the *axial* free valences on the connecting carbons opposite ( $180^\circ$ ) to each other. (This bond should be encircled by a red marker. If it is not, *click* on the bond.) To do so, hold down both the **Alt** key and the left button and *drag* the mouse up and down. *Equatorial* free valences on the cyclohexane carbons adjacent to those that connect the two rings are now positioned to form the **B** ring (in a chair geometry).


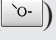


Select  $sp^3$  carbon () and *click* on the two opposing *equatorial* free valences on the rings you have just joined. The resulting structure should appear as follows.



Click on  and, one after another, *click* on the (*equatorial*) free valences that nearly overlap to create the **B** ring. *Click* on . Check to see that all rings are chairs and that the stereochemistry is as in the figure in the previous page. If not, you can start over (**Clear** from the **Edit** menu), or undo the last operation (**Undo** from the **Edit** menu).

Next, form the **D** ring starting from the carbon that will become a carbonyl group. Select **Carbonyl** from the **Groups** menu and *click* on the appropriate free valence on the C ring. Select  $sp^3$  carbon () from the model kit and *click* first on the single free valence on the carbonyl carbon and then on a free valence of the  $sp^3$  carbon just added. *Click* on  and join the free valence from the second  $sp^3$  carbon added in the previous step and the appropriate *equatorial* free valence on the C ring.

Click on . You have made one of the forms of perhydrophenanthrene and its name will appear at the bottom right of the screen.

Add two methyl groups. Select  $sp^3$  carbon () from the model kit. Click on the *axial* free valence at the appropriate juncture of the **C** and **D** rings, and then on the *axial* free valence at the proper juncture of the **A** and **B** rings. Add the hydroxy group. Select  $sp^3$  oxygen () from the model kit, and click on the appropriate *axial* position on the **A** ring. Click on  and then on . The name androsterone should appear at the bottom right of the screen.

Androsterone incorporates several chiral centers. To assign them as R or S, select **Configure...** (**Model** menu) and check **R/S** under **Atom** in the dialog that appears. Click on **OK**. R/S labels will appear next to each of the chiral centers. (To remove the labels, select **Labels** from the **Model** menu.)

2. Select **Calculations...** from the **Setup** menu and specify calculation of an energy using the **Hartree-Fock 3-21G** model. Select **AM1** from the menu to the right of **StartFrom**. Click on **Submit** and accept the name *androsterone*.
3. When completed (several minutes), select **Properties (Display menu)**. Write down the surface area of a space-filling (CPK) model given in the **Molecule Properties** dialog.
4. Select **Surfaces (Setup or Display menu)**. Click on **Add** and select **Density** from the menu. Leave the dialog on screen.
5. When completed, click on **Density** in the **Surfaces** dialog. Click on the density surface to bring up the **Surface Properties** dialog. This provides both the surface area and the volume of the electron density surface. The second box to the right of **IsoVal** contains the percentage of the total number of electrons that are enclosed by the density surface. Adjust this value by typing over the existing value and pressing the **Enter** key (**return** on Mac) such that the surface area of the density surface is nearly the same as that for the CPK model. What is the percentage enclosure?
6. Remove *androsterone* from the screen.

# Chapter 5

## Groups of Organic Molecules

---

*This chapter introduces and illustrates a number of basic operations involved in processing groups of molecules, as well as the associated spreadsheet for organizing and fitting data and facilities for making plots.*

Computational investigations, like experimental investigations, are rarely restricted to single molecules, but typically involve a series of related molecules. Here, it may be of interest to compare geometries, energies or other calculated properties, or to compare trends in calculated and measured properties. **Spartan** provides facilities for these purposes. In particular, it allows molecules to be grouped, either manually, or automatically as a result of a conformational search, from following a particular vibrational motion, or from a scan of one or more geometric variables. Once put into a group, molecules may be aligned based either on their structure or chemical functionality. Calculations may be performed either on individual molecules or, just as simply, on the complete group of molecules. The results of calculations on a group can be examined and analyzed individually or altogether to seek out trends.

Associated with a group (or an individual molecule) is a spreadsheet. This allows convenient access to virtually any calculated quantity that can be given a value. Additionally, data may be entered manually into the spreadsheet or transferred from another application such as Excel. Data in the spreadsheet may be manipulated, linear regression analyses performed and both 2D and 3D plots constructed. Alternatively, the data in a **Spartan** spreadsheet may be transferred to Excel for further analysis. The Excel worksheet can then be embedded in the **Spartan** document (see **Embedded Data** under the **File** menu; **Chapter 12**).

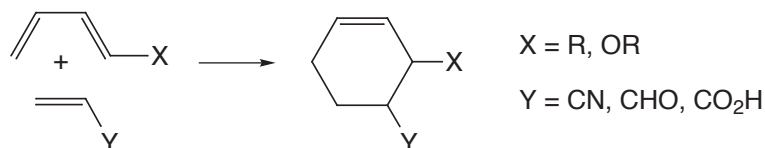
The tutorials in this chapter introduce and illustrate a number of the basic group operations available in **Spartan**. These include building a




group from scratch, processing groups resulting from a conformational search and from varying a torsion angle. Also provided is an example of fitting an experimental observable to one or more calculated properties by way of linear regression. The use of the **Reactions** dialog is illustrated both as it operates on molecules in a group and on molecules derived from substitution of grouped molecules.

## Dienophiles in Diels-Alder Cycloadditions



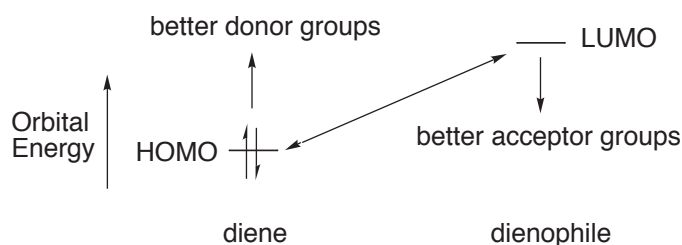The most common Diels-Alder reactions involve electron-rich dienes and electron-deficient dienophiles.






The rate of these reactions generally increases with increasing  $\pi$ -donor ability of the diene substituent, and with increasing  $\pi$ -acceptor ability of the dienophile substituent. This can be rationalized by noting that donor groups raise the energy of the highest-occupied molecular orbital (HOMO) on the diene, while acceptor groups lower the energy of the lowest-unoccupied molecular orbital (LUMO) on the dienophile. Thus, the HOMO-LUMO gap is reduced, leading to enhanced stabilization.






To test such an hypothesis, you will examine the extent to which experimental relative rates of Diels-Alder cycloadditions involving cyclopentadiene and a variety of cyanoethylenes correlate with dienophile LUMO energies.


1. Build acrylonitrile,  $\text{H}_2\text{C}=\text{C}(\text{H})\text{CN}$ , (select **New** from the **File** menu or *click* on ) , minimize () and place it on the clipboard. *Right click* on the background and choose **Copy** from the resulting menu.

Copy	Ctrl+C
Paste	Ctrl+V
Delete Selected	
Properties	

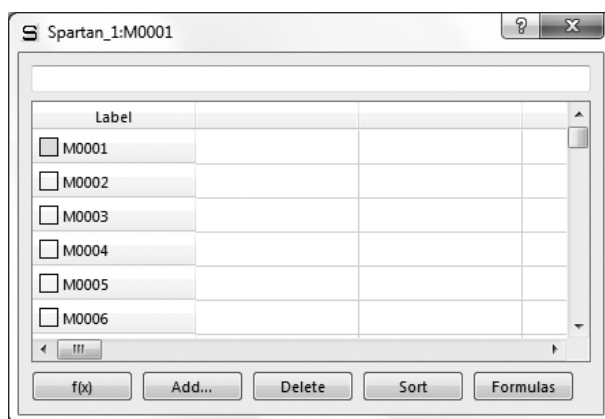
2. Select **New Molecule** (not **New**) from the **File** menu. This requests that a molecule that has yet to be built is to be appended onto the end of a group of molecules in which the currently selected molecule (acrylonitrile) is a member. The screen will be cleared. *Click* on **Clipboard** at the bottom of the model kit and *click* on screen. Acrylonitrile will appear. *Click* on **Groups** in the model kit, select **Cyano** and add to the appropriate free valence on acrylonitrile to make 1,1-dicyanoethylene. *Click* on .
3. Repeat this procedure (**New Molecule**, followed by **Clipboard**, followed by **Groups**, followed by ) four more times to build *cis* and *trans*-1,2-dicyanoethylene, tricyanoethylene and tetracyanoethylene. When you are all done (six molecules in total), *click* on  to remove the model kit.
4. The molecules have been put into a single group, allowing calculated properties to be accessed via a spreadsheet. Select **Spreadsheet** from the **Display** menu.



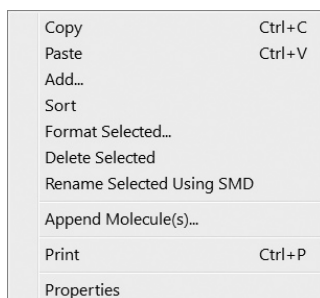To select individual molecules, *click* on their labels (**M0001**, ...) in the left hand column, or use the  and  keys at the bottom left of the screen. Position the cursor in the header field above the left-hand column and *left click*, then *right click* and select **Rename Selected Using SMD** from the menu that results.



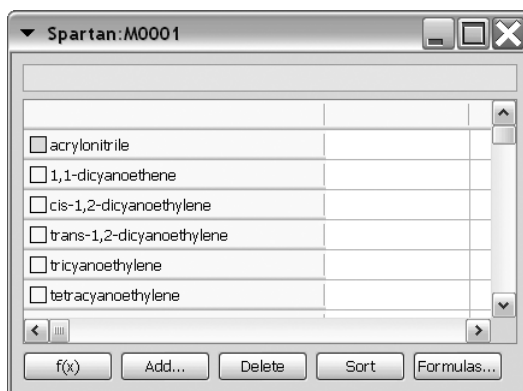
This results in replacement of the default labels by proper names from the Spartan Molecular Database.




5. Enter the **Calculations** dialog (**Calculations** under the **Setup** menu), select **Energy** using the **Hartree-Fock 3-21G** model, and **AM1** from the **Start from** menu. A 3-21G energy calculation will be preceded by calculation of geometry using the AM1 model. Make certain that **Global Calculations** is *checked* to ensure that the calculations are to be applied to all the molecules.
6. Submit the job. Name it ***Diels-Alder dienophiles***. After it completes, enter experimental relative rates into the spreadsheet.

Experimental relative rates for Diels-Alder cycloadditions of cyclopentadiene*			
dienophile	log <sub>10</sub> (relative rate)	dienophile	log <sub>10</sub> (relative rate)
acrylonitrile	0	1,1-dicyanoethylene	4.64
<i>trans</i> -1,2-dicyanoethylene	1.89	tricyanoethylene	5.66
<i>cis</i> -1,2-dicyanoethylene	1.94	tetracyanoethylene	7.61

\* J. Sauer, H. Weist and A. Mielert, *Chem. Ber.*, **97**, 3183 (1964).

Double click on the header cell for a blank column in the spreadsheet, type **Log(rate)=** and press the **Enter** (return on Mac) key. You need to press the **Enter** key (return on Mac) following each entry or use the  key. Sort by relative rate. Click on the column header **Log(rate)**, and then click on **Sort** at the bottom of the spreadsheet.

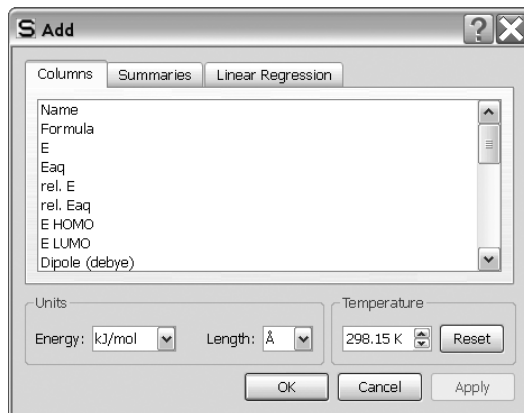- Click on a header cell for a blank column and click on **Add...** at the bottom of the spreadsheet. Alternatively, right click inside the header cell for a blank column in the spreadsheet, and select **Add...** from the menu that results. Select **E LUMO** from the list in the dialog that results and **eV** from the **Energy** menu. Click on **OK**.




Two data columns are added to the spreadsheet.

	Log(rate)	E LUMO (eV)
<input type="checkbox"/> acrylonitrile	0.00000000	2.66409914
<input type="checkbox"/> trans-1,2-dicyanoethylene	1.89000000	0.86282668
<input type="checkbox"/> cis-1,2-dicyanoethylene	1.94000000	0.935615905
<input type="checkbox"/> 1,1,1,2-dicyanoethene	4.46000000	1.20423400
<input type="checkbox"/> tricyanoethylene	5.66000000	-0.351754439
<input type="checkbox"/> tetracyanoethylene	7.61000000	-1.47917760


The spreadsheet has now served its purpose by collecting both calculated (LUMO energy) and experimental (relative rate) data for all six dienophiles. You can, if you wish, now remove it from

the screen (*click* on  at the top).

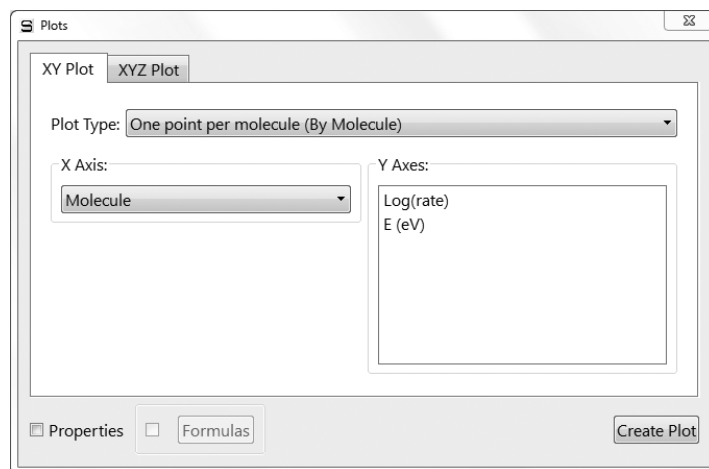
8. Select **Plots...** from the **Display** menu. This leads to the **Plots** dialog.



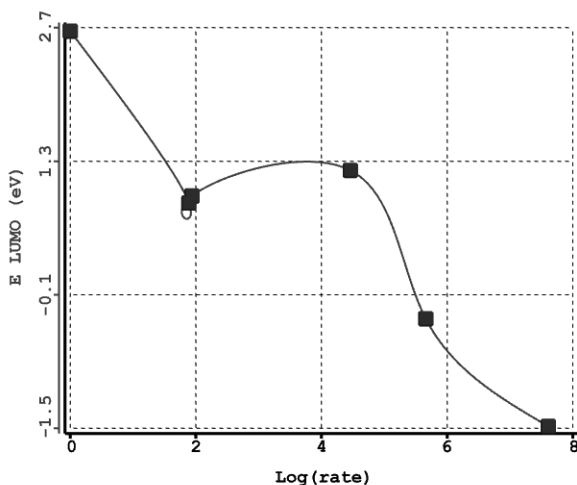
Make certain that the **XY Plot** tab at the top of the dialog is selected. Select **Log(rate)** from the list of items in the **X Axis** menu and **E LUMO(eV)** from the **Y Axes** list. *Click* on **Create Plot** to remove the dialog and show the plot.



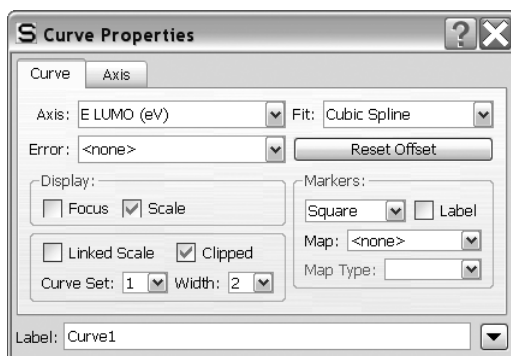
The default curve is a so-called cubic spline that smoothly connects the data points. To establish a correlation, use a simple linear least-squares fit. Select **Properties** from the **Display** menu and *click* anywhere on the curve that you have just drawn. The **Curve Properties** dialog appears.



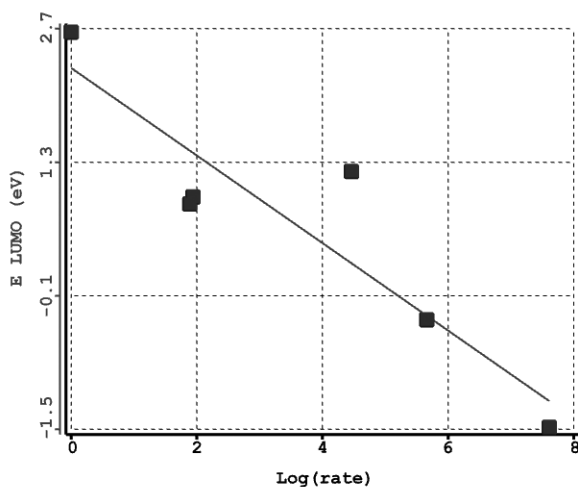Select **Linear LSQ** from the **Fit** menu. A new plot appears.




Note the correlation between relative rates and LUMO energy.



9. If the **Database Preview** dialog is not already on screen, *click* on  at the bottom of the screen to bring it up. *Check* the box to the left of B3LYP/6-31G\* in the dialog and *click* on **Replace All**. The LUMO energies you have just calculated will be replaced by values from the Spartan Molecular Database, and the plot will be updated.
10. Remove all molecules and any remaining dialogs from the screen.

## Addition vs. Substitution

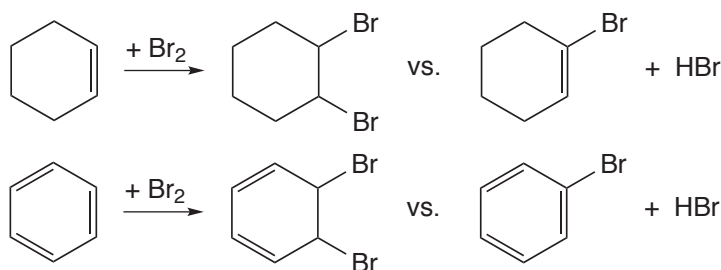
Alkenes normally undergo addition reactions whereas aromatic compounds normally undergo substitution reactions. For example, bromine reacts with cyclohexene to give *trans*-1,2-dibromocyclohexane (the addition product) not 1-bromocyclohexene, whereas it reacts with benzene to give bromobenzene (the substitution product) not *trans*-5,6-dibromo-1,3-cyclohexadiene.



In this tutorial you will use results from the Spartan Molecular Database as accessed from *Spartan's* reaction energy calculator in an effort to establish the change in preferred product. Specifically, heats of formation from the T1 thermochemical recipe will be employed. No calculations are involved.

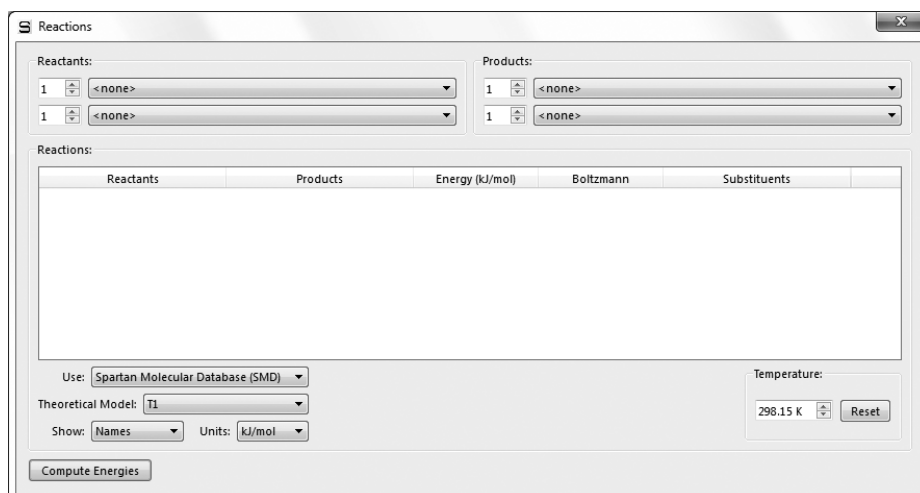1. One after another, build cyclohexene, *trans*-1,2-dibromocyclohexane, 1-bromocyclohexene, benzene, *trans*-5,6-dibromo-1,3-cyclohexadiene, bromobenzene, bromine (Br<sub>2</sub>) and hydrogen bromide (eight molecules in total). Put all in the same document. Use **New** from the **File** menu () to enter the builder for the first molecule, and **New Molecule (File menu)** for each successive molecule.
2. Select **Spreadsheet** from the **Display** menu. *Left click* inside the header cell for the leftmost column, then *right click* and select **Rename Selected Using SMD** in the menu that appears. *Click* on **OK** in the dialog that results. The molecules can now be referenced using chemical names. Close the spreadsheet (*click* on  at the top).
3. Select **Reactions** from the **Display** menu.






Choose **SMD Database** from the menu to the right of **Use** and **T1** from the menu to the right of **Theoretical Model** at the bottom of the dialog. Compute the energy for Br<sub>2</sub> addition to cyclohexene: select *cyclohexene* and *bromine* as **Reactants** and *trans-1,2-dibromocyclohexane* (and *<none>*) as **Products** and *click* on the **Compute Energies** button at the bottom left of the dialog. Record the reaction energy. Repeat for the corresponding substitution energy (same reactants but the products are 1-bromocyclohexene and hydrogen bromide) and for both addition and substitution reactions of benzene (reactants are *benzene* and *bromine* and products are *trans-5,6-dibromomo-1-3-cyclohexadiene* and *<none>* for addition and *bromobenzene* and *hydrogen bromide* for substitution).


Are all reactions thermodynamically favorable (*exothermic*)? Identify any reactions that are not and provide a rationale as to why they are not. Why is there a change in preferred reaction in moving from the alkene to the arene?

4. Close the (unnamed) document and any open dialogs.

## Allyl Vinyl Ether

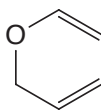

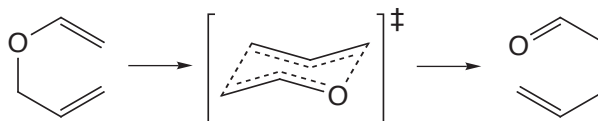Allyl vinyl ether undergoes Claisen rearrangement, the mechanism presumes a chair arrangement of the reactant.




Is this the lowest-energy conformer (global minimum) or is additional energy required to properly orient the molecule for reaction? To find out, you need to locate all the conformers of allyl vinyl ether, identify the best chair structure and evaluate its energy relative to that of the actual global minimum. You will first carry out a conformational search using the MMFF molecular mechanics model, and then obtain relative conformer energies using single-point 6-31G\* Hartree-Fock calculations based on 3-21G Hartree-Fock equilibrium geometries.

1. Bring up the organic model kit and build allyl vinyl ether from a sequence of  $sp^2$  and  $sp^3$  carbons and an  $sp^3$  oxygen. Click on to produce a refined geometry. Click on .
2. Bring up the **Calculations** dialog and select **Conformer Distribution** from the top menu to the right of **Calculate** and **Molecular Mechanics** and **MMFF** from the two bottom menus. Click on **OK**.
3. Submit the job. Accept the name *allyl vinyl ether*. When completed, it will give rise to a series of low-energy conformers\* placed in a new file *allyl vinyl ether.Conf.M0001*.\*\* A prompt will ask you if you want to open this file. Click on **OK**.\*\*\* Select **Spreadsheet** from the **Display** menu.


\* By default, only conformers within 40 kJ/mol of the global minimum will be kept. This can be changed (see **Conformational Search** in **Appendix C**).


\*\* **M0001** is the default label of the molecule you built. You can change it by altering the **Labels** field in the **Molecule Properties** dialog (**Properties** under the **Display** menu; **Chapter 18**).

\*\*\* To avoid confusion, it is a good idea to close the original file *allyl vinyl ether*.

Size the spreadsheet such that all rows (corresponding to different conformers) are visible at one time. *Click* on **Add...** at the bottom of the spreadsheet. Select **E** from the list at the top of the dialog which appears, **kJ/mol** from the **Energy** menu and *click* on **OK**. Energies for each of the different conformers will be added to the spreadsheet. Examine the lowest-energy conformer (the top entry). Is it in a chair conformation suitable for Claisen rearrangement? If not, identify the lowest-energy conformer that is suitable. You can keep two or more conformers on screen at the same time by *checking* the boxes immediately to the left of the molecule labels (the leftmost column) in the spreadsheet. To get a clearer idea of structural similarities and differences, align the conformers. Select **Align** from the **Geometry** menu (or *click* on the  icon at the top of the screen) and select **Structure** from the **Align by** menu at the bottom right of the screen. A message will appear at the bottom left of the screen.


Select atoms.

*Click* on oxygen and on the two carbons of the vinyl group bonded to oxygen to designate them as alignment centers. Each will be marked by a red circle. If you make a mistake, *click* on the circle and it will disappear. When you are done, *click* on the **Align by** button at the bottom right of the screen. *Click* on .

- To obtain a better estimate of the energy required to twist allyl vinyl ether into a conformer suitable for Claisen rearrangement, perform 6-31G\* energy calculations using 3-21G geometries. Make a copy of *allyl vinyl ether.M0001* () and name it *allyl vinyl ether Hartree-Fock*. Using the copy, delete all conformers except the global minimum and the lowest-energy conformer poised for Claisen rearrangement. Select each conformer to be discarded, and then select **Delete Molecule** from the **File** menu. Alternatively, *click* on the cell in the spreadsheet containing the label for the molecule to be deleted, and then *click* on **Delete** at the bottom of the spreadsheet. When you are done, the spreadsheet should contain only two entries.

5. Enter the **Calculations** dialog, and specify calculation of an energy using the **Hartree-Fock 6-31G\*** model. Also, specify **3-21G** from the **Start from** menu. Make certain that **Global Calculations** at the bottom of the dialog is checked to specify that the dialog settings apply to both conformers. *Click* on **Submit**.
6. When completed, examine the conformer energies. How much energy is needed to go from the global minimum to a conformer poised to undergo Claisen rearrangement?
7. Remove all molecules and any remaining dialogs from the screen.

## Internal Rotation in Dimethylperoxide

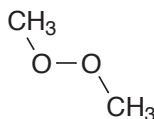


Quantum chemical calculations, in particular, Hartree-Fock molecular orbital calculations, density functional calculations and MP2 calculations, may be called on to furnish data to parameterize empirical energy functions for use in molecular mechanics and/or molecular dynamics calculations. Most important are data relating to torsional motions, for it is here that experimental data are most scarce.




Note that the function chosen to represent the energy of rotation about a single bond needs to reflect the inherent periodicity, that is, it must repeat itself in  $360^\circ$ . One suitable choice is a truncated Fourier series.

$$E^{\text{torsion}}(\omega) = k^{\text{torsion1}} (1 - \cos(\omega - \omega^{\text{eq}})) + k^{\text{torsion2}} (1 - \cos 2(\omega - \omega^{\text{eq}})) \\ + k^{\text{torsion3}} (1 - \cos 3(\omega - \omega^{\text{eq}}))$$

Here,  $\omega^{\text{eq}}$  is the ideal dihedral angle and  $k^{\text{torsion1}}$ ,  $k^{\text{torsion2}}$  and  $k^{\text{torsion3}}$  are parameters. The first (one-fold) term accounts for the difference in energy between *syn* and *anti* conformers, the second (two-fold) term for the difference in energy between planar and perpendicular conformers, and the third (three-fold) term for the difference in energy between eclipsed and staggered conformers.



In this tutorial, you will generate a potential energy function for rotation about the oxygen-oxygen bond in dimethylperoxide using 6-31G\* Hartree-Fock calculations (and optionally from B3LYP/6-31G\* density functional calculations) and then fit this potential to a truncated Fourier series.

1. Build dimethylperoxide. If the molecule is not already in an *anti* conformation, *click* on  and set the COOC dihedral angle to **180** ( $180^\circ$ ) by *typing* **180** in the box at the lower right of the screen and *pressing* the **Enter** key (**return** key on Mac). ***Do not minimize.***

2. Select **Constrain Dihedral** from the **Geometry** menu (or *click* on the  icon at the top of the screen). Select the COOC torsion, and then *click* on  at the bottom right of the screen. The icon will change to  indicating that a dihedral constraint is to be applied.



*Check* the box to the left of **Profile** at the bottom right of the screen. This will result in two additional text boxes.



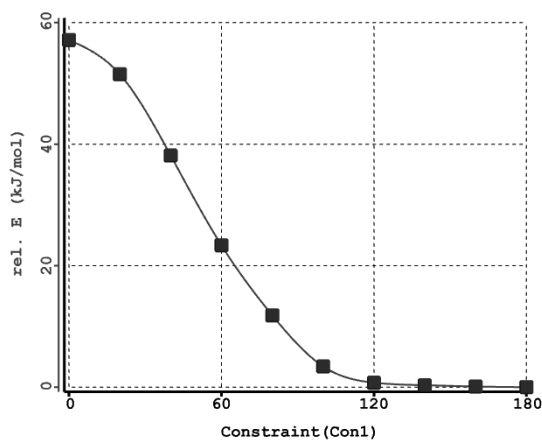
- Leave **180** ( $180^\circ$ ) in the original (leftmost) box alone, but change the contents of the box to the right of **to** to **0** ( $0^\circ$ ). You need to *press* the **Enter (return)** key after you type in the value. **Steps** should be **10**. If it is not, *type 10* and *press* the **Enter (return)** key. What you have specified is that the dihedral angle will be constrained first to  $180^\circ$ , then to  $160^\circ$ \*, etc. and finally to  $0^\circ$ . *Click* on .
4. Bring up the **Calculations** dialog and select **Energy Profile** from the top menu to the right of **Calculate**, and **Hartree-Fock** and **6-31G\*** from the two bottom menus. *Click* on **Submit** and accept the name *dimethylperoxide*.
  5. When the calculations have completed, they will go into a file named *dimethylperoxide.Prof.M0001*. A prompt will ask you if you want to open this file. *Click* on **OK**. Align the conformers. *Click* on , select **Structure** from the **Align by** menu at the bottom right of the screen and, one after the other, *click* on both oxygens and on one of the carbons. Then *click* on the **Align by** button at the bottom right of the screen. Bring up the spreadsheet (**Display** menu), and enter both the energies relative to the  $180^\circ$  conformer, and the COOC dihedral angles. First *click* on the label (**M0001**) for the top entry in the spreadsheet (the  $180^\circ$  conformer), then *click* on the header cell for the leftmost blank column, and finally, *click* on **Add...** at the bottom of the spreadsheet. Select **rel. E** from the quantities in the dialog that results, **kJ/mol** from


---


\* The difference between constraint values is given by: (final-initial)/(steps-1).

the **Energy** menu and *click* on **OK**. To enter the dihedral angle constraints, *click* on , *click* on the constraint marker attached to dimethylperoxide and *click* on **P** at the bottom of the screen (to the right of the value of the dihedral angle). *Click* on .

6. Select **Plots...** (**Display** menu). Select **Constraint(Con1)** from the **X Axis** menu and **rel. E(kJ/mol)** from the **Y Axes** list. *Click* on **Create Plot**.



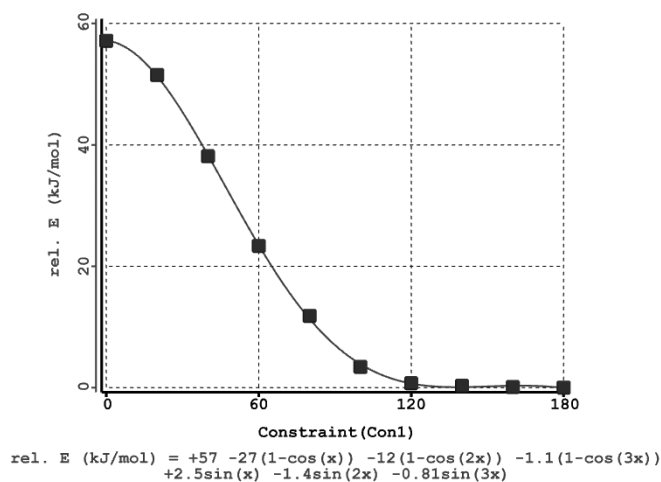A smooth curve (called a cubic spline) connects the data points. To fit the points to a Fourier series, bring up the **Properties** dialog (**Display** menu) and *click* anywhere on the curve. Select **Fourier LSQ** from the **Fit** menu and *click* on  to remove the dialog. A new plot appears.



The actual fit expression appears at the bottom of the plot.


### 7 and 8 optional

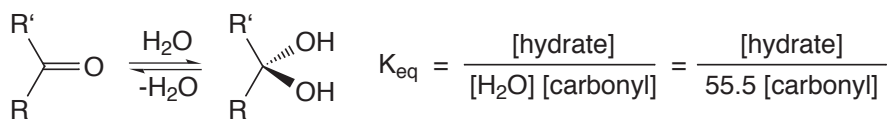7. Use energies from B3LYP/6-31G\* density functional calculations, together with the Hartree-Fock geometries to provide a better fitting function. First, make a copy of *dimethylperoxide.Prof.M0001*. Name it *dimethylperoxide density functional*. Next, enter the **Calculations** dialog and specify calculation of **Energy** using the **B3LYP/6-31G\* Density Functional** model. Make certain that **Global Calculations** (at the bottom of the dialog) is *checked* to signify that energy calculations are to be performed on all conformers. *Click* on **Submit**.
8. It involves energy calculations for all ten conformers and will require 10-20 minutes to complete. When it is done, draw a new energy plot and compare it to the energy plot produced earlier.
9. Remove all molecules and dialogs from the screen.



## Hydration of Carbonyl Compounds

The hydration of carbonyl compounds has been extensively studied primarily because it serves as a model for a number of important reactions, nucleophilic addition to carbonyl compounds among them.



Experimental $K_{\text{eq}}$ for hydration of carbonyl compounds*			
	$\log(k_{\text{eq}}/55.5)$		$\log(k_{\text{eq}}/55.5)$
PhCOCH <sub>3</sub>	-6.8	CF <sub>3</sub> COCH <sub>3</sub>	-0.2
CH <sub>3</sub> COCH <sub>3</sub>	-4.6	PhCOCF <sub>3</sub>	0.1
PhCHO	-3.8	H <sub>2</sub> CO	1.6
<i>t</i> -BuCHO	-2.4	CF <sub>3</sub> CHO	2.7
CH <sub>3</sub> CHO	-1.7	CF <sub>3</sub> COCF <sub>3</sub>	4.3

\* J.P. Guthrie, Can. J. Chem., **53**, 898 (1975); **56**, 962 (1978).



In this tutorial, you will use *Spartan*'s linear regression analysis tool to correlate calculated properties of carbonyl compounds with measured equilibrium constants for their hydration.

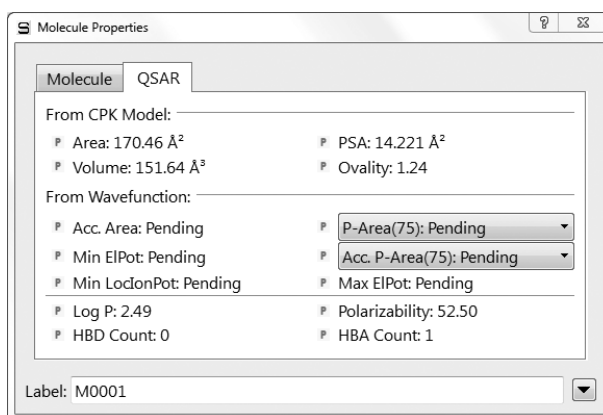
1. Build the compounds\* listed above. Use **Carbonyl (Groups menu)** as a starting point for each to ensure consistent atom numbering. Put in one list; select **New Molecule** following construction of each molecule, that is, build and minimize acetophenone, select, **New Molecule**, build and minimize acetone, select **New Molecule**, etc. Click on .
2. Enter the **Calculations** dialog and specify calculation of **Equilibrium Geometry** using the **Hartree-Fock 3-21G** model.\*\*

\* Keep the phenyl ring and the carbonyl group coplanar for benzaldehyde, acetophenone and trifluoroacetophenone.

\*\* You can save computer time by specifying calculation of an **Energy** using the **Hartree-Fock 3-21G** model, and then selecting **AM1** from the **Start from** menu. This requests a semi-empirical AM1 geometry in place of the 3-21G geometry. As all of these molecules are in SMD, you can save even more time by bringing up the **Preview** dialog (click on  at the bottom right of the screen), clicking on the yellow square to the left of **3-21G** and then clicking on **Replace All** at the bottom of the dialog. (You need to be certain that SMD and not SSPD is selected in the menu at the bottom of *Spartan*'s screen.) In this case, skip the next step as all the results you require are already available.

Make certain that **Global Calculations** is *checked*. Click on **Submit**. Name the job *carbonyl compounds*.

3. After the calculations complete, bring up the spreadsheet. *Double click* inside the header cell of an empty column, *type* **Log(Keq)** and *press* the **Enter** key (**return** key on Mac). Then enter the experimental equilibrium constants from the table on the previous page. You need to *press* the **Enter** (**return**) key following each data entry (or use the  key). Sort the list according to the value of Log(Keq). *Click* inside the header cell **LogKeq** and then *click* on **Sort** at the bottom of the spreadsheet. Also, replace the molecule identifiers in the spreadsheet (*M0001*, ...) by correct names (*acetophenone*, ...). Position the cursor inside the header cell of the leftmost column of the spreadsheet, *left click*, then *right click* and choose **Rename Selected Using SMD** from the menu that appears.
4. Select **Properties** (**Display** menu) to bring up the **Molecule Properties** dialog. *Click* on  to the left of the following entries: **E HOMO**, **E LUMO**, **Dipole Moment** and **Weight**.
5. *Click* on the **QSAR** tab at the top of the **Molecule Properties** dialog.




*Click* on  to the left of **Volume** under **From CPK Model**.


6. *Click* on the oxygen atom for whatever compound is displayed. *Click* on  to the left of **Electrostatic** (under **Charges**) in the **Atom**

**Properties** dialog (that has replaced the **Molecule Properties** dialog), to place oxygen charges into the spreadsheet.

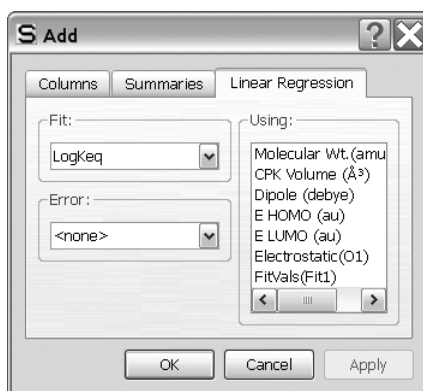7. Click on **Add...** at the bottom of the spreadsheet, and then *click* on the **Linear Regression** tab at the top of the dialog that results. This brings up the **Add Linear Regression** dialog.



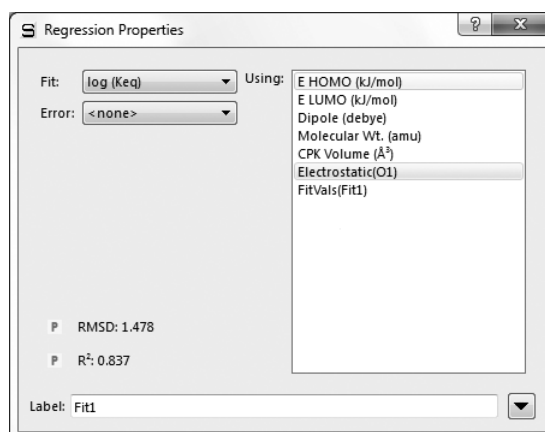Select **Log(Keq)** from the **Fit** menu and **Electrostatic (Ox)**\* from the **Using** list. Click on **OK**. Two new rows will be added at the bottom of the spreadsheet. With a **Properties** dialog on screen, *click* anywhere on the row in the spreadsheet labelled **Fit1** (the first fit performed on this data set). The **Regression Properties** dialog appears.



---

\* Numbering depends on the order that atoms were introduced during building, but the reference is to the carbonyl oxygen.

This provides information about the fit of  $\text{Log}(\text{Keq})$  to the charge on oxygen, in particular, the value of  $R^2$ . This will tend to unity as the quality of the fit improves.


8. Inside the **Regression Properties** dialog, try fitting to different (single) properties in the **Using** list. (The entries toggle on and off with repeated *clicking*.) Maximize the value of  $R^2$ . Next, try all combinations of two properties. The combination of **E LUMO** and **Electrostatic (Ox)** should give the highest value of  $R^2$ . With this combination selected, bring up the **Plots** dialog (**Display** menu) and select **Log(Keq)** from the **X Axis** menu, and **FitVals (Fit1)** from the **Y Axes** list and *click* on **Create Plot**. Next, *click* on the curve (not on the axes) to bring up the **Curve Properties** dialog (which replaces the **Regression Properties** dialog). Select **Linear LSQ** (to replace **Cubic Spline**) from the **Fit** menu. The resulting plot should show a reasonable correlation.
9. Close *carbonyl compounds* and remove any remaining dialogs from the screen.

## Acidities of Carboxylic Acids

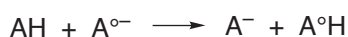
Acid strength is among the most important molecular properties. It is readily available from calculation, either in terms of absolute deprotonation energy,



or, more commonly, as the deprotonation energy relative to that of some standard acid ( $\text{A}^\circ\text{H}$ ).



The value of the electrostatic potential in the vicinity of the acidic hydrogen in the neutral acid provides a possible additional measure of (relative) acid strength. Electrostatic potential maps are known to uncover gross trends in acidity, for example, the acidic hydrogen in a strong acid, such as nitric acid, is more positive than that in a weak acid, such as acetic acid, which in turn is more positive than that in a very weak acid, such as ethanol.

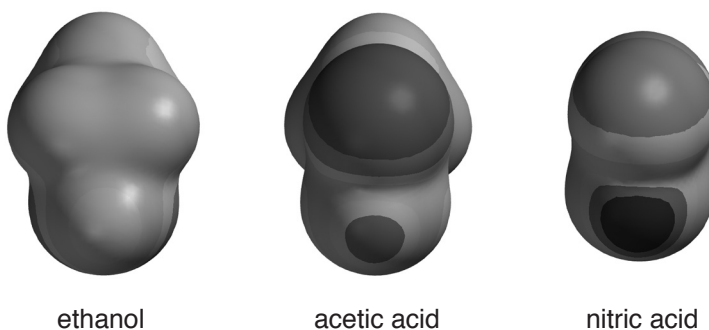

In this tutorial, you will explore the use of electrostatic potential maps to quantify changes in acid strength due to subtle variations in structure. You will use the Spartan Molecular Database of calculated structures, energies and properties to shorten the time required for an investigation.

1. One after the other, build trichloroacetic, dichloroacetic, chloroacetic, formic, benzoic, acetic and pivalic acids. Put all molecules into the same document (use **New Molecule** from the **File** menu following the first molecule). *Click* on  when you are done.

2. All of the molecules that you have built are available in the Spartan Molecular Database. *Click* on  to the left of the name of the presently selected molecule at the bottom of the screen. The **Preview** dialog will appear. (You need to be certain that **SMD** and not **SSPD** is selected in the menu at the bottom of *Spartan*'s screen.) This displays the structure of the selected molecule on the left and a list of theoretical models for which calculated results are available on the right. *Check* the yellow square to the left of **3-21G**, *click* on **Replace All** at the bottom of the dialog and respond **OK** to the message that follows. Structures obtained from Hartree-Fock 3-21G calculations will replace those you have built. However, you need to calculate a wavefunction for each of the molecules in order to produce electrostatic potential maps.
3. Enter the **Calculations** dialog (**Calculations** from the **Setup** menu). Request a **Hartree-Fock 3-21G Energy** calculation, and make certain that **Global Calculations** is *checked*. *Click* on **OK**. Enter the **Surfaces** dialog (**Surfaces** from the **Setup** menu). Make certain that **Global Surfaces** at the bottom of the dialog is checked. *Click* on **Add** at the bottom of the dialog and select **Electrostatic Potential Map** from the menu. Leave the **Surfaces** dialog on screen. Submit the job (**Submit** from the **Setup** menu). Name it *carboxylic acids*.
4. After the calculations have completed, bring up the spreadsheet (**Spreadsheet** from the **Display** menu). Expand it so that you can see all seven molecules, and that two data columns are available. Position the cursor in the leftmost header cell of the spreadsheet, *right click* and select **Rename Selected Using SMD** in the menu that appears. *Double click* inside the header cell of the next available data column, *type* **pKa** and *press* the **Enter** key (**return** key on Mac). Enter the experimental pK<sub>a</sub>'s (given on the next page) into the appropriate cells under this column. You need to *press* the **Enter** (**return**) key following each entry. Alternatively, use the  key to move to the next cell.


acid	pK <sub>a</sub>	acid	pK <sub>a</sub>
trichloroacetic (Cl <sub>3</sub> CCO <sub>2</sub> H)	0.7	benzoic (C <sub>6</sub> H <sub>5</sub> CO <sub>2</sub> H)	4.19
dichloroacetic (Cl <sub>2</sub> CHCO <sub>2</sub> H)	1.48	acetic (CH <sub>3</sub> CO <sub>2</sub> H)	4.75
chloroacetic (ClCH <sub>2</sub> CO <sub>2</sub> H)	2.85	pivalic ((CH <sub>3</sub> ) <sub>3</sub> CCO <sub>2</sub> H)	5.03
formic (HCO <sub>2</sub> H)	3.75		

Experimental data from: E.P. Sargeant and B. Dempsey, *Ionization Constants of Organic Acids in Aqueous Solution*, IUPAC no. 23, Pergamon Press, 1979.


5. Arrange the seven molecules on screen such that you can clearly see the acidic hydrogen on each. To display all molecules at once, *check* the box on the right hand side of the **Label** column (in the spreadsheet) for each entry. To manipulate them independently of one another, *click* to deselect **Coupled** from the **Model** menu.
6. Select **Electrostatic Potential Map** inside the **Surfaces** dialog to turn on the electrostatic potential map for each molecule. Bring up the **Properties** dialog (**Properties** from the **Display** menu). If the **Surface Properties** dialog is not already displayed, *click* on the electrostatic potential map to display it. *Click* on  to the right of **Max** inside the **Surface Properties** dialog. The maximum value of the electrostatic potential (corresponding to the acidic proton) will be pasted to the spreadsheet. Leave the dialog on screen.
7. Plot experimental pK<sub>a</sub> vs. maximum in the electrostatic potential. Bring up the **Plots** dialog (**Plots** from the **Display** menu), select **pKa** under the **X Axis** menu and **Property Max (Surface)** from the **Y Axes** list, and *click* on **Create Plot**. The data points are connected by a smooth curve (a so-called cubic spline). To get a least squares fit, *click* on the curve, and select **Linear LSQ** from the **Fit** menu in the **Curve Properties** dialog. Does there appear to be a correlation between pK<sub>a</sub> and the maximum value of the electrostatic potential?
8. Close **carboxylic acids** and any remaining dialogs from the screen.

## Positional Selectivity of Substituted Naphthalenes

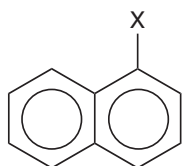Thermochemical stabilities of positional isomers may depend on several factors, including a tendency to minimize unfavorable non-bonded intramolecular interactions (*sterics*) and to minimize overall dipole moment (*electrostatics*). Naphthalene offers a good example with only two different positions.




1-substituted naphthalene

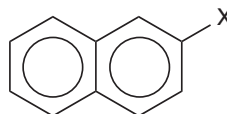


2-substituted naphthalene


In this tutorial you will use results from the Spartan Molecular Database as accessed from *Spartan's* reaction energy calculator to assign in an effort to establish the change in preferred product. Specifically, heats of formation from the T1 thermochemical recipe will be employed. Instead of building a series of substituted naphthalenes, you will use *Spartan's* substituent model kit to construct a pair of molecules, one marked for substitution on the 1-position and the other on the 2-position. No calculations are involved.

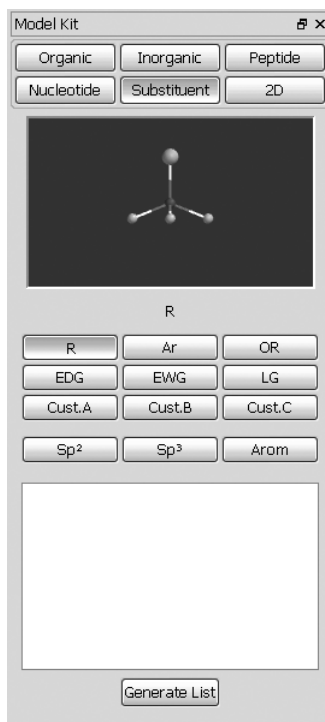1. One after the other build methane, ammonia, water hydrogen fluoride, hydrogen cyanide and formic acid. Put all in the same list (use **New** to enter the builder for the first molecule and **New Molecule** for each successive molecule). Bring up the spreadsheet (**Spreadsheet** under the **Display** menu), *double click* inside the leftmost cell for each molecule and replace the identifier (*M0001*, ...) by proper functional group names (*methyl*, *amino*, *hydroxy*, *fluoro*, *cyano* and *carboxylic acid*). Put the contents on the clipboard. *Right click* inside the header cell for the leftmost column of the spreadsheet and select **Copy** for the menu that results. Close the document without saving.
2. Select **New** from the **File** menu (or *click* on ) and build naphthalene. *Click* on the **Substituent** tab at the top of the organic model kit to bring up the substituent model kit.






Click on **Cust. A** near the middle of the model kit, then *right click* inside the box near the bottom and select **Paste** from the menu that results. Names of the six functional groups that you built in the previous step will appear in the box. As you *click* through the list, ball-and-wire models of their structures will appear in the small screen at the top of the model kit. Note that one of the hydrogens (white balls) for each is highlighted in gold. This hydrogen will be removed to make a free valence for attachment. The carboxylic acid group has two different hydrogens, one on carbon (leading to a carboxylic acid) and one on oxygen (leading to a formate). If the hydrogen on the carbon is not already highlighted, *click* on it and it will be highlighted. *Click* on the free valence on naphthalene that will lead to 1-substituted naphthalenes. The molecule on screen will now appear with a marker to indicate substitution at the 1-position.

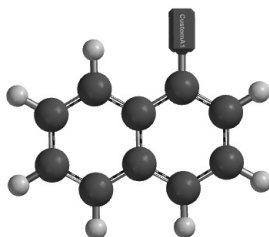





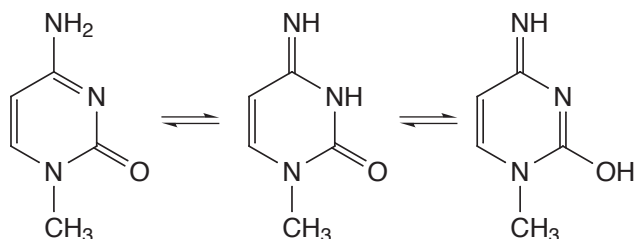

3. Select **New Molecule** (**File** menu). Bring up the organic model kit and again build naphthalene. Reenter the substituent model kit, make certain that **Cust. A** is selected (if not, *click* on it) and *click* on the free valence on naphthalene that will lead to 2-substituted naphthalenes. *Click* on  to remove the model kit.
4. Bring up the spreadsheet (**Spreadsheet** from the **Display** menu). *Double click* inside the leftmost cell for the first molecule and change the default identifier (**M0001**) with ***1-substituted naphthalenes***. *Double click* inside the corresponding cell for the second molecule and replace the default name by ***2-substituted naphthalenes***.
5. Bring up the **Reactions** dialog (**Reactions** from the **Display** menu). Choose **SMD Database** from the menu to the right of **Use** and **T1** from the menu to the right of **Theoretical Model**. Select ***1-substituted naphthalenes*** as **Reactants** and ***2-substituted naphthalenes*** as **Products**, and *click* on **Compute Energies** at the bottom of the dialog. Reactions will be written (and reaction energies computed) for all six substituents.
6. Close the (unnamed) document and any open dialogs.

## Tautomers of Nucleotide Bases

Protons bound to heteroatoms in heterocyclic compounds are likely to be very mobile in solution and, where two or more heteroatoms are present in a structure, different isomers (tautomers) may be in equilibrium. As a case in point, consider the nucleotide base cytosine (where a methyl group has replaced the sugar-phosphate backbone at the 1-position).



The existence of a low-energy tautomer could have far-reaching consequences, given that the valence structure of cytosine is key to hydrogen bonding in DNA. In this tutorial, you will examine the possible tautomers of 1-methylcytosine for evidence of low-energy structures.\*

1. Build 1-methylcytosine. Start with **Amide** from the **Groups** menu. Select planar trigonal nitrogen for the external amino group. Minimize (*click* on ) and *click* on .
2. Note that the word **Tautomer** appears at the bottom right of the screen, indicating that tautomers exist. Select **Identify Tautomers** from the **Search** menu (or *click* on the  icon at the top of the screen). Step through the tautomers using the  and  keys that appear at the bottom right of the screen. To put the tautomers in a list, *click* on  at the bottom right of the screen and *click* on **OK** in the dialog that results.
3. Enter the **Calculations** dialog with the list of tautomers (close the molecule you built to get the tautomers). Specify calculation of

\* Note, however, that were the energy of an alternative tautomer only 10 kJ/mol higher than that for the normal structure, this would translate into a relative abundance of only about 1% at room temperature. Thus, any alternative tautomers would need to be very close in energy to the lowest-energy tautomer to have noticeable effect.

**Equilibrium Geometry** using the **3-21G Hartree-Fock** model. Submit the job with the name *cytosine tautomers*.

4. After the calculations have completed, bring up the spreadsheet and *click* on the line corresponding to the structure for cytosine that you originally built. *Click* on **Add** at the bottom of the spreadsheet, select **rel.E** and **Boltzmann Distribution** from the available quantities and **kJ/mol** from the **Energy** menu and *click* on **OK**. Are either of the alternatives close in energy to the normal form of cytosine?

### **5 and 6 optional**

To get a better estimate of relative tautomer energies, perform energy calculations using the T1 model.

5. Make a copy of *cytosine tautomers*. Name it *cytosine tautomers T1*. Enter the **Calculations** dialog (**Setup** menu) with this copy. Specify **Energy** from the top menu to the right of **Calculate** and **Thermochemical Recipes** and **T1** from the two bottom menus. Submit the job.
6. When the calculations have completed (several minutes), bring up the spreadsheet and again examine relative tautomer energies and Boltzmann populations. Are either of the alternative structures competitive?
7. Remove all molecules and remaining dialogs from the screen.



# Chapter 6

## Spectra of Organic Molecules

---

In addition to equilibrium geometries and conformations, reaction energies and conformational energy differences and diverse molecular properties, calculations are able to account for molecular spectra. These include, infrared and Raman spectra, NMR spectra and UV/visible spectra. The first two arise from the transitions between ground and excited vibrational states, the third from transitions between nuclear spin states and the fourth from transitions between ground and excited electronic states. Most commonly, experimental spectra are used to provide clues to the structure of an unknown molecule, that is, features in the spectrum are taken as evidence for features in the molecular structure. On the other hand, a calculated spectrum starts with a known structure. A high degree of similarity with a measured spectrum may be taken as evidence that the molecule on which the calculation was made is the same as that for which a spectrum was measured. Lack of similarity suggests that the two molecules are not the same.

The first tutorial in this chapter deals with the infrared spectrum of methyl formate. It details the steps needed to calculate a spectrum and to relate the spectrum to the underlying molecular structure. The second tutorial illustrates the use of the Spartan Infrared Database to identify a molecule based on its infrared spectrum. The remaining tutorials in the chapter deal with NMR spectroscopy. The first details the steps involved in calculating and displaying a proton spectrum and a 2D COSY spectrum for 2-norbornene. The second details the steps involved in calculating and interpreting both an  $^{13}\text{C}$  spectrum and a DEPT spectrum for a more complex organic molecule. The two final tutorials also focus on  $^{13}\text{C}$  NMR. The first of these compares NMR spectra for “similar” flexible and rigid molecules while the second examines the dependence of carbon chemical shifts on stereochemistry.

All of the tutorials make use of the EDF2/6-31G\* density functional model. Structures and (infrared and NMR) spectra for all the molecules considered actually exist in the Spartan Spectra and Properties Database (SSPD). As the calculations involved in the last three tutorials involve significant computer time, you might want to obtain some or all of the data from SSPD. On the other hand, running the calculations will give you a more accurate impression of what is actually needed to calculate spectra. If you do decide to use SSPD, you need to change the selection at the bottom of *Spartan*'s interface from **SMD** to **SSPD**.

## Infrared Spectrum of Methyl Formate

In the harmonic approximation, the frequency at which a diatomic molecule vibrates is proportional to the square root of the ratio of the force constant (the second derivative of the energy with respect to change in bond length) and the reduced mass (the product of the masses of the two atoms divided by their sum). This is in the infrared region of the electromagnetic spectrum. Frequency increases with increasing force constant or “stiffness” of the bond and decreases with increasing masses of the atoms involved in the bond.

$$\text{frequency} \propto \sqrt{\frac{\text{force constant}}{\text{mass}}}$$

This expression may be generalized to a polyatomic molecule and requires change into a coordinate system (so-called *normal coordinates*) that leads to a diagonal matrix of second energy derivatives. The normal coordinates for a polyatomic molecule will typically involve the motions of several (and perhaps all) of the atoms.



The degree (or intensity) of absorption of infrared radiation by a diatomic molecule is proportional to the change in the dipole moment with change in bond length. Note that there is no change in the dipole moment of a homonuclear diatomic molecule with change in bond length, meaning that the intensity of absorption is zero. Molecules such as N<sub>2</sub> and O<sub>2</sub> are transparent in the infrared.

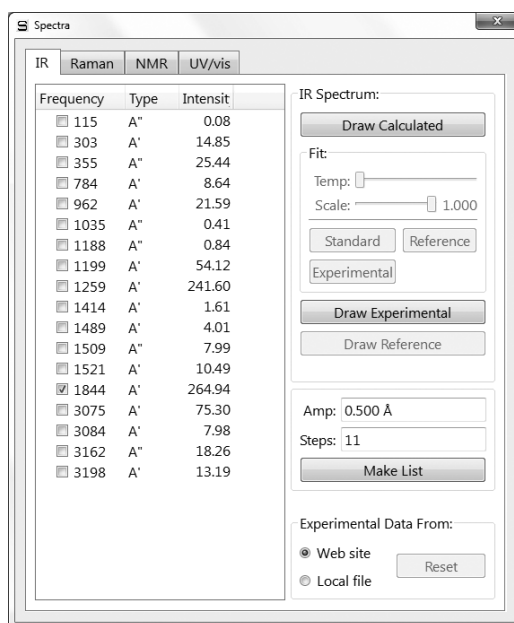
Were this not the case, the earth would be far hotter as radiation of the energy due to the sun would be trapped in the atmosphere (the greenhouse effect). The atmospheric concentration of carbon dioxide, the major atmospheric component that absorbs in the infrared, is < 400 ppm.

The intensity of each of the individual lines in an infrared spectrum of a polyatomic molecule follows from the change in dipole moment along the associated normal coordinate. Note that some of the normal coordinates in a polyatomic molecule may lead to no change in dipole moment, for example, the symmetric stretch in carbon dioxide where both CO bonds are simultaneously moved.



This tutorial illustrates the steps required to calculate and display an infrared spectrum and to compare it with an experimental spectrum. You will see how the ease or difficulty of molecular motion and atomic masses affect frequency.


1. Build methyl formate and *click* on  . Bring up the **Calculations** dialog (**Setup** menu), and request an equilibrium geometry using the EDF2/6-31G\* density functional model. *Check IR* to the right of **Compute** and *click* on **Submit**. Accept the name.
2. After the calculation has completed (several minutes), select **Spectra** from the **Display** menu and *click* on the **IR** tab.




*Click* on both **Draw Calculated** and **Draw Experimental** at the right of the **IR Spectra** dialog. If you are online, select **Web Site** under **Experimental Data from** at the bottom of the dialog before *clicking* on **Draw Calculated**. If not or if you are having trouble downloading the experimental spectrum, select **Local File** at the bottom right of the dialog. Navigate to the *spectra of organic molecules* subdirectory in the *tutorials* directory\* and

\* For Windows, this is found in *Program Files/Wavefunction/Spartan10* or *Program Files(x86)/Wavefunction/Spartan10*. For Linux, this is found in the directory where *Spartan* was installed. For Macintosh, this is located in the *Spartan.10.full* folder on the installation CDROM.

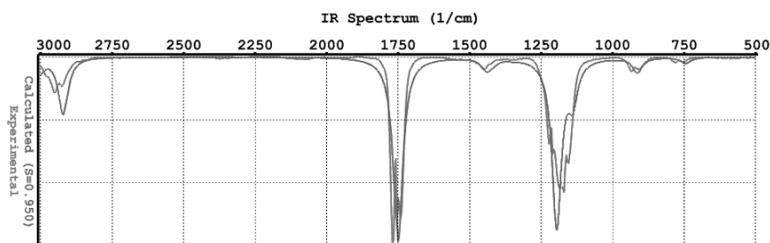*click* on **methyl formate**. You need to select **All Files** from the menu to the right of the file name in the **Open** directory. The calculated frequencies have been scaled to take account of a systematic error due primarily to the harmonic approximation. You can alter the default value (0.962) with the **Scale** slider bar. You can also alter the width of the lines in the calculated spectrum (which can be loosely associated with temperature) with the **Temp** slider bar. To get a best fit to the experimental spectrum (shown below), *click* on the **Experimental** button below the two slider bars. The scale factor has been reduced slightly (to 0.950) and the “temperature” has increased. You can, if you wish, return to the default values by *clicking* on the **Standard** button.



3. Calculated frequencies scaled by the value on the slider bar are listed at the left of the dialog. *Check* the box to the left of the smallest frequency. Describe the motion. Repeat for the largest frequency. Which motion should be more difficult? Why? Finally, *check* the box corresponding to the most intense line, and describe the motion. Why is this line intense?
4. Make a copy of methyl formate (); name it **methyl formate d3**. Select **Properties** from the **Display** menu and *click* on one of the three hydrogen atoms on the methyl group to bring up the **Atom Properties** dialog. Change **Mass Number** from **1** to **2 Deuterium**. Repeat for the other two methyl group hydrogen atoms. Resubmit the calculation (it will require only a few seconds). Compare the frequencies of the undeuterated and deuterated forms of methyl formate, and identify which change the most and which change the least. Explain your results.
5. Close all molecules and dialogs that are open on screen.

## Searching *Spartan*'s Infrared Spectral Database



Pattern matching (“fingerprinting”) a measured infrared spectrum to experimental spectra contained in a library (database) has been common practice for many years. This has not previously been extended to comparisons of measured spectra to databases of calculated spectra. The likely reason is that the results of calculations (a set of vibrational frequencies and intensities to go along with these frequencies) do not look like experimental infrared spectra, at least, they do not look like spectra obtained at normal temperatures. However, a fit of the calculated data to a Lorentzian function in which peak width at half height is treated as a parameter (one value for all peaks) makes the two visually quite similar. One could argue that this parameter loosely corresponds to temperature, which is consistent with the fact that infrared spectra measured at low temperature comprise a series of sharp lines (rather than bands).

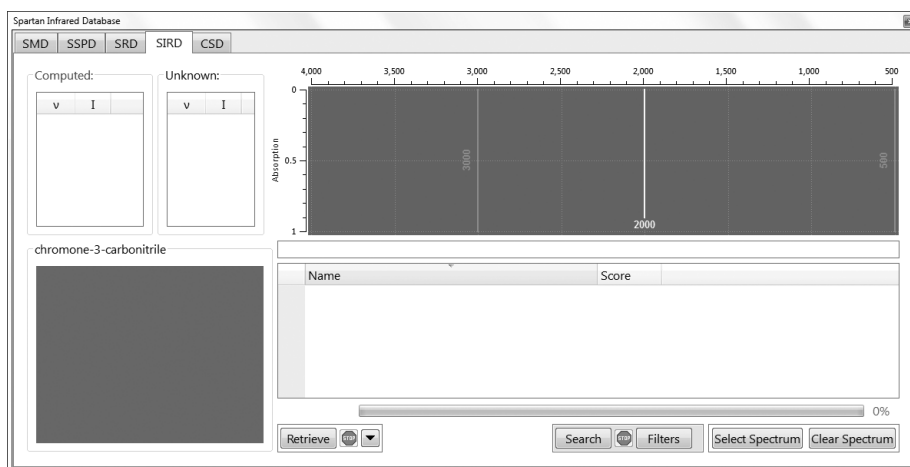
The frequency of an infrared absorption is proportional to (the square root of) the second derivative of the energy at a minimum in the potential surface for a particular vibrational coordinate, that is, the curvature of the surface. Practical frequency calculations replace the real potential at the energy minimum by a quadratic function. This leads to a surface that is too steep (and an energy second derivative that is too large), meaning that calculated frequencies will be larger than measured frequencies. For the most part, the error is systematic, with calculated frequencies being between 3% and 15% larger than measured frequencies, depending on the theoretical model. Hartree-Fock models show an error near the top of the range, while density functional models show an error near the bottom of the range. Semi-empirical models do not show a consistent error pattern. Such a systematic error can be compensated for, at least in part, by incorporating a single linear scaling parameter into the fitting function.

In summary, *Spartan* fits calculated and measured infrared using a Lorentzian function that incorporates two parameters, a non-linear parameter accounting for peak width and a linear scaling parameter.

This tutorial illustrates the way in which searches of the *Spartan*


Infrared Database (SIRD) are setup and performed and the results examined. You may choose from a small selection of measured infrared spectra (or supply your own spectra as a .dx file).

1. Click on  and then select **Databases** from the **Search** menu (or click on ). Click on the **SIRD** tab to bring up the **Spartan Infrared Database** dialog.




Click on **Select Spectrum** at the bottom right of the dialog. Navigate to the *spectra of organic molecules* subdirectory under the *tutorials* directory.\* You need to select **All Files** from the menu to the right of the file name inside the **Open** dialog.

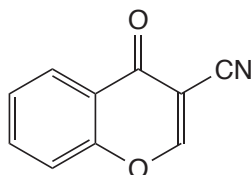
Select one of the following molecules and click on **Open**.



chromone-3-carbonitrile

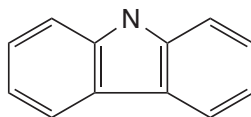

dibenzopyrrole

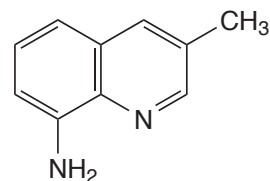


8-amino-2-methylquinoline


The experimental infrared spectrum of the selected molecule will appear in the window at the top right of the dialog with the name underneath. Immediately to the left is a scroll box containing

\* For Windows, this is found in *Program Files/Wavefunction/Spartan10* or *Program Files(x86)/Wavefunction/Spartan10*. For Linux, this is found in the directory where *Spartan* was installed. For Macintosh, this is located in the *Spartan.10.full* folder on the installation CDROM.

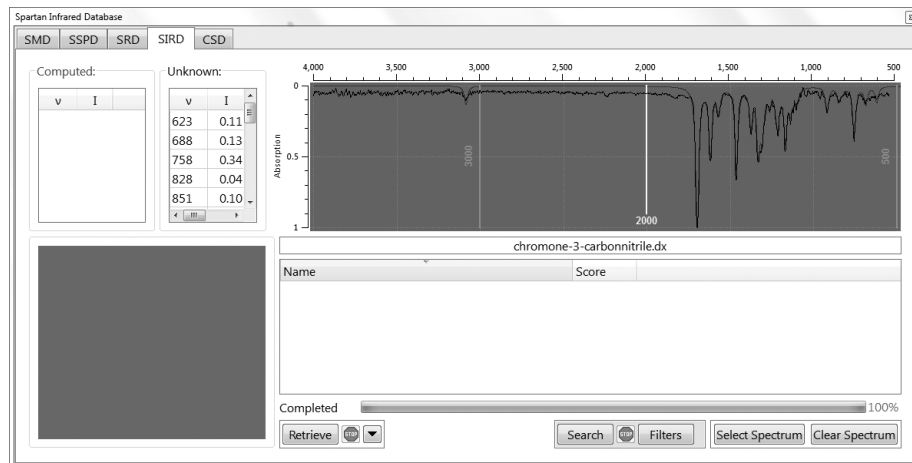the most intense lines in the spectrum (obtained from a fit of the experimental spectrum to a Lorentzian function).



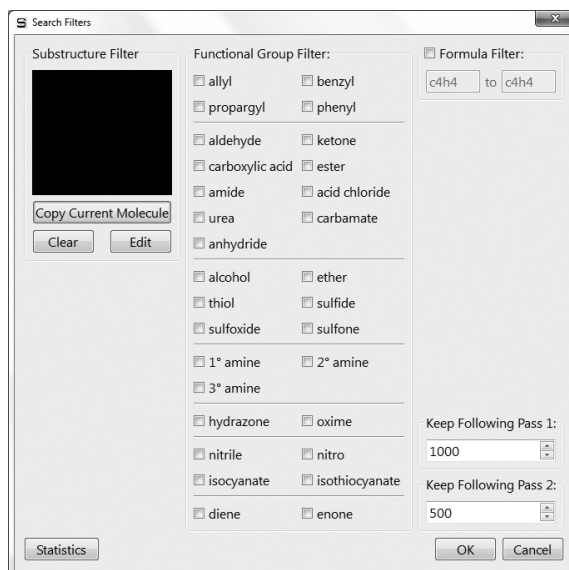2. Click on **Filters** at the bottom of the dialog to bring up the **Search Filters** dialog.



This provides for substructure and formula filters (neither of which is used in this tutorial) as well as functional group filters. You can if you wish *check* an appropriate entry, for example, *nitrile*, if you have selected chromone-3-carbonitrile. Click on **OK** to exit the dialog.

3. *Click* on **Search** at the bottom of the (**Spartan Infrared Database**) dialog. The search may require a minute or two. When it has completed, scan the list of hits at the lower right of the dialog for the name of your query. It will be at or near the top of the list.\* *Click* on the name. The calculated infrared spectrum will be superimposed onto the experimental spectrum (in the window at the top right of the dialog).
4. Close all molecules and dialogs that are open on screen.

---

\* This will not always be the case. The obvious exception occurs where the query is not in the database. Problems will also occur when the (query) infrared spectrum lacks sufficient structure to distinguish it from similar molecules.

## Proton NMR Spectrum of 2-Norbornene

Proton NMR spectroscopy was the first tool available to chemists that allowed definitive assignment of the molecular structures of complex organic molecules. By the 1970's, it had largely replaced infrared spectroscopy and to a large extent chemical proofs of structure.  $^{13}\text{C}$  NMR is now more dominant, but proton NMR remains an essential tool in the chemist's arsenal.

NMR is based on the fact that nuclei possess spins that can either align parallel or antiparallel to an applied magnetic field, giving rise to different nuclear spin states. The relative energy of these states ( $\Delta E$ ) depends on the nucleus and on the strength of the applied magnetic field, by way of a simple relationship.

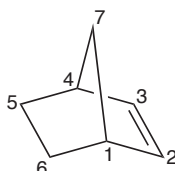
$$\Delta E = \gamma \hbar B_0$$

$\gamma$  is the gyromagnetic ratio (a constant for a given type of nucleus),  $\hbar$  is Planck's constant divided by  $2\pi$  and  $B_0$  is the strength of the magnetic field *at the nucleus*. While the two nuclear spin states are normally in equilibrium, this equilibrium can be upset by applying a second magnetic field. The absorption of energy as a function of field strength (a *resonance*) between the states can then be detected.

The key to the utility of the magnetic resonance experiment is that the energy at which a nucleus resonates depends on its location in the molecule, and is different for each (chemically) distinct nucleus. The reason is that the applied magnetic field is weakened by electrons around the nucleus. Nuclei that are well shielded by the electron cloud will feel a lesser magnetic field than those that are poorly shielded, and will show a smaller energy splitting. The difference, given relative to a standard, is termed a *chemical shift*. By convention, both proton and  $^{13}\text{C}$  chemical shifts (treated later in this chapter) are reported relative to tetramethylsilane (TMS) as a standard.






While each *unique* proton in a molecule gives rise to a single line (resonance) in the spectrum, the spins on nearby nuclei add and subtract to the external magnetic field. This leads to a splitting of lines, the splitting pattern depending on the number of neighboring protons and their geometry. Discounting splitting, the intensity of

the lines is approximately proportional to the number of equivalent protons that contribute. For example, the proton NMR spectrum of 2-norbornene shows six lines, two with unit intensity corresponding to the two different protons on the methylene bridge (C<sub>7</sub>), and four with twice the intensity corresponding to protons at C<sub>1</sub> (C<sub>4</sub>), C<sub>2</sub> (C<sub>3</sub>) and the two different (*endo* and *exo*) protons on C<sub>5</sub> (C<sub>6</sub>).



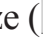In this tutorial, you will calculate the proton NMR spectrum of 2-norbornene and compare it with that contained in the Spectral Database of Organic Compounds (SDBS), freely-accessible on the web from the National Institute of Advanced Industrial Science and Technology (AIST) in Japan.\*

You will also make a COSY plot which identifies pairs of protons that are separated by three bonds. This requires HH coupling constants that are not calculated by *Spartan* but evaluated empirically based on the three-dimensional geometry of the molecule.

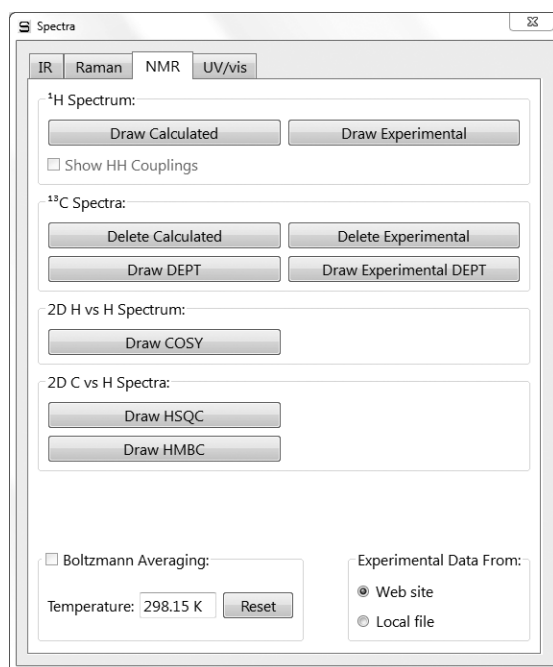
1. Build 2-norbornene. Start from cyclohexane (**Rings** menu), add an sp<sup>3</sup> carbon () to an *axial* ring position and bond to () to a free valence on the opposite side of the ring. Make a double bond between the appropriate ring carbons (). Minimize () and *click* on (). The name *2-norbornene* will appear at the bottom right of the screen.
2. Enter the **Calculations** dialog (**Calculations** from the **Setup** menu), and specify an **Equilibrium Geometry** calculation using the **EDF2 Density Functional** model with the **6-31G\*** basis set. *Check* **NMR** to the right of **Compute**. *Click* on **Submit**. Accept the name *2-norbornene*.
3. The calculation will require several minutes. When it has completed, proton shifts can be found in **Atom Properties**

\* Web address: [http://riodb1.ibase.aist.go.jp/sbds/cgi-bin/cre\\_index.cgi](http://riodb1.ibase.aist.go.jp/sbds/cgi-bin/cre_index.cgi). Unfortunately, electronic access to SDBS is prohibited by AIST.




dialog. Select **Properties** from the **Display** menu and one after another *click* on the hydrogens. Propose assignments for the six resonances observed in the experimental spectrum: 5.99, 2.84, 1.61, 1.31, 1.08 and 0.95 ppm.


4. Select **Spectra** from the **Display** menu and click on the **NMR** tab. The **NMR Spectra** dialog appears.



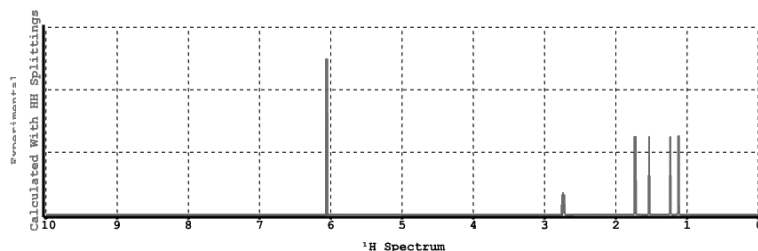
Click on **Draw Calculated** and *check* the box to the left of **Show HH Couplings** (under  $^1\text{H}$  Spectrum) at the top left of the dialog. The following proton spectrum appears.



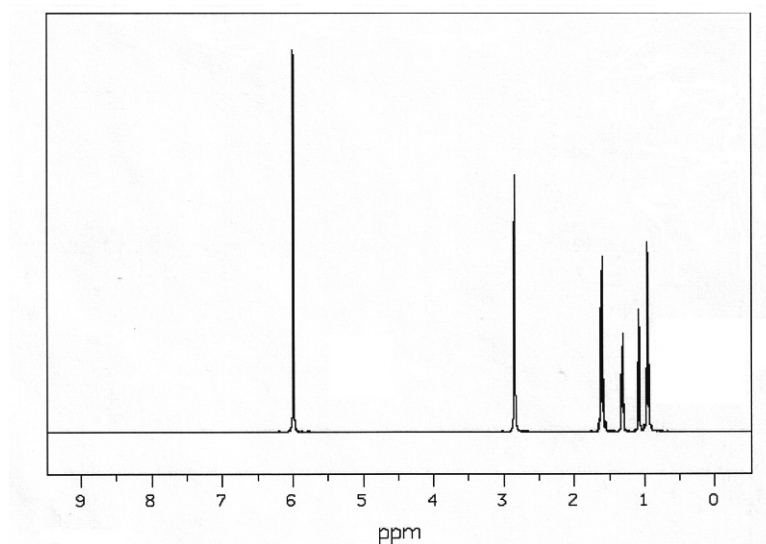
Compare the calculated proton spectrum with that in SDBS (see below or examine it on-line).



5. You may have a difficult time seeing that the lines are split (in both calculated and observed spectra). This is because the splittings are very small relative to differences in the magnitudes of the chemical shifts themselves. *Spartan* provides a magnifier tool to assist. With the **Properties** dialog on screen, click on the horizontal axis of the spectra plot. *Click* on the **Magnifier** button at the bottom right of the **Plot Properties** dialog that results. *Click* on the border of the green magnifier box that results to select it. Move it (as you would any graphical object) to a position directly above the resonance around 6 ppm. The single line (due to protons on the double bond) will be split in two due to interaction (coupling) with the magnetic spins of the protons on the bridgehead carbons.

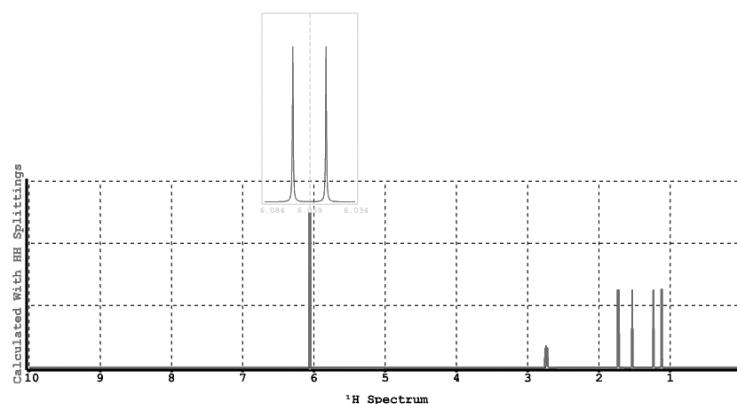

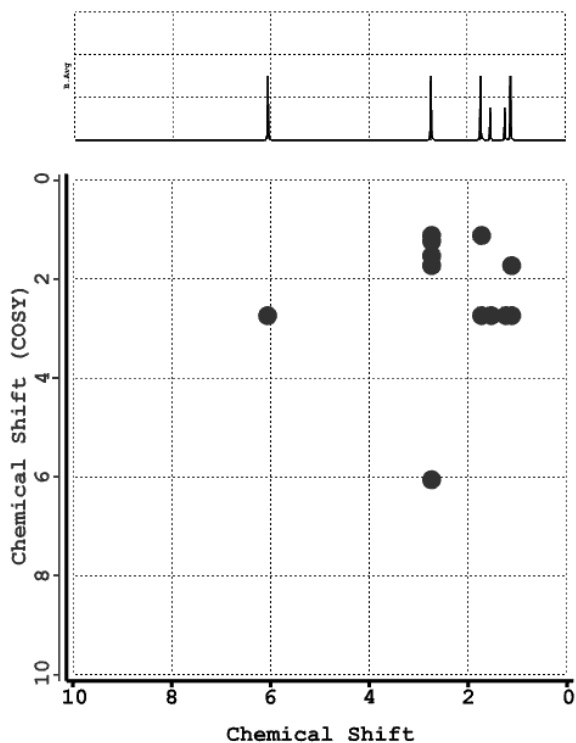6. Draw the COSY spectrum for 2-norbornene (in addition to the proton spectrum). *Click* on **Draw COSY** under **2D H vs. H Spectra**.




The circles indicate interactions between neighboring protons. Note that coupling constants are not obtained from quantum chemical calculations, but rather are based on geometrical relationship between protons in the calculated structure.


7. Close *2-norbornene* and any remaining dialogs when you are finished.

## $^{13}\text{C}$ Spectrum of Caulophylline

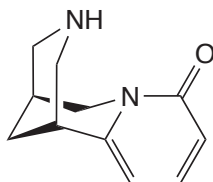There are several reasons why NMR spectroscopy, in particular,  $^{13}\text{C}$  NMR, is one of the most important routine analytical techniques available for characterizing organic molecules. The experiment, which is straightforward and can be carried out rapidly, requires relatively small samples and is non-destructive. The resulting (proton decoupled) spectrum is quite simple, comprising but a single line for each and every unique carbon. However, assigning  $^{13}\text{C}$  spectra is by no means trivial, even for molecules that might appear to be quite “simple”. The problem is that the positions of the lines in the spectrum (the chemical shifts) are very sensitive to the environment in which the carbons find themselves. Where two or more carbons in a molecule reside in what appears to be “similar” environments, it may be very difficult to distinguish them.

This tutorial uses the alkaloid caulophylline (cytisine) to illustrate the use of calculated  $^{13}\text{C}$  spectra to assist in assigning the measured spectrum of the molecule.



1. Build caulophylline and minimize () . Bring up the **Calculations** dialog (**Setup** menu) and specify calculation of equilibrium geometry with the EDF2/6-31G\* density functional model. *Check NMR* to the right of **Compute**. *Click* on **Submit** and specify a name. The job will require upwards of one hour to complete\*.
2. When the calculation is done (or after you have accessed information from SSPD), bring up the **NMR Spectra** dialog (**Spectra** under the **Display** menu) and *click* on the **NMR** tab. Make certain that the box to the left of **Boltzmann Averaging**

---

\* An optimized structure and NMR spectrum for caulophylline are available in the Spartan Spectra and Properties Database (SSPD). If you decide to use this instead of doing the calculations, select SSPD from the menu at the bottom of *Spartan*'s main screen, and *click* on **Replace** inside the **Preview** dialog.

at the bottom of the dialog is not checked (if it is, *click* inside the box to remove the checkmark). *Click* on **Draw Calculated** in the section designated  $^{13}\text{C}$  **Spectra**. Next, *check* **Local File** in the section designated **Experimental Data From**, *click* on **Draw Experimental** (in the  $^{13}\text{C}$  **Spectra** section), move to the *spectra of organic molecules* subdirectory under the *tutorials* directory\* and *click* on *caulophylline carbon NMR*. You need to select **All Files** from the menu to the right of the file name inside the **Open** dialog. The reason that two of the resonances are missing from both calculated and experimental spectra is because they are outside the default scale (0 to 150 ppm). To change the scale, select **Properties** from the **Display** menu and *click* on the horizontal axis of the plot, to bring up the **Plot Properties** dialog. *Click* on the **X-Axis** tab at the top of the dialog. Change “150” at the left below **Range** to “175”, increase the number of horizontal increments (**Tics**) from “6” to “8” (top right of dialog) and close the dialog. The plot now shows all ten resonances.


3. Use the calculated spectrum to associate the individual lines in the experimental  $^{13}\text{C}$  spectrum with specific carbons in the structure of caulophylline. Predict the DEPT spectrum, where CH groups point up,  $\text{CH}_2$  groups point down and quaternary carbons are missing. Leave the  $^{13}\text{C}$  plots on screen and *click* on **DEPT Draw** and on **Draw Experimental DEPT** (in the  $^{13}\text{C}$  section of the **NMR Spectra** dialog). The experimental spectrum is *caulophylline DEPT NMR* in the *spectra of organic molecules* subdirectory under the *tutorials* directory. Are your peak assignments consistent with the experimental DEPT spectrum?
4. Close all molecules and dialogs that are open on screen.

---

\* For Windows, this is found in *Program Files/Wavefunction/Spartan10* or *Program Files(x86)/Wavefunction/Spartan10*. For Linux, this is found in the directory where *Spartan* was installed. For Macintosh, this is located in the *Spartan.10.full* folder on the installation CDROM.

## $^{13}\text{C}$ Spectrum of *cis*-1,2-Dimethylcyclohexane

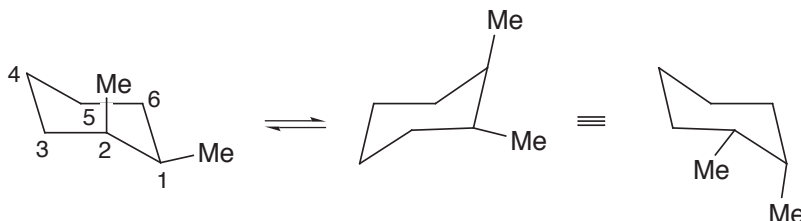At normal temperatures, the NMR spectrum of a conformationally-flexible molecule represents a (Boltzmann-weighted) average of the NMR spectra of all accessible conformers. Only when the temperature is lowered will the spectrum reveal its components, and its complexity. A simple example is *cis*-1,2-dimethylcyclohexane, a molecule with two equivalent conformers.



The room temperature  $^{13}\text{C}$  NMR spectrum contains only four lines at 34.9, 31.9, 24.2 and 16.4 ppm relative to TMS, corresponding to an equal weighting of  $\text{C}_1$  and  $\text{C}_2$ ,  $\text{C}_3$  and  $\text{C}_6$ , and  $\text{C}_4$  and  $\text{C}_5$  and the two methyl carbons, respectively. Only when the sample is cooled does the spectrum reveal eight distinct lines.

In this tutorial, you will calculate the  $^{13}\text{C}$  NMR spectrum of *cis*-1,2-dimethylcyclohexane, and then average chemical shifts for equivalent carbons to yield a “room temperature” spectrum. You will then obtain a  $^{13}\text{C}$  spectrum for *cis*-decalin to see if it serves as a model for static *cis*-1,2-dimethylcyclohexane.


1. Build and minimize *cis*-1,2-dimethylcyclohexane. Inside the **Calculations** dialog, specify calculation of **Equilibrium Geometry** using the **Density Functional EDF2 6-31G\*** model. *Check NMR* to the right of **Compute**, and submit the job. Accept the name *cis-1,2-dimethylcyclohexane*.
2. After the calculation has completed (several minutes), record the  $^{13}\text{C}$  chemical shifts\* for all ring carbons in *cis*-1,2-dimethylcyclohexane, either from the **Atom Properties** dialog (select **Properties** from the **Display** menu and *click* on a

\*  $^{13}\text{C}$  chemical shifts from B3LYP/6-31G\* and EDF2/6-31G\* calculations (only) have been corrected to account for local environment. The two correction schemes are different with that for B3LYP/6-31G\* employing fewer empirical parameters.

carbon) or by labeling the model (select **Configure...** from the **Model** menu and *check Chem Shift* under **Atom** in the dialog that results). In the latter case, you may want to remove the hydrogens (*click* to deselect **Hydrogens** from the **Model** menu) to simplify the display. Finally, average the calculated chemical shifts in order to mimic what is seen at room temperature. Are the resulting  $^{13}\text{C}$  shifts consistent with what is actually observed? Elaborate.

3. *cis*-Decalin is a plausible model for static *cis*-1,2-dimethylcyclohexane, in that it is unable to undergo (rapid) ring inversion.









Build *cis*-decalin. Start with cyclohexane, add  $\text{sp}^3$  carbons at adjacent *cis* free valences and connect these carbons by two additional  $\text{sp}^3$  groups. Instead of calculating its equilibrium geometry and NMR spectrum, get them from the Spartan Spectra and Properties Database (SSPD). *Click* on  to the right of the name (*cis-decalin*) at the bottom of the screen to bring up the **Preview** dialog. *Click* on the **Replace** button at the bottom of the dialog. The entry in SSPD will replace the molecule you built.

4. Select **Spectra** from the **Display** menu and *click* on the **NMR** tab. *Click* on **Draw Calculated** under  $^{13}\text{C}$  **Spectrum**. Because of the ( $\text{C}_2$ ) symmetry of *cis*-decalin, your spectrum will contain only five lines, all of equal intensity. Relate these to the six lines corresponding to the ring carbons in static *cis*-1,2-dimethylcyclohexane. Pay particular attention to the differences in chemical shifts between carbons which average in the equilibrating molecule. Is *cis*-decalin a good model for *cis*-1,2-dimethylcyclohexane?
5. Close all molecules and any open dialogs.

## Stereochemical Assignments from $^{13}\text{C}$ Spectra

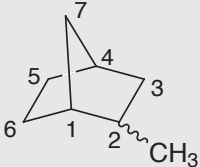
NMR spectroscopy, in particular  $^{13}\text{C}$  spectroscopy, is without doubt the method of choice to establish the three-dimensional structure of organic molecules. Only X-ray diffraction provides more definitive results, although the requirement of a crystalline sample severely limits its application. It is now practical to routinely calculate the NMR spectra of organic molecules. In time, the availability of a “virtual NMR spectrometer” will offer organic chemists an entirely new paradigm for structure determination, that is direct comparison of a measured spectrum with calculated spectra for one or more chemically reasonable candidates.

Here, you will obtain  $^{13}\text{C}$  chemical shifts for *endo* and *exo* stereoisomers of 2-methylnorbornane, and compare these to experimental shifts. The EDF2/6-31G\* density functional model will be employed. You will establish the extent to which the calculations are able to reproduce differences in chemical shifts as a function of stereochemistry.

1. Build norbornane. Start with **Cyclohexane** from the **Rings** menu. Add  $\text{sp}^3$  carbon  to one of the *axial* free valences (leaving you with *axial*-3-methylcyclohexane). Click on  and then click on one of the free valences of the methyl group and the *equatorial* free valence on cyclohexane that is directly across. Click on . The name norbornane should appear at the bottom right of the screen. Put norbornane onto the clipboard. *Right click* and select **Copy** from the menu that appears.
2. Select  $\text{sp}^3$  carbon  from the model kit and click on the *exo* free valence of the carbon in the 2-position. The name ***exo*-2-methylnorbornane** should appear at the bottom right of the screen.
3. Select **New Molecule** from the **File** menu. Click on **Clipboard** at the bottom of the model kit and click on screen. Norbornane appears. Select  $\text{sp}^3$  carbon  and click on the *endo* free valence of the carbon in the 2-position. The name ***endo*-2-methylnorbornane** should appear at the bottom right of the screen. Click on .



- Name your molecules. Bring up the spreadsheet (**Spreadsheet** under the **Display** menu), *right click* on the header cell of the leftmost column and select **Rename Selected Using SSPD** from the menu that results.
- Inside the **Calculations** dialog, specify calculation of **Equilibrium Geometry** using the **EDF2 6-31G\* Density Functional** model. *Check NMR* (to the right of **Compute**). Submit the job with the name **2-methylnorbornane**.
- When completed (several minutes), compare the  $^{13}\text{C}$  chemical shifts with experimental values with particular attention to differences between *endo* and *exo* stereoisomers.

	<i>endo</i>	<i>exo</i>	$\Delta$
			
C <sub>1</sub>	43.5	42.2	-1.3
C <sub>2</sub>	36.8	34.6	-2.2
C <sub>3</sub>	40.2	40.7	0.5
C <sub>4</sub>	37.3	38.2	0.9
C <sub>5</sub>	30.3	30.6	0.3
C <sub>6</sub>	29.0	22.4	-6.6
C <sub>7</sub>	35.0	38.9	3.9
CH <sub>3</sub>	22.3	17.4	-4.9

Can the calculations be used to distinguish between the two isomers?

- Close **2-methylnorbornane**.

## <sup>13</sup>C Chemical Shifts Depend on Conformation

The time for nuclear spin relaxation is very long relative to time required for equilibration among conformers, at least at normal temperatures. This means that the each of the carbon chemical shifts in the NMR spectrum of a flexible molecule, <sup>13</sup>C, will be a weighted average of the shifts NMR of the individual conformers, <sup>13</sup>C<sub>i</sub>.

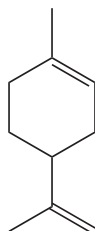
$$^{13}\text{C} = \sum_i \omega_i ^{13}\text{C}_i$$

The weight,  $\omega_i$ , is given by the Boltzmann equation, and depends on its energy,  $\epsilon_i$  (relative to that of the lowest-energy conformer,  $\epsilon_0$ ) and on temperature,  $T$ . Summation is over all conformers (including the lowest-energy conformer),  $g_i$  is the number of times that conformer  $i$  appears in the overall distribution and  $k$  is the Boltzmann constant.

$$\omega_i = g_i \exp - [\epsilon_i - \epsilon_0/kT] / \sum_j \{ \exp - [\epsilon_j - \epsilon_0/kT] \}$$

In practice, conformers that are 10 kJ/mol or more above the lowest-energy conformer make a negligible contribution to the total (at room temperature).


1. Open *limonene* from the *spectra of organic molecules* subdirectory under the *tutorials* directory. This molecule with one freely-rotatable bond has three distinct conformers.



2. Select **Spreadsheet** under the **Display** menu. *Click* on **Add...** at the bottom of the spreadsheet, select both *rel. E* and **Boltzmann Distribution** from the menu and *click* on **OK**. You will see that all three conformers are within 2 kJ/mol, meaning that all should contribute significantly to the NMR spectrum.
3. Select **Spectra** under the **Display** menu and *click* on the **NMR** tab. *Click* on both **Draw Calculated** and **Draw Experimental** under

<sup>13</sup>C Spectra (the experimental spectrum has been attached to the .spartan file). Move among the three conformers and examine the changes to the calculated <sup>13</sup>C spectrum and the “quality of fit” to the experimental spectrum. Do any of the chemical shifts differ by more than a few ppm from one conformer to another?

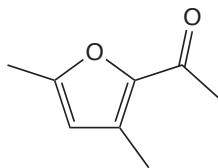4. Click on **Boltzmann Averaging** at the bottom of the **NMR Spectra** dialog. The calculated spectrum will now be based on the Boltzmann average (it will not change as you move from one conformer to another). Close *limonene* when you are done.
5. Repeat the above steps for *3,5-dimethyl-2-furyl methyl ketone* (also available in the *spectra of organic molecules* subdirectory under the *tutorials* directory).



The two conformers of this molecule are different in energy by 17 kJ/mol and the minor conformer contributes very little to the overall average NMR spectrum. Does the large difference in (conformer) energy anticipate a large difference in <sup>13</sup>C shifts (at least for those carbons that are directly involved in the change)?

6. Close *3,5-dimethyl-2-furyl methyl ketone* and any open dialogs.

# Chapter 7

## Organic Reactions

---

*This chapter outlines and illustrates strategies for locating and verifying transition states for organic reactions as well as exploring changes in product distributions as a function of substituents and reactant stereochemistry.*

The treatment of chemical reactions adds an entirely new dimension to the application of quantum chemical models. Whereas unique valence structures may generally be written for most molecules and, based on these structures, reasonable guesses at bond lengths and angles may be made, it is often difficult to designate appropriate valence structures for transition states, let alone specify detailed geometries. While there is a complete absence of experimental data for the structures of transition states, calculated transition-state geometries are now commonplace. **Spartan** provides both an extensive library of calculated transition-state geometries and a facility for automatically matching an entry in this library with the reaction at hand.\* This library is also available as the Spartan Reaction Database (SRD), and may be searched by substructure to yield all available transition states of reactions related to the one of interest.

**Spartan** also provides a procedure for driving user-defined coordinates. Aside from conformational analysis (see discussion in previous chapter), the major application of this is to force reactions, thereby permitting identification of transition states.

The first four tutorials in this chapter illustrate **Spartan's** automatic procedure for guessing transition-state geometries based on its library of organic reactions. The last two of these show how a transition state for one reaction may be used to provide a good guess at the transition state of a closely-related reaction. These tutorials also illustrate the

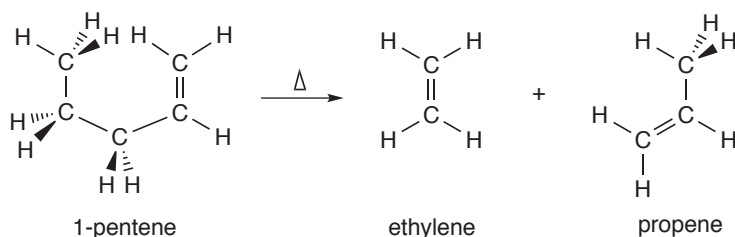
---

\* Where a reaction is unknown to **Spartan's** library, a fallback technique which averages reactant and product geometries (similar to the so-called linear synchronous transit method) is invoked.




use of vibrational analysis (infrared spectroscopy) to verify that a particular structure corresponds to a transition state and to show the motion connecting it to reactants and products. The fifth tutorial illustrates how a reaction may be “driven” through a transition state. This tutorial, along with the fourth tutorial, draw the connection between relative activation energies and kinetic product distributions. The last tutorial in this chapter illustrates how transition states may be extracted from the Spartan Reaction Database.


An example of a transition state calculation for an organometallic reaction is provided in **Chapter 10**.

## Ene Reaction of 1-Pentene





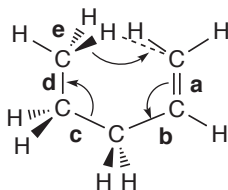
The proposed mechanism of the ene reaction involves simultaneous transfer of a hydrogen atom and CC bond cleavage. In this tutorial, you will obtain the transition state for the ene reaction of 1-pentene from 3-21G Hartree Fock calculations, and examine the reaction coordinate for evidence of concerted motion. Optionally, you will perform B3LYP/6-31G\* density functional energy calculations to estimate the activation energy for the reaction.

1. Bring up the organic model kit and build 1-pentene in a conformation in which one of the terminal hydrogens on the ethyl group is poised to transfer to the terminal methylene group. To rotate about a (single) bond, first *click* on it to select (it will be marked by a yellow cylinder), hold down the left button and **Alt** key (**option** key on Mac) and move the mouse up and down. Minimize and *click* on .
2. First, save 1-pentene as *1-pentene density functional* for optional use later (). Also make a copy for immediate use; name it *ene reaction 1-pentene*. In both cases you will need to replace the suggested name (*1-pentene*).
3. Select **Guess Transition State** from the **Search** menu (or *click* on the  icon at the top of the screen). *Click* on bond **a** in the figure on the following page and then *click* on bond **b**. A curved arrow from double bond **a** to single bond **b** will be drawn.



Next, *click* on bond **c** and then on bond **d**. A second curved arrow from bonds **c** to **d** will be drawn. Finally, *click* on bond **e** and then, while holding down the **Shift** key, *click* on the (methyl) hydrogen to be transferred and on the terminal (methylene) carbon to receive this hydrogen. A third curved arrow from bond **e** to the center of a dotted line that has been drawn between the hydrogen and oxygen, will appear.

If you make a mistake, you can remove an arrow by selecting **Delete** from the **Build** menu (*click* on ) and then *clicking* on the arrow. (You will need to select  to continue.) Alternatively, hold down the **Delete** key as you *click* on an arrow. With all three arrows in place, *click* on  at the bottom right of the screen. Your structure will be replaced by a guess at the ene transition state. If the resulting structure is unreasonable, then you have probably made an error in the placement of the arrows. In this case, select **Undo** from the **Edit** menu to return to the model with the arrows and modify accordingly.

4. Enter the **Calculations** dialog (**Setup** menu), and specify calculation of transition-state geometry using the 3-21G Hartree-Fock model. Select **Transition State Geometry** from the top menu to the right of **Calculate**, and choose **Hartree-Fock** and **3-21G** from the two bottom menus. Finally, *check* **IR** to the right of **Compute**. This will allow you to confirm that you have found a transition state, and that it smoothly connects reactant and product. *Click* on **Submit**.
5. When the job completes, animate the motion of atoms along the reaction coordinate. Select **Spectra** from the **Display** menu and *click* on the **IR** tab. *Click* on the top entry in the list in the **IR** dialog that results. It corresponds to an imaginary frequency, and will be designated with an **i** in front of the number.

A vibrational frequency is proportional to the square root of the ratio of the force constant (reflecting the curvature of the potential surface along a particular coordinate) divided by a combination of the masses of atoms involved in motion along that coordinate. At a transition state (a maximum in the reaction coordinate), the curvature is negative. Since the mass term is positive, the quantity inside the square root is negative and the frequency is an imaginary number.




Is the vibrational motion consistent with an ene reaction of interest and not with some other process?

6. Controls at the bottom right of the **IR** dialog allow for changing both the amplitude of vibration (**Amp**) and the number of steps that make up the motion (**Steps**). The latter serves as a speed control. Change the amplitude to 0.3. *Type 0.3* in the box to the right of **Amp** and *press* the **Enter** key (**return** key on Mac). Next, *click* on **Make List** at the bottom of the dialog. This will give rise to a group of structures that follow the reaction coordinate down from the transition state both toward reactant and product. Remove the original transition state from the screen. *Click* on **ene reaction 1-pentene** (the vibrating molecule) and close it, along with the **IR** dialog.
7. Enter the **Calculations** dialog and specify calculation of **Energy** using the **3-21G Hartree-Fock** model (the same level of calculation used to obtain the transition state and calculate the frequencies). Make certain that **Global Calculations** is checked. Next, enter the **Surfaces** dialog and specify evaluation of two surfaces: a bond density surface and a bond density surface onto which the electrostatic potential has been mapped. Before you do so, make certain that **Global Surfaces** is checked. *Click* on **More Surfaces...\***, select **density (bond)** for **Surface** and **none** for **Property** and *click* on **Apply**. Select **density (bond)** for **surface** and **potential** for **Property** and *click* on **OK**.

---




\* You need to use **More Surfaces** rather than **Add** as the electrostatic potential is mapped on the bond density surface.




- Submit for calculation. Name it *ene reaction 1-pentene sequence*. Once the job has completed, enter the **Surfaces** dialog and examine the surfaces that you have calculated. For each, step through the sequence of structures ( and ) keys at the bottom of the screen) or animate the reaction (). Note, in particular, the changes in bonding revealed by the bond density surface. Also pay attention to the value of the potential on the migrating atom. This reflects its charge. Is it best described as a proton (blue), hydrogen atom (green) or hydride anion (red)?
- Close *ene reaction 1-pentene sequence* as well as any remaining dialogs.





### 10 to 13 optional

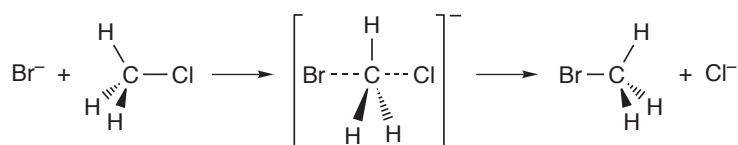
Methods that account for electron correlation are generally needed to furnish accurate estimates of absolute activation energies. Perform B3LYP/6-31G\* density functional energy calculations on both 1-pentene and on the ene reaction transition state (using 3-21G geometries).


- Open *ene reaction 1-pentene* () and make a copy (). Name it *ene reaction 1-pentene density functional*.
- Enter the **Calculations** dialog, and specify calculation of **Energy** using the **B3LYP/6-31G\* Density Functional** model. Remove the checkmark from **IR**. Submit the job.
- Open *1-pentene density functional* (). Specify calculation of **Energy** using the **B3LYP/6-31G\* Density Functional** model. Also select **3-21G** from the **Start from** menu. This designates that a 3-21G equilibrium structure is to be used (that is, the same as employed for the transition state geometry). Submit the job.
- Obtain the activation energy (difference in total energies between 1-pentene and the ene reaction transition state).
- Remove any molecules and dialogs from the screen.





## S<sub>N</sub>2 Reaction of Bromide and Methyl Chloride

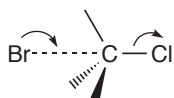


1. The S<sub>N</sub>2 reaction passes through a transition state in which carbon is in a trigonal bipyramid geometry and the entering and leaving groups are colinear. To build it, first construct methyl chloride. Then select bromine from the palette of icons in the model kit, *click* on the **Insert** button at the bottom right of the screen or hold down the **Insert** key (**option** key on Mac) and *click* anywhere on screen. Alternatively, *double click* in a blank area of the screen following selection of bromine. Two detached fragments, methyl chloride and hydrogen bromide, appear on screen. *Click* on  and then *click* on the free valence on bromine. *Click* on . Alternatively, hold down the **Delete** key and *click* on the free valence on bromine. You are left with methyl chloride and bromine atom (bromide). Manipulate the two such that bromide is poised to attack methyl chloride from the backside (as in the transition state above). (Recall that translations and rotations normally refer to both fragments, but can be made to refer to a single fragment by first *clicking* on the fragment and then holding down on the **Ctrl** key while carrying out the manipulations.) **Do not minimize** (If you do so by accident, select **Undo** from the **Edit** menu). *Click* on .
2. Select **Guess Transition State** from the **Search** menu or *click* on . *Click* on bromide and, while holding down the **Shift** key, *click* again on bromide and then on carbon. A dotted line will be drawn from bromine to carbon, together with an arrow from bromine to the center of this line. Next, *click* on the CCl bond and then *click* on the chlorine. A second arrow from the carbon-chlorine bond to the chlorine will be drawn.






Click on  at the bottom right of the screen. Your structure will be replaced by a guess at the transition state.

- Click on  and then on the CBr bond. Replace the current CBr distance in the box at the bottom right of the screen by **3.8** (3.8Å) and *press* the **Enter** key (**return** key on Mac). You have now made a complex representing the reactant.
- Select **Constrain Distance** from the **Geometry** menu (or *click* on the  icon at the top of the screen). *Click* on the CBr bond, and then *click* on  at the bottom right of the screen. The icon will change to  indicating a constraint is to be applied to this distance. *Check* the box to the left of **Profile** at the bottom of right of the screen. This will result in two additional text boxes.

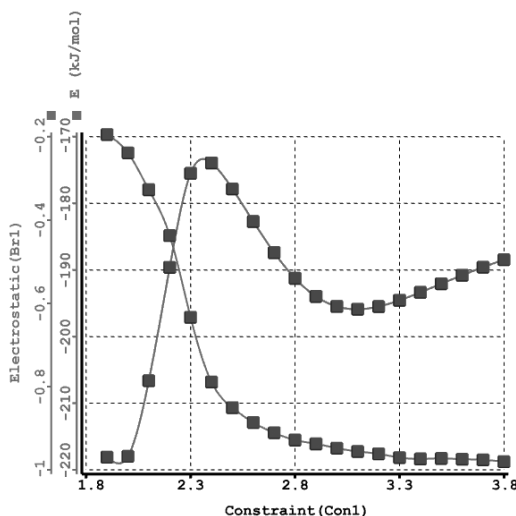


Leave the value **3.8** (3.8Å) in the original (leftmost) box alone, but change the number in the box to the right of **to** to **1.9** (1.9Å) and *press* the **Enter** (**return**) key. Change the number in the box to the right of **Steps** from **10** (the default) to **20**. 20 Calculations with CBr bond lengths constrained from 3.8Å (the starting point) to 1.9Å (the ending point) will be performed. The transition state should have a CBr distance in between these values. *Click* on .

- Enter the **Calculations** dialog, and select **Energy Profile**, **Semi-Empirical** and **AM1** from the appropriate menus to the right of **Calculate**. You need to change **Total Charge** to **Anion**.
- Submit the job. Name it *bromide+methyl chloride*. When completed, it will give rise to a sequence of calculations placed in *bromide+methyl chloride.Prof.M0001*. You will be prompted as to whether you want to open this file. *Click* on **OK**. Align the molecules. *Click* on , select **Structure** from the **Align by** menu at the bottom right of the screen and, one after the other, *click* on the three hydrogens. Finally *press* the **Align by** button at the bottom right of the screen. *Click* on .
- Bring up the **Spreadsheet** (**Spreadsheet** from the **Display** menu)

and *click* on **Add...** Select **E** from among the quantities listed at the top of the dialog, **kJ/mol** from the **Energy** menu, and *click* on **OK**. Next, enter the (constrained) CBr distances and bromine charges in the spreadsheet. *Click* on , *click* on the constraint marker and *click* on **P** at the bottom right of the screen. *Click* on **V**. Bring up the **Properties** dialog. *Click* on bromine and *click* on **P** to the left of **Electrostatic** under **Charges** in the **Atom Properties** dialog. Finally, bring up the **Plots** dialog, and select **Constraint (Con1)** (the distance at which the CBr bond has been constrained) from the **X Axis** menu, and both **E (kJ/mol)** and **Electrostatic (Br1)** from the **Y Axes** list. *Click* on **Create Plot**.

One plot gives the energy as the reaction proceeds and the other gives charge on bromine. Are the two related? Explain.






$S_N2$  reactions involving charged species normally need to be carried out in highly-polar media, for example, water. *Spartan* provides two solvation models in which the solvent is represented by a reaction field (**Appendix A**). The SM5.4 model only alters the energy and not the charges. The more recent SM8 model does affect charges but is available only for Hartree-Fock and density functional models.


8. Add aqueous phase data (based on the SM5.4 model) to the spreadsheet. *Click* on an empty column header, *click* on **Add...**, select **Eaq** from the list of available quantities (**kJ/mol** from the

Energy menu), and *click* on **OK**. Bring up the **Plots** dialog and select **Constraint (Con1)** from the **X Axis** and **Eaq(kJ/mol)** from the **Y Axes** list. *Click* on **Create Plot**.


### 9 to 12 optional

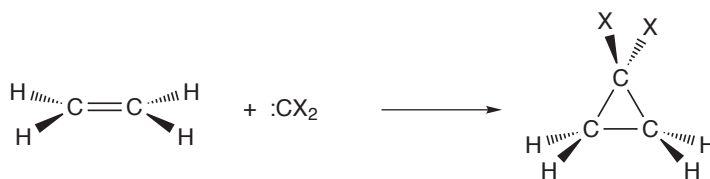9. With *bromide+methyl chloride.Prof.M0001* selected, enter the **Calculations** dialog. Specify **Energy (Semi-Empirical and AM1)** should be selected) and *click* on **OK**.
10. Enter the **Surfaces** dialog. *Click* on **More Surfaces....** Select **density (bond)** from the **Surface** menu and **none** from the **Property** menu and *click* on **Apply**. Select **density (bond)** from the **Surface** menu, but this time **potential** from the **Property** menu. *Click* on **OK**.
11. Submit the job. When completed, enter the **Surfaces** dialog and examine the bond density surface that you generated. *Click* on  at the bottom left of the screen to animate the display. Note, that bonds are smoothly broken and formed during the course of reaction. *Click* on  at the bottom of the screen when you are done.
12. Reenter the **Surfaces** dialog. Turn display of the bond density off, and then turn display of the electrostatic potential mapped onto the bond density on. *Click* on . Relate the migration of negative charge during reaction as indicated by colors in the electrostatic potential map to the plot you constructed earlier in step 7. Recall, that colors near red indicate maximum negative potential.
13. Remove all molecules and any remaining dialogs from the screen.


## Carbene Additions to Alkenes




Singlet carbenes add to alkenes to yield cyclopropanes. Since a singlet carbene possesses both a high-energy occupied molecular orbital in the plane of the molecule, and a low-energy, out-of-plane unoccupied molecular orbital, this reaction presents an interesting dilemma. Clearly it would be more advantageous for the low-lying vacant orbital on the carbene, and not the high-lying filled orbital, to interact with the olefin  $\pi$  system during its approach.

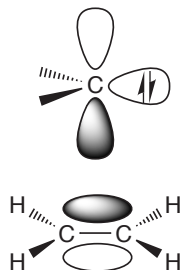

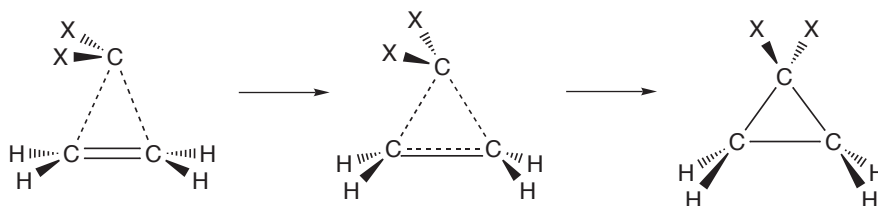
However, this leads to a product with an incorrect geometry. The carbene must twist by  $90^\circ$  during the course of reaction.



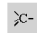In this tutorial, you will use the Hartree-Fock 3-21G model to find the transition state and to analyze the motion of the fragments.


1. Bring up the organic model kit and build ethylene.
2. Select  from the model kit. *Click* on the **Insert** button at the bottom right of the screen or hold down the **Insert** key (option key on Mac) and then *click* anywhere on screen. Alternatively, *double click* on a blank area of the screen following selection of  $sp^3$

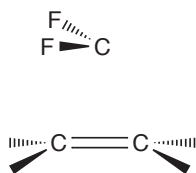carbon. Next, select  $\ominus\text{F}$  from the model kit and *click* on two of the free valences on the  $\text{sp}^3$  carbon. Next, *click* on  $\times$  and, one after the other, *click* on the remaining two free valences on the  $\text{sp}^3$  carbon. Finally, *click* on  $\oplus$ . You are left with two fragments, ethylene and difluorocarbene. Orient the two as to be poised for reaction.\*




Translations and rotations normally refer to the complete set of fragments, but if you *click* on a fragment (not on a free valence) to select it, and then hold down the **Ctrl** key they will refer to an individual fragment.

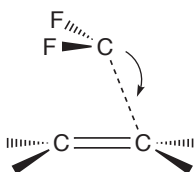3. Select **Guess Transition State** from the **Search** menu or *click* on . *Click* on the carbon on the  $\text{CF}_2$  fragment and then, while holding down the **Shift** key, *click* on the  $\text{CF}_2$  carbon and on one of the carbons on the ethylene fragment. A dotted line is drawn between the two carbons that are to be bonded along with an arrow from the  $\text{CF}_2$  carbon to the center of this line.



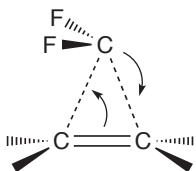
*Click* on the CC double bond and then, while holding down the **Shift** key, *click* on the other ethylene carbon and on the  $\text{CF}_2$  carbon. A second dotted line and arrow will be drawn.







---

\* Proper orientation of the two fragments is not crucial in this case, but is primarily to allow you to associate the arrows with the intended reaction. Proper orientation is, however, essential where different stereochemical outcomes are possible.




Click on  at the bottom right of the screen. Your structure will be replaced by a guess at the transition state.



4. Enter the **Calculations** dialog. Specify calculation of a **Transition State Geometry** using the **3-21G Hartree-Fock** model. *Check IR* to the right of **Compute**, and *click* on **Submit**. Name the job *difluorocarbene+ethylene*.
5. When the job is complete, examine the geometry of the transition state. In light of the previous discussion, would you describe your structure as corresponding to an early or late transition state? Animate the vibration corresponding to the reaction coordinate. Bring up the **IR** dialog (**Spectra** under the **Display** menu and *click* on the **IR** tab). *Click* on the imaginary frequency at the top of the list of frequencies. Does the animation show that the carbene reorients as it approaches the double bond? Turn the animation off by again *clicking* on the imaginary frequency.
6. Select **Properties** (**Display** menu) and, in turn, *click* on each of the four hydrogens in the transition state. Change the value in the **Mass Number** menu in the **Atom Properties** dialog (**Properties** under the **Display** menu) from **1** to **2 Deuterium**. Resubmit the job. (No additional quantum chemical calculations are involved, but the vibrational analysis needs to be repeated.) When complete, examine the new set of vibrational frequencies. Note that they are uniformly smaller than those for the undeuterated system, and that the largest changes are for vibrational motions where hydrogens are involved.

### 7 to 10 optional


7. Make a copy of *difluorocarbene+ethylene* () name it *difluorocarbene+cyclohexene*. The original file will be closed.
8. Bring up the organic model kit with the copy (*click* on ). Select **Freeze Center** from the **Geometry** menu (or *click* on the  icon at the top of the screen). One after the other, *click* on all atoms in the transition state except two of the hydrogens on one



side of the ethylene fragment.\* Each of the selected atoms will be given a magenta colored marker, indicating that it is to be frozen (not moved during optimization). Build a chain of four  $sp^3$  carbons  off one of the free valences that you have not frozen. Select  and *click* on the free valence on the terminal  $sp^3$  carbon and on the remaining free valence that you have not frozen to make a cyclohexene ring. *Click* on . Note that the only four methylene groups move and not the underlying transition state. You have produced an excellent guess at the transition state for difluorocarbene addition to cyclohexene based on your previous transition state for the simpler process.\*\*

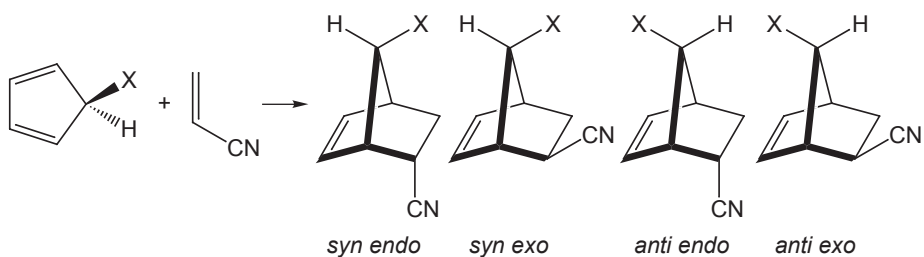
9. Enter the **Calculations** dialog. It should already specify calculation of **Transition State Geometry** using the **3-21G Hartree-Fock** model. Remove the checkmark on **IR** (to the right of **Compute**). **Do not** check **Frozen Atoms** (to the right of **Subject To**). *Click* on **Submit**.
10. To compare the geometries of the two transition states, you need to put them into the same document. Select **Append Molecule(s)...** from the **File** menu, *click* on *difluorocarbene+ethylene* in the browser that results and *click* on **Open**. Select **Spreadsheet** from the **Display** menu and *check* each of the boxes to the right of their identifiers. This allows the two transition states to be displayed simultaneously on screen. *Click* on **Add** at the bottom of the spreadsheet and select **Alignment Scores** from the menu that results. *Click* on , select **structure** from the **Align by** menu at the bottom right of the screen, *click* on the three carbons that make up the cyclopropane ring and *click* on **Align by** at the bottom right of the screen. *Click* on . Are the two transition states very similar as expected (an alignment score of 1.0 is perfect alignment)?
11. Remove all molecules and dialogs from the screen.

---

\* You could do this more quickly by first *clicking* on , then drawing a selection box around all atoms to be frozen, and then *clicking* inside this box.


\*\* It is not necessary to replace the deuteriums by hydrogens. You will not request an IR spectrum and calculated geometry does not depend on atomic mass.


## Stereospecific Diels-Alder Reactions



Diels-Alder cycloaddition of 5-substituted cyclopentadienes with acrylonitrile can lead to four stereoproducts, in which the substituent, X, at the 5 position is *syn* or *anti* to the dienophile, and the nitrile is *endo* or *exo*. *Anti* products are preferred when X is alkyl (consistent with sterics), while *syn* products are favored when X is halogen or alkoxy. In general, *endo* adducts are kinetically favored over *exo* adducts (see following tutorial). In this tutorial, you will use AM1 calculations to obtain both *syn* and *anti* transition states for *endo* addition of both 5-methylcyclopentadiene and 5-fluorocyclopentadiene and acrylonitrile, and then use Hartree-Fock 3-21G calculations to estimate relative activation energies. All four transition states are likely to be very similar to the transition state for the parent cycloaddition (cyclopentadiene and acrylonitrile). You will take advantage of this, and obtain a transition state for the parent reaction, and then to use it as a starting point for locating transition states for the substituted systems.

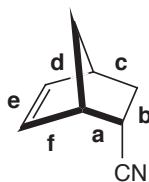

1. Construct the following substituted norbornene (the product of *endo* addition of cyclopentadiene and acrylonitrile).




2. Select **Guess Transition State** from the **Search** menu or *click* on . *Click* on bond **a** (see figure above) and then on bond **b**. A curved arrow will be drawn from **a** to **b**. Next, *click* on bond **c** and then on bond **d**; a second curved arrow from **c** to **d** be drawn. Finally, *click* on bond **e** and then on bond **f**, leading to a third curved arrow. The model on screen should now appear as follows.




Click on  at the bottom right to produce a guess at the transition state.

3. Inside the **Calculations** dialog, specify calculation of **Transition State Geometry** using the **Semi-Empirical AM1** model. *Check IR* under **Compute** to specify calculation of vibrational frequencies. This will allow you to verify that you have indeed located a transition state that is consistent with a Diels-Alder reaction.
4. Submit the job. Name it *cyclopentadiene+acrylonitrile*. When completed, examine the resulting structure and calculated frequencies. Verify that there is one and only one imaginary frequency. Animate the motion associated with the imaginary frequency (**Spectra** under the **Display** menu, *click* on the **IR** tab and *check* the box to the left of the imaginary frequency in the dialog that results) to verify that it corresponds to the expected motion along the reaction coordinate.
5. Place the transition-state structure onto the clipboard. Right *click* on the background and select **Copy** from the resulting menu. Close *cyclopentadiene+ acrylonitrile* (.
6. Bring up the organic model kit (). *Click* on **Clipboard** (at the bottom of the model kit) and then *click* anywhere on screen. The transition-state structure will appear. Add a fluorine to the methylene group of the cyclopentadiene fragment of the transition state either *syn* or *anti*. **Do not minimize**. Your starting structure should provide an excellent guess at the transition state for the substituted system.
7. Select **New Molecule** (not **New**) from the **File** menu. The screen will clear. *Click* on **Clipboard** and then *click* anywhere on screen. Add fluorine to the other methylene group position. Repeat the

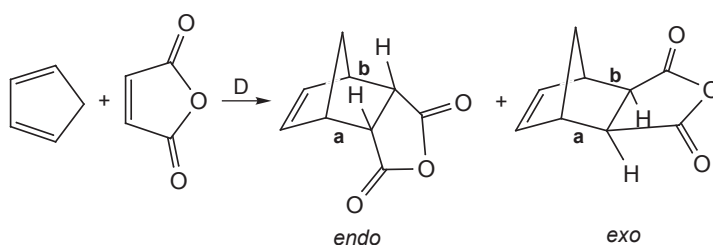
process two more times (starting with **New Molecule**), to add methyl both *syn* and *anti* on the methylene group. When you are done (four substituted transition states in total), *click* on .

Bring up the spreadsheet and replace the molecule identifiers in the leftmost column (*M0001*, etc.) but more descriptive names (*F syn*, etc.).






8. Enter the **Calculations** dialog and specify calculation of **Transition State Geometry** using the **Semi-Empirical AM1** model. *Check IR* to the right of **Compute**. Make certain that **Global Calculations** is checked. Submit the job (name it ***Diels-Alder stereochemistry***). When it has completed, verify that all four structures correspond to transition states (using **Spectra** from the **Display** menu).
9. Make a copy of ***Diels-Alder stereochemistry*** (). Name it ***Diels-Alder stereochemistry Hartree-Fock***. Enter the **Calculations** dialog with this copy and specify a **3-21G Hartree-Fock Energy** calculation. Make certain that you remove the checkmark on **IR** and that **Global Calculations** is checked before you exit the dialog.
10. Submit the job. When it completes, bring up the **Reactions** dialog (**Reactions** from the **Display** menu). Select *syn* and *anti* F compounds (and then *syn* and *anti* Me compounds) for **Reactants** and **Products**, respectively, select **Current Document** from the **Use** menu, **kJ/mol** from the **Units** menu and *click* on the **Compute Energies** button at the bottom of the dialog. Use the results to identify the more stable product in each case.
11. Remove all molecules and any remaining dialogs from the screen.

## Thermodynamic vs. Kinetic Control

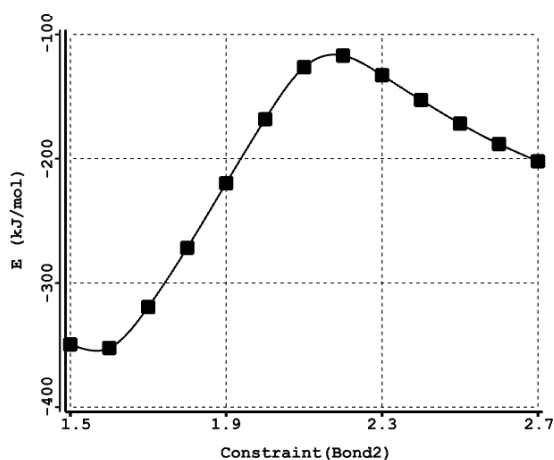
Chemical reactions may yield different products depending on the conditions under which they are carried out. High temperatures and long reaction times favor the most stable (thermodynamic) products, whereas low temperatures and short reaction times favor the most easily formed (kinetic) products. For example, the kinetic product in Diels-Alder cycloaddition of cyclopentadiene and maleic anhydride is the *endo* adduct, whereas it seems likely that the less-crowded *exo* adduct is the thermodynamic product.



Here you will first obtain pathways for reactions leading to both *endo* and *exo* adducts using the PM3 semi-empirical model and then follow these by 3-21G Hartree-Fock energy calculations to get a better estimate of the difference in activation energies.

1. Build the *endo* Diels-Alder adduct of cyclopentadiene and maleic anhydride and minimize ()
2. Select **Constrain Distance** from the **Geometry** menu (or *click* on ) and *click* on bond **a** in the figure above. *Click* on the icon  at the bottom right of the screen (it will then turn to ). *Check* the box to the left of **Profile** at the bottom right of the screen. This will give rise to two additional check boxes. Change the value in the leftmost box from 1.54Å to 1.5Å by *typing* **1.5** inside the box and then *pressing* the **Enter** key (**return** key on Mac). *Type* **2.7** (2.7Å) into the box to the right of **to** and *press* the **Enter** (**return**) key. Set **Steps** to **13**. *Type* **13** and *press* the **Enter** (**return**) key. You have specified that bond **a** will be constrained first to 1.5Å, then to 1.6Å, then to 1.7Å, etc. and finally to 2.7Å. *Click* on .

- Repeat the process for bond **b**. When you are done, both bonds **a** and **b** will be constrained from 1.5Å to 2.7Å in 13 equal steps. Click on .
- Bring up the **Calculations** dialog and select **Energy Profile** from the top menu to the right of **Calculate**, and **Semi-Empirical** and **PM3** from the two bottom menus. Click on **Submit** at the bottom of the dialog. Name it *cyclopentadiene+maleic anhydride endo*.
- When completed the job will give rise to a new document *cyclopentadiene+maleic anhydride endo.Prof.M0001*. This contains 13 calculations corresponding to the 13 steps that make up the energy profile. You will be prompted as to whether you want to open this file. Click on **OK**.\*
- Bring up the **Spreadsheet**, and click on **Add...** at the bottom. Select **E** from among the entries, **kJ/mol** from the **Energy** menu, and click on **OK**. Next, click on  and click on one of the two CC bonds varied in the energy profile. Click on  at the bottom right of the screen. Click on . Finally, bring up the **Plots** dialog (**Plots** from the **Display** menu), and select **Constraint (BondX)\*\*** from among the items in the **X Axis** menu and **E (kJ/mol)** from the **Y Axes** list. Click on **Create Plot**.





\* To avoid confusion, it is a good idea to close the original file *cyclopentadiene+maleic anhydride endo*.

\*\* Bonds are numbered in the order they were formed upon initial construction of the molecule.

Identify both the reactant and transition state from the plot and estimate the activation energy for the cycloaddition reaction.



7. Repeat steps 1 to 6 for the *exo* adduct. Compare the activation energy for *exo* addition to that for *endo* addition (above). What is the kinetic product?

### 8 and 9 optional

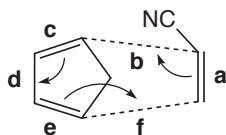
8. Open *cyclopentadiene+maleic anhydride endo.Prof.M0001* () and make a copy (). Name it *cyclopentadiene +maleic anhydride endo Hartree-Fock*. Enter the **Calculations** dialog and specify a **Hartree-Fock 3-21G Energy** calculation. Submit. When completed, perform the same spreadsheet and plot operations you did for the PM3 calculations.
9. Repeat the above procedure for the *exo* adduct and compare the two activation energies. What is the kinetic product?
10. Remove all molecules and dialogs from the screen.


## Activation Energies of Diels Alder Reactions



In an earlier tutorial *Dienophiles in Diels-Alder Reactions* (Chapter 5), you examined the extent to which LUMO energies for a series of related dienophiles correlated with relative rates of Diels-Alder cycloadditions involving cyclopentadiene. In this tutorial, you will compare calculated activation energies for this same set of reactions with the experimental rates. You will use the previous set of calculations on the dienophiles and obtain transition states from the Spartan Reaction Database (SRD). The only quantum chemical calculation that is required is for cyclopentadiene.


1. Build cyclopentadiene and minimize (). Before you proceed with the next step, copy cyclopentadiene onto the clipboard. *Right click* and select **Copy** from the menu that appears.
2. Enter the **Calculations** dialog (**Setup** menu) and specify an **Equilibrium Geometry** calculation using the **Hartree-Fock 3-21G** model. Submit the job with the name *cyclopentadiene*. Move it off to the side of the screen for use later on.
3. Enter the builder (). *Click* on **Clipboard**. The structure of cyclopentadiene will appear at the top of the model kit. *Click* anywhere on screen and cyclopentadiene will appear.
4. Select **Alkenyl** from the **Groups** menu, *press* the **Insert** button (at the bottom right of the model kit) or hold down the **Insert** key (**alt** key on Mac) and *click* on screen, or *double click* in a blank area on screen. Both cyclopentadiene and ethylene will appear on screen, but they will not be connected. Select **Cyano** from the **Groups** menu and *click* on one of the free valences on ethylene. Both cyclopentadiene and acrylonitrile will appear on screen.
5. Orient the two molecules such that they are poised for a Diels-Alder reaction leading to an *endo* product.





Select **Define Reaction Query** from the **Search** menu (*click on \**). *Click* on the bond marked **a** in the above figure and, while holding down the **Shift** key, *click* on the two carbons that when connected will lead to the bond marked **b**. A dotted line will be drawn between these carbons and an arrow drawn from the center of bond **a** to the center of bond **b**. Next, *click* on bond **c** and then on bond **d**. A second arrow will be drawn. Finally, *click* on bond **e** and, while holding down the **Shift** key, *click* on the two carbons that when connected will lead to the bond marked **f**. A second dotted line and a third arrow will be drawn.

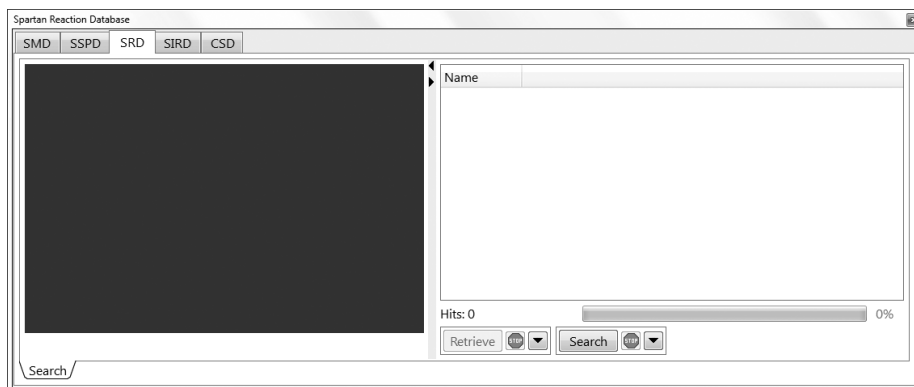
6. Select **Define Structure Query** from the **Search** menu (*click on *). *Click* on all three free valences on the acrylonitrile fragment. They will be replaced by orange cones, meaning that anything (including hydrogen) may be attached to these positions.
7. Select **Open Databases** from the **Search** menu (*click on *) and *click* on the **SRD** tab at the top of the dialog that results. This leads to the **Spartan Reaction Database (SRD)** dialog.



Before you begin the search, *click* on  to the right the **Search** button at the bottom of the dialog. This brings up the **Search**


\* **Define Search Query** is equivalent to **Guess Transition State** for the purpose of defining “curly arrows”. It lacks the ability to automatically guess a transition state.

## Options dialog.

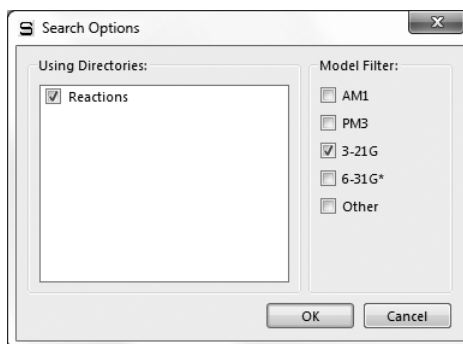

Check **3-21G** under **Method Filters** and *click* on **OK**. *Click* on the **Search** button to search the Spartan Reaction Database for Diels-Alder transition states between cyclopentadiene and substituted acrylonitriles.

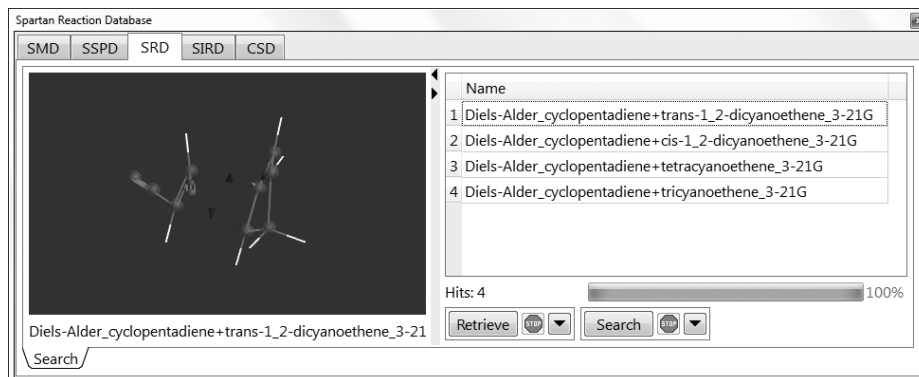- When the search completes, a listing of hits appears at the right of the dialog.



This should contain four Diels-Alder transition states for *endo* additions of cyclopentadiene and *trans*-1,2-dicyanoethylene, *cis*-1,2-dicyanoethylene, tricyanoethylene and tetracyanoethylene. All need to be retrieved into a single document. *Click* on  to the right of the **Retrieve** button at the bottom of the dialog to bring up the **Retrieve Options** dialog. Make sure that **New Document** is selected in this dialog and *click* on **OK**. Hold down the **Ctrl** key, select (*click* on) all six Diels-Alder transition states and *click* on **Retrieve** at the bottom of the dialog. A new

- (unnamed) document will be created. Open a spreadsheet for this document.
9. Open *Diels-Alder dienophiles* that you created in an earlier tutorial\* and bring up a spreadsheet for this document. Select all entries in this document. *Click* inside the lefthand cell of the first entry, hold down the **Shift** key and *click* on the last entry. *Right click* and select **Copy** from the menu that results. Move the cursor over the leftmost cell of the first empty column of the new (unnamed) document, *right click* and select **Paste** from the menu that results. The contents of *Diels-Alder dienophiles* (the dienophiles and cyclopentadiene) will be copied into the new document. Close *Diels-Alder dienophiles*.
  10. Bring up the **Reactions** dialog (**Reactions** from the **Display** menu). One after the other, calculate activation energies for the six Diels-Alder reactions. Use cyclopentadiene and one of the dienophile for **Reactants** and the appropriate transition state and *none* for **Products**. Select **Current Document** under **Use** at the bottom of the dialog and *click* on **Compute Energies**. Is there a correlation between calculated activation energies and relative rates (experimental relative rates are available in the tutorial **Diels-Alder Dienophiles in Chapter 5**)?
  11. Close all molecules and any open dialogs.

---

\* If you have not completed this tutorial, a copy of the document *Diels-Alder dienophiles* may be found in the *organic reactions* sub-directory (*tutorials* directory). For Windows, this is found in *Program Files/Wavefunction/Spartan10* or *Program Files(x86)/Wavefunction/Spartan10*. For Linux, this is found in the directory where *Spartan* was installed. For Macintosh, this is located in the *Spartan.10.full* folder on the installation CDROM.

# Chapter 8

## Medicinal Chemistry

---

*This chapter illustrates applications of **Spartan** to problems of relevance to medicinal chemists. The first involves an investigation of the relationship between polar surface area (defined both by atom counts and by electrostatic potential maps) and the rate of transport between blood and brain. The second uses a published pharmacophore to decide whether the antihistamine terfenadine should be considered a potential potassium channel blocker. The third seeks to find molecules with similar analgesic activity as morphine, based first on similar structure and then on similar chemical functionality.*

Many molecular properties may now be calculated using quantum chemical methods, some more accurately than they can be measured. Still, some quantities of key interest to medicinal chemists are not yet subject to calculation. A good example is the rate at which a molecule passes through a membrane. Here, overall molecular polarity plays an obvious role. Too polar and the lipid will not take up the molecule; not sufficiently polar and it will be unable to dissolve in the aqueous environment on either side of the membrane. One descriptor commonly used to account for polarity is the so-called polar surface area, defined as the area of space-filling model due to nitrogen and oxygen together with any bonded hydrogens. The first tutorial in this chapter examines the correlation between polar surface area and blood-brain transport and then proposes an alternative based on electrostatic potential maps.

Many (if not most) drugs act by binding to a protein or nucleotide, either promoting a particular process (agonists) or inhibiting one (antagonists). New drugs can be found by identifying molecules that are structurally and/or functionally similar to compounds that are known or presumed to bind. The second tutorial examines whether the antihistamine terfenadine is likely to show serious side effects because of a binding affinity in potassium ion channels. Here, the

similarity analysis is not carried out between molecules but between a molecule and an environment. The third tutorial makes use of similarity analysis to suggest compounds that look like morphine and hence might exhibit similar action.

## Anticipating Blood-Brain Transport

A drug must be able to be transported throughout the body. This includes intestinal absorption and passage across the blood/brain barrier. In this tutorial, you will first examine how well the polar surface area, the area of space-filling model due to nitrogen and oxygen atoms together with any attached hydrogens, correlates with  $\log(C_{\text{brain}}/C_{\text{blood}})$ \*.

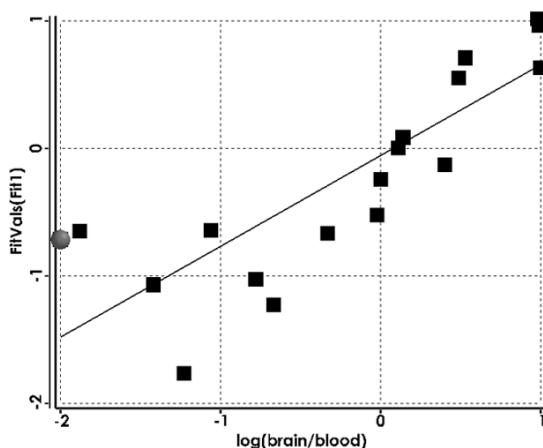
1. Open **blood brain transport** from the **medicinal chemistry** sub-directory (**tutorials** directory).\*\* This comprises a list of drugs/drug candidates for which experimental data relating to the ratio of concentrations in the brain and in the blood ( $C_{\text{brain}}/C_{\text{blood}}$ ) are available. These span a range of nearly 5 log units.
2. Bring up a spreadsheet (**Spreadsheet** under the **Display** menu) and size it to show all compounds, several extra rows and three columns, in addition to the column of experimental data **log(brain/blood)** that is already displayed. Bring up the **Molecule Properties** dialog and *click* on the **QSAR** tab. *Click* on the **P** to the left of **PSA**. Polar surface area values will fill one column in the spreadsheet.
3. *Click* on the **Add** button at the bottom of the spreadsheet and select (*click* on) the **Linear Regression** tab at the top of the dialog that appears (most likely the **Add Columns** dialog). This leads to the **Add Linear Regression** dialog. Select **log(brain/blood)** from the **Fit** menu and *click* on **PSA** from the quantities in the box below **Using**. *Click* on **OK**. Bring up the **Regression Properties** dialog by first selecting **Properties** from the **Display** menu and then *clicking* on the cell labelled **Fit1** under **Label** in the spreadsheet. RMSD and  $R^2$  values are reported. The smaller RMSD and the closer to unity  $R^2$ , the better the fit. Leave the **Regression Properties** dialog on screen.

---

\* J. Kelder *et al.*, *Pharmaceutical Res.*, **16**, 1514 (1999). The 19 named compounds considered in this paper have been included.

\*\* For Windows, this is found in **Program Files/Wavefunction/Spartan10** or **Program Files(x86)/Wavefunction/Spartan10**. For Linux, this is found in the directory where **Spartan** was installed. For Macintosh, this is located in the **Spartan.10.full** folder on the installation CDROM.


4. Make a plot of  $\log(\text{brain/blood})$  vs. the regression fit. Select **Plots** from the **Display** menu to bring up the **Plots** dialog. Select **log (brain/blood)** from the **X Axis** menu, **Fit Vals (Fit 1)** from the **Y Axis** list and *click* on **OK**. Initial presentation uses a smooth curve (cubic spline) to connect the data points. *Click* on the curve (not on the axes) to bring up the **Curve Properties** dialog. Select **Linear LSQ** from the **Fit** menu to give a least squares line.



Next, you will consider an alternative definition of polar surface area, based on electrostatic potential maps.

5. Select **Surfaces** from the **Display** menu and select **Electrostatic Potential Map** in the dialog that results. An electrostatic potential map for the selected drug appears (maps for the remaining drugs may be seen by stepping through the list). The overall size and shape is that of the electron density and corresponds roughly to a conventional space-filling model. The colors indicate the value of the electrostatic potential. Colors toward red designate areas of negative potential (where a positive charge is most likely to be attracted), while colors toward blue designate areas of positive potential (where a positive charge is least likely to be attracted). To see the molecular skeleton underneath the electrostatic potential map, change to a transparent or mesh display. *Click* on the map to select it and select **Transparent** or **Mesh** from the **Style** menu that appears at the bottom right of the screen.

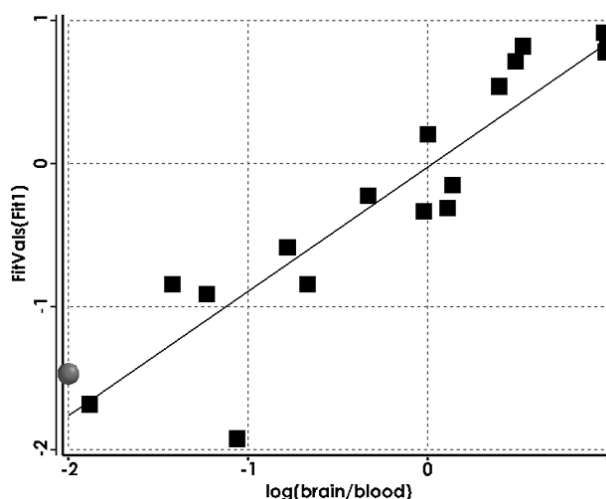6. Select **Properties** from the **Display** menu. *Click* on the map. *Click* on the  button to the left of **P-Area** inside the **Surfaces Properties** dialog that results. Polar areas, defined as the area for which the absolute value of the electrostatic potential is  $> 100$  kJ/mol\*, will be added to the spreadsheet.
7. *Click* on the **Add** button at the bottom of the spreadsheet and *click* on the **Linear Regression** tab at the top of the dialog that results. Select **log(brain/blood)** from the **Fit** menu and *click* on **Polar Area** from the quantities in the box below **Using**. *Click* on **OK**. Bring up the **Regression Properties** dialog by *clicking* on the cell labelled **Fit2** under **Label** in the spreadsheet. Note that the  $R^2$  value is better (closer to unity) for this fit than for the previous fit (to PSA defined in the usual manner). Make a plot of  $\log(\text{brain/blood})$  vs. the new regression fit. Select **Plots** from the **Display** menu, choose **log (brain/blood)** from the **X-Axis** menu and **FitVals (Fit 2)** (instead of **Fit Vals (Fit 1)**) from the **Y Axis** list and *click* on **OK**. Switch to a least square fit by *clicking* on the curve that results and selecting **Linear LSQ** from the **Fit** menu in the **Curve Properties** dialog. Which property provides better correlation with transport across the blood/brain barrier, PSA or polar area?




8. Close *blood brain transport* when you are finished.


\* This value may be changed in the **Settings Preferences** dialog (**Preferences** under the **Options** menu; **Chapter 20**).



## Terfenadine. A Potassium Channel Blocker?

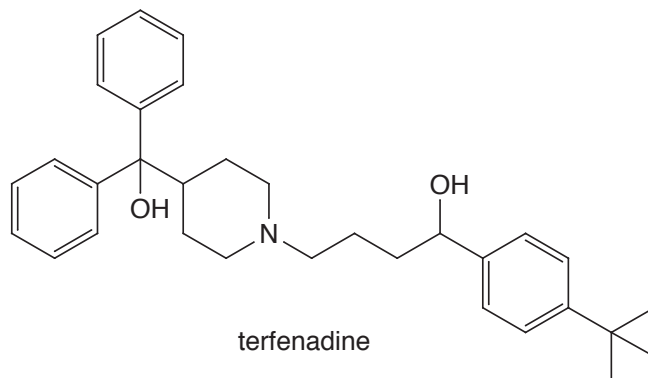Cardiovascular toxicity due to blocking of potassium ion channels is commonly screened early in the drug development process. One way this is done is to compare drug candidates to a pharmacophore deduced from 3D QSAR studies. This tutorial uses a published pharmacophore\* to see whether the antihistamine, terfenadine should be considered a potential channel blocker.




1. Open *terfenadine conformer library* from the *medicinal chemistry* sub-directory (*tutorials* directory).\*\* This contains a selection of several hundred conformers of terfenadine (the number is reported in the **Molecule Properties** dialog) obtained from a conformer library calculation. It has been supplied in order to save computer time. Close *terfenadine conformer library*.
2. Open *potassium channel blocker pharmacophore* from the *medicinal chemistry* sub-directory (*tutorials* directory). This is a five-point pharmacophore comprising four hydrophobes (blue spheres) and one positive ionizable center (red sphere).
3. Select **Set Similarity Centers** from the **Geometry** menu (or *click* on the  icon at the top of the screen) and then select **CFD** from the menu at the bottom right of the screen. In turn, *click* on each of the (five) pharmacophore elements. In response,

\* S. Ekins *et al*, J. Pharmacology and Experimental Therapeutics, **301**, 427 (2002).

\*\* For Windows, this is found in *Program Files/Wavefunction/Spartan10* or *Program Files(x86)/Wavefunction/Spartan10*. For Linux, this is found in the directory where *Spartan* was installed. For Macintosh, this is located in the *Spartan.10.full* folder on the installation CDROM.


each will be surrounded by a violet circle indicating that it is to be used in the similarity analysis.

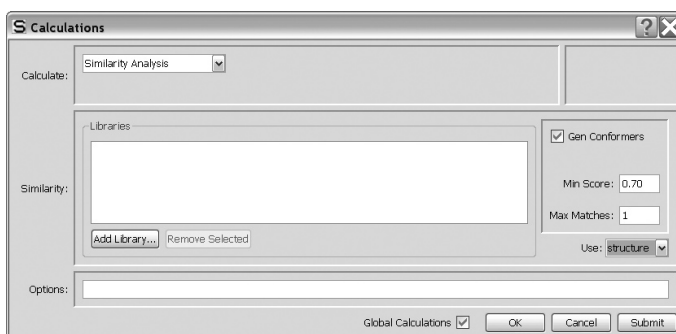4. Select **Calculations...** from the **Setup** menu. Select **Similarity Analysis** from the top left menu to the right of **Calculate** inside the **Calculations** dialog. This leads to a new form of the **Calculations** dialog.



*Click* on **Add Library** to bring up a file browser. Locate **terfenadine conformer library** in the **medicinal chemistry** sub-directory (**tutorials** directory) and *click* on **OK**. This indicates that the library of terfenadine conformers will be searched for a match to the pharmacophore. *Click* on **Submit**.

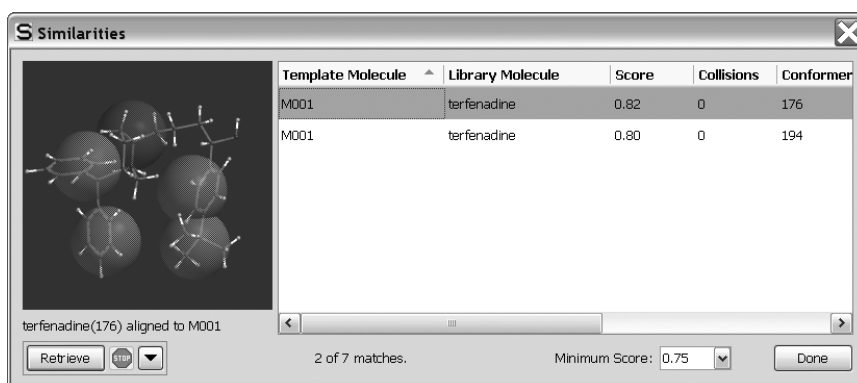5. When the analysis has completed, select **Similarities...** from the **Display** menu to bring up the **Similarities** dialog.




One or more hits will appear in the box at the right of the dialog, together with a similarity score (limiting on 1.0 which means a perfect match). If there is more than one hit, *click* on that

with the largest score. A composite graphic of terfenadine and the pharmacophore will appear in a window at the left of the dialog. You can manipulate the two (as a single object) in the usual way (you need to position the cursor inside the window). Note that the positive ionizable center is positioned above the nitrogen, while the four hydrophobes are above the phenyl rings and the *tert*-butyl group.

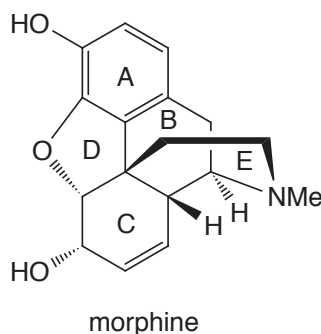The scoring algorithm not only accounts for position but also for direction of any nitrogen or oxygen centers that overlap with (non-hydrophobe) pharmacophore elements.





6. You can, if you wish, retrieve your hit (without retrieving the pharmacophore), by *clicking* on the **Retrieve** button at the bottom left of the **Similarities** dialog. This can be used for further analysis.
9. Close all open documents on screen.

## Morphine. Structure vs. Pharmacophore




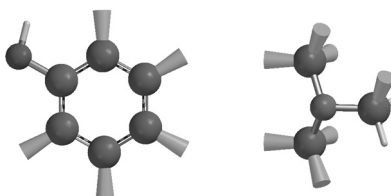
Three elements of morphine appear to be required in order for the molecule to act as an analgesic. These are: i) the nitrogen center (assumed to be protonated in the protein-bound complex), ii) the aromatic ring and iii) the hydroxyl group attached to the aromatic ring. Loss of any of these results in significant loss of activity. This knowledge may either be used directly, to identify other likely analgesics based on structure, or indirectly, to construct a simple 3-point pharmacophore that in turn may be employed to find compounds with potentially comparable analgesic activity.

This tutorial comprises three parts. In the first part, you will identify molecules in SMD that incorporate the three required elements. In practice, you will look for molecules that incorporate a nitrogen center substituted by a methyl group and two other ( $sp^3$ ) carbons and a substituted phenol. You will then select a single hit and generate a conformer library for similarity analysis. (In a comprehensive investigation, you would likely use several hits or all hits; the restriction here is simply in order to save computer time.)


1. Bring up the organic model kit () . Select **Benzene** from the **Rings** menu and *click* anywhere on screen. Select  $sp^3$  oxygen () and *click* on one of the free valences of benzene. You have made phenol. Select  $sp^3$  nitrogen () , *click* on the **Insert** button at the bottom of the screen or hold down the **Insert** key (**option** key on Mac) and *click* anywhere on screen. Alternatively, *double click* on a blank region on screen following selection of  $sp^3$  nitrogen. Two molecules (phenol and ammonia) appear on screen. *Click* on  $sp^3$  carbon () and *click* on all three free

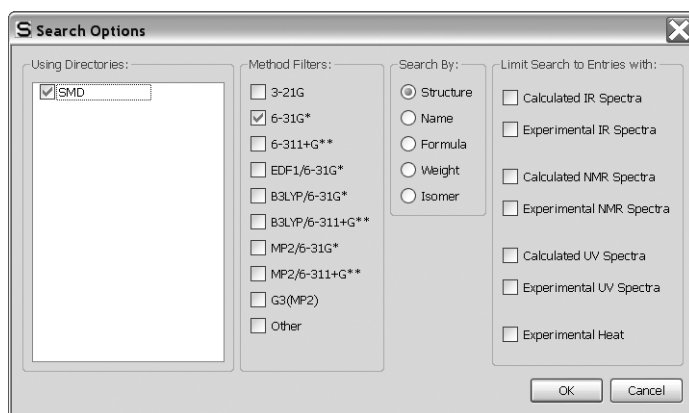
valences on the nitrogen. Phenol and trimethylamine now appear on screen. *Click* on  remove the model kit.

2. Select **Define Structure Query** () from the **Search** menu. *Click* on all five free valences on the carbons in the phenol fragment and on all three free valences on two of the three methyl groups and on one free valence for the remaining methyl group in the trimethylamine fragment. Orange cones will appear at all (twelve) selected positions.




Anything can be grown from these positions (including hydrogen) but any hit must contain phenol and a nitrogen bonded to three  $sp^3$  carbon centers.



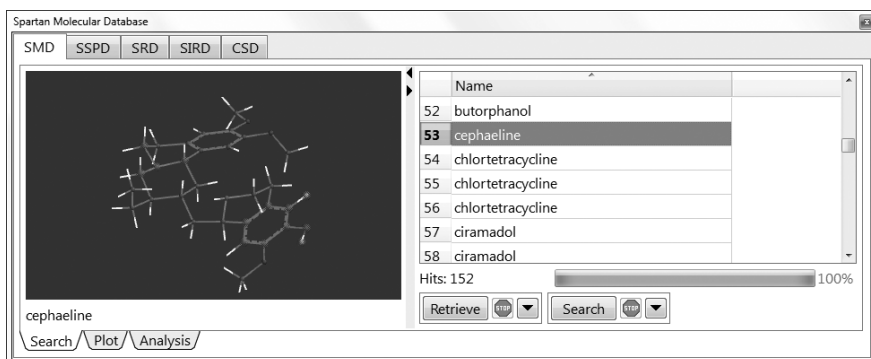3. Select **Open Databases** from the **Search** menu (). *Click* on the **SMD** tab in the dialog that results. *Click* on  to the right of the **Search** button in the **SMD** dialog to bring up the **Search Options** dialog.






*Check 6-31G\** under **Method Filters** and **Structure** under **Search by** in the **Search Options** dialog. *Click* on **OK**. Finally, *click* on the **Search** button at the bottom of the **SMD** dialog.

4. The search will result in several hits. Find and *click* on **cephaeline** from among those hits listed inside in the scroll box at the right of the **SMD** dialog.

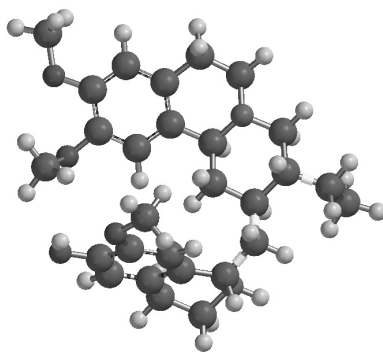


Its structure will be displayed in the window at the left of the dialog\*. You can manipulate the model in the usual way by positioning the cursor inside this window. *Click* on  to the right of the **Retrieve** button at the bottom of the **SMD** dialog to bring up the **Retrieve Options** dialog. *Check New Document* under **Retrieve Options** and *click* on **OK**. Finally, *click* on the **Retrieve** button. Dismiss the **SMD** dialog by *clicking* on  at the top right-hand corner.

Close the molecule used for the search (phenol and trimethylamine) without saving (). Save the molecule you retrieved under the name **cephaeline** (.

5. To shorten computer time, remove some of the degrees of conformational freedom from **cephaeline**. Select **Set Torsions** from the **Geometry** menu (or *click* on ). Your model will be augmented with yellow cylinders to designate flexible bonds. Remove the markers from the bonds connecting the three methoxy groups to the aromatic rings by *double clicking* on each in turn. These bonds will no longer be included in the conformational search. Your model should now appear as follows with only three single bonds selected.

\* *Right click* on **Name** at the top of the list to sort in alphabetical order.




If more are selected, *double click* on the appropriate torsion markers to remove them.

- Bring up the **Calculations** dialog (**Calculations...** from the **Setup** menu). Select **Conformer Library** from the top left menu to the right of **Calculate**. No further information is required. *Click* on **Submit** at the bottom of the dialog. Name it *cephaeline library*. You can continue with the tutorial while you are waiting for the job to complete. When *cephaeline library* completes, close this file (🗑️).


In the second part of the tutorial, you will identify the key structural components in morphine and perform a similarity analysis based on structure on the conformer library for cephaeline.


- You can either build morphine or (to save time) open it: *morphine* from the *medicinal chemistry* sub-directory (*tutorials* directory).

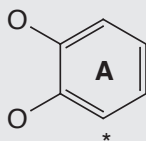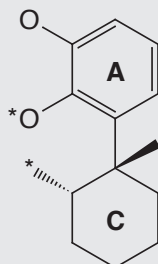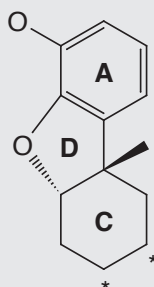
To build morphine, start with **Benzene** from the **Rings** menu. Add two (*ortho*)  $sp^3$  oxygens (🗑️) to make catechol and orient as below (free valences not drawn).

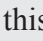



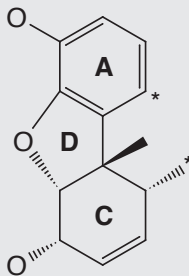
Select **Cyclohexane** from the **Rings** menu and *click* on the free valence on the \*'ed carbon in catechol. Orient as below, with the free valence on the cyclohexane carbon used to attach to benzene pointing straight up.





Select **Make Bond** from the **Build** menu (or *click* on ) . *Click* on the \*'ed free valences indicated in the structure above, one on oxygen and one on carbon pointing down (relative to the free valence on the carbon used to attach the two rings). A bond will be drawn to make the **D** ring.

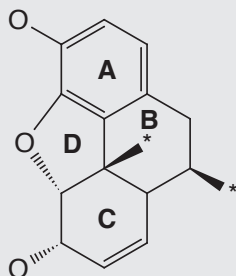


**Make Bond** is still selected. *Click* on *axial* free valences on the \*'ed carbons in the cyclohexane ring to introduce a double bond. Also, add a hydroxyl group to this ring () adjacent to the double bond. (Pay attention to stereochemistry.) Select **Minimize** from the **Build** menu (or *click* on ) to refine your structure.




Add  $sp^3$  carbons () to free valences at the \*'ed positions in the figure above (pay attention to stereochemistry on the cyclohexene ring). *Click* on  and then *click* on free valences on each of the two  $sp^3$  carbons that you have just added. You have made the **B** ring. Check that the stereochemistry is as it should be.





Connect the two \*'ed free valences in the figure above with a chain of two  $sp^3$  carbons ( $\text{>C-}$ ) and an  $sp^3$  nitrogen ( $\text{>N}$ ). Add an  $sp^3$  carbon to the nitrogen. You have made morphine (its name should appear at the bottom of the screen). *Click* on  $\text{E}$  to produce a refined structure and on  $\text{V}$  to remove the model kit.

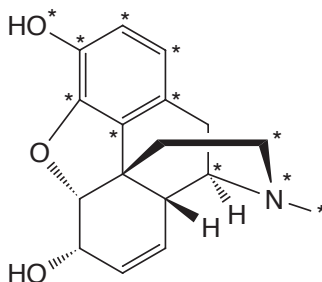
8. Select **Set Similarity Centers** from the **Geometry** menu (or *click* on  $\text{G}$ ) and select **structure** from the menu at the bottom right of the screen. *Click* on the eleven \*'ed atoms in the figure below. A red circle will be drawn around each. If you make a mistake and select the wrong atom, *click* on the circle to deselect it.



9. Select **Calculations...** from the **Setup** menu. Select **Similarity Analysis** from the top left menu to the right of **Calculate** inside the **Calculations** dialog. *Click* on **Add Library** to bring up a file browser, locate and select *cephaline library*, and *click* on **OK**. Make certain that **Structure** is selected from the **Use** menu at the lower right of the dialog. You have requested that the library of cephaline conformers be searched for a match to the key structural elements of morphine. *Click* on **Submit**.<sup>\*</sup> Supply the name *morphine*.


<sup>\*</sup> Make sure *cephaline library* is closed before carrying out the similarity analysis.

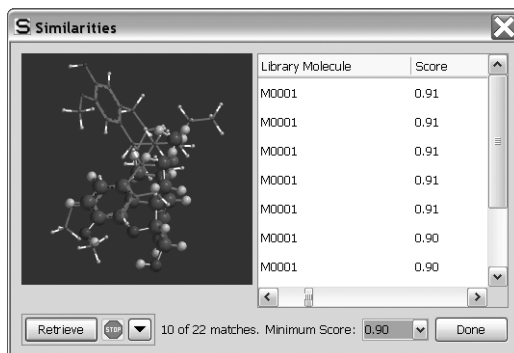
10. The similarity analysis will require a few seconds. When it completes, select **Similarities...** from the **Display** menu.



11. One or more hits appear in the box at the right of the dialog together with a score. Sort according to score (best on top) by *clicking* on **Score** at the top of the box. Select (*click* on) one of the best matches. A composite graphic consisting of morphine and the matching cephaeline conformer from the library will appear in a window at the left. You can manipulate the two structures (as a single object) in the usual way.

Repeat the similarity analysis using CFD's rather than structure.

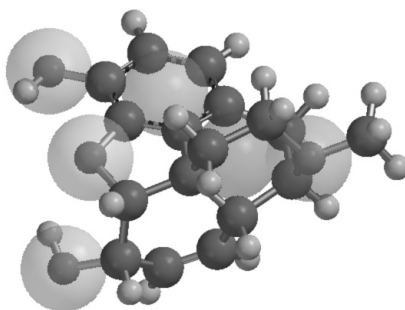12. Select **Set Similarity Centers** from the **Geometry** menu () , but this time select **CFD** from the menu at the bottom right of the screen. CFD's will augment the morphine structure.



13. Designate three of the CFD's as similarity centers. Select (*click* on) the CFD over the phenolic oxygen. A violet circle will surround this center indicating that it is to be used in the similarity analysis. (If you make a mistake and select the wrong

- CFD, *click* on the circle and it will disappear.) Select (*click* on) the CFD at nitrogen and the CFD at the middle of the phenol ring.
14. Enter the **Calculations** dialog. This should already designate **Similarity Analysis** and the appropriate library (*cephaeline*). Change the entry in the **Use** menu to **CFD**. *Click* on **Submit**.
  15. When the analysis is complete, bring up the **Similarities** dialog (**Similarities...** under the **Display** menu). One or more hits should appear in the box at the right of the dialog. Select a hit to get a composite graphic (CFD and structure from the library). You can easily see the extent to which the two are matched.
  16. Close any molecules and dialogs that are open on screen.

# Chapter 9

## Polypeptides to Proteins

---

*This chapter introduces the peptide model kit and shows how a structure from the Protein Data Bank (PDB) may be brought into **Spartan** for further analysis and a Ramachandran plot displayed.*

At present, calculations on biopolymers (proteins and RNA/DNA strands) are limited to molecular mechanics models. However, quantum chemical calculations are now routinely applicable to molecules bound into biopolymers and to the investigation of the electronic and steric factors that are responsible.

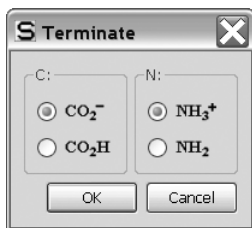
This chapter comprises three tutorials. The first illustrates **Spartan**'s peptide model kit to construct triglycine and decaglycine in idealized helical arrangements. Comparisons are made between zwitterion and non-zwitterion forms in and out of water and electrostatic potential maps are employed to investigate charge separation.

The second tutorial illustrates **Spartan**'s on-line access to PDB. A small molecule is extracted from a protein structure and the question is posed as to the energy of its (bound) conformer compared to that of the free conformer.

The third tutorial also involves extraction from a PDB file, but this time the quantity of interest is not the molecule itself, but rather its environment. This is turned into a so-called pharmacophore (a simplified representation of the environment), and a search for other molecules that may also be consistent with this environment is carried out.












Choose **CO<sub>2</sub>H** for the C terminating group and **NH<sub>2</sub>** for the N terminating group. *Click* on **OK**.

Amino acids replace atomic fragments, functional groups, rings and ligands as the basic building blocks in the peptide model kit, and these other building blocks are missing. Therefore, most modifications of peptides, aside from modifications in sequence of amino acids and in overall conformation, need to be carried out using either the organic or inorganic model kits.

5. Select **Freeze Center** from the **Geometry** menu (or *click* on ). *Click* on **Freeze Heavy** at the bottom right of the screen. This indicates that all heavy (non-hydrogen) atoms are not to be moved during molecular mechanics minimization.\* *Click* on  to produce a refined structure (subject to the restriction that heavy atoms are kept in place).
6. Select **New Molecule** from the **File** menu. The triglycine fragment is already specified. *Click* anywhere on screen. *Click* on **Terminate** inside the model kit and this time select **CO<sub>2</sub><sup>-</sup>** and **NH<sub>3</sub><sup>+</sup>** for C and N terminating groups, respectively. *Click* on **OK**. *Click* on  and then *click* on **Freeze Heavy** at the bottom of the screen. Finally, *click* on .
7. *Click* on . Both forms of triglycine are now in the same document. Enter the **Calculations** dialog and specify an **Energy** calculation using the **Hartree-Fock 3-21G** model. Do not change **Total Charge** (from **Neutral**) for the zwitterionic form of triglycine. Even though its Lewis structure incorporates formal + and - charges, the molecule is neutral. Also make certain that **Global Calculations** is checked. *Click* on **Submit** at the bottom

---

\* This is necessary to maintain the idealized helix structure.

of the dialog. Name the job *triglycine*.






8. After the calculations complete, bring up the spreadsheet, *click* on the header cell of an empty column, and *click* on **Add...** at the bottom of the spreadsheet. Select both **E** and **Eaq** from the list of available properties at the top of the dialog which appears, and **kJ/mol** from the **Energy** menu, and then *click* on **OK**.

Which form of triglycine, zwitterion or non-zwitterion, is favored in the gas phase? Why? Does solvent better stabilize the neutral or zwitterionic form? Why?



9. Remove *triglycine* and remaining dialogs from the screen.

### 10 to 14 optional

An electrostatic potential map for the zwitterionic form of a larger polypeptide, decaglycine, clearly shows charge separation, and explains why such a structure is not favored in the gas phase.

10. Select **New** from the **File** menu () and bring up the peptide model kit. *Click* seven times on  (glycine) from the selection of amino acid codes. \* A sequence of ten glycines has now been requested. Make certain that  **$\alpha$  Helix** is selected. *Click* anywhere on screen. *Click* on **Terminate** in the model kit and select **CO<sub>2</sub><sup>-</sup>** and **NH<sub>3</sub><sup>+</sup>** for C and N terminating groups, respectively, and then *click* on **OK**. *Click* on  and then *click* on **Freeze Heavy** at the bottom of the screen. *Click* on , and then *click* on .
11. Enter the **Calculations** dialog and specify **Energy** using the **AM1 Semi-Empirical** model. *Check* **Converge** (near the bottom right of the dialog). Semi-empirical calculations on molecules of this complexity often have trouble converging and this option invokes a number of special convergence techniques. *Click* on **OK**.
12. Bring up the **Surfaces** dialog and request calculation of an electrostatic potential map. *Click* on **Add** and select **Electrostatic Potential Map** from the menu.

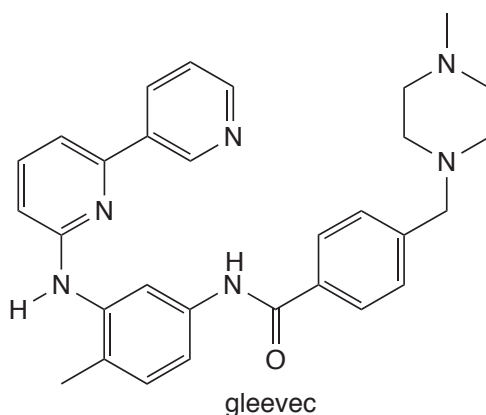
---

\* This assumes that the box at the top of the peptide model kit already contains a sequence of three glycines. If you are starting from scratch, first make certain that  is selected and **Sequence** is checked, and then *click* nine times on .



13. Submit the job. Name it *decaglycine zwitterion*. When completed, again enter the **Surfaces** dialog. *Click* on the line *Electrostatic Potential Map*. Examine the electrostatic potential map. Choose **Properties** from the **Display** menu and then *click* on the electrostatic potential map. *Click* on the selection box next to **Legend** to provide a graphical legend for the value of electrostatic potential. Recall that colors near red correspond to areas on the surface that are negatively charged (attracted to a point-positive charge), while colors toward blue designate regions that are positively charged. Intermediate colors (greens) represent neutral regions. What do you conclude about the environment in which decaglycine wishes to exist?
14. Remove decaglycine and any remaining dialogs from the screen.

## Gleevec. Protein Bound vs. Free Conformer



Gleevec (Glivec, Imatinib) is a protein kinase inhibitor used in anticancer therapy. It specifically targets a protein kinase coded for in the rogue gene *bcr-abl*. Several crystal structures for gleevec docked in protein kinases have appeared in the literature and are available in the PDB. In this tutorial, you will examine one of them\* to establish whether or not the conformation of gleevec in the protein is identical (or similar) to that of the free (gas-phase) molecule, and if not, what energy penalty is to be paid to adopt the protein-bound conformation. The next tutorial will return to this same structure and extract a pharmacophore (gleevec's footprint).

1. Retrieve a protein structure from the PDB. Select **Access PDB Online...** from the **File** menu. Type *1opj* into the box to the right of PDB ID and *click* on **Open**.

If you are not online, you may skip this step as the structure is available as *1opj* in the *polypeptides to proteins* sub-directory (*tutorials* directory). Open this document. For Windows, it is found in *Program Files/Wavefunction/Spartan10 or Program Files(x86)/Wavefunction/Spartan10*. For Linux, it is found in the directory where *Spartan* was installed. For Macintosh, it is located in the *Spartan.10.full* folder on the installation CDROM.


The protein will be represented by a ribbon model. Red sections correspond to parts of the chain that are  $\alpha$  helices, while blue

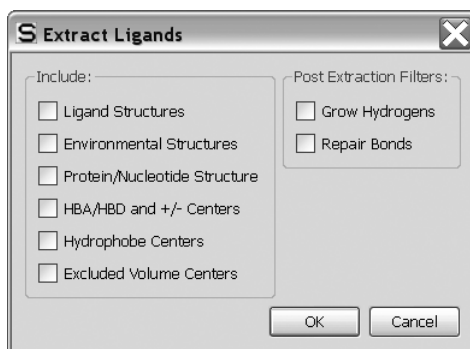
\* B. Nager *et al*, Cell, **112**, 859 (2003): PDB identification 1opj.

regions correspond to parts that are  $\beta$  sheets (the remaining regions are green). Note that these assignments are taken from the PDB file. A Ramachandran plot (the  $\phi$  and  $\psi$  torsional angles connecting the individual amino acids) allows you to see the clustering of  $\alpha$  helices and  $\beta$  sheets. To display, select **Ramachandran Plot** from the **Model** menu. The colors of the dots in the resulting plot map one-to-one with the colors of the ribbon model. When you are done, go back to the **Model** menu and deselect **Ramachandran Plot** to remove the plot.

The colors on the ribbon display may be modified, by selecting **Configure...** from the **Model** menu and *clicking* on the **Ribbons** tab. Coloring **By Residue** will give each amino acid its own unique color.

The protein you have brought in from PDB incorporates two sets of bound molecules, depicted in the model as two sets of translucent spheres. Gleevec is the larger of the two incorporated molecules and is designated by its **PDB HET** code **STI**.

2. Select **Extract Ligands** from the **Search** menu (or *click* on the  icon at the top of the screen). Select (*click* on) one of the gleevec molecules inside the protein structure and *click* on the **Extract Ligands** button at the bottom right of the screen. This leads to the **Extract Ligands** dialog.




**Check Ligand Structures** under **Include** and **Grow Hydrogens** under **Extraction Filters** inside the dialog and *click* on **OK**.

3. A ball-and-spoke model for gleevec will appear on screen. (You no longer need the protein; select and close it.) Put a copy of gleevec on the clipboard. With the molecule selected, select **Copy** from the **Edit** menu.
4. Select **New Molecule** (not **New**) from the **File** menu. Select **Clipboard** from the organic model kit and *click* anywhere on screen. *Click* on . You now have two copies of gleevec in a single document. You can step between them using the step buttons at the bottom left of the screen.
5. Select the first copy and then select **Calculations...** from the **Setup** menu. **Equilibrium Geometry** from the top left menu to the right of **Calculate** and **Molecular Mechanics** and **MMFF** from the two bottom menus. The conformation will remain unchanged, but bond lengths and angles will be optimized.

This is necessary in order to establish the energy difference between protein-bound and free conformations of gleevec. The resolution of protein X-ray crystallography is not adequate to establish bond lengths and angles to chemical accuracy.


Remove the checkmark from **Global Calculations** at the bottom of the dialog to signify that your choice (geometry optimization) only applies to this copy of gleevec. *Click* on **OK** (not on **Submit**).

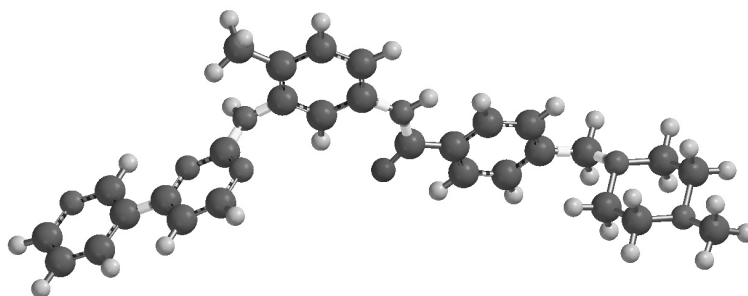6. Select the second copy of gleevec which will be used as a starting point for a conformational search. First, remove some degrees of conformational freedom in order to shorten the computer time for the conformational search. *Click* on . Your model will be augmented with yellow cylinders and circles to designate flexible bonds and rings, respectively. Remove the marker from the bond connecting the phenyl ring to the carbonyl end of the amide group, by *double clicking* on it. Also, remove the (two) circles on the piperazine ring by *double clicking* on each in turn. The conformation of the piperazine ring will not change during the search. The model should now appear as below.




7. Re-enter the **Calculations** dialog. Select **Equilibrium Conformer** from the top left menu to the right of **Calculate** and (as before) **Molecular Mechanics** and **MMFF** from the two bottom menus. *Click* on **Submit**, and provide the name *gleevec protein bound vs. free conformer*.
8. When the calculations have completed, bring up the spreadsheet (**Spreadsheet** under the **Display** menu), *click* inside leftmost the cell for the second molecule (the one on which a conformational search was performed) then *click* inside the header cell for a blank data column. *Click* on **Add** at the bottom of the dialog and select **Rel. E** from the scroll box in the dialog that appears and **kJ/mol** from the **Energy** menu. The relative energy of the protein bound conformer is provided in the spreadsheet. Is it close to zero meaning that the protein-bound conformer will be present in significant amount in a sample of free gleevec?
9. Close *gleevec protein bound vs. free conformer* and any dialogs remaining on screen.

## Gleevec. Making a Pharmacophore from PDB

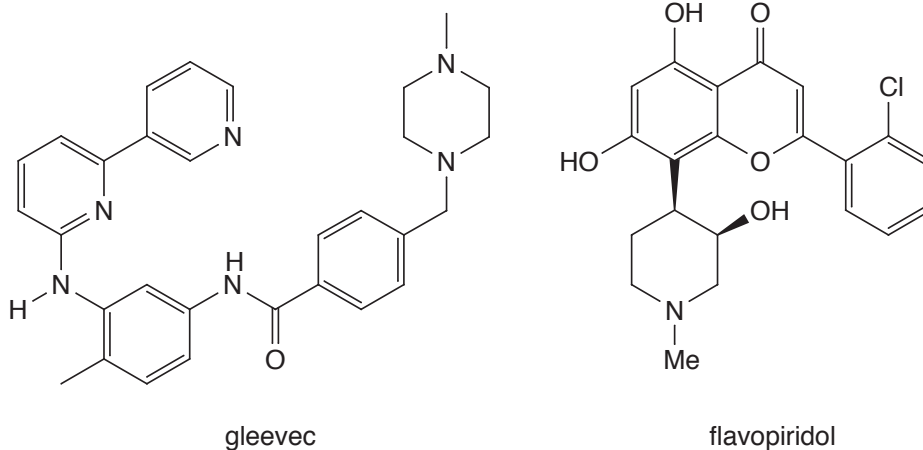


In the previous tutorial, you extracted gleevec from a PDB structure of a protein kinase in order to see if the bound conformer was the same as the free conformer. Here, you will start from the same protein structure\*, but instead of abstracting the molecular structure of gleevec, you will abstract its “footprint”, that is, the locations from which it interacts with its protein host in terms of hydrogen-bond or charge-charge interactions (non-steric contacts). This information, together with knowledge of either the steric requirements of the guest or the space occupied by the protein (excluded volumes), constitutes a structural pharmacophore, that is, a template for scrutiny of other possible guests. In this tutorial, you will use a subset of the non-steric contacts together with excluded volume elements to see to what extent this pharmacophore fits another protein kinase inhibitor, flavopiridol.

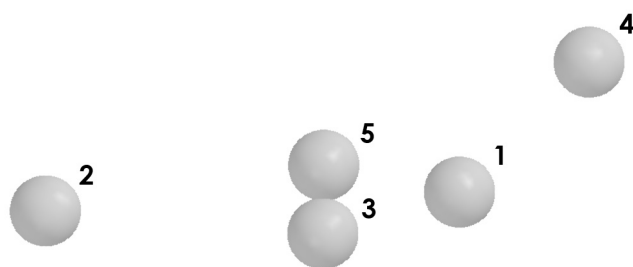
1. Either obtain the PDB file *1opj* from the PDB database (see previous tutorial) or open *1opj* from the *polypeptides to proteins* sub-directory (*tutorials* directory)\*\*. Select **Ligands** from the **Search** menu (or *click* on ) , *click* on the model for gleevec inside the protein structure (displayed as a set of transparent spheres and look much like a space-filling model) to select

\* B. Nager *et al.*, *Cell*, **112**, 859 (2003): PDB identification 1opj.



\*\* For Windows, this is found in *Program Files/Wavefunction/Spartan10* or *Program Files(x86)/Wavefunction/Spartan10*. For Linux, this is found in the directory where *Spartan* was installed. For Macintosh, this is located in the *Spartan.10.full* folder on the installation CDROM.

it, and *click* on the **Exact Ligands** button. Inside the **Extract Ligands** dialog, *check* both **HBA/HBD** and **+/- Centers** and **Excluded Volume Centers**. If any other items are selected, make certain to *deselect* them before you *click* on **OK**.

2. You are finished with *Iopj* and can close it () . Also, simplify the display of the extracted information by removing the excluded volume elements from view. Select **Configure...** from the **Model** menu and *click* on the **CFDs** tab. *Uncheck* **Excluded Volumes** and *click* on **OK**.
3. The model that now appears comprises five blue spheres representing non-steric contacts between gleevec and its protein host. Turn on labels. Select **Configure...** from the **Model** menu, *click* on the **Labels** tab, *check* **CFD Labels** and *click* on **OK**.



These have been assigned based on the observation that each is within a predefined distance of (one or more) complementary groups from amino acids in the protein host. Because hydrogen positions are not established (in the experimental X-ray structure) it may be impossible to say whether a particular site on the guest is acting as a (hydrogen-bond) donor or acceptor. In fact, it may not even be possible to say whether a particular site is protonated or deprotonated (although pH may dictate this). This being the case, all non-steric contacts have initially been given four CFD definitions: hydrogen-bond acceptor (**HBA**), hydrogen-bond donor (**HBD**), positive ionizable (+) and negative ionizable (-). Examination of the structure of gleevec suggests more focused assignments: **1** (**HBA**, +), **2** (**HBA**), **3** (**HBD**), **4** (**HBA**) and **5** (**HBA**). (See above figure for numbering.)

4. You will use only three elements of this pharmacophore (together with all the excluded volume elements) to assess whether or not flavopiridol is likely to fit into the same host environment as gleevec. Select **Properties** from the **Display** menu. Then select **Set Similarity Centers** from the **Geometry** menu (or *click* on ). Make certain that **CFD** (not **Structure**) is selected in the lower right hand corner of *Spartan*'s interface. *Click* on the pharmacophore element designated as **1** (see figure on previous page). A violet circle will surround the element designating that it is to be used in the similarity analysis. If you make a mistake and select the wrong element, *click* on the circle and the designation will be removed. Inside the **CFD Properties** dialog, remove checkmarks from **HBD** and – **Ionizable** (you are left with **HBA** and + **Ionizable**). Repeat the process for the pharmacophore element **3**, designating it as **HBD** only and **5**, designating it as **HBA** only. When you are done, only three of the five pharmacophore elements should be marked. *Click* on .
5. Select **Calculations...** from the **Setup** menu and then **Similarity Analysis** from the top left menu to the right of **Calculate**. *Click* on **Add Library...** and then locate and select *flavopiridol* in the *polypeptides to proteins* sub-directory (*tutorials* directory) and *click* on **OK**. Make certain that **CFD** is selected next to **Use:** in the lower right of the **Calculations** dialog. *Click* on **Submit**. Name the document *flavopiridol fit to gleevec pharmacophore*.
6. When the similarity analysis has completed, bring up the **Similarities** dialog (**Similarities...** from the **Display** menu). Sort the hits according to similarity score (*click* on **Score** at the top of the dialog). One after the other, *click* on the top scoring entries in the box to the right of the dialog and examine the fits in the window at the left.
7. Close any open documents and dialogs on screen.





# Chapter 10

## Inorganic and Organometallic Molecules

---

*This chapter shows how to construct inorganic and organometallic molecules using **Spartan**'s inorganic model kit. It also identifies quantum chemical models that appear to be suitable for calculations on molecules incorporating transition metals.*

Many molecules are made up of a relatively few elements and obey conventional valence rules. They may be easily built using the organic model kit. However, others cannot be assembled with this model kit either because they incorporate other elements, or do not conform to normal valence rules or they involve ligands. Important among these are inorganic and organometallic compounds involving transition metals. They may be constructed using the inorganic model kit.

Transition-metal inorganic and organometallic compounds may also require different quantum chemical methods from those that have proven to be satisfactory for organic molecules. In particular, Hartree-Fock and MP2 models are known to produce poor results where transition metals are involved. The PM3 and PM6 semi-empirical models have been parameterized for most transition metals, and generally provide a good account of equilibrium geometries. Density functional models are also successful with regard to the calculation of geometries of molecules incorporating transition metals. There is no reason to believe that density functional models should not provide a solid account of the energies of reactions involving transition-metal systems, although there are very few experimental thermochemical data with which to compare.


Finally, note that there is almost no experience with molecules incorporating lanthanides and actinides, and few practical quantum chemical methods are available.

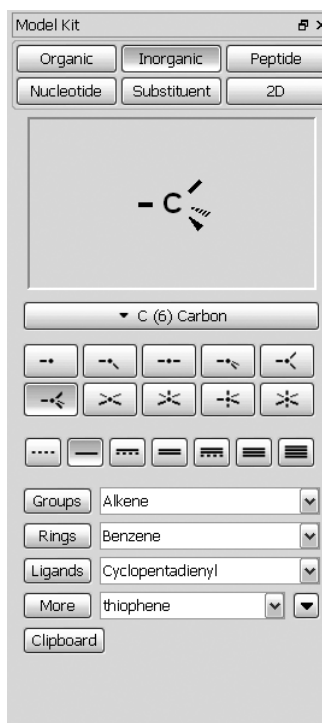
The tutorials in this chapter illustrate construction of inorganic and organometallic molecules, as well as the selection of suitable quantum chemical models.

## Sulfur Tetrafluoride



Sulfur tetrafluoride cannot be constructed using *Spartan's* organic model kit. This is because sulfur is not in its normal bent dicoordinate geometry, but rather in a trigonal bipyramid geometry with one of the *equatorial* positions vacant. However, the molecule can easily be made using the inorganic model kit.

1. Bring up the inorganic model kit by *clicking* on  and then *clicking* on the **Inorganic** tab at the top of the model kit.





The inorganic model kit comprises an atom bar (*clicking* on which bring up the *Periodic Table*<sup>\*</sup>) followed by a selection of atomic hybrids, then bond types, and finally **Rings**, **Groups**,





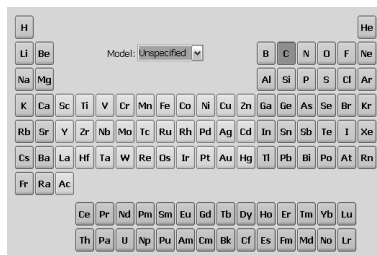\* Not all methods are available for all elements listed. Elements for which a specific method (selected in the **Calculations** dialog) are available will be highlighted following selection of a theoretical model from the **Model** menu that appears in the center of the *Periodic Table*.

**Ligands**, **More** and **Clipboard** menus (all except for **Ligands** are the same as found in the organic model kit).

2. Click on the atom bar to bring up the *Periodic Table*.



Select (click on) **S** in the *Periodic Table* and the five coordinate trigonal bipyramid structure  from the list of atomic hybrids. Click on screen. A trigonal bipyramid sulfur will appear at the top of the model kit.

3. Again, click on the atom bar, select **F** in the *Periodic Table* and the one-coordinate entry  from the list of atomic hybrids. One after the other, click on both *axial* free valences of sulfur, and two of the three *equatorial* free valences.
4. It is necessary to delete the remaining free valence (on an *equatorial* position); otherwise it will become a hydrogen. Click on  and then click on the remaining *equatorial* free valence.
5. Click on . Click on  to remove the model kit.
6. Select **Calculations...** from the **Setup** menu. Specify calculation of **Equilibrium Geometry\*** using the **Hartree-Fock 3-21G** model\*\* and click on **OK**.
7. Select **Surfaces** from the **Setup** menu, click on **Add** and select **HOMO** from the menu that results. Leave the dialog on screen.

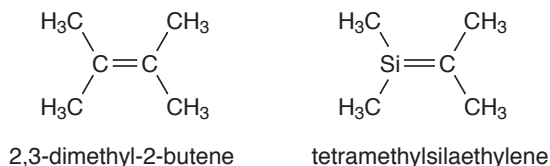
---

\* It should be noted that were an incorrect geometry specified at the outset, optimization would lead to the correct structure, as long as the starting geometry possessed no symmetry ( $C_1$  point group). Thus, square planar  $SF_4$  in  $D_{4h}$  symmetry would remain square planar, while an almost square planar structure (distorted only slightly from  $D_{4h}$  symmetry to  $C_1$  symmetry) would collapse to the proper structure.





\*\* The 3-21G basis set is extended by a set of d-type valence functions for sulfur (and other second-row elements).

8. Submit the job. Accept the name *sulfur tetrafluoride*. When completed, select **Properties (Display menu)** and *click* on an atom, for example, sulfur. Three different atomic charges will appear in the (**Atom Properties**) dialog (corresponding to different methods for establishing atomic charge). Of these, the procedure based on fitting the electrostatic potential is generally considered to be the best. Are the electrostatic charges consistent with covalent or ionic bonding?
9. *Check* the box to the left of **HOMO** inside the **Surfaces** dialog. Is the highest-occupied molecular orbital in sulfur tetrafluoride consistent with the notion that sulfur is surrounded by six electron pairs? Elaborate.
10. Remove *sulfur tetrafluoride* and any remaining dialogs from the screen.

## Reactivity of Silicon-Carbon Double Bonds




With the exception of so-called phosphorous ylides, compounds incorporating a double bond between carbon and a second-row element are quite rare. Most curious perhaps is the absence of stable compounds incorporating a carbon-silicon double bond. This can be rationalized by using local ionization potential and LUMO maps to compare the reactivities of olefins and silaolefins.

- Build 2,3-dimethyl-2-butene with the organic model kit, minimize () and copy onto the clipboard. (*Right click* on the molecule and select **Copy** from the resulting menu.) Select **New Molecule** (not **New**) from the **File** menu, *click* on **Clipboard** inside the organic model kit and *click* anywhere on screen. *Click* on the **Inorganic** tab (at the top of the organic model kit) to bring up the inorganic model kit. *Click* on the atom bar and select **Si** from the *Periodic Table*. *Double click* on one of the  $sp^2$  carbons to substitute a silicon atom for a carbon atom and to make tetramethylsilaethylene. *Click* on  and then on . You are left with 2,3-dimethyl-2-butene and tetramethylsilaethylene in the same document.
- Use the structures in SMD to save computer time. *Click* on  to the left of the molecule name (either 2,3-dimethyl-2-butene or tetramethylsilaethylene) at the bottom of the screen, select **3-21G** from the entries in the **Preview** dialog and *click* on **Replace All**.
- Select **Calculations** from the **Setup** menu. Specify calculation of **Energy** using the **Hartree-Fock 3-21G** model. *Click* on **OK**.
- Select **Surfaces** from the **Setup** menu. First request a local ionization potential map. *Click* on **Add** and select **Local Ionization Potential Map** from the menu. Then request a LUMO map. *Click* on **Add** and select **|LUMO|map** from the menu.


5. Submit the job. Name it *tetramethylsilaethylene*. (Even though both molecules are in SMD, *Spartan* does not attempt to assign names to documents that contain more than one molecule.) When completed, bring up the spreadsheet (**Spreadsheet** under the **Display** menu) and *check* the box at the far left for each entry in the spreadsheet. This allows simultaneous display of both molecules. Select **Coupled** (**Model** menu) to uncouple the motions of the molecules. Position the two molecules side by side on screen.
6. Inside the **Surfaces** dialog, select **Local Ionization Potential Map**. Compare local ionization potential maps for the olefin and silaolefin, recognizing that the more red the color, the lower the ionization potential and the more susceptible toward electrophilic attack. Which molecule do you conclude is likely to be more reactive? Turn the local ionization potential maps off by again selecting it. Select **|LUMO|map** to turn the LUMO maps on. Here, the more blue the color, the greater the concentration of the LUMO and the more susceptible toward nucleophilic attack. Which molecule do you conclude is likely to be more reactive? Close *tetramethylsilaethylene* when you are finished.


### 7 - 10 Requires Access to the Cambridge Structural Database

Identify molecules in the Cambridge Structural Database (CSD) that incorporate a carbon-silicon double bond.

7. Build silaethylene,  $\text{H}_2\text{Si}=\text{CH}_2$ , by starting from ethylene and (using the inorganic model kit) substituting a silicon for one of the carbons.
8. Select **Define Structure Query** from the **Search** menu (). *Click* on all four free valences in silaethylene. They will be marked by orange cones.
9. Bring up the **Cambridge Structural Database (CSD)** dialog by first selecting **Open Databases** from the **Search** menu and then *clicking* on the **CSD** tab. *Click* on the **Search** button at the bottom of the dialog. In a few seconds, hits will begin to appear



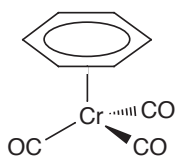
in the box at the right of the dialog. Examine each in turn (in the window at the left of the dialog) by *clicking* on its CSD reference code (**REFCODE**). Identify a structure that appears to be representative for closer scrutiny. Before you retrieve it, *click* on  to the right of the **Retrieve** button at the bottom of the dialog. Make certain that all five boxes under **Post Retrieval Filters** in the **Retrieve Options** dialog are *checked*. *Click* on **OK**. *Click* on the selected **REFCODE** and then *click* on the **Retrieve** button. Dismiss the **CSD** dialog.

10. *Click* on . Except for hydrogens, all atoms in the structure you have retrieved will be fixed and energy minimization will only change hydrogen positions. This is to compensate for the fact that hydrogens are generally not well located in X-ray crystal structure determinations.

Change the display to a space-filling model. Pay particular attention to the incorporated silicon-carbon double bond. Is it accessible in all of the compounds or is it encumbered by the surrounding groups? What, if anything, do your observations suggest about why this particular compound can be (has been) made?

11. Close silaethylene and any open dialogs.


## Benzene Chromium Tricarbonyl



benzene chromium tricarbonyl


Comparison of electrostatic potential maps for benzene chromium tricarbonyl and free benzene will allow you to classify  $\text{Cr}(\text{CO})_3$  as an electron-donor or an electron-acceptor substituent.

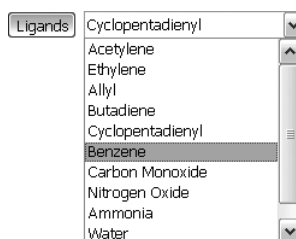1. Click on and bring up the inorganic model kit. Click on the atom bar and select **Cr** from the *Periodic Table*. Select the four-coordinate tetrahedral structure from the list of atomic hybrids. Click anywhere on screen.
2. Click on **Ligands** in the model kit, select **Benzene** from the menu of available ligands.



Click on one of the free valences on the four-coordinate chromium center.

3. Select **Carbon Monoxide** from the **Ligands** menu, and click on the remaining (three) free valences on chromium. Click on to produce a refined structure.
4. Select **New Molecule** from the **File** menu. Click on **Rings**, select **Benzene** and click on the screen. Click on and then on . The document now contains both benzene chromium tricarbonyl and benzene.
5. Select **Calculations...** (**Setup** menu). Specify calculation of **Equilibrium Geometry** with the **Semi-Empirical PM3** model. Make certain that **Global Calculations** (at the bottom of the

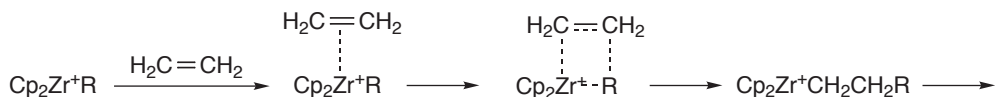
- dialog) is checked. *Click* on **OK**.
6. Select **Surfaces** (**Setup** or **Display** menu). *Click* on **Add** and choose **Electrostatic Potential Map** from the menu. Make certain that **Global Surfaces** is checked.
  7. Submit the job. Name it *benzene chromium tricarbonyl*. When completed bring up the spreadsheet (**Spreadsheet** from the **Display** menu), and *check* the box to the right of the label for both entries. This allows the two molecules to be displayed simultaneously on screen. If **Coupled** (**Model** menu) is *checked*, remove the checkmark by selecting it. The two molecules may now be moved independently. Orient each molecule so that you can clearly see the benzene face (exposed face in the case of the organometallic).
  8. *Check **Electrostatic Potential Map** in the **Surfaces** dialog. Compare electrostatic potential maps for both free and complexed benzene, with attention to the exposed benzene face.\* Does the  $\text{Cr}(\text{CO})_3$  group donate or withdraw electrons from the ring? Would you expect the aromatic ring in benzene chromium tricarbonyl to be more or less susceptible to electrophilic attack than free benzene? More or less susceptible to nucleophilic attack?*
  9. Both benzene and benzene chromium tricarbonyl are available from the Spartan Molecular Database. Make a copy of *benzene chromium tricarbonyl* (); name it *benzene chromium tricarbonyl density functional*. With either molecule in the copy selected, *check **B3LYP/6-31G\**** from the list of available models in the **Preview** dialog and *click* on **Replace All**. Inside the **Calculations** dialog, specify an **Energy** calculation using the **B3LYP/6-31G\* Density Functional** model. Submit. When completed, examine the electrostatic potential maps. Are they qualitatively similar to those from the PM3 calculations?
  10. Remove all molecules and any remaining dialogs from the screen.

---


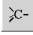
\* Electrostatic potential maps (as well as other maps) for molecules in a group will be put onto the same (color) scale. This allows comparisons to be made among different members.


## Ziegler-Natta Polymerization of Ethylene

Ziegler-Natta polymerization involves a metallocene. This first complexes an olefin, which then inserts into the metal-alkyl bond.





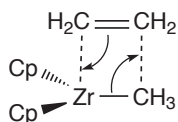In this tutorial, you will use PM3 calculations to obtain a transition state for insertion of ethylene into  $\text{Cp}_2\text{ZrCH}_3^+$  and, (optionally), estimate the activation energy using density functional calculations.

1. Bring up the inorganic model kit. *Click* on the atom bar and select **Zr** from the *Periodic Table*. Select  from the list of hybrids and *click* on screen. Select **Cyclopentadienyl** from the **Ligands** menu and *click* on two of the free valences on zirconium.
2. Move to the organic model kit. Select  $\text{sp}^3$  carbon  and *click* on the remaining free valence on zirconium. Select **Alkenyl** from the **Groups** menu, *click* on the **Insert** button at the bottom of the model kit or hold down the **Insert** key (**option** key on Mac) and *click* anywhere on screen or *double click* in a blank area on screen.
3. Orient the two fragments ( $\text{Cp}_2\text{ZrCH}_3$  and ethylene) as shown below:




(To move the two fragments independently, hold down the **Ctrl** key.)

4. Select **Guess Transition State** from the **Search** menu (). *Click* on the ZrC (methyl) bond, hold down the **Shift** key and, one after another, *click* on the methyl carbon and on one of the ethylene carbons. Next, *click* on the ethylene double bond, hold down the **Shift** key and, one after another, *click* on the other ethylene carbon and on zirconium. *Click* on  at the bottom right of the screen. In a few seconds, a guess at the transition state appears.

5. Bring up the **Calculations** dialog. Specify a **Transition State Geometry** calculation using the **PM3 Semi-Empirical** model. Change **Total Charge** to **Cation** and *check IR* to the right of **Compute**. *Click* on **Submit** and name the job **Cp2ZrMe cation + ethylene**.
6. When the job has completed, bring up the **IR** dialog (**Spectra** from the **Properties** menu and *click* on the **IR** tab) and *click* on the imaginary frequency at the top left of the dialog. Would you describe the process as concerted or occurring in discrete steps?

### 7 to 9 optional

7. Perform B3LYP/6-31G\* density functional energy calculations using PM3 geometries to obtain a better estimate for the energy barrier for ethylene insertion. *Spartan* will automatically make use of a pseudopotential for Zr instead of the all-electron 6-31G\* basis set. Make a copy of **Cp2ZrMe cation+ethylene** (); name it **Cp2ZrMe cation+ethylene density functional**. Enter the **Calculations** dialog with this copy, and specify calculation of **Energy** using the **B3LYP/6-31G\* Density Functional** model. Remove the checkmark on **IR** (to the right of **Compute**). **Total Charge** should still be set to **Cation**. *Click* on **Submit**.
8. Build both ethylene and  $\text{Cp}_2\text{ZrCH}_3^+$  (name them *ethylene density functional* and *Cp2ZrMe cation*, respectively). For  $\text{Cp}_2\text{ZrCH}_3^+$ , start with three-coordinate trigonal Zr, and then add two cyclopentadienyl ligands and a four-coordinate tetrahedral carbon. For each, enter the **Calculations** dialog, and specify calculation of **Energy** using the **B3LYP/6-31G\* Density Functional** model. Select **PM3** from the **Start from** menu to designate use of a PM3 geometry. For *Cp2ZrMe cation density functional* (only), set **Total Charge** to **Cation**.
9. Submit both jobs. When they have completed, calculate an activation energy for the insertion reaction.
10. Remove all molecules and any remaining dialogs from the screen.

# Chapter 11

## “Dry Labs”: Using the Spartan Spectra and Properties Database

---

*The tutorials in this chapter comprise examples where structures, energies, spectra and other properties can be extracted from the Spartan Spectra and Properties Database doing away with the need for actually having to carry out new quantum chemical calculations.*

The Spartan Spectra and Properties Database (SSPD) is a collection of molecular structures, spectra and diverse properties obtained from a single theoretical model, specifically the EDF2/6-31G\* density functional model. SSPD includes many of the molecules found in the Spartan Molecular Database (SMD) as well as numerous molecules not found in SMD, but differs from (and we believe, improves on) this previous collection in several important ways. First, for the vast majority of entries, the choice of conformer is based on the T1 thermochemical recipe rather than on the MMFF molecular mechanics model. T1 is far more likely to correctly identify the lowest-energy conformer (or at least a conformer that is very close in energy to the lowest-energy conformer) than is MMFF. Second, SSPD entries provide a much wider variety of molecular properties than available in SMD, including a selection of QSAR descriptors. They also include both the infrared and NMR spectra (neither of which is available in SMD). Third, all SSPD entries include the wavefunction, allowing “on-the-fly” generation of molecular orbitals, electron densities and other graphical surfaces and the electrostatic potential and other property maps. In short, SSPD provides a rich source of diverse information for a wide selection of molecules.

The limitation of SSPD to a single theoretical model is a deliberate attempt to focus attention away from the characteristics of a particular model (and the question “which model is best”) and onto the chemistry at hand. While the choice EDF2/6-31G\* can be seen as a compromise between accuracy and practicality, in fact, EDF2/6-31G\* provides excellent results for a variety of properties that can be accurately measured, including molecular geometry and infrared and NMR spectra. Furthermore, the T1 heat of formation, a byproduct of using the T1 recipe to establish best conformer, is included in the SSPD entry. We believe that over the wide variety of molecules considered, T1 provides a reliable an indicator of internal energy as any (calculated or experimental) measure presently available.

Collections of calculated results such as the Spartan Spectra and Properties Database can serve different roles. At one end, information in a database may be used to supplement the calculation capabilities of *Spartan* and eliminate the need to perform calculations that have already been done. At the other end, is to use the database as the exclusive source of information, doing away altogether with the need for quantum chemical calculations. From the instructional perspective (teaching about molecular models and molecular modeling) this allows use of more complex examples than might otherwise be practical. The tutorials in this chapter fall into the latter category, and can be thought as the ultimate in “dry labs”.

To complete the tutorials in this chapter, SSPD needs to be installed on *Spartan* (the sample version included with the program will suffice or the full version can be used). In addition, database selection menu at the bottom of *Spartan*'s main screen needs to be set to **SSPD**.

## Isomeric C<sub>5</sub>H<sub>8</sub> Dienes

In addition to 2-methyl-1,3-pentadiene, *cis* and *trans*-1,3-pentadiene (conjugated dienes) and 1,4-pentadiene (a non-conjugated diene), there are three other dienes with the same formula: 3-methyl-1,3-butadiene (1,1-dimethylallene) 1,2-pentadiene (ethyl allene) and 2,3-pentadiene (1,2-dimethylallene).

1. Build all seven C<sub>5</sub>H<sub>8</sub> dienes (place in a single document) and bring up the spreadsheet (**Spreadsheet** from the **Display** menu). With any of the dienes selected, *click* on the **Replace All** button in the **Preview** dialog at the bottom of the screen. This will substitute entries in SSPD for all seven C<sub>5</sub>H<sub>8</sub> dienes. Molecule names are available but are not used unless requested explicitly. To use the names in SSPD, *click* on the top left header cell (**Label**) in the spreadsheet, hold down the right mouse button and select **Rename Selected Using SSPD** from the menu. Bring up **Molecule Properties** dialog (**Display** menu) and *click* on the **P** buttons to the left of both **Expt. Heat** and **T1 Heat**. This enters the experimental heat of formation (from the NIST database) as well as the value obtained from the T1 thermochemical recipe into the spreadsheet. The T1 heat of formation should provide a more accurate account of relative stability than available from the EDF2/6-31G\* model itself.

The T1 model is presently defined only for uncharged, closed-shell molecules comprising H, C, N, O, F, Si, P, S, Cl and Br only. T1 heats of formation are not available for the other molecules in SSPD (~2% of the total).





2. First verify that the T1 model correctly reproduces the ordering of experimental heats. Then use the T1 heats to answer the following questions:

What is the difference in stability (heat) between the best of the three conjugated dienes and 1,4-pentadiene? What is the approximate energy gain due to conjugation? What is the difference in stabilities between the best of the substituted allenes



and the best of the conjugated dienes?

### 3 Optional


3. There are several other  $C_5H_8$  isomers in addition to the dienes discussed above, a number of which will be in SSPD (many more in the full version than in the sample). With any of the dienes selected, bring up the **Databases** window (**Open Databases** from the **Search** menu or *click* on  and *click* on the **SSPD** tab). *Click* on  to the right of the **Search** button in the **SSPD** dialog and select **Isomer** under **Search By** in the **Search Options** menu that results. *Click* on the **Search** button. Select several (or all) of the non-diene isomers, by *clicking* on them in the list at the right of the dialog while holding down the **Ctrl** key. *Click* on  to the right of the **Retrieve** button, then *click* on **Current Document** in the dialog that results. Remove the **Databases** window by selecting **Close Databases** (that has replaced **Open Databases**) from the **Search** menu or *clicking* on . Bring up the spreadsheet. Add both experimental and T1 heats of formation to the spreadsheet. Note that not all the experimental heats of formation are available. In fact, experimental heats of formation exist for only around 2,000-3,000 molecules. According to the T1 model, are any of the non-diene isomers more stable than the best diene?
4. Close any molecules and dialogs that are open on screen.

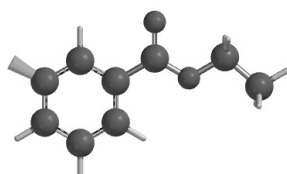
## Using Infrared Spectroscopy to Identify an Unknown Ethyl Benzoate Derivative


Unlike NMR spectra (at least carbon spectra) infrared spectra are rich in detail. Some chemists would view this as a curse, in that the complexity of an infrared spectrum masks its connection with the three-dimensional molecular structure. Others might see complexity as a blessing, providing a detailed pattern or fingerprint with which to identify a molecule. In fact, this difference reflects the different ways that chemists actually use NMR and infrared spectroscopy. NMR is most commonly employed to assign structures of new molecules whereas IR is most commonly used to identify previously characterized (known) molecules.

Extensive databases of experimental infrared spectra exist, and sophisticated programs have been developed to search these databases for matches to unknown spectra. It turns out that infrared spectra calculated from density functional models and adjusted to account for known systematic errors in frequencies and for finite temperature, closely match experimental spectra. This suggests that they might be used in lieu of experimental spectra to identify unknowns. This tutorial illustrates the point, attempting to identify which of a series of *meta*-substituted ethyl benzoates contained in the Spartan Infrared Database (SIRD) best matches the infrared spectrum of an unknown molecule.

SIRD uses the infrared spectra contained in the Spartan Spectra and Properties Database (SSPD). The fact that there are two different tabs simply reflects the different search paradigms (spectra matching for SIRD vs. substructure matching for SSPD).

1. Build ethyl benzoate. Then select **Define Structure Query** from the **Search** menu (  ) and *click* on one of the free valences on the phenyl ring that is *meta* to the ester group.



2. Select **Open Databases** from the **Search** menu (  ) and *click* on the **SIRD** tab at the top of the dialog that results. Size the dialog to occupy the greater part of the screen. *Click* on **Select Spectrum**, move to the *using SSPD* subdirectory under the *tutorials* directory and *click* on *unknown ethyl benzoate*. Make certain that **All Files** is selected in the menu to the right of the file name in the **Open** dialog. In a few seconds, the experimental infrared spectrum of the unknown will appear at the top right of the dialog. A list of the frequencies for all peaks may be found under **Unknown:** immediately to the left is the spectrum. This has been obtained by fitting the experimental spectrum to a Lorentzian function. You can obtain the frequency of a line in the spectrum using the yellow measurement bar. Position the cursor on top of the bar, hold the left button and move the mouse left of right until the bar is over the peak.
3. You could carry out your search with nothing other than the unknown spectrum (see the optional part of this tutorial), or you could restrict it using structural or other information that you have. Let's assume that we know that the unknown is a *meta*-substituted ethyl benzoate. *Click* on the **Filters** button and then *click* on **Copy Current Molecule** at the top left of the **Search Filters** dialog that results. *Click* on **OK** to exit the dialog. What you have done is to restrict the search to *meta*-substituted ethyl benzoates, that is, examine only molecules that contain this particular substructure (ethyl benzoate substituted in the *meta* position). *Click* on the **Search** button and wait for the search to complete. Hits resulting from the search appear at the bottom right of the **SIRD** dialog, in numerical order starting from the best. As you *click* on each, the corresponding calculated infrared spectrum will be superimposed onto the experimental spectrum (of the unknown) at the top left of the display. Suggest a structure for the unknown.

#### 4 Optional

4. Remove the substructure filter and repeat the search. *Click* on the **Filters** button, *click* on the **Clear** button at the top left of

the **Search Filters** dialog and *click* on **OK**. *Click* on the **Search** button. Be patient; the search will require a minute or two. Unlike searches of experimental databases where spectra are directly compared, each of the calculated spectra in the database are adjusted for frequency scale and peak width (roughly associated with temperature) as it is compared to the unknown. Removing the substructure filter greatly increases the number of spectra that need to be examined. Does the *meta* substituted ethyl benzoate that came in at the top in the previous search, also come in at or very near the top in the full search?

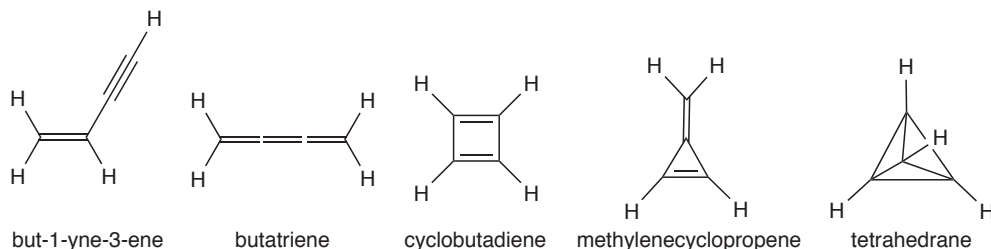
A search of the full database will require more computer time than the search of the 5,000 molecule sample database supplied with *Spartan*, simply because it contains an order of magnitude more molecules. However, the search is actually divided into preliminary and final parts, only the first of which depends on the overall size of the database.




5. Close any molecules and dialogs that are open on screen.


## Infrared Spectra of Short-Lived Molecules

While infrared spectroscopy is by no means as formidable a tool as NMR for elucidating the structures of molecules, it does offer significant advantages which can make it the method of choice in some instances. It is more sensitive, meaning that it can be used to detect transient species, for example, molecules trapped in an inert matrix at very low temperature. On the practical side, an infrared spectrometer is inherently simpler (more robust) and can be put into a much smaller package than an NMR spectrometer. This means that the spectrometer can be brought to the sample rather than the other way around. Imagine, putting an NMR spectrometer on the Mars rover! The problem is that IR spectroscopy unlike NMR typically does not provide an unambiguous structure. However, a close match to a calculated spectrum may offer convincing evidence for the veracity of a particular structure (or a mismatch convincing evidence that a structure is incorrect).

The infrared spectrum of an unknown  $C_4H_4$  isomer shows strong absorptions at 854, 1608, 2994 and 3080  $cm^{-1}$ . Possible structures include but-1-yne-3-ene, butatriene, singlet-1,3-butadiene, methylenecyclopropene and tetrahedrane.



1. Build any one of these molecules. Select **Open Databases** from the **Search** menu () and *click* on the **SSPD** tab. *Click* on  to the left of the **Search** button, select **Isomer** under **Search By:** inside the **Search Options** dialog and *click* on **OK**. *Click* on the **Search** button. Structures for all  $C_4H_4$  isomers will appear in a list at the right of the **Databases** dialog. Select all of them (*click* on each in turn while holding down the **Ctrl** key). Retrieve to a new document. *Click* on  to the left of the **Retrieve** button, select **New Document** and *click* on **OK**. *Click* on the **Retrieve**

button. Remove the **Databases** dialog (select **Close Databases** from the **Search** menu or *click* on ).

2. Give proper names to your compounds. Open the spreadsheet (**Spreadsheet** under the **Display** menu), *click* on **Label** (the leftmost header cell), then *right click* and select **Rename Selected Using SSPD** from the menu that appears.
3. Select **Spectra** from the **Display** menu and *click* on the **IR** tab at the top of the resulting dialog. Numerical values of the frequencies and infrared intensities of the selected molecule in the document will appear at the left of the dialog. Note that the frequencies have been scaled to take account of systematic errors in the calculations (primarily due to the harmonic approximation). To see unscaled frequencies, move the slider bar to the right of **Scale:** from 0.962 (the default) to 1.000. To return to scaled frequencies, *click* on the **Standard** button. Display the infrared spectrum of the selected isomer. *Click* on **Display Calculated** at the right of the dialog. Move through the list to bring up the spectra for the other isomers.

The calculated infrared spectrum of which if any of the  $C_4H_4$  isomers best fits the unknown infrared spectrum?

4. Close any molecules and dialogs that are open on screen.

## Using $^{13}\text{C}$ NMR to Distinguish Structural Isomers

NMR spectroscopy, in particular  $^{13}\text{C}$  NMR spectroscopy, is arguably the single most important tool available to an organic chemist to establish molecular structure. Obtaining an NMR spectrum is straightforward and the spectrum itself is very simple, comprising but a single line for each unique carbon the molecule. Because a mass spectrum is normally also available, the molecular formula will be known, and assignment of the NMR reduces to deciding among possible isomers. Chemical evidence (how the molecule was made) can usually be counted on to eliminate some choices and to strengthen the case for others. Still, pinning an NMR spectrum to a particular molecule can be difficult and fraught with error.

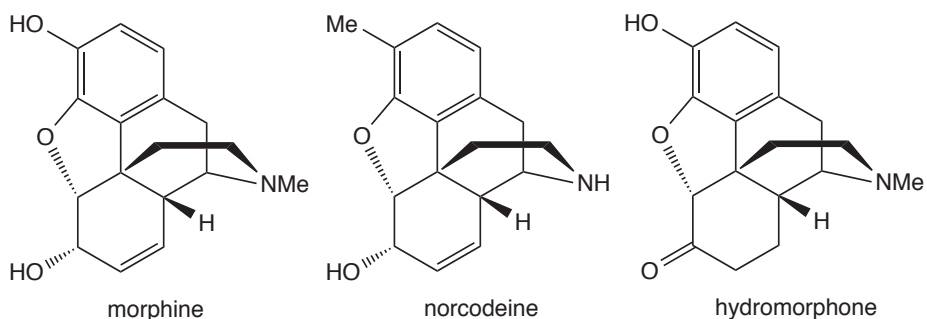
NMR assignment problems can be reduced but not altogether eliminated by requiring that both proton and  $^{13}\text{C}$  NMR spectra are consistent with a particular structure. This is routine practice. Further help comes from taking a DEPT spectrum, allowing methylene carbons to be distinguished from methyne and methyl carbons and (in conjunction with the  $^{13}\text{C}$  spectrum) from quaternary carbons. This is less common but becoming more and more routine. In principal, all ambiguity can be eliminated by cross-correlating the results of the  $^{13}\text{C}$  NMR with those of proton NMR using a variety of so-called 2D spectra, most important, the HMBC spectrum. However, 2D NMR experiments are not yet the norm, if for no other reason because they require significant instrument time on a spectrometer that may cost upward of \$1,000,000.

The fact that NMR spectra for reasonable size organic molecules may now be routinely calculated on a personal computer costing \$1000 or less raises the possibility for another tool to assist with spectral assignments. Whether this turns into reality ultimately depends on the ability of the calculations to obtain NMR spectra that are sufficiently accurate to distinguish among the different isomers.

In this tutorial we first see which of three structural isomers with molecular formula  $\text{C}_{17}\text{H}_{19}\text{NO}_2$  (morphine, norcodeine and hydromorphone) best fits the NMR spectrum of an unknown isomer, or if clear assignment is not possible. Then, we draw comparisons between calculated and experimental DEPT spectra (the former

just a different presentation of the  $^{13}\text{C}$  spectra) to see if this further clarifies the assignment.

1. Build morphine, norcodeine and hydromorphone and put in the same document.



If the name for each fails to appear at the bottom of the screen (all three isomers are in SSPD), then you have made a mistake. When you are done, *click* on the **Replace All** button at the bottom of the **Preview** dialog. This will substitute structures from SSPD for what you have built. Give proper names to the three isomers. Bring up a spreadsheet (**Spreadsheet** under the **Display** menu), *click* on **Label** (the leftmost header cell), then *right click* and choose **Rename Selected Using SSPD** from the menu.

2. Bring up the **NMR Spectra** dialog (**Spectra** under the **Display** menu) and *click* on the **NMR** tab. First make certain that the box to the left of **Boltzmann Averaging** at the bottom of the dialog is not checked (if it is, *click* inside the box to remove the checkmark). *Click* on **Draw Calculated** in the section designated  $^{13}\text{C}$  **Spectra**. Next, *check* **Local File** in the section designated **Experimental Data From**, *click* on **Draw Experimental** (in the  $^{13}\text{C}$  **Spectra** section), move to the **using SSPD** subdirectory under the **tutorials** directory and *click* on **unknown C17H19NO2 carbon NMR**. Make certain that **All Files** is selected in the menu to the right of the file name in the **Open** dialog. You will need to repeat this two more times (providing the same experimental spectrum for each of the three candidates). Calculated and experimental  $^{13}\text{C}$  spectra now appear on the same plot. Step through the three isomers. Which,




if any, of the calculated  $^{13}\text{C}$  spectra best fits the experimental spectrum? To associate lines with specific carbons in each of the calculated spectra, bring up the **Atom Properties** dialog. Select **Properties** from the **Display** menu and *click* on an atom. Alternatively, you can attach chemical shifts as labels on your structures. Select **Configure** from the **Model** menu and *click* on **Chem Shift** under **Atom** in the dialog that results. Simplify your models by changing to a tube model (**Tube** under the **Model** menu) and by removing hydrogens from the model (**Hydrogens** under the **Model** menu).

3. *Click* on both **Delete Calculated** and **Delete Experimental** buttons in the dialog (which have replaced the corresponding **Draw** buttons) to remove the  $^{13}\text{C}$  spectra. In their place, show the corresponding calculated and experimental DEPT spectra. *Click* on the **Draw DEPT**. *Click* on **Draw Experimental**, then move to the *using SSPD* subdirectory under the *tutorials* directory and *click* on *unknown C17H19NO2 morphine DEPT NMR*. Again, you will need to do this three times. Recall that lines in a DEPT spectrum that point up correspond to carbons with either one or three attached hydrogens, lines that point down to carbons with two attached hydrogens and that missing lines correspond to quaternary carbons. How many methyne and methyl carbons are in the structure of the unknown? How many methylene carbons are there? How many quaternary carbons are there? Which if any of the three isomers examined provide the same count? What is the structure?
4. Close any molecules and dialogs that are open on screen.

## Using $^{13}\text{C}$ NMR Spectra to Distinguish Stereoisomers

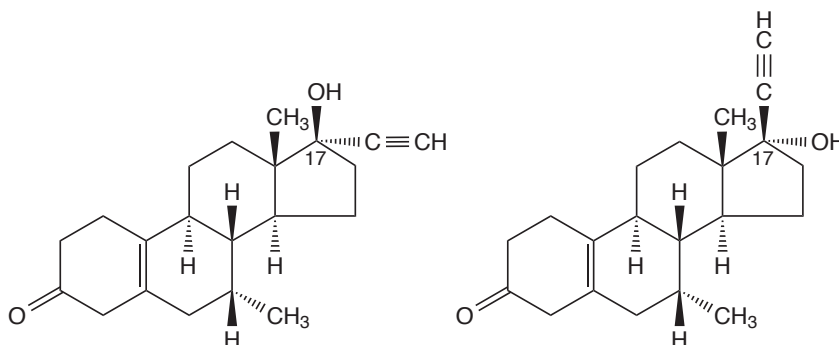As commented in the previous tutorial, the combination of 1D and 2D NMR spectra is almost always able to provide unambiguous assignment of molecular structure, including assignment of stereochemistry. However, as also noted, practical considerations (time on an expensive instrument) often make 2D experiments (in particular HMBC) the exception rather than the routine. Without the additional experiments, distinguishing stereoisomers based on their  $^{13}\text{C}$  NMR spectra may be difficult experimentally, simply because the alternatives are likely to be structurally very similar. For the same reason, it is likely to be an ideal case for calculations as comparisons between molecules with similar structures should benefit from cancellation of errors. The bottom-line question is whether or not calculated  $^{13}\text{C}$  NMR spectra good enough to be able to clearly be able to distinguish between stereoisomers? This tutorial provides an illustration.

The hydroxyl CO bond in tribolone, shown below on the left, is *gauche* to the methyl group at  $\text{C}_{17}$  (the alkyne is *anti*), whereas it is *anti* (and the alkyne is *gauche*) in the stereoisomer shown on the right. The experimental  $^{13}\text{C}$  spectrum of a molecule that has been assigned as tribolone shows a line at 12.7 ppm attributed to the methyl group at  $\text{C}_{17}$  and lines at 88.9 and 74.8 ppm attributed to the internal and external alkyne carbons, respectively. The full  $^{13}\text{C}$  spectrum may be found in the *using SSPD* subdirectory under the *tutorials* directory as *tribolone proposed carbon NMR*. Make certain that **All Files** is selected in the menu to the right of the file name in the **Open** dialog.




1. Build both stereoisomers and replace your structures with those in SSPD. Select **Properties** from the **Display** menu and

for each stereoisomer, record the calculated  $^{13}\text{C}$  shifts for the methyl group carbon and for the two alkyne carbons. Do the calculations support or refute the experimental assignment, or are they ambiguous?

2. Which stereoisomer is more stable according to the T1 heats of formation accessible from SSPD? Offer a reasonable explanation for the observed preference.
3. Close any molecules and dialogs that are open on screen.




**Formula Editor** dialog). A formula for the query appears in the box at the center of the **Formula Editor** dialog, preceded by the name of the query *F0*. Replace this name by *enol to keto*, click on the **Enter** button (to the right of the query name and formula). You have now defined the query as the difference between the T1 heats of formation of the keto and enol forms of (substituted) propen-2-ol. Leave the **Formula Editor** dialog open.

4. Select **Open Databases** from the **Search** menu () and click on the **SSPD** tab at the top of the dialog. Drag the query name *enol to keto* from the formula editor onto an empty column header at the top right of the **SSPD** dialog, and click on the **Search** button. In a few seconds, hits (molecules in SSPD formed from propen-2-ol by single or double substitution at C<sub>1</sub>) will begin to appear at the right of the dialog. Reaction energies for any entries where T1 data exists both for actual search query (the substituted enol) and for its associated keto form will appear in the column designated *enol to keto*. *Pending* signifies that an entry exists for the search query (the enol) but not for the keto form.
5. After the search has completed, position the cursor over the *enol to keto* column, hold down the right mouse button and select **Sort enol to keto** from the menu that appears. Molecules where the equilibrium most favors the keto form now appear at the top. Scroll down the list to identify molecules where the enol form is preferred.
6. Close any molecules and dialogs that are open on screen.

## Electrophilic Reactivity of Polycyclic Aromatics

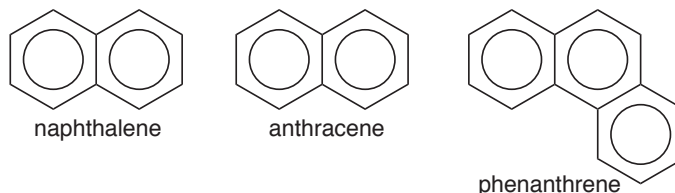
Naphthalene, anthracene and phenanthrene undergo electrophilic substitution just as does benzene. Are they likely to be more or less reactive? Local ionization potential maps should be able to provide an answer. The lower the potential on the aromatic ring, the more likely it is to give up an electron, that is, the more susceptible it will be to electrophilic attack.




1. Build benzene, naphthalene, anthracene and phenanthrene and put into a single document. With any of the four molecules selected, *click* on **Replace All** from the **Preview** dialog. Give the molecules proper names. Bring up a spreadsheet (**Spreadsheet** under the **Display** menu), *click* on **Label** (leftmost header cell), then *right click* on an empty header cell and choose **Rename Selected Using SSPD** from the menu.
2. Bring up the **Surfaces** dialog from the **Setup** or **Display** menu and make certain that **Global Surfaces** is *checked*. *Click* on **More Surfaces...** at the bottom of the dialog to bring up the **Add Surfaces** dialog. (You need to use this route to specify fixed surfaces and save computer time.) Select **density** from the **Surfaces** menu and **ionization** from the **Properties** menu (**medium** from the **Resolution** menu), *check* **Fixed** to the right of **IsoValue** and *click* on **OK**. When the graphics calculations have completed, bring up the spreadsheet (**Spreadsheet** under the **Display** menu) and *check* the box to the left of each of the four entries. The four molecules will now be simultaneously displayed, but they will be coupled. To move them independently, select **Coupled** from the **Geometry** menu (this will remove the checkmark to the left). Position the four molecules on screen such that the  $\pi$  system is visible, and *check* **Local Ionization Potential Map** inside the **Surfaces** dialog. Which of the four molecules is indicated by its local ionization

potential map to be most susceptible to electrophilic attack? Which positions on naphthalene, anthracene and phenanthrene are most susceptible?

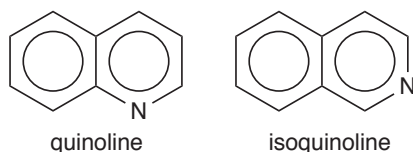3. Compare local ionization potential maps for naphthalene, quinoline and isoquinoline. Focusing on the carbon ring, which of the three is likely to be most susceptible to electrophilic attack? Which is likely to be least susceptible?



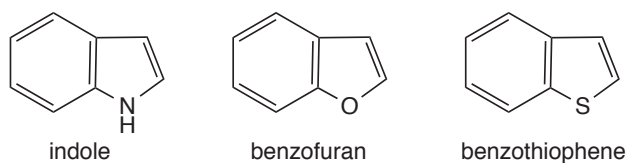4. Compare the local ionization potential map for naphthalene with those of indole, benzofuran and benzothiophene. Focusing on the carbon ring, which of the four compounds is most susceptible to electrophilic attack? Which is the least susceptible?



5. Close any molecules and dialogs that are open on screen.

# Section III

## Features and Functions

---

This section describes the functions available under the menus incorporated into *Spartan'10*, and is intended to serve as a general reference to the program. The coverage follows the order of the menus presented in *Spartan's* user interface: **File**, **Edit**, **Model**, **Geometry**, **Build**, **Setup**, **Display**, **Search**, **Options**, **Activities** and **Help**. In turn, it describes the functions and usage of each of the menu entries. Entries under the **File**, **Model**, **Geometry** and **Build** menus deal primarily with the input and construction of 3D structures and with their display and query. Entries under the **Setup** menu specify molecular mechanics or quantum chemical models as well as similarity analysis tools, designate what properties and spectra are to be obtained, request one or more graphical models and submit jobs for calculation in the background. Entries under the **Display** menu access calculated quantities while those under the **Search** menu access databases of calculated information. Entries under the **Options** menu set program defaults, designate file paths and control job queues.

This chapter provides only limited commentary about the performance and requirements of different computational methods and the utility of the different graphical models. Additional information is provided under *Topics* accessed from the *Activities* menu and to *A Guide to Molecular Mechanics and Quantum Chemical Calculations*, available as a PDF under the **Help** menu.



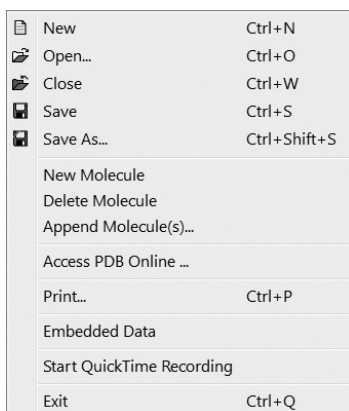


# Chapter 12

## The File Menu

---

*Operations under the **File** menu access model kits to build, edit and substitute molecules and the file system to read and write both native and non-native files, print text and on-screen graphics, access the online PDB database of protein and nucleotide structures, embed text and other files into **Spartan** documents and create QuickTime movies.*



### **New** (📄)

Brings up a model kit and clears the screen. Model kits are discussed in **Chapter 16**.

### **Open...** (📂)

Opens a file that contains all information associated with a particular molecule or pharmacophore (or list of molecules and/or pharmacophores). In addition to native (.spartan) files (documents), supported are files containing 2D drawings: CambridgeSoft

(ChemDraw) CDX (.cdx), and MDL SDF (.sdf)\*, SKC (.skc), and TGF (.tgf), files containing 3D structures: MacroModel (.mac), SYBYL Mol (.mol), SYBYL Mol2 (.mol2), and PDB (.pdb), and files containing 1D “strings”: SMILES (.smi). Finally, infrared, NMR and UV/visible spectra may be input (discussion of file formats is provided in **Appendix I**).

We use the word “file” to refer to any information contained in the “file system” (the disk on your computer) and the word “document” to refer to files that are native to *Spartan*. Thus, an MDL SDF *file* may be read into *Spartan* and then written as a *document*. Documents may contain information on a single molecule or single pharmacophore or on a collection of molecules and/or pharmacophores. This collection is loosely referred to as a “list” or as a “spreadsheet” (in the context that it is presented in a spreadsheet).

Cartesian coordinates from any of these files (except spectra files which do not necessarily contain coordinates) may be replaced by coordinates generated based on atomic connectivity using **Replace Coords.** in the **Molecule Utilities** dialog (accessible from **Properties** under the **Display** menu; **Chapter 18**). Additionally, hydrogen atoms that are missing in the input structure may be provided using **Grow Hydrogens** in the **Molecule Utilities** dialog. Non-native files are normally hidden from view, but may be seen by selecting **All Files** from the **Files of type** menu at the bottom of the dialog.

### Close (🗑️)

Closes the document containing the selected molecule, as well as the spreadsheet, graphics, plots, spectra and document specific dialogs. If the document has not previously been saved, a name is requested.

---

\* SDF files may contain information in addition to geometry, for example, the name of the molecule or its CAS number. Each quantity appears in a line following <xyz> where xyz is a character string denoting the quantity, for example, CAS to denote the CAS number.

<CAS>  
22090-26-2

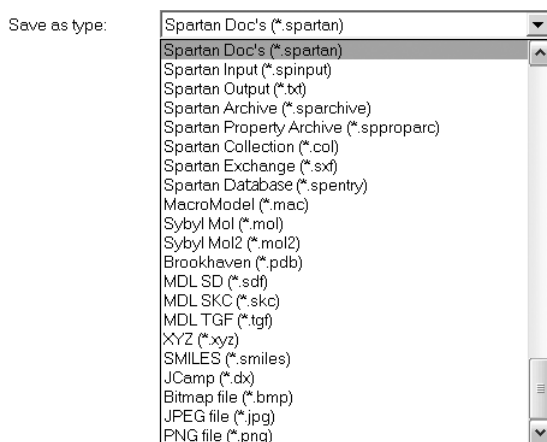
This information may be accessed in the spreadsheet by typing the character string into a header cell, for example, typing CAS will result in the CAS number.

If a previously-saved document has been altered, verification that the changes are to be saved is requested.

## Save

### Save As... (📁)

Saves the document containing the selected molecule *exactly as it appears on screen*. Opening the document will bring it on screen exactly as it was last saved. If the document has not previously been named, **Save** behaves as **Save As...** Documents may be either be saved in native format or in one of the formats listed under **Open**. In addition, bitmap (.bmp), JPEG (.jpg) and PNG (.png) graphics file formats are supported as is the Spartan Database (.spentry) format for creating custom databases (see **Appendix G**). Support is also provided for writing QuickTime (movie) files (see discussion later in this chapter). Selection is made under the **Save as type** menu in the **Save As** dialog.



## New Molecule

Brings up a model kit and clears the screen. This differs from **New**, in that the resulting molecule is appended to the end of the list associated with the molecule that is presently selected.\*

---

\* Examples are provided in the tutorials *Dienophiles in Diels-Alder Cycloadditions and Hydration of Carbonyl Compounds* in Chapter 5, *Stereospecific Diels-Alder Reactions* in Chapter 7, *Polyglycine* in Chapter 8 and *Benzene Chromium Tricarbonyl* in Chapter 10.

## Delete Molecule

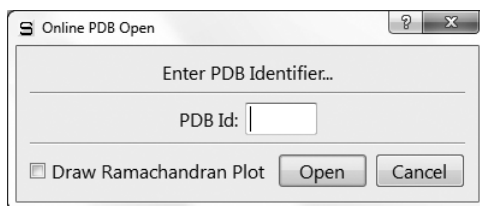
Deletes the selected molecule(s) from a document. Deleting the last molecule leads to an empty document.

## Append Molecule(s)...

Appends one or more documents onto the end of the document that contains the selected molecule. \* **Append Molecule(s)...** leads to a file browser from which one or more documents need to be selected. \*\*

## Access PDB Online...

Provides access to the online **Protein Data Bank (PDB)**\*\*\* comprising nearly 70,000 protein and nucleotide structures. Selection results in a dialog.



To access a PDB structure, enter the four character identification code in the box to the right of **By PDB ID** at the bottom of the dialog and *click* on the **Search** button. If the PDB entry contains more than one structure (typical with protein structures obtained from NMR

---

\* An example is provided in the tutorial *Carbene Additions to Alkenes* in **Chapter 7**.

\*\* Alternatively, molecules may be appended onto an existing document either by copy/paste operations using the clipboard or by *dragging* from an external window. Both require that the spreadsheet associated with the destination document be open on screen.

To copy a molecule open on screen onto the clipboard, first select (*click* on) it, and then select **Copy** from the **Edit** menu. Alternatively, *click* on its label in its spreadsheet (in the leftmost column), and then select **Copy** from the **Edit** menu. The latter permits several molecules to be selected (and copied) at once using the **Shift** and **Ctrl** keys in the usual manner. Once on the clipboard, the molecule or molecules may be moved to the destination list by *clicking* on an empty row header in the spreadsheet (for the destination document), and then selecting **Paste** from the **Edit** menu.

To copy a document from an external window, *drag* it onto the open spreadsheet (associated with the destination document) inside of **Spartan**. Several documents can be *dragged* at once using the **Shift** and **Ctrl** keys in the usual manner.


\*\*\* The web address is <http://www.rcsb.org>

spectroscopy) and/or the PDB ID yields more than one entry, all structures will be returned in a list.

PDB access will typically require a few seconds. The PDB ID will appear at the right and a ribbon model of the protein or nucleotide will appear on screen. A Ramachandran plot associated with a protein structure may either be upon initial retrieval of the PDB file by *checking* the box to the left of **Draw Ramachandran Plot** or later from **Ramachandran Plot** under the **Model** menu (**Chapter 14**).

## Print/Print Output/Print Spreadsheet

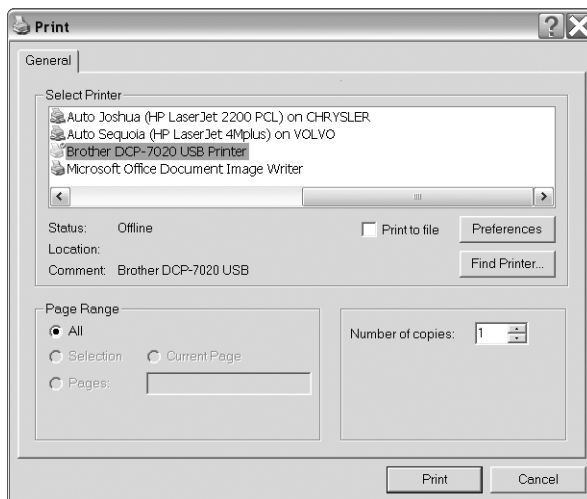
Selection leads to a dialog in order to designate a printer, specify print layout and number of copies. It also allows printing to a file.



By default, prints whatever is presently displayed on screen. To print text output, bring up and select an output window (**Output** under the **Display** menu) and select **Print Output...** (that has replaced **Print...**) from the **File** menu. To print the contents of the spreadsheet, bring up a spreadsheet (**Spreadsheet** under the **Display** menu) and select **Print Spreadsheet...** (that has replaced **Print...**) from the **File** menu. The results of a reaction energy calculation (**Reactions...** under the **Display** menu; **Chapter 18**), may be printed using **Print** from the contextual menu.

## Embedded Data

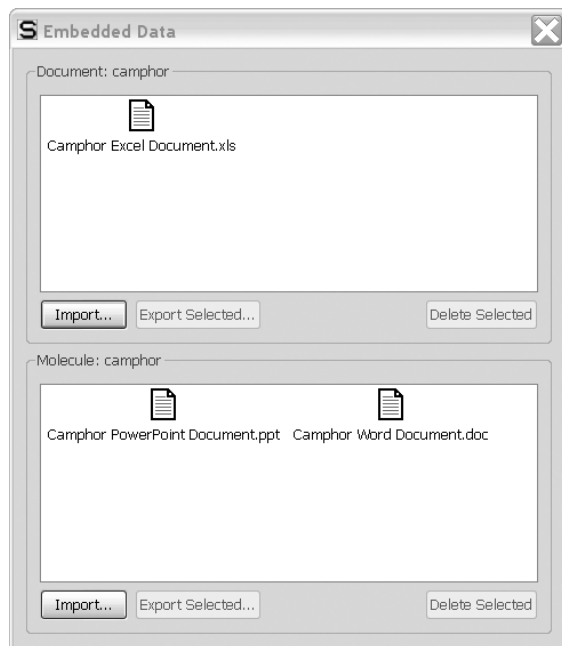This allows non-native files, for example, Word, Excel or PowerPoint files, to be embedded in a *Spartan* document. These may be associated either with an individual molecule or with all the molecules (or both). Selection results in the **Embedded Data** dialog, the upper half of which relates to document level files and the lower half of which relates to molecule level files.



A file may be imported (at the document or molecule level) by *clicking* on the appropriate **Import** button. This leads to a file browser from which a file may be selected. Once embedded in a *Spartan* document, a non-native file may be exported by first selecting (*clicking* on) it and then *clicking* on **Export Selected**. It may be deleted by first selecting it and then *clicking* on **Delete Selected**.

Note that *applications cannot be launched* from the **Embedded Data** dialog. They need to be copied outside prior to launching.

## Start/Stop QuickTime Recording

This allows QuickTime movies to be created. To start making a movie, select **Start QuickTime Recording**. Any motions of all molecules in *Spartan*'s main screen will be captured. Dialogs (including the builders) will not be captured. Note the use of "tumbling" (see **Settings Preferences** under **Preferences** in the **Options** menu; **Chapter 20**) in making QuickTime movies. To stop making a movie, select **Stop QuickTime Recording** (which has replaced **Start QuickTime Recording** in the **File** menu) and supply the requested file name.

## Exit

Exits *Spartan*, that is, clears the screen and closes all open documents. A prompt for a name is provided for each document that has not previously been saved.





# Chapter 13

## The Edit Menu

---

Operations under the **Edit** menu provide for undoing commands, copying items to and from the clipboard, finding text and graphics, centering molecules on screen and clearing the selected molecule.

Undo	Ctrl+Z
Cut	Ctrl+X
Copy	Ctrl+C
Paste	Ctrl+V
Find...	
Find Next	F3
Center	
Clear	

### Undo

Undoes the last operation from the **Build** and **Edit** menus and from the **Molecule Utilities** dialog. Undoes transition-state formation and retrieval from *Spartan*'s Databases.

### Cut, Copy, Paste

**Cut** moves the selected item to the clipboard and removes it from the document. **Copy** copies the item to the clipboard. The item is unaffected. **Paste** transfers the contents of the clipboard to the selected location. The contents of the clipboard are unaffected.

Among the important uses of the clipboard within *Spartan* are:

- (i) Transferring on-screen graphics into other applications such as Microsoft Word<sup>®</sup> and PowerPoint<sup>®</sup>.
- (ii) Temporary storage of a molecular structure for use in molecule building.
- (iii) Transferring data between *Spartan* spreadsheets and between a *Spartan* spreadsheet and other applications such as Microsoft Excel<sup>®</sup>.

- (iv) Making multi-molecule documents and/or transferring molecules between documents.

**Cut** operations for (i) and (ii) require drawing a selection box. First position the cursor slightly above and slightly to the left of the item to be transferred. Then, while holding down both buttons, *drag* the mouse to a location slightly below and slightly to the right of the item to be transferred. Finally, release both buttons. **Copy** operations for (i) and (ii) also refer to the contents of a selection box *if one has been drawn*, but to the selected molecule if a box has not been drawn. **Copy** operations from a *Spartan* spreadsheet refer to all information associated with a molecule if selection is made on the header cell of the leftmost column, but only to the selected (text) information if selection is made on any other column.

Further discussion relating to use of the clipboard in molecule building is provided in **Chapter 16** and for moving data in and out of the spreadsheet is provided in **Chapter 18**.

### **Find..., Find Next**

**Find** locates a text string defined in the **Find** dialog if an output window or a spreadsheet is selected, or a structure sequence defined on the clipboard if an on-screen model is selected. **Find Next** locates the next occurrence of a text string or a structure sequence.

### **Center**

Centers on screen all molecules in the document for which the selected molecule is a member (only the selected molecule is displayed).

### **Clear**

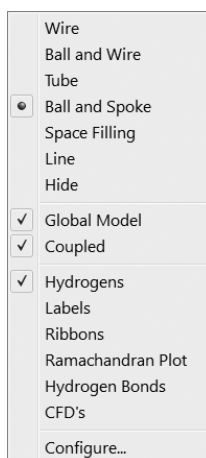
Clears (deletes) the structure and other information for the selected molecule, and brings up a model kit. No information is actually removed from the file system until the document is saved.

# Chapter 14

## The Model Menu

---

*Structure models available under the **Model** menu include wire, ball-and-wire, tube, ball-and-spoke, space-filling (CPK) and line models, with or without hydrogens, with or without hydrogen bonds indicated, with or without chemical function descriptors (CFD's) shown and with or without atom labels, as well as ribbon displays for polypeptides and polynucleotides, with or without labels and with or without hydrogen bonds indicated. It allows drawing a Ramachandran plot for a protein structure retrieved from the Protein Data Bank (PDB). The menu also provides for configuring atom labels to display element name, R/S chirality, mass number, charge or chemical shift, and for specifying color coding and display style for ribbon labels, as well as turning a variety of other labels on and off. Finally, it allows model style to be applied globally (to all molecules in a document) and models to be manipulated in concert.*




Only one model style **Wire**, **Ball and Wire**, **Tube**, **Ball and Spoke**, **Space Filling** or **Hide**) may be selected. The selected model is designated by a check mark ✓ in front of its entry in the menu. **Global Model**, **Coupled**, **Hydrogens**, **Labels**, **Ribbons**, **Hydrogen Bonds** and **CFD's** operate as toggle switches. A ✓ in front of the entry in the menu indicates that it is turned on.

All structure models, and graphics may be displayed either in orthogonal or perspective projections. The latter may be valuable in helping to visualize large molecules. Selection is done in the **Settings Preferences** dialog (**Preferences...** under the **Options** menu; **Chapter 20**). Both structure models and graphics may be presented in 3D stereo. This is also controlled from the **Settings Preferences** dialog as well as from the **3** key. Stereographic displays require perspective projections.

## Wire

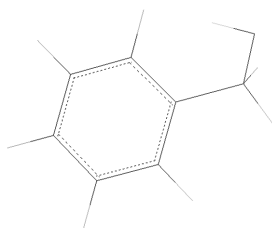This represents the molecule as a wire model where the vertices represent the atoms.

Wire Model




The bonds are drawn in two colors, one for each of the atoms making up the bond. Default atom colors are given in **Table 14-1**.

Atom colors apply globally (to all atoms of given type), and may be changed using the **Set Colors** dialog (**Colors** under the **Options** menu; **Chapter 20**). Colors of individually selected atoms may be set using the **Atom Style** dialog (**Properties** under the **Display** menu; **Chapter 18**). All models use the same color scheme for atoms, and provide for the same mechanism of changing colors globally or individually.

## Ball and Wire

This represents atoms by small balls and bonds by wires.

Ball-and-Wire Model

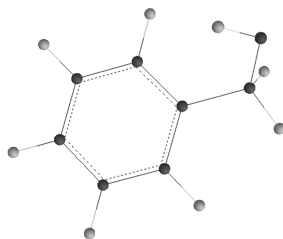

**Table 14-1: Default Atom Colors**

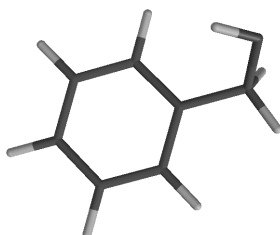
<b>main group</b>		<b>main group (con't)</b>	
Hydrogen	white	Indium	orange
Lithium	orange	Tin	gray
Beryllium	green	Antimony	orange
Boron	orange	Tellurium	orange
Carbon	gray	Iodine	violet
Nitrogen	blue-gray	<b>transition metals</b>	
Oxygen	red	Scandium-Zinc	green
Fluorine	pale yellow	Yttrium-Cadmium	green
Sodium	yellow	Lanthanum-Mercury	green
Magnesium	purple	<b>lanthanides</b>	
Aluminum	magenta	Cerium-Lutetium	violet
Silicon	gray	<b>actinides</b>	
Phosphorus	orange	Thorium-Lawrencium	dark blue
Sulfur	yellow	<b>noble gases</b>	
Chlorine	green	Helium	dark orange
Potassium	red	Neon	dark orange
Calcium	red	Argon	dark orange
Gallium	orange	Krypton	dark orange
Germanium	gray	Xenon	orange
Arsenic	dark orange		
Selenium	red orange		
Bromine	dark red		
Rubidium	red		
Strontium	red		

The balls are color coded according to atom type, and the wires representing bonds are drawn in two colors (as in wire models).

## Tube

This is similar to the wire model, except that tubes instead of wires are used to represent bonds.

Tube Model




Tubes may either be solid cylinders or be split to represent multiple bonds depending on whether **Split Tubes** in the **Settings Preferences** dialog (**Preferences...** under the **Options** menu; **Chapter 20**) is turned off or on. As with wire models, bonds are drawn in two colors.

## Ball and Spoke

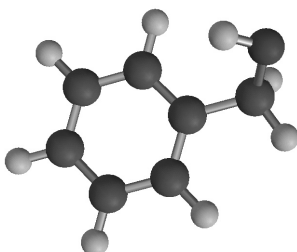This represents atoms by balls (the size and color of which depend on atom type), and bonds by spokes.

Ball-and-Spoke Model




Spokes may either be cylinders or be split to represent multiple bonds depending on whether **Split Tubes** in the **Settings Preferences** dialog (**Preferences...** under the **Options** menu; **Chapter 20**) is turned off or on. Bond (spoke) color is gray by default but it may be changed using the **Set Colors** dialog (**Colors** under the **Options** menu; **Chapter 20**). Colors of individually selected bonds may be set using the **Bond Style** dialog (accessible from **Properties** under the **Display** menu; **Chapter 18**).

## Space Filling

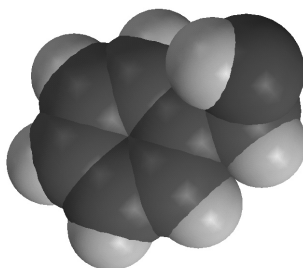This represents the molecule as a composite of spheres, the radii of which have been chosen to approximate van der Waals contact distances.\* Also known as CPK models, space-filling models are intended to portray overall molecular size and shape.


Space-Filling Model



Volume, surface area and polar surface area (PSA)\*\* displayed under the **QSAR** tab in the **Molecule Properties** dialog (accessible from **Properties** under the **Display** menu; **Chapter 18**) correspond to a space-filling model.

## Line

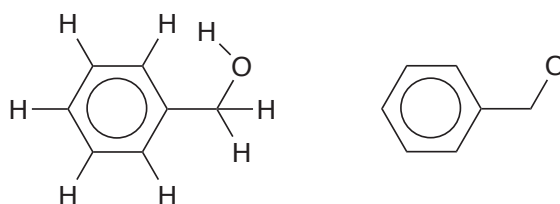The **Line** model allows *Spartan* to display native structures from 2D drawing programs such as ChemDraw and IsisDraw. It designates carbons as the termini or the intersection of lines and non-carbons by their elemental symbols without any use of color.



These models are perhaps best shown minus hydrogens (see **Hydrogens** later in this chapter), and include stereochemical cues (— and ..... ) if these have been provided in the 2D drawing or if

---

\* Default values for van der Waals radii may be changed from the **VDW Radii** dialog accessible from **Preferences** under the **Options** menu (**Chapter 20**). Settings apply to all atoms of given atomic number.

\*\* By default, polar surface area is defined as the area due to nitrogen and oxygen and any hydrogens attached to nitrogen and oxygen. Other combinations of atoms may be employed to define a polar surface area. See **Formulas** under the **Display** menu (**Chapter 18**).



they are specified from the **Bond Properties** dialog (accessible from **Properties** under the **Display** menu; **Chapter 18**).

## Hide

This removes the structure model from the screen. This may be desirable where its display may lead to unnecessary crowding, for example, in proteins where ribbon displays are more appropriate. A structure model may be restored by selecting it from the **Model** menu.

Different parts of a molecule may be rendered in terms of different model styles and colors. This allows focus on specific interesting regions of a molecule. Regions may be individual atoms and/or bonds or any collection of atoms and/or bonds contiguous or not. They are selected either by *clicking* on an individual atom or bond or, with the aid of the **Ctrl** (select multiple) or **command** key (Macintosh), **Shift** (select over a range), and **Alt** (select all bonded) or **option** key (Macintosh), by *clicking* on a set of atoms and/or bonds, or by defining a selection box. Discussion has already been provided in **Chapter 2**. Style and color of selected model components is by way of the “**Selected**” **Style** dialog accessible from **Properties** under the **Display** menu (**Chapter 18**).

## Global Model

If *checked* (turned on), this signifies that all molecules in a document will share attributes. These include presentation of hydrogens, atom and other labels, hydrogen bonds, CFD's and ribbon displays. Global model style is controlled from the **Molecule Preferences** dialog (**Preferences...** under the **Options** menu; **Chapter 20**). **Global Model** acts in a toggle manner, switching between global and local display. **Global Model** is normally on.

## Coupled

If *checked* (turned on), this signifies that all molecules in a document selected for simultaneous display will be moved together. **Coupled** is turned on following molecule alignment (see **Align** under the **Geometry** menu; **Chapter 15**). **Coupled** acts in a toggle manner, that is, repeated selection couples and decouples the molecules.

## Hydrogens

If *checked*, this signifies that hydrogens are to be included in the model. Note that structures retrieved from the Cambridge Structural Database or the Protein Data Bank may lack hydrogens. These need to be added to the model before they can be displayed, before CFD's (Chemical Function Descriptors) can be assigned (see discussion later in this chapter) and before calculations may be performed. This may be done either upon retrieval (see **Open/Close Databases** and **Extract Ligands** under the **Search** menu; **Chapter 19**) or using **Grow Hydrogens** in the **Molecule Utilities** dialog (accessible from **Properties** under the **Display** menu; **Chapter 18**). **Hydrogens** acts in a toggle manner, that is, repeated selection turns the display of hydrogens on and off.

## Labels

If *checked*, this signifies that labels associated with atoms, ribbons and bonds as well as with other attributes specified in **Configure...** (see discussion later in this chapter) are to be displayed in the model. **Labels** acts in a toggle manner, that is, repeated selection turns display of labels on and off. **Labels** is automatically turned on following selection of **Apply** or **OK** in the **Configure** dialog.

## Ribbons

If *checked*, this signifies that ribbons are to be displayed along with the selected model. (If only ribbons are desired, for example, in proteins, select **Hide** for the model.) **Ribbons** acts in a toggle manner, that is, repeated selection turns display of ribbons on and off.

## Ramachandran Plot

If *checked*, this draws a Ramachandran plot for a protein input from the Protein Data Bank (see **Access PDB Online** under the **File** menu; **Chapter 12**). **Ramachandran Plot** acts in a toggle manner, that is, repeated selection turns the plot on and off. Coloring of the points on the plot (red for  $\alpha$ -helices, blue for  $\beta$ -sheets, green otherwise) is based on assignments in the PDB file.

## Hydrogen Bonds

If *checked*, this signifies that hydrogen bonds are to be drawn as part of the model. Hydrogen bonds are defined as non-bonded contacts between a nitrogen or oxygen and a hydrogen attached to nitrogen or oxygen separated by a distance ranging from 1.6 to 2.1 Å and making an X–H--Y (X,Y = N, O) angle of >120°. **Hydrogen Bonds** acts in a toggle manner, that is, repeated selection turns display of hydrogen bonds on and off.

## CFD's

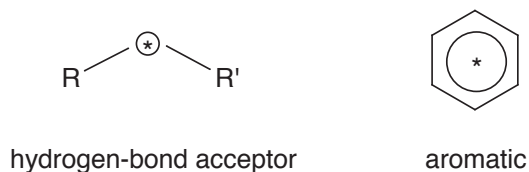
If *checked*, this signifies that **Chemical Function Descriptors (CFD's)** are to be displayed along with the structure model. CFD's are descriptors given to a molecule in order to characterize or anticipate its chemical behavior or to identify commonality among molecules with different structures. They parallel terms in a chemist's vocabulary such as *lone pair* (to suggest the role of a hydrogen-bond acceptor) and *sterically crowded* (to suggest that getting close may be difficult).

*Spartan* uses seven different kinds of CFD, the first six of which may be thought of as attributes of a molecule.


CFD	chemical language
1. hydrophobe	sterically-crowded region
2. aromatic	aromatic $\pi$ system
3. hydrogen-bond donor	acidic hydrogen
4. hydrogen-bond acceptor	lone pair
5. positive ionizable site	basic site
6. negative ionizable site	acidic site

Default CFD assignments follow from atomic connectivity and rules covering common organic functional groups. They may be modified on a case-by-case basis using the **CFD Properties** dialog (accessible from **Properties** under the **Display** menu; **Chapter 18**). CFD's (represented by \*'s in the figures below) may either be atom centered, or they may be centered in between atoms. For example, a hydrogen-bond acceptor CFD is placed on an ether oxygen, while

an aromatic CFD is placed at the center of a benzene ring.



CFD's may also be associated with a number of common functional groups. For example, a single CFD is provided for the carboxylic acid group, sited in-plane and equidistant from the carbon and two oxygens. The CFD's for this and other functional groups may depend on pH. At neutral (or basic) pH this CFD should be designated both as a hydrogen-bond donor and a hydrogen-bond acceptor, whereas at acidic pH it should be designated as a negative ionizable site.

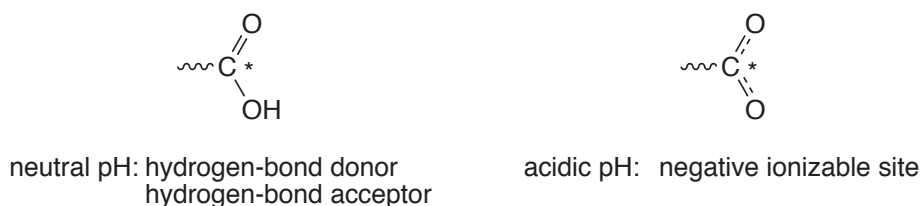

Finally, note that default assignments for hydrogen-bond donors and acceptors depend not only on atomic type (N, O, . . . ) and (in the case of hydrogen-bond donors) on the availability of hydrogens, but also on known (or presumed) chemical behavior. For example, because an amide nitrogen would generally not be considered to be a hydrogen-bond acceptor, no (hydrogen-bond acceptor) CFD is provided. It is assigned to be a hydrogen-bond donor (if hydrogens are available) and the amide oxygen is assigned as a hydrogen-bond acceptor.

Note that, hydrogen-bond acceptor and donor CFD's may be related to the electrostatic potential obtained from quantum chemical calculations. For example, a negative potential associated with an oxygen center suggests that it might serve as a hydrogen-bond acceptor, while a positive potential associated with a hydrogen attached to such a center suggests its role as a hydrogen-bond donor. These are, of course, the same conclusions that a chemist would draw simply by looking at the molecule. However, the calculations may be able to say which acceptor or donor sites are likely

to be strong and which are likely to be weak.


Calculated electrostatic potentials may offer significant advantage over CFD's and other qualitative descriptions when there are no heteroatoms and hydrogen bonding is not possible. The  $\pi$  systems of benzene and other aromatics provide the most common examples. The same conditions that make the ring susceptible to electrophilic aromatic substitution mean that it repels other electron-rich regions (and attracts electron-poor regions).

A seventh CFD type, excluded volume, derives from knowledge of a molecule incorporated into a host. Excluded volumes may be obtained from PDB files using the **Ligands** dialog (**Extract Ligands** under the **Search** menu; **Chapter 19**).

### Configure...

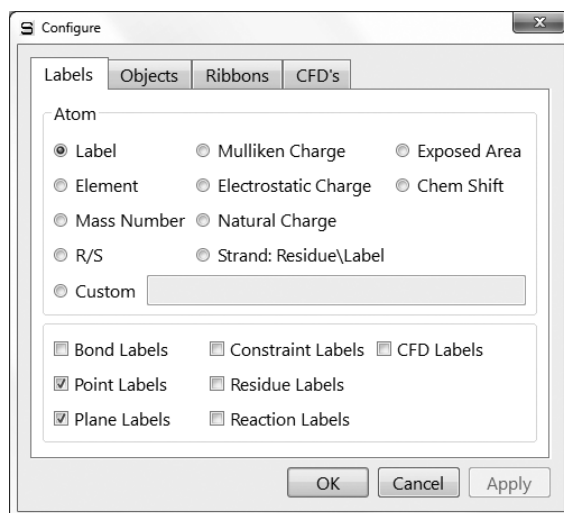This selects the types of labels attached to atoms, ribbons and CFD's.

### Configure Labels




Atom labels may be selected from among the following: **Labels**, a unique element/number combination that may be changed from the **Atom Properties** dialog (accessible from **Properties** under the **Display** menu; **Chapter 18**), **Element**, **Mass Number**, **R/S** (chirality), **Mulliken Charge**, **Electrostatic Charge**, **Natural Charge**, **Strand: Res\Label** (polypeptides and polynucleotides),

**Exposed Area** (of an atom in a space-filling model), **Chem Shift**. Selection of **Custom** leads to a text box into which a formula defining an atomic label may be entered. Details are provided in **Formulas** under the **Display** menu (**Chapter 18**). In addition, **Bond Labels**, **Point Labels**, **Plane Labels**, **Constraint Labels**, **Residue Labels**, **Reaction Labels** and/or **CFD Labels** may be provided. Default settings (for a new molecule) are made in the **Molecule Preferences** dialog (**Preferences** under the **Options** menu; **Chapter 20**).

## Configure Objects

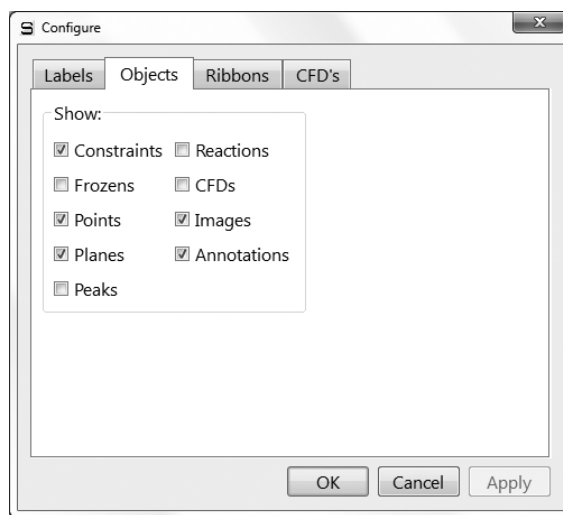*Clicking on the **Objects** tab leads to the **Configure Objects** dialog.*



If checked, **Constraint** and **Frozen** markers, **Points** and **Planes**, infrared and Raman **Peaks**, **Reaction** arrows and **CFD's** attach to the model. If not checked, these are shown only in the respective modes, for example, **Frozen** markers are shown only if **Freeze Center** is selected.

## Configure Ribbons

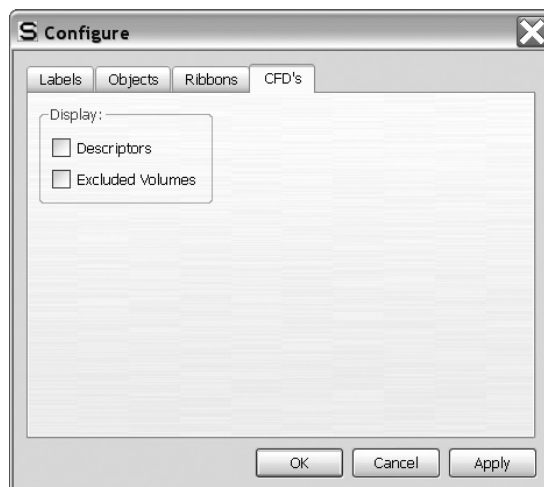
*Clicking on the **Ribbons** tab leads to the **Configure Ribbons** dialog.*



Ribbon coloring may be selected from among the following: **Monochrome**, **By Secondary Structure**, **By Strand** or **By Residue**. Ribbon style may be selected from among the following: **Ribbons**, **Beads** or **Lines**.


### Configure CFD's

*Clicking on the **CFD's** tab leads to the **Configure CFDs** dialog.*



This allows for turning CFD descriptors on (identifying the type of CFD) identifiers and CFD's corresponding to excluded volumes.

The **Configure** dialog is removed from the screen with all

selections maintained by *clicking* on **OK**. *Clicking* on **Cancel** or on  removes the dialog but selections are lost. *Clicking* on **Apply** maintains the selections but leaves the **Configure** dialog on screen. Note, that **Labels** (from the **Model** menu) will be turned on following either *clicking* on **OK** or on **Apply**.






# Chapter 15
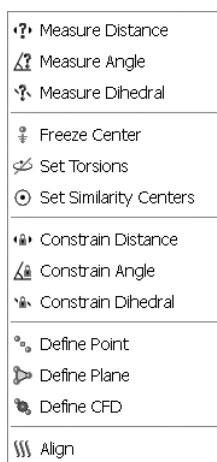
## The Geometry Menu

*Functions available under the **Geometry** menu allow querying and changing bond lengths, angles and dihedral angles, defining points and ligand points, planes and chemical function descriptor (CFD) locations, setting geometrical constraints or ranges of constraints, freezing atomic centers, altering default bond and ring conformer assignments, introducing NOE's in conformational searching, selecting atoms or CFD's for similarity analysis and aligning molecules in a list.*



**Measure Distance** ()


**Measure Angle** ()

**Measure Dihedral** ()


**Measure Distance** displays the distance (in Ångstroms) between two atoms, whether or not they are bonded. Selection results in a message at the bottom left of the screen.


Select two atoms, a bond, ...

*Clicking* on two atoms displays the distance at the bottom right of the screen.

Distance(C1,C2) = 1.500 Å 

Alternatively, *clicking* on a bond displays its length.

Length(Bond1) = 1.531 Å 

**Measure Distance** may also be used to alter the distance between atoms (as long as both are not incorporated into the same ring), by altering the contents of the box to the right of **Distance (A,B) =** or **Length (A)=**, and then *pressing* the **Enter** key (**return** key on Mac). The distance (length) may be entered into the spreadsheet by *clicking* on  to the right of its display (see **Spreadsheet** under the **Display** menu; **Chapter 18**). Alternatively, the label “**Distance (A,B)=**” or “**Length (A)=**” may be *dragged* into the spreadsheet, the formula editor (see **Formulas** under the **Display** menu; **Chapter 18**) or into any of the dialogs associated with the Spartan Molecular Database or Spartan Spectra and Properties Database (see **Open Databases** under the **Search** menu; **Chapter 19**).

Angle and dihedral angle queries are handled in a similar manner. Angles require that three atoms or two bonds be identified in the proper order while dihedral angles require that four atoms or three bonds be identified in the proper order.


### **Freeze Center**

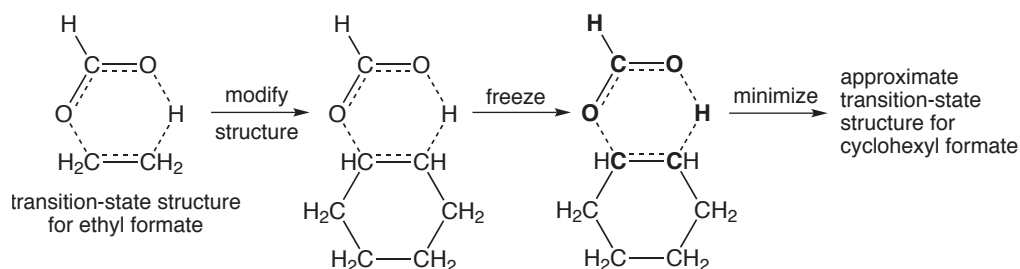This allows one or more atoms to be held in place during minimization (in build mode) or (optionally) during equilibrium or transition-state geometry optimization, conformational searching, or energy profile generation. Outside of build mode, use of frozen atoms needs to be explicitly indicated in the **Calculations** dialog (**Chapter 17**).

Atom freezing may be useful in a number of situations, among them guessing a transition-state geometry for a reaction that is closely related to one for which a transition state is available. For example, a good guess at the transition state for pyrolysis of cyclohexyl formate will be obtained by modifying the transition state for pyrolysis of ethyl formate, freezing all but the modified sections (designated in bold in the figure below) and then minimizing.\*

---

\* An example of this is presented in the tutorial *Carbene Additions to Alkenes* in **Chapter 7**.



Selection of **Freeze Center** leads to a message at the bottom left of the screen.

Select atom to freeze.

*Clicking* on an atom or free valence\*, freezes it; *clicking* again thaws it. Buttons at the bottom right of the screen allow for freezing all atoms (**Freeze All**), freezing all heavy (non-hydrogen) atoms (**Freeze Heavy**) and for thawing all atoms (**Thaw All**).

Another important use of frozen atoms is to “correct” XH bond lengths in structures resulting from a search of the Cambridge Structural Database (see **Open/Close Databases** under the **Search** menu; **Chapter 19**). Hydrogen positions are more often than not poorly located in small molecule X-ray structures, and XH bond lengths are commonly as much as 0.1 to 0.2Å shorter than they should be. Structures incorporating such bond lengths are clearly inappropriate for energy and property calculations and may also be problematic as starting geometries in quantum chemical calculations. Protein crystal structures always lack hydrogens, and hydrogen positions of bound ligands extracted from PDB files may be poorly defined. A reasonable solution to both problems is to freeze all heavy (non-hydrogen) atoms (**Freeze Heavy**) and then to carry out molecular mechanics minimization using either **Minimize** from the **Build** menu (see discussion in **Chapter 16**) or from the **Calculations** dialog under the **Setup** menu (see discussion in **Chapter 17**).

Frozen atoms are indicated by magenta colored markers (☞). Whether or not these are included with the model (outside of freeze center mode) is controlled from the **Molecule Preferences** dialog under **Preferences...** in the **Options** menu (**Chapter 20**).

\* The bond distance in this case is that appropriate for hydrogen being added to the free valence.

## Set Torsions (🔄)

*Spartan* automatically identifies bonds and rings for conformational searching and specifies default step sizes.

Two different rule sets have been provided. **Normal** is intended for use in determining the lowest-energy conformation or in establishing a distribution of low-energy conformers, whereas **Trimmed** is intended for use in building conformer libraries for use in similarity analyses. (Different modes of conformational analysis are discussed in **Chapter 17**.) Selection of rule set is made in the **Settings Preferences** dialog (**Preferences...** under the **Options** menu; **Chapter 20**).


**Set Torsions** allows these defaults to be altered. Selection results in rotatable bonds each being marked by a gold cylinder, and flexible rings each being marked by a gold circle around one or more atoms, and a message at the bottom left of the screen.

Select bonds or ring atoms.

*Clicking* on a bond or an atom contained in a ring selects it for rotation. In the case of a ring, rotation means that the atom is to be puckered up and down (restricted rotation). The default rotation is provided in a box to the right of **Fold** at the bottom right of the screen. This is typically 2 or 3 for a single bond (step size of 180° and 120°, respectively) and 3 for a flexible ring. A value of either 0 or 1 indicates that the bond (ring) is not to be rotated (rotation by 360°). Other integer values may be entered into the box, followed by *pressing* the **Enter** key (**return** key on Mac). The original settings may be retrieved by *clicking* on **Defaults** at the bottom right of the screen. *Double clicking* on a marked bond or ring atom removes the marker and eliminates the conformational degree of freedom.

Saturated six-membered rings for example, cyclohexane, are treated as special cases with only chair conformers examined. This is indicated by two circles on opposite atoms. Other conformers, for example, twist-boat conformers, can be generated by selecting additional ring atoms.



◀ and ▶ buttons at the bottom right of the screen are available to step

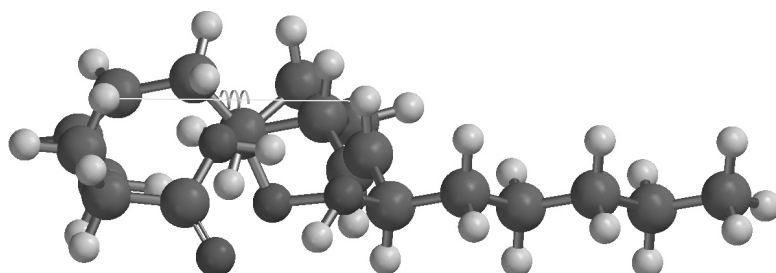
through the possible single-bond conformers\* (ring conformers are not provided). Geometries have not been optimized and some conformers may be severely crowded. Any conformer can be selected in lieu of the initial structure. The full set of single-bond conformers may be placed in an unnamed document by *clicking* on  (to the right of the step buttons), and then *clicking* **OK** in the dialog that results. Note, that duplicate conformers have not been removed.

The total number of conformers (excluding ring conformers) that will be examined in a systematic search (see discussion in **Chapter 17**) is indicated at the bottom right of the screen “**Conformers = xx**”. This is an upper bound to the number of conformers that will actually be kept as duplicates and high-energy conformers are eliminated. This can be *dragged* into the spreadsheet (see discussion in **Chapter 18**), into the **Formulas** dialog (see discussion in **Chapter 19**), or into any of the dialogs associated with the Spartan Molecular Database or Spartan Spectra and Properties Database (see discussion in **Chapter 19**).

**Set Torsions** is also used to specify non-bonded distances that need to be kept below a threshold value. These follow from NOE (Nuclear Overhauser Effect) measurements and may be referred to as NOE conditions or simply NOEs. NOEs are specified by *clicking* on two (non-bonded) atoms while holding down the **Shift** key. In response, a message appears at the bottom of the screen.

NOE(H2,H7) =  

*Clicking* on  changes it to  and enters the default value for the NOE threshold into the box. This value can be changed. A line is drawn between atoms that are to be kept within the threshold value.



---

\* A limit of 1000 conformers is enforced.

Once set, NOEs are used in conformational searching (conformer distribution only) without further user intervention. They act as a post-processor filter to eliminate conformers that do not satisfy the constraints and do not affect the speed of the search.

### Set Similarity Centers

This specifies the set of atomic centers, or **Chemical Function Descriptors (CFD's)** or pharmacophore elements that are to be used in similarity analysis. Choice of whether the similarity analysis is to be based on structure or CFD's (pharmacophore elements) is made from a menu that appears at the bottom right of the screen following selection of **Set Similarity Centers**. This menu is also accessible from the **Similarity Analysis** task in the **Calculations** dialog (**Calculations...** under the **Setup** menu; **Chapter 17**).

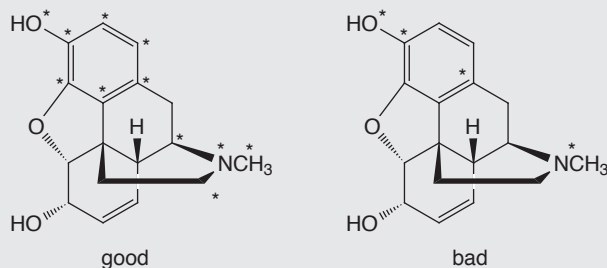
#### Similarity Based on Structure

Selection of **structure** from the **Similarity by** menu at the bottom right of the screen results in a message at the bottom left of the screen.

Select atoms.

*Clicking* on an atom designates it as an alignment center and it marks it with a violet circle. *Clicking* on the circle removes the designation (and the circle). Hydrogens are not used in similarity analysis and may not be selected.

The computer time required for similarity analysis based on structure depends on the number of molecules in the query, the number of molecules in the library to which comparisons are being made and the choice of similarity centers. Single atoms should never be selected as similarity centers unless they are unique. Rather, contiguous groups of atoms that are likely to be unique should be chosen. For example, analysis carried out on morphine based on selection of its phenol and trialkylamine components (left) will be much faster and much more informative than that based on selection of the phenol oxygen, two carbons on benzene and the nitrogen (right).



This is because there will be many few occurrences of phenol and trialkylamine substructures in the library entries than there will be aromatic carbons and nitrogen atoms and fewer permutations that need to be examined and fewer false positive results.

### Similarity Based on CFD's

Similarity analysis may also be based on CFD's. These indicate whether a particular center is likely to act as a hydrogen-bond donor or hydrogen-bond acceptor, is likely to take on positive or negative charge or is a hydrophobe. This is accomplished by selecting **CFD** from the **Similarity by** menu at the bottom right of the screen. In response, the selected molecule is augmented with CFD's, and a message appears at the bottom left of the screen.

Select CFD's.

*Clicking* on a CFD designates it as a similarity center and marks it with a violet circle. *Clicking* on the circle removes the designation (and the circle).

The computer time required for similarity analysis based on CFD's depends on the number of molecules and/or pharmacophores in the query, the number of molecules and/or pharmacophores in the library and the nature of the selected CFD's in the query. For each query/library pair, where the query is represented by  $i$  CFD's of type A,  $j$  CFD's of type B,  $k$  CFD's of type C, etc., the number of different alignment combinations that need to be considered is given by:

$$i!j!k! \dots$$

(! is the factorial symbol.) This rapidly becomes unmanageable. For example, if there are 6 CFD's of type A, 3 of type B and 2 of type C, the



number of combinations is  $6!3!2!$  or 8640. Considerable savings can be achieved by limiting the number of CFD's of a given type in the overall representation.

## Similarity Based on a Pharmacophore

Finally, it is possible to base similarity on a pharmacophore. This is set up in exactly the same way as similarity based on CFD's insofar as the elements of a pharmacophore are CFD's (see previous section).

Selection of similarity centers either needs to be repeated for the remaining molecules in the list or the selections may be designated as global, that is, applying to all molecules in the list. The latter is applicable where the molecules in the list are conformers or are closely related at least insofar as the alignment centers are concerned. Global designation is accomplished by *checking* the **Global Similarity Settings** at the bottom of the screen.

**Constrain Distance** ()

**Constrain Angle** ()

**Constrain Dihedral** ()

These introduce one or more geometrical constraints during structure minimization (in build mode), and (if requested) during equilibrium or transition-state geometry optimization or conformational searching using methods under the **Calculations** dialog (**Chapter 17**). They also allow for setting a range of constraints needed for generation of energy profiles. Constraints may be useful in a number of situations, among them:


- (i) constructing conformational energy profiles where one or more dihedral angles need to be fixed while other geometrical variables are optimized,
- (ii) optimizing molecular structures where the values of certain key parameters are known, for example, optimizing the geometry of a molecule with an intramolecular hydrogen bond or a disulfide linkage, and

- (iii) building molecules with unusual geometries, for example, molecules with very long bonds, as might be required in the construction of transition states and intermolecular complexes.


Selecting **Constrain Distance** results in a message at the bottom left of the screen.





Select two atoms, a bond, ...

*Clicking* on two atoms, or a bond results in a message at the bottom right of the screen.

Constraint(C1,C2) =  

*Clicking* on  changes it to  and shows the current distance.


Constraint(C1,C2) =  

This (constraint) distance can now be changed by altering the contents of the box and then *pressing* the **Enter** key (**return** key on Mac). Alternatively, the existing distance may be used as the constraint distance. If the selected distance had previously been constrained, the icon  would have been initially displayed. In this case, *clicking* on  turns the constraint off and returns the icon to the . Finally, the value of the constraint, that may be different from the value of the current distance\*, may be entered into the spreadsheet by *clicking* on  to its right. Alternatively, the label “**Constraint (A,B)=**” may be *dragged* into the spreadsheet, the formula editor (see **Formulas** under the **Display** menu; **Chapter 18**) or into any of the dialogs associated with the Spartan Molecular Database or Spartan Spectra and Properties Database (see **Open Databases** under the **Search** menu; **Chapter 19**).

This sequence of operations (bond identification followed by turning the constraint on and off) may be repeated as many times as necessary. Any bonds or non-bonded distances on which constraints are to be imposed are indicated by magenta colored markers. Any constraints introduced are automatically enforced only upon energy minimization

---

\* Note, however, that it may be problematic to carry out a constrained geometry optimization starting from a structure that is very different from that satisfying one or more constraints.

in add fragment mode () , but need to be requested for methods under the **Calculations** dialog (**Chapter 17**).

Angle and dihedral angle constraints are handled in a similar manner. Note that *points and planes may not be used to define constraints*.

Constraints may be dragged into the spreadsheet (see discussion in **Chapter 18**), into the **Formulas** dialog (see discussion in **Chapter 18**), or into the databases window (see discussion in **Chapter 19**). Constraints may also be modified as well as posted to the spreadsheet from the **Constraint Properties** dialog (accessible from **Properties** in the **Display** menu; **Chapter 18**). In addition, a sequence of constraints may be defined (from some initial value to some final value in a given number of steps).

Locking in a constraint leads to two additional text boxes at the bottom right of the screen. This allows generation of an energy profile along a predefined set of coordinates\* (see discussion under **Calculations** from the **Setup** menu; **Chapter 17**).



Whereas the leftmost box sets the initial value of the constraint, the middle box to the right of **to** sets the final value and the rightmost box to the right of **Steps:** sets the number of steps. For example, were the initial value set to 0°, the final value to 180° and the number of steps to 10, then a series of ten constraints (0°, 20°, 40°, ... 180°) would be specified. This can also be accomplished using the **Constraint Properties** dialog.

Whether or not constraint markers are included as part of the model (outside of constrain distance, constrain angle or constrain dihedral mode) is controlled from the **Molecule Preferences** dialog (**Preferences...** under the **Options** menu; **Chapter 20**).

---


\* Examples of this are provided in the tutorial *Internal Rotation in Dimethylperoxide* in **Chapter 6**, and the tutorial *Thermodynamic vs. Kinetic Control* in **Chapter 7**.

## Define Point ()

This defines a point as the geometric (unweighted) center of selected atoms (or points) previously defined. Selection results in display of a message at the bottom left of the screen.

Select atoms. Repeat to terminate.

*Clicking* on atoms (or points) in any order, and *clicking* a second time on any one of the atoms (or points) defines a point (depicted as a small sphere). As many points as desired can be defined and these are treated in the same way as an atom in defining distances, angles, etc. Points move with the molecule as its geometry is altered.

Selecting **Define Point** (or *clicking* on ) while holding down on the **Shift** key, followed by *clicking* on the appropriate atoms, leads to a **ligand point**. This is a point of attachment directed perpendicular to the geometric center of the plane defined by three atoms (or best plane in the case of four or more atoms). A ligand point shares all the characteristics of a normal point, but may also be used to bond to atomic fragments, functional groups, etc. See **Make Bond** under the **Build** menu (**Chapter 16**) for a discussion. Ligand points move with the molecule as geometry is altered.

Delete from the **Build** menu () or the **Delete** key may be used to remove a point or ligand point.

Whether or not points and ligand points are shown as part of the model is controlled from the **Molecule Preferences** dialog (**Preferences...** under the **Options** menu; **Chapter 20**).

## Define Plane ()

This defines and displays a reference plane. Selection results in display of a message at the bottom left of the screen.

Select three atoms.

*Clicking* on three atoms or points defines a plane. As many planes as desired may be defined, and these may be used in defining

distances, angles, etc. Planes move with the molecule as its geometry changes.

**Delete** from the **Build** menu () or the **Delete** key may be used to remove a plane.

Whether or not planes are included as part of the model for an is controlled from the **Molecule Preferences** dialog (**Preferences...** under the **Options** menu; **Chapter 20**).

### **Define CFD** ()

This defines the position of a new CFD. The **CFD Properties** dialog (accessible from **Properties** under the **Display** menu; **Chapter 18**) allows CFD size and type to be selected. Selection results in display of a message at the bottom left of the screen.

Select atoms. Repeat to terminate.

*Clicking* on atoms (or points) in any order, and *clicking* a second time on any one of the atoms (or points) defines the position of a CFD (depicted as a transparent sphere labelled Null). *Clicking* twice on only one atom locates the CFD on the atom. As many CFD's as desired can be defined. These move with the molecule as its geometry is altered.

Whether or not CFD's are included as part of the model is controlled from the **Molecule Preferences** dialog (**Preferences** under the **Options** menu; **Chapter 20**).

### **Align** ()

This aligns the selected molecule to all other molecules in the same document based either on structure or on CFD's. If the selected molecule is a pharmacophore, aligns the pharmacophore to all molecules in the same document (based on CFD's). Choice of whether alignment is based on structure or CFD's is made from the **Align by** menu that appears at the bottom right of the screen.

The terms “alignment” and “similarity analysis” have subtly different meanings inside of *Spartan*. Alignment refers to the operation in which the molecules and/or pharmacophores in a document are reoriented to best coincide with the selected molecule (pharmacophore). Where a molecule is selected, alignment may either be based on a set of atoms or on a set of CFD’s. Where a pharmacophore is selected, alignment is based on a set of CFD’s. A single reorientation is provided for each molecule in the document that can be aligned to the selected molecule (pharmacophore).

Similarity analysis refers to the operation in which all molecules in a user-specified library comprising one or more documents are compared one-on-one to all molecules and/or pharmacophores in the selected document (in the graphical user interface). As with alignment, either structural elements or CFD’s may be employed. The result is a similarity score from 0 to 1 (1 is perfect) and a manipulatable model of the two molecules (or molecule and pharmacophore). It is likely that most comparisons will not be successful, and also that some comparisons may lead to more than one score (corresponding to different orientations of the two components). Similarity analysis both by structure and CFD as well as to a pharmacophore is discussed in detail in **Chapter 17**.

In short, the role of alignment is to orient molecules and/or pharmacophores in order to highlight their common features, while the role of similarity analysis is to rank pairs of molecules or molecule/pharmacophore pairs according to the extent of their common features.


### **Align Based on Structure**

Selection of **Structure** from the **Align by** menu results in a message at the bottom left of the screen.

Select atoms .

*Clicking* on an atom designates it as an alignment center, and marks it with a red circle. *Clicking* on the circle removes the designation (and the circle). Following selection of alignment centers, *clicking* on the **Align by** button at the bottom right of the screen aligns the molecules. If no atoms are selected prior to *clicking* on **Align by**, then alignment is based on all (non-hydrogen) atoms.

Following alignment, two or more molecules may be displayed at once using spreadsheet functions (see **Spreadsheet** under the

**Display** menu; **Chapter 18**). Their motions (coordinates) will be coupled following alignment, but may be uncoupled allowing the aligned molecules to move independently (see **Coupled** under the **Model** menu; **Chapter 14**). Note that alignment center selections are kept and molecules can be realigned by again selecting **Align** from the **Geometry** menu (or *clicking* on ) followed by *clicking* on the **Align by** button.

### **Align Based on CFD's**

It is also possible to base molecule alignment on CFD's. These indicate whether a particular center is likely to act as a hydrogen-bond donor or acceptor, is likely to take on positive or negative charge or is a hydrophobe. This is accomplished by choosing **CFD** from the **Align by** menu at the bottom right of the screen followed by *clicking* on the **Align by** button. In response, the selected molecule is augmented with CFD's, and a message appears at the bottom left of the screen.


Select CFD's.

*Clicking* on a CFD designates it as an alignment center and marks it with a red circle. *Clicking* on the circle removes the designation (and the circle).

The **CFD Properties** dialog (accessible from **Properties** under the **Display** menu; **Chapter 18**) allows definitions of the individual CFD's to be altered.

Following selection of CFD's, *clicking* on the **Align by** button at the bottom right of the screen aligns the molecules. If no CFD's are selected prior to *clicking* on **Align by**, then alignment is based on all CFD's.

Following alignment, one or more molecules may be displayed at once using spreadsheet functions (see **Spreadsheet** under the **Display** menu; **Chapter 18**). Their motions (coordinates) will be coupled following alignment, but may be uncoupled allowing the aligned molecules to move independently (see **Coupled** under the **Model** menu; **Chapter 14**). Note that alignment center selections

are kept and molecules can be realigned by again selecting **Align** from the **Geometry** menu (or *clicking* on ) followed by *clicking* on the **Align by** button.

### **Align to a Pharmacophore**

Finally, it is possible to align molecules to a pharmacophore contained in the same list. This is set up in exactly the same way as alignment based on CFD's insofar as a pharmacophore element is a CFD (see previous section).

For both structure and CFD alignment (and for alignment to a pharmacophore), an alignment score from 0 to 1 (perfect alignment), is available in the spreadsheet. This is accessed by *clicking* on the **Add** button at the bottom of the spreadsheet, and selecting **Alignment Score** from the list of available properties (see **Spreadsheet** under the **Display** menu; **Chapter 18**). A score of 0 is assigned to molecules that cannot be aligned to the selected molecule.





# Chapter 16

## The Build Menu

---

The **Build** menu provides model kits and tools for building and editing organic, inorganic and organometallic molecules, as well as polypeptides and polynucleotides, for generating lists of substituted molecules, and a molecular mechanics procedure for structure refinement. This menu also provides seamless access to ChemDraw.

### Model Kits

**Spartan** provides five different model kits for assembling a variety of molecular systems: an organic model kit for most organic molecules, an inorganic model kit for organic molecules not well represented in terms of an uncharged (non-zwitterionic) valence structure, as well as inorganic and organometallic molecules, model kits for polypeptides and polynucleotides and a model kit for attaching substituents (lists of chemical groups) to previously-built molecules. The organic and inorganic model kits utilize atomic fragments, functional groups and rings (and ligands in the inorganic model kit), while the peptide model kit uses the set of natural amino acids as building blocks, and the nucleotide model kit the set of nucleotide bases.

Two-dimensional molecular structures (drawings) produced using ChemDraw\* can also be brought into **Spartan** and converted to three-dimensional structures.

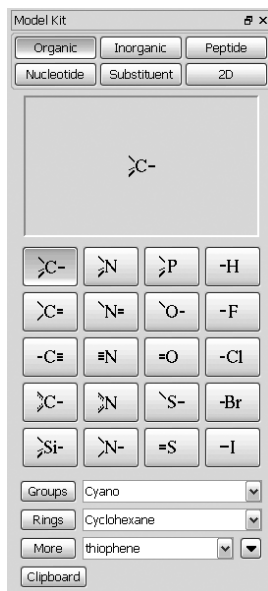
Molecule construction in **Spartan** proceeds much in the same manner as a chemist would assemble a structure from a model kit, that is, pieces are taken from the kit one at a time and added sequentially to the molecule under construction.

---

\* Seamless access to ChemDraw is available for the Windows version only. All versions are able to read ChemDraw files. ChemDraw is not included with **Spartan** and must be obtained from CambridgeSoft at <http://www.cambridgesoft.com>.

## Organic Model Kit

The organic model kit contains a suite of molecule building/editing tools specifically designed to construct organic molecules.



In the center of the model kit are a selection of atomic fragments, which from left to right and then top to bottom, correspond to:

C(sp <sup>3</sup> )	N(sp <sup>3</sup> )	P(sp <sup>3</sup> )	H
C(sp <sup>2</sup> )	N(sp <sup>2</sup> )	O(sp <sup>3</sup> )	F
C(sp)	N(sp)	O(sp <sup>2</sup> )	Cl
C(aromatic)	N(aromatic)	S(sp <sup>3</sup> )	Br
Si(sp <sup>3</sup> )	N(planar)	S(sp <sup>2</sup> )	I

A fragment is chosen by *clicking* on its icon, which is then displayed at the top of the model kit. Once selected, the fragment may be used to initiate building, to add alongside of an existing structure or appended onto an existing structure. To initiate building, *click* anywhere on screen. To add alongside of an existing structure, first *click* on the **Insert** button at the bottom right of the screen or hold down the **Insert** key (**option** key on Mac), and then *click* anywhere on screen or *double click* in a blank area on screen. To bond to an existing structure, *click* on a free valence (**not an atom**). (Free valences are colored yellow on the selected molecule.) Bond type in the case of atomic fragments


with multiple bond types, for example,  $sp^2$  carbon, depends on the nature of the free valence selected.

While only H, C, N, O, F, Si, P, S, Cl, Br and I are available from the organic model kit, other elements may be substituted using atom replacement feature available in the inorganic model kit (see **General Molecule Building Functionality** later in this chapter). Similarly, bond types may be altered in the inorganic model kit. The latter may be necessary in order to alter bonding in structures retrieved from the Cambridge Structural Database or ligands obtained from PDB files. Atom and bond types may also be altered using the **Atom** and **Bond Properties** dialogs, respectively (accessible from **Properties** under the **Display** menu; **Chapter 18**). Note that bond type (or the fact that bonds may be present) is not considered in quantum chemical calculations, and only atom identities and coordinates are required.

Menus inside the model kit provide access to a number of pre-built fragments corresponding to functional groups (**Groups**) and rings (**Rings**), and to additional libraries of rings (as well as any user-defined structures) stored in *Spartan*'s file system (**More**). The model kit also accesses the clipboard (**Clipboard**).

## Groups

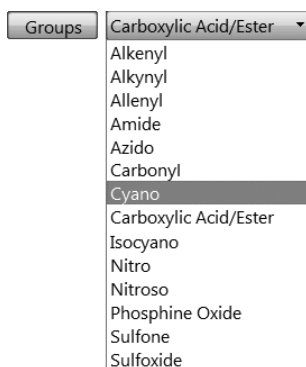*Clicking* on **Groups** brings up a menu of groups, and displays an icon of one group at the top of the model kit.




Once selected from the menu, a group may be used to initiate building, to add alongside of an existing structure on screen, or to add to an existing structure.

The amide and carboxylic acid/ester groups have more than one different free valence. The free valence that is to be used is marked with a gold • (in the icon at the top of the model kit). The marked position circulates among the possible positions with repeated *clicking* on the icon.

## Rings

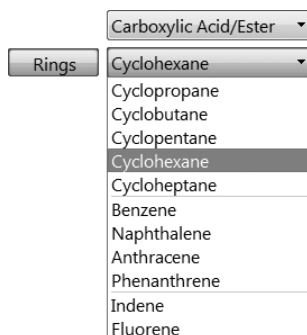
*Clicking* on **Rings** brings up a menu of hydrocarbon rings, and displays an icon of one ring at the top of the model kit.



Once selected from the menu, a ring may be used to initiate building, to add alongside of an existing structure on screen, or to add to an existing structure.

Cyclohexane, cycloheptane, naphthalene, anthracene, phenanthrene, indene and fluorene have more than one different free valence. The one that is to be used is marked with a gold • (in the icon). The marked position circulates among the available positions with repeated *clicking* on the icon. Selection of an *axial* or *equatorial* free valence in cyclohexane and cycloheptane is indicated by the label **ax** or **eq** appearing alongside the icon. All rings in this menu are hydrocarbons, but heteroatoms may be substituted (see **General Molecule Building Functionality** later in this chapter).

## More

This provides access to a broader selection of rings as well as to access user-defined entities (rings, groups, ligands, etc.). Upon initial entry, the menu to the right of **More** will be empty. It can be populated, by *clicking* on  to the far right. This brings up a file browser that has been set to point toward a directory containing

several *Spartan* files of common rings.

nitrogen heterocycles	saturated nitrogen rings
oxygen heterocycles	saturated oxygen rings
sulfur heterocycles	saturated sulfur rings
mixed heterocycles	saturated mixed rings

*Clicking* on a file followed by *clicking* on **Open** or *double clicking* on a file fills the menu to the right of **More**. Menu entries are selected in the usual way. In response, a ball-and-wire model of the selected ring will appear at the top of the model kit. This may be manipulated (rotated, translated, zoomed) using the usual mouse/keyboard commands (you need to position the cursor inside the box). The ring may be used to initiate building, to add alongside of an existing structure, or to add to an existing structure. In the latter case, the attachment point (on the ring in the window) needs to be identified by *clicking* on the appropriate free valence.

Three other collections (ligands, chelates and high-coordination fragments) appear in the selected directory. They are intended for use with the inorganic model kit as discussed in the next section. User-defined libraries may also be accessed (they are normal *Spartan* documents).

## Clipboard

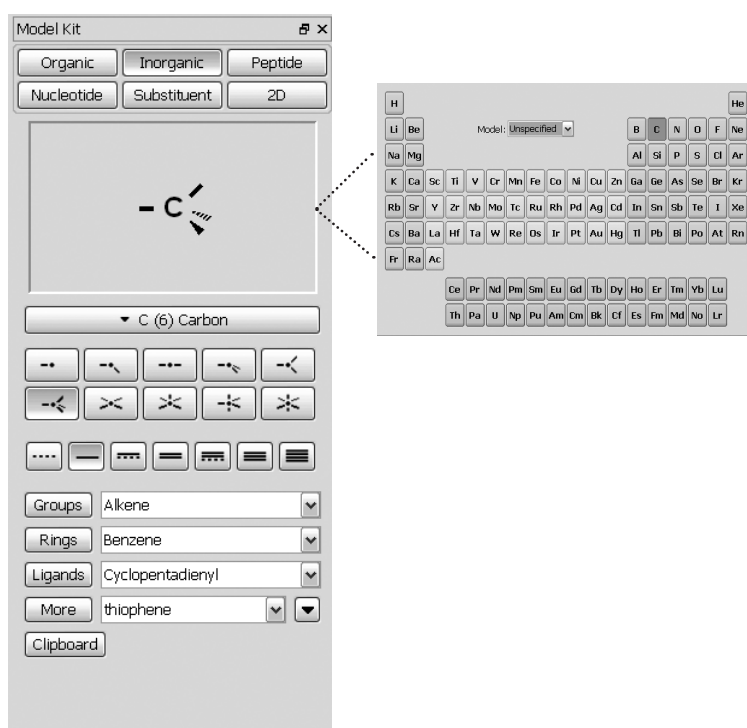
*Clicking* on **Clipboard** accesses the clipboard. A ball-and-wire model of whatever is on the clipboard is displayed at the top of the model kit. This may be manipulated using the usual mouse/keyboard commands (you need to position the cursor inside the box). Once selected, the molecule may be used to initiate building, to add alongside of an existing structure, or to add to an existing structure. In the latter case, the attachment point needs to be identified by *clicking* on the appropriate free valence in the clipboard.

An empty clipboard will be signaled by:




## Inorganic Model Kit

*Spartan*'s inorganic model kit allows construction of a much wider class of molecules (including inorganic and organometallic species) than possible with the organic model kit. Structures that violate conventional bonding rules may be constructed, as this model kit purposefully provides no checking. The inorganic model kit is reached by *clicking* on the **Inorganic** tab, located at the top of the organic (peptide, nucleotide or substituent) model kit. (Return to the organic, peptide, nucleotide or substituent model kit or to ChemDraw is accomplished by *clicking* on the **Organic**, **Peptide**, **Nucleotide**, **Substituent** or **2D** tab, respectively, in the inorganic model kit).



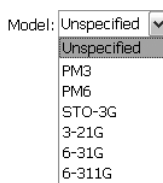Atoms may be selected by *clicking* on the center atom button. This leads to a full *Periodic Table*. Main-group elements are colored red, transition metals are colored green and lanthanides and actinides are colored blue. The **Model** menu inside the *Periodic Table* contains a selection of theoretical models (basis sets except for semi-empirical models).



Selecting an entry from this menu leads to recoloring of the *Periodic Table*. A light green color is used to indicate elements for which the selected model may be used.\* Immediately below is a selection of atomic hybrids.

Selection of atom type is effected by *clicking* on the appropriate element in the *Periodic Table*. The entry will be highlighted. Selection of an atomic hybrid follows by *clicking* on the appropriate icon which will then be highlighted.\*\* This combination (atom type + atomic hybrid) may be used to initiate building, to add alongside of an existing structure or to append onto an existing molecular fragment. To initiate building, *click* anywhere on screen. To add alongside of an existing structure, first *click* on **Insert** at the bottom right of the screen or *press* the **Insert** key (**option** key on Mac), and then *click* anywhere on screen, or *double click* in a blank area on screen. To bond to an existing fragment, *click* on the appropriate free valence.

Two of the hybrids (trigonal bipyramidal and square-based pyramidal) may bond either *axially* or *equatorially*. Selection of the appropriate bonding point, marked by a •, is effected by repeatedly *clicking* on the icon; the bonding point alternates between the two sites.

Atoms are connected with whatever bond type (partial single, single, aromatic, double, triple or quadruple) is selected from a menu below atom types in the model kit. A bond type may be changed by first selecting a type and then *double clicking* on the bond. Bond types have no impact on quantum chemical calculations, but do affect molecular

---

\* The existence of either an all electron basis set or a pseudopotential is sufficient to classify a model as available for a particular element. Note that while molecular mechanics models are available for all elements, they have been carefully parameterized for only a relatively few of these. Use caution in interpreting molecular mechanics geometries for molecules incorporating heavy main-group elements, transition metals, lanthanides and actinides.

\*\* Hybrids for a number of high-coordination centers are available as a library reachable from **More** (see discussion under organic model kit).




mechanics calculations which reference a Lewis structure (including minimization in the builder; see discussion later in this chapter).

No valence checking is performed in the inorganic model kit, and the user is free to construct any arrangement of atoms.

Menus under **Groups**, **Rings** and **More** are the same as those described for the organic model kit as is **Clipboard**. One additional fragment collection is provided:

### Ligands

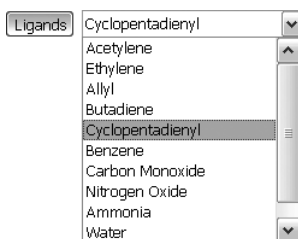This provides access to a number of pre-built ligands, useful in the construction of inorganic and organometallic molecules. Its operation is analogous to that for the **Groups** and **Rings** menus. *Clicking* on **Ligands** brings up a menu of available ligands, and results in an icon of one ligand from this menu being displayed at the top of the model kit.



A ligand may be used to initiate building or to add alongside or to an existing structure. Additional ligands are accessible from **More** (see previous discussion). Ligands may also be built with the aid of ligand points (**Define Point** in the **Geometry** menu; **Chapter 15**).



inorganic model kits.

There are two different modes of operation: single amino acid mode and polypeptide mode. The former is used to initiate building with a single amino acid, to add a single amino acid alongside of an existing structure or to add a single amino acid to an existing structure, while the latter is used to construct amino acid sequences (polypeptides). **Sequence** off (unchecked) corresponds to single amino acid mode, and on (checked) corresponds to polypeptide mode.

Peptide construction (**Sequence** on) is accomplished in three steps:

### **Specification of Amino Acid Sequence**

This is accomplished by *clicking* in the desired order on the amino acid codes. Building occurs from the N end to the C end of the peptide. In response to each selection, the three-letter code is appended to the sequence of codes in the box at the top of the model kit. The stereochemical configuration of the amino acid is by default the l configuration; this may be changed to the d configuration prior to selection of the amino acid, by *checking d* to the right of **stereoisomer** in the model kit. (It may be changed back to l by *checking l*). d amino acids are indicated by **.d** following the code in the box.

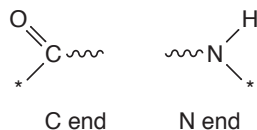
The sequence may be altered by changing the text in the box. Existing amino acid codes may be deleted or changed or new codes can be added. The entire sequence may be specified in this way if desired. Specification of a non-existent code will result in an error message. The sequence can be cleared by *clicking* on **Clear**.

### **Specification of Macroscopic Structure**


Once sequencing is complete, macroscopic structure ( $\psi$  and  $\phi$  angles), is specified by *clicking* on one of  $\alpha$  **Helix**,  $\beta$  **Sheet** or **Other**. In the case of the first two, preset angle values are displayed on the right. In the case of specification of **Other**, boxes appear, into which the desired dihedral angles need to be entered.


## Termination

The peptide is not yet terminated, and the two ends are still set up for addition of further amino acids.



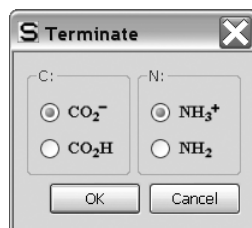where the \* indicates a free valence. *Clicking* on **Terminate** at the bottom of the model kit leads to the **Terminate** dialog.




C and N terminating groups may be selected by repeated *clicking* on the C and N icons, respectively. Selection will rotate among the available terminating groups. *Clicking* on **OK** removes the dialog and terminates the polypeptide. *Clicking* on **Cancel** or  removes the dialog but does not terminate the polypeptide.

The peptide (or single amino acid) may now be used either to initiate building, by *clicking* anywhere on screen or added alongside of an existing structure, by first *clicking* on **Insert** at the bottom right of the screen or *pressing* the **Insert** key (**option** key on Mac), followed by *clicking* anywhere on screen, or *double clicking* in a blank area on screen. If unterminated, it may also be joined onto an existing structure by *clicking* on a free valence. In the latter case, attachment is made from the N end, unless the free valence corresponds to an unterminated peptide fragment, in which case the appropriate end required to make an amide bond is used.

## Nucleotide Model Kit

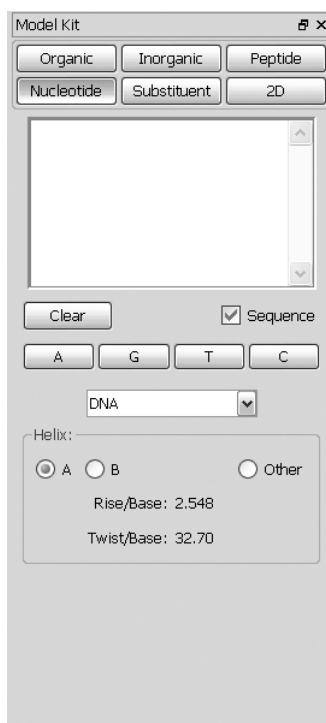Finally, *Spartan* provides a model kit for construction of polynucleotides. It is reached by *clicking* on the **Nucleotide** tab which is located at the top of the organic, inorganic, peptide and substituent model kits. (Return to the organic, inorganic, peptide or substituent model kit or to ChemDraw is accomplished by *clicking* on the **Organic, Inorganic, Peptide, Substituent** or **2D** tab, respectively, in the nucleotide model kit.)



At the middle of the model kit is a menu designating the type of polynucleotide.

DNA  
DNA (single strand)  
DNA-RNA  
RNA  
RNA (double strand)  
RNA-DNA

Immediately below this menu are icons, designating the nucleotide

bases. Selection of DNA, DNA (single strand) or DNA-RNA from the menu leads to one set of icons.



Selection of RNA, RNA (double strand) or RNA-DNA leads to a second set, the only difference is that uracil (U) has been substituted for thymine (T).



A nucleotide base is selected by *clicking* on its letter, following which either an icon of the base is displayed in the box at the top of the model kit, or the letter for the base is appended to the sequence of letters in the box. Nucleotide bases replace atomic fragments, functional groups, rings and ligands as the building blocks in the nucleotide model kit. Because these other building blocks are missing, modifications of nucleotides, aside from modifications in sequence and helical structure, need to be carried out using either the organic or inorganic model kits.

There are two different modes of operation: single base mode and polynucleotide mode. The former is used to place a single base or base pair on screen, to add a single base or base pair alongside of an existing structure, or to add a single base or base pair to an existing structure, while the latter is used to construct strands of DNA or RNA (or mixed strands). **Sequence** off (unchecked) corresponds to single base (base pair) mode and on (checked) corresponds to polynucleotide mode.

Polynucleotide construction (**Sequence** on) is accomplished in two steps:


### **Specification of Base Sequence**

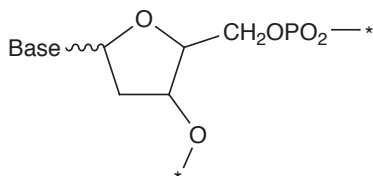This is accomplished by *clicking* in order on the base codes. In response to each selection, the letter code is appended to the sequence of codes in the box at the top of the model kit. The sequence may be altered by editing the contents of the box. Existing base codes may be deleted or changed or new codes added. The entire sequence can be specified in this way if desired. The sequence may be cleared by *clicking* on **Clear**.

## Specification of Helical Structure

Once sequencing is complete, a helical structure may be specified by *clicking* on **A** or **B**. These correspond to A and B helices, respectively. Selecting **Other** allows user modification of the rise (in Å) per base (**Rise/Base**) and twist (in degrees) per base (**Twist/Base**).

Note that the polynucleotide is not yet terminated, and the two ends are still set up for addition of further bases or base pairs.



\* indicates a free valence. Hydrogens occupy all free valences (except the \*'ed positions at the two ends of the chain).


The polynucleotide (or single base pair) may now be used to either initiate building, by *clicking* anywhere on screen, added alongside of an existing structure, by first *clicking* on **Insert** at the bottom right of the screen or *pressing* the **Insert** key (**option** key on Mac) followed by *clicking* on screen, or *double clicking* on a blank area on screen, or joined onto an existing structure by *clicking* on a free valence. In the latter case, attachment is made from the phosphate end.

## Substituent Model Kit

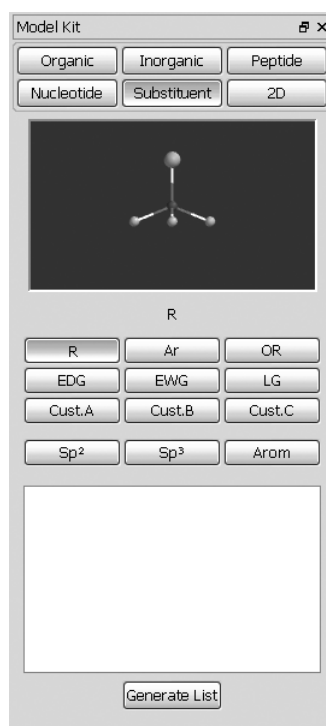The substituent model kit designates one or more locations for libraries of substituents to be added to a molecule. Additionally, it designates one or more carbon hybrids for replacement by their heteroatom “equivalents”, for example, sp<sup>3</sup> carbon by sp<sup>3</sup> nitrogen. There are two uses for the **Substituent** model kit. The first is to facilitate assembling lists of molecules that differ only in substitution and/or heteroatom. The second application is to use “molecules” incorporating one or more substituent libraries as queries in searches of the Spartan Molecular Database and Spartan Spectra and Properties Database (see **Open/Close Databases** under the **Search** menu; **Chapter 19**)

or as elements of reaction energy calculations (see **Reactions** under the **Display** menu; **Chapter 18**). Here, the search or reaction energy calculation is not confined to a single molecule but rather is open to a set of molecules that are related by substitution. Further discussion is provided in the appropriate sections of **Chapters 18** and **19**.

The substituent model kit is accessed by *clicking* on the **Substituent** tab located at the top of the organic, inorganic, peptide and nucleotide model kits. (Return to the organic, inorganic, peptide or nucleotide model kit or to ChemDraw is accomplished by *clicking* on the **Organic**, **Inorganic**, **Peptide**, **Nucleotide** or **2D** tab, respectively, in the substituent model kit.)



At the center of the model kit are four rows of buttons. Buttons in the first two rows provide access to predefined substituent libraries: **R** (alkyl groups), **Ar** (aromatic groups), **OR** (alkoxy groups), **EDG** (electron donor groups), **EWG** (electron withdrawing groups) and **LG** (leaving groups). Buttons in the third row allow access to user-defined substituent libraries (**Cust. A**, **Cust. B**, **Cust. C**), while those in the last row permit substitution of heteroatom hybrids for  $sp^3$ ,  $sp^2$



and aromatic carbons incorporated in a molecule. Immediately below the buttons, a dialog lists the entries in the individual libraries (the three custom libraries should be blank upon initial entry). As entries are selected, they will appear in the window at the top of the dialog. Predefined libraries may not be altered and individual entries may not be renamed, but existing library entries in custom libraries may be renamed or deleted and new entries added. Modifications to existing substituent libraries need to be performed *following attachment* using the **Substituent Properties** dialog (accessible from **Properties** under the **Display** menu; **Chapter 18**).

Creating lists of substituted molecules starting from an unsubstituted molecule occurs in three stages: library preparation, attachment and (optionally) list generation.

### Library Preparation

Six substituent libraries have been provided.

**R** (alkyl group): methyl, ethyl, *n*-propyl, isopropyl, *n*-butyl, *sec*-butyl, isobutyl, *t*-butyl, *t*-amyl

**Ar** (aromatic ring): phenyl, *p*-tolyl, *p*-nitrophenyl, 1-naphthyl, 2-naphthyl, 2-pyridyl, 3-pyridyl, 4-pyridyl

**OR** (alkoxy group): methoxy, ethoxy, *n*-propoxy, isopropoxy, *t*-butoxy


**EDG** (electron donor group): hydroxy, methoxy, ethoxy, phenoxy, amino, dimethylamino, mercapto, methylthio, phenylthio

**EWG** (electron withdrawing group): cyano, nitro, trifluoromethyl, carboxy, formyl, methoxycarbonyl, ethoxycarbonyl, methyl sulfonyl, methyl sulfoxyl

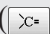
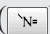





**LG** (leaving group): bromo, iodo, tosylate, triflate, mesylate

The contents of an existing library may be displayed in the scroll box at the bottom of the model kit by *clicking* the appropriate button. A ball-and-wire model of a substituent may be seen at the top of the model kit by *clicking* on its name in the scroll box. The model may be manipulated in the usual way. The attachment point is indicated by a gold ball surrounding a free valence.

To create a custom substituent library, first *click* on the **Cust. A**, **Cust. B** or **Cust. C** buttons inside the model kit. The scroll box is empty. If a substituent (set of substituents) has previously been copied to the clipboard, it may be pasted into the library by *right clicking* inside the scroll box and selecting **Paste** from the menu that appears. Alternatively, a substituent (set of substituents) may be appended to the document (including an empty document) by *right clicking* inside the scroll box, selecting **Append** from the menu that results and selecting a *Spartan* document comprising one or more molecules from the browser. Finally, any *Spartan* document may be *dragged* onto the scroll box from the file system. As many substituents (sets of substituents) as desired may be pasted, appended and/or *dragged* into a library. To remove a substituent from the library, *right click* on its name in the scroll box and select **Delete** from the menu that appears. **Ctrl** (command for Mac) and **Shift** keys may be used in the usual way.

Substituents need to be specified in terms of both their structure and their connection point. One way to do this is to specify the latter is to replace all free valences *except the one that will serve as the connection point* by hydrogens. For example, to designate a 2-pyridyl substituent, add hydrogens () from the organic model kit) to four of the five free valences in pyridine, leaving only one free valence in the 2-position available. Alternatively, a connection point can be specified by *clicking* on the appropriate position as the substituent displayed at the top of the model kit.

The three buttons in the last row designate heteroatom equivalents of carbon-based hybrids.

**sp<sup>2</sup>**: carbon () , nitrogen ()  
**sp<sup>3</sup>**: carbon () , nitrogen () , oxygen ()  
**Arom.**: carbon () , nitrogen ()


## Attachment

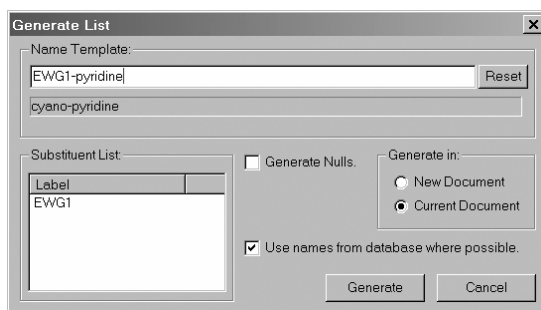Following selection, a substituent library may be used to initiate building, to add alongside of an existing structure or to attach to an existing structure. The model is augmented by an icon containing the library name from the button, for example **EWG**, followed by a number. As many substituent libraries as desired may be attached to a molecule. Each library is given a unique label, that can be changed prior to list generation using the **Substituent Properties** dialog. A substituent library may be removed from the molecule by selecting **Delete** from the **Build** menu and then *clicking* on the library icon.

A document containing one or more molecules with attached substituent libraries may be saved prior to list generation just as any other *Spartan* document. Any queries (bond distances, angles, etc.) and calculations refer to the unsubstituted molecule, that is, hydrogens replacing substituent markers. Database searches (on the Spartan Molecular Database or Spartan Spectra and Properties Database) and reaction energy calculations based on entries in the Spartan Molecular Database or Spartan Spectra and Properties Database pertain to all possible substituted molecules.

## List Generation

Generation of a list of substituted molecules follows specification and (optionally) editing of one or more substituent libraries and attachment of each of these libraries at one or more positions in a molecule. Duplicates are eliminated. *Clicking* on **Generate List** at the bottom of the model kit, leads to the **Generate List** dialog.



A box at the top of the dialog allows the user to specify (or modify)

a **Name Template**, from which the names of substituted molecules will be derived. Upon entry, this takes the form.

Ax-By-Cz-molecule name

A, B, C are the names of the substituent libraries (**R**, **EWG**, etc.), x, y, z are integers and molecule name is the name of the unsubstituted molecule if and as it appears in the Spartan Molecular Database (SMD) or in the Spartan Spectra and Properties Database (SSPD). If the unsubstituted molecule is not in the database, then the box below **Name-Template** will be blank. If the substituted molecule is in SMD or SSPD, its name in the database may be used (instead of the name generated from the template) by *checking* the box to the left of **Use names from SMD or SSPD where possible**.

Nulls are molecules with less than the full complement of substituents, for example, monosubstituted molecules in the case where two different positions have been substituted. These will not be included in the list of substituted molecules unless **Generate Nulls** is checked.

Note, that the conformations of the substituted molecules are arbitrary and they have not been subjected to energy minimization to relieve unfavorable crowding interactions. Therefore, the list entries should be examined, conformations adjusted as appropriate and minimized with molecular mechanics prior to using them for quantum chemical calculations.

The list may either be written to a **New Document** or appended to the end of the **Current Document**.

### Accessing ChemDraw (NAL, NAM)\*

The ChemDraw program may be seamlessly accessed inside of *Spartan*\*\* allowing chemical drawings to be produced in a familiar environment and then brought over as 3D structures\*\*\*. The conversion

---

\* ChemDraw files (.cdx) may be read with all versions of *Spartan*.

\*\* ChemDraw is not provided with *Spartan* but must be obtained from CambridgeSoft at <http://www.cambridgesoft.com>.

\*\*\* Transfer is one directional only. 3D structures that have been altered may not be transferred back to ChemDraw.

is unambiguous as long as all stereochemical cues are in place in the 2D drawing. Note that the conformation implied by the 2D drawing may be ambiguous and need to be altered.

To access ChemDraw, enter the builder, *click* on the **2D** tab and *click* on **Edit** at the bottom of the panel that results. ChemDraw will appear. To return to *Spartan*, close ChemDraw. The 2D drawing will appear at the center of the panel and manipulatable 3D structure will appear at the top of the panel. *Clicking* on screen will move the 3D structure into *Spartan's* main window.

## General Molecule Building Functionality

Additional capabilities are available with **Add Fragment** selected:

### Multiple Fragments

Multiple fragments may result either from bond breakage (see **Break Bond** later in this chapter) or from use of the **Insert** button or **Insert** key (**option** key on Mac), or *double clicking* in a blank area on screen. A fragment is selected by *clicking* on it, following which the associated free valences are colored gold (free valences for any non-selected fragments are colored white). All builder functions apply to the selected fragment only. Rotation and translation also apply to the selected fragment, but may be made to apply to the entire set of fragments by holding down the **Ctrl** key while carrying out these operations.

Fragments may be attached using **Make Bond** (see discussion later in this chapter).

### Bond Rotation/Bond Stretching

In addition to molecule rotation, translation and scaling, the mouse is used to rotate about and stretch bonds not incorporated into rings. This is accomplished via the following sequence of operations:

- (i) *Clicking* on the bond, which is then marked by a red cylindrical arrow. (The bond connecting the last atom, group,

ring or substituent added to the molecule is automatically selected.)

- (ii) Simultaneously holding down the **Alt** key (**option** key on Mac) and the left mouse button while *dragging* the mouse up and down, for bond rotation, or the **Alt** (**option**) key and the right mouse button for bond stretching.

### Atom/Fragment Replacement

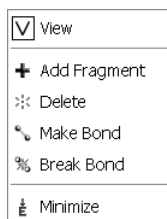
Another function of the mouse is atom replacement. *Double clicking* on an atom (not a free valence) while an atomic fragment in the organic model is highlighted, replaces the atom by selected fragment. Free valences are adjusted to accommodate the replacement, for example, replacement of  $sp^3$  carbon by  $sp^3$  oxygen results in two free valences being removed. Atom replacements that violate valence rules or that would disrupt substituents are not permitted. *Double clicking* on an atom (**not a free valence**) while an element in the *Periodic Table* from the inorganic model kit is selected, replaces the atom by the selected element. The latter merely changes the atomic number. No changes in the number or arrangement of free valences is made, and no checking is done. Atom replacement is not available in the peptide, nucleotide and substituent model kits.

### Chirality Inversion

In the **Add Fragment** mode, *double clicking* on a chiral atom with the **Ctrl** key (**command** key on Mac) depressed inverts the chirality of the atom (R→S or S→R). This is not available in fused ring systems. *Double clicking* on any atom with both **Ctrl** (**command** key on Mac) and **Shift** keys depressed inverts the absolute configuration of the molecule.

## Building/Editing Menu Functions

Molecule building/editing functions are found under the **Build** menu.



### View (V)

This exits build mode, and removes the model kit from the screen.

Initial entry into build mode is by way of **New** or **New Molecule** under the **File** menu (Chapter 12). **Add Fragment**, **Delete**, **Make Bond**, **Break Bond** and **Minimize** are for modifying existing structures.

### Add Fragment (+)

In addition to the capabilities discussed under **General Molecule Building Functionality**, this allows access to the libraries of atomic fragments, groups, rings, ligands and substituents, as well as the file system and the clipboard. *Clicking* on any buttons or menus in the organic, inorganic, peptide, nucleotide or substituent model kits, leads to **Add Fragment**. (If a model kit is not on screen, selection brings up the last-accessed model kit.) A fragment may be used to initiate building by *clicking* anywhere on screen, to add alongside an existing structure on screen by *clicking* on the **Insert** button at the bottom right of the screen or *pressing* the **Insert** key (**option** key on Mac) followed by *clicking* anywhere on screen, or by *double clicking* in a blank area on screen, or be added to an existing structure by *clicking* on the appropriate free valence. Fragment addition can be terminated by selection of any other function.

## Delete (✖)

This allows atom, free valence and substituent library removal from a structure. Selection leads to a message at the bottom left of the screen.

Select object to delete.

Subsequent *clicking* on an atom, free valence or substituent library results in its deletion. Deletion of an atom results in deletion of all of its associated free valences. Free valences for any atoms to which the deleted atom was previously connected are restored. Note that atom deletion may result in one or more detached fragments. *Clicking* inside a selection box results in deletion of everything inside the box. Selection of **Delete** does not bring up a model kit nor does it remove a model kit that is present on screen. **Delete** may be terminated by selection of any other function.

**Delete** may also be used to delete points and planes, CFD's and pharmacophore elements as well as spectra and plots (or individual curves that make up plots).

Deletion may also be accomplished by holding down on the **Delete** key while *clicking* on the item to be deleted.

## Make Bond (°)

This allows bonds to be drawn between free valences and/or atoms. Selection leads to a message at the bottom left of the screen.

Select two free valences.

*Clicking* on two free valences (on different atoms) will cause these atoms to be linked by a single bond. Alternatively, *double clicking* on each of two atoms will bond them, and *clicking* on a free valence on one atom and *double clicking* on a different atom will bond the two atoms. Note that available free valences are consumed as a result of bond formation, irrespective of whether free valences or atoms are selected.\* If the selected atoms are

---

\* Free valences can be “protected” without altering the molecule by adding hydrogens to them (H) from the organic model kit). One use of this is in conjunction with custom substituent libraries (see discussion of the substituent model kit earlier in this chapter).



already bonded, this will result in the bond order being increased by one, that is, single → double, double → triple. Selection of **Make Bond** does not bring up a model kit nor does it remove a model kit that is already present on screen. **Make Bond** may be terminated by selection of any other function.

### **Break Bond** ()

This allows breaking an existing bond resulting in free valences. Selection leads to a message at the bottom left of the screen.


Select bond to break.

*Clicking* on a bond breaks it and restores free valences. Note that bond breaking may result in detached fragments. Selection of **Break Bond** does not bring up a model kit nor does it remove a model kit that is present on screen. **Break Bond** may be terminated by selection of any other function.

### **Minimize** ()

This uses MMFF molecular mechanics to refine the geometry. Selection leads to a message at the bottom left of the screen.

Minimizer is active.

The molecular mechanics energy\* in kJ/mol, displayed at the bottom right of the screen, is continually updated during the minimization process. Minimization may be stopped at any time by *clicking* on the  icon at the bottom right of the screen. Any geometrical constraints imposed on the structure (see **Constrain Distance**, **Constrain Angle**, **Constrain Dihedral** in **Chapter 15**) are enforced during minimization. NOE's are not taken into account. Also, any frozen atoms in the structure (see **Freeze Center** under the **Geometry** menu; **Chapter 15**) remain frozen. Following completion, **Minimize** reverts back to **Add Fragment** only if a model kit is on screen.

---

\* The mechanics energy is a combination of the strain energy which is necessarily positive and the non-bonded or intramolecular interaction energy which can be either positive or negative. It will most commonly be a positive quantity, although it can be slightly negative.

# Chapter 17

## The Setup Menu

---

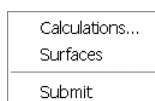
*This chapter describes functions available under the **Setup** menu. **Calculations** is used to specify molecular mechanics calculations, semi-empirical calculations, Hartree-Fock molecular orbital calculations and correlated calculations, including density functional calculations, Møller-Plesset calculations, coupled cluster calculations and quadratic configuration interaction calculations for ground-state species, and configuration interaction calculations and (time dependent) density functional calculations for excited-state species. Tasks include calculation of energy, equilibrium structure and conformation, transition-state structure and constructing energy profiles, although not all tasks are available for some methods. In addition, Hartree-Fock and density functional calculations may be extended to include an empirical solvation term, allowing “solvated” equilibrium and transition-state structures, infrared spectra and molecular properties to be obtained. (Solvated NMR, Raman and UV/visible spectra are not available.) A wide variety of all-electron Gaussian basis sets are supported for Hartree-Fock and correlated calculations as are pseudopotentials for calculations on molecules incorporating elements for which all-electron basis sets are not available. Also provided are a number of thermochemical recipes, including G3 and G3(MP2), as well as the T1 recipe that has been parameterized to closely reproduce G3(MP2) heats of formation with 2-3 orders of magnitude less computation time. **Calculations** also requests IR, Raman, NMR and UV/visible spectra, and calculation and printing of a variety of molecular properties and QSAR descriptors. Finally, **Calculations** specifies conditions for identifying similar molecules, based either on molecular structure or chemical functionality, as well as for identifying molecules that are compatible with a particular molecular environment (a pharmacophore).*

***Surfaces** is used to designate graphical surfaces, including electron*

*and spin densities, electrostatic potentials, hydride potentials, local ionization potentials and molecular orbitals, for later display as surfaces, property maps and contour plots. Accessible regions on electron density surfaces and property maps based on these surfaces may be demarked.*

***Submit** is used to direct jobs for calculation, either locally or on a remote Linux server. Dual and quad-core processors may either be used to simultaneously process different jobs or (for most tasks involving Hartree-Fock, density functional and RI-MP2 models) may be directed at a single job.*

The **Setup** menu provides access to dialogs for specifying molecular mechanics and quantum chemical calculations, for specifying conditions for similarity analysis as well as specifying surfaces and property maps and for submitting jobs for calculation.



## Calculations...

Calculations modules perform: molecular mechanics calculations using the SYBYL and MMFF force fields; MNDO (MNDO/d for second-row elements), AM1, RM1, PM3 and PM6 semi-empirical calculations, including PM3 calculations on transition-metal inorganic and organometallic systems; Hartree-Fock molecular orbital calculations; BP, BLYP, EDF1, EDF2 and B3LYP among other density functional calculations, including calculations based on a user-specified mix of existing exchange and correlation functionals; MP2 and RI-MP2 Møller-Plesset calculations; configuration interaction calculations (for excited-state species). All are applicable to energy calculation, equilibrium and transition-state geometry determination, conformational searching and making an energy profile, although some methods may not be practical for some tasks. Available for energy calculations only, are a series of high-order correlated models, including MP3 and MP4 Møller-Plesset models, coupled cluster models and quadratic configuration interaction models.

A wide selection of basis sets is available, ranging from the STO-3G minimal basis set (recommended for Hartree-Fock models only), to the 3-21G split-valence basis set (recommended for Hartree-Fock models only), to 6-31G\*, 6-311G\*, and cc-pVDZ, cc-pVTZ, cc-pVQZ and cc-pV5Z polarization basis sets. In addition, pseudopotentials for elements which all-electron basis sets do not exist are available. All basis sets, except STO-3G and 3-21G, may be supplemented with polarization functions on hydrogen, additional and/or higher-order polarization functions on heavy atoms and/or diffuse functions on hydrogens and/or heavy atoms.

Finally, a series of recipes for estimating heats of formation are available. These include the G3 and G3 (MP2) recipes as well as the T1 recipe that is based on G3(MP2) and provides nearly identical heats of formation, but is two to three orders of magnitude less costly in terms of computation.

Quantum chemical calculations also result in a variety of atomic and molecular properties and QSAR descriptors as well as IR, Raman, NMR and UV/visible spectra. IR spectra are available for molecular mechanics, semi-empirical, Hartree-Fock, density functional and MP2 and RI-MP2 models. Other spectra are available only for Hartree-Fock and density functional models.

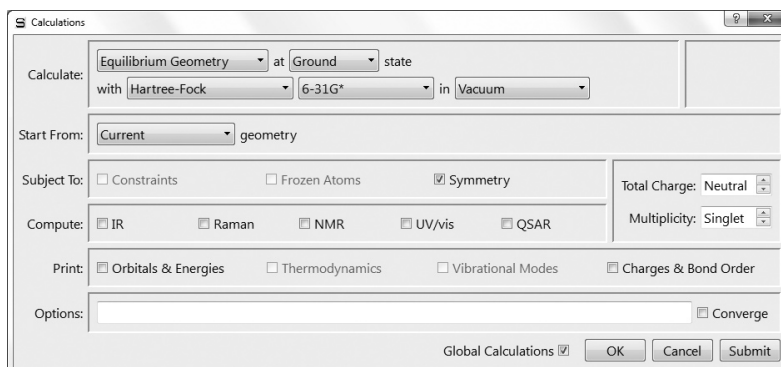
Solvent may be introduced in one of two ways, both of which involve empirical formulations. The simpler and less costly in terms of computation is to add the aqueous solvation energy obtained by the SM5.4 model to the energy of any molecular mechanics or quantum chemical calculation. Equilibrium and transition-state geometries as well as spectra and any properties derived from the wavefunction are unaffected. A second and computationally more costly method is based on the more recent SM8 model. Here, the correction affects the wavefunction in addition to the energy. This means that it can be used to obtain equilibrium and transition-state geometries as well as infrared spectra and most properties that depend on the wavefunction. The SM8 solvation correction is limited to Hartree-Fock and density functional models, and applies to a variety of solvents (including water). A listing of available solvents is provided in **Appendix C**. Solvent

Corrected Raman, NMR and UV/visible spectra are not available.

Discussion of molecular mechanics and quantum chemical methods in general, focusing on the specific methods available in *Spartan*, is provided both in **Topics** under the **Activities** menu and in *A Guide to Molecular Mechanics and Quantum Chemical Calculations* which is provided as a PDF under the **Help** menu. The guide also discusses how different methods may be combined, as well as illustrates graphical modeling techniques.

Finally, the **Calculations** dialog accesses a procedure to quantify the extent to which two molecules or a molecule and a pharmacophore are similar. This does not involve either molecular mechanics or quantum chemical calculations but depends only on molecule (pharmacophore) geometry.

Selection of **Calculations...** results in the **Calculations** dialog\*.



This contains a number of pull-down menus, buttons, check boxes and a space to enter additional information (keywords):

## Calculate

This section contains a series of menus and check boxes that specify the task to be accomplished, the electronic state (ground or first excited) the type of calculation, method, basis set and solvent to be employed, as well as details specific to each calculation type.

---

\* The dialogs associated with similarity analysis and with conformer library generation are different from that associated with the other entries under the **Calculate** menu, and will be presented later in this chapter.

## Task

Specification of **Task** is by way of a pull-down menu:

Energy
Equilibrium Geometry
Transition State Geometry
Equilibrium Conformer
Conformer Distribution
Conformer Library
Energy Profile
Similarity Analysis

**Energy** specifies calculation of energy (and in the case of quantum chemical methods, a wavefunction) at a single geometry.

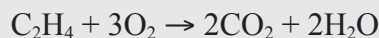
*Spartan* reports energies from molecular mechanics calculations in kJ/mol, from semi-empirical calculations and from T1, G3(MP2) and other Gx thermochemical recipes as heats of formation in kJ/mol, and from Hartree-Fock and correlated calculations as total energies in atomic units (hartrees).

The molecular mechanics energy comprises two parts: the strain energy and the non-bonded energy. The strain energy is the difference in energy between a molecule and its strain free analog. It is nearly always positive and less than a few hundred kJ/mols in magnitude. The non-bonded energy accounts for attraction or repulsion between atomic centers that are not connected due to van der Waals and Coulombic interactions. Because the strain energy of every molecule relates to a different standard, molecular mechanics energies cannot be used to obtain reaction energies (unless there are no changes in bonding between reactants and products).

The heat of formation is to the enthalpy at 298K of a balanced chemical reaction in which a molecule is converted to a set of standard products. For example, the heat of formation of ethylene is given by reaction,




where graphite and hydrogen molecule are the carbon and hydrogen standards, respectively. In practice, the actual measurement is typically carried out for a combustion reaction, for example, for ethylene:



Heats of formation may either be positive or negative quantities and generally span a range of only a few thousand kJ/mol.

Heats of formation are not suitable for presenting energy data from most quantum chemical calculations (except the Gx and T1 recipes), simply

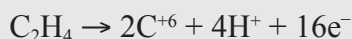because the standards for several important elements (most notably, carbon) are not well-defined isolated species. In its place is the energy of a reaction that splits a molecule into isolated nuclei and electrons, for example, for ethylene:



Total energies, as the energies of such reactions are termed, are always negative and may be very large (tens of thousands of kJ/mol). They are most commonly given in atomic units (hartrees).

$$1 \text{ atomic unit} = 2625 \text{ kJ/mol}$$

Aside from a difference in units (see **Appendix D**), it makes no difference whatsoever which standard is employed to investigate thermochemistry.

To summarize, the heat of formation differs from the total energy both with regard to the standard reaction and with regard to units. Either provides a suitable basis for thermochemical calculations.

**Equilibrium Geometry** specifies that the nearest energy minimum will be located, and **Transition State Geometry** that the nearest *first-order saddle point* (energy maximum in one dimension and energy minima in all other dimensions or commonly known as a *transition state*) will be located. **Equilibrium Conformer**, **Conformer Distribution** and **Conformer Library** attempt to characterize the conformers available to a molecule based on different criteria (see discussion following). **Energy Profile** steps along a user-defined coordinate or set of coordinates. **Similarity** quantifies the likeness among molecules based either on structure or chemical functionality or between molecules and a template (a pharmacophore).

Three different conformational analysis procedures are available in *Spartan*. The first two (**Equilibrium Conformer** and **Conformer Distribution**) may be employed with a variety of molecular mechanics and quantum chemical models (although practical concerns may limit application to molecular mechanics), while the third (**Conformer Library**) is restricted to the MMFF molecular mechanics model.

**Equilibrium Conformer** replaces the initial conformer of a molecule by the lowest-energy conformer. This procedure is typically used to get a good guess at the best (lowest-energy) conformer for calculation of properties using quantum chemical methods.

**Conformer Distribution** provides a selection of low-energy conformers. This procedure is commonly used to identify a set of energetically accessible conformers and construct a *Boltzmann distribution* for calculation of average molecular properties.

Except for molecules with very few degrees of conformational freedom full *systematic searches* (looking everywhere) are not likely to be practical. A *Monte-Carlo search* follows a path that biases in favor of low-energy conformers (but that does not completely exclude high-energy conformers). While there is no guarantee that the lowest-energy conformer (the *global minimum*) will actually be located, it can be shown that the set of conformers examined (and kept) approach a Boltzmann distribution. (This is strictly true only for a large number of search steps, but is closely approximated for searches with reasonable numbers of steps.)

**Conformer Library** replaces the initial conformer of a molecule by the lowest-energy conformer and attaches the coordinates of a set of conformers spanning all shapes accessible to the molecule without regard to energy. Used to build libraries for similarity analysis.

The following recipe is used to generate a conformer library:

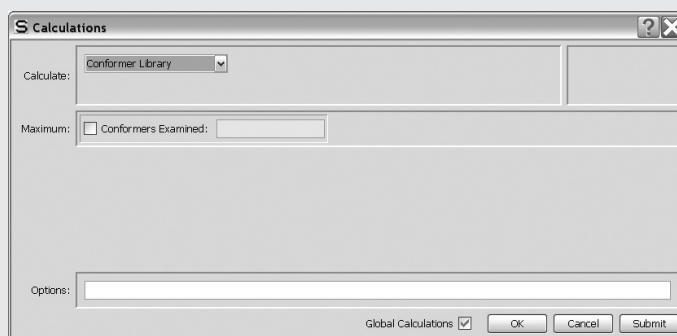
1. Step through all possible conformers in a systematic manner and calculate equilibrium structures and energies for each using the MMFF force field. Eliminate duplicates and very high-energy conformers.
2. Place spheres on each nitrogen, and oxygen atom and at the center of each hydrophobic region for each conformer.
3. Find the smallest set of conformers for which spheres do not overlap, that is, eliminate any conformer for which all of its spheres overlap with the corresponding spheres on another conformer.

Because this procedure involves a systematic search of all possible conformers, it is desirable to use a set of conformer assignment rules that eliminate degrees of freedom, but are unlikely to produce significant shape changes. See **Set Torsions** under the **Geometry** menu (**Chapter 15**) and **Settings Preferences** under the **Options** menu (**Chapter 20**). A limit may be placed on the total number of conformers actually considered. This is known as a *restricted systematic search*, a procedure that walks



through all possible conformers but randomly discards some without further minimization in order to enforce a fixed number of conformers. **Conformers Examined** sets this number and **Conformers Kept** ( $\leq$  to **Conformers Examined**) sets the number that are actually kept. This is different from the Monte-Carlo method, in that it strives to span all possible shapes open to a flexible molecule without consideration of energy.



Selection of **Conformer Library** results in a different dialog.



**Conformers Examined**, sets the maximum number of conformers.

**Energy** applies to all methods. **Equilibrium Geometry**, **Equilibrium Conformer**, **Conformer Distribution** and **Energy Profile** does not apply to advanced correlated methods and thermochemical recipes. **Transition State Geometry** does not apply to molecular mechanics methods, advanced correlated methods and thermochemical recipes. **Conformer Library** applies only to MMFF molecular mechanics. **Similarity Analysis** does not involve any molecular mechanics or quantum chemical methods.

## State

Specification of **State** (electronic state) is by way of a pull-down menu. The default setting **Ground** may be changed to **First Excited** by *clicking* on  to the right of the box and back again to **Ground** by *clicking* on . Different methods are available to handle ground and excited-state species.\* Hartree-Fock, Møller-Plesset

---

\* Note that gradients for density functional methods for excited states need to be calculated numerically. This means that determining an equilibrium geometry for an excited state of a molecule is significantly slower than determining an equilibrium geometry for the ground state. Analytical gradients for CIS and CIS(D) methods for excited states are available.

as well as a variety of advanced correlated methods are available only for ground-state species, while configuration interaction methods are available only for excited-state species. Density functional methods are available for both ground and excited-state species. **State** is ignored for molecular mechanics calculations, and semi-empirical calculations are limited to the ground state.

## Type of Calculation

Specification of **Type of Calculation** is by way of a menu, the selections in which depend on electronic state. For ground states:

Molecular Mechanics
Semi-Empirical
Hartree-Fock
Density Functional
Møller Plesset
Advanced Correlated
Thermochemical Recipes

For excited states:

Density Functional
Configuration Interaction

Some calculation types require additional information:

Molecular Mechanics	<table border="1"><tr><td>MMFF</td></tr><tr><td>MMFFaq</td></tr><tr><td>SYBYL</td></tr></table>	MMFF	MMFFaq	SYBYL	Møller-Plesset	<table border="1"><tr><td>MP2</td></tr><tr><td>RH-MP2</td></tr></table>	MP2	RH-MP2				
MMFF												
MMFFaq												
SYBYL												
MP2												
RH-MP2												
Semi-Empirical	<table border="1"><tr><td>AM1</td></tr><tr><td>RM1</td></tr><tr><td>PM3</td></tr><tr><td>PM6</td></tr><tr><td>MNDO</td></tr></table>	AM1	RM1	PM3	PM6	MNDO	Advanced Correlated**	<table border="1"><tr><td>MP3</td></tr><tr><td>MP4</td></tr><tr><td>CCSD</td></tr><tr><td>CCSD(T)</td></tr></table>	MP3	MP4	CCSD	CCSD(T)
AM1												
RM1												
PM3												
PM6												
MNDO												
MP3												
MP4												
CCSD												
CCSD(T)												
Density Functional*	<table border="1"><tr><td>BP</td></tr><tr><td>BLYP</td></tr><tr><td>EDF1</td></tr><tr><td>EDF2</td></tr><tr><td>B3LYP</td></tr><tr><td>MD6</td></tr><tr><td>wB97X-D</td></tr></table>	BP	BLYP	EDF1	EDF2	B3LYP	MD6	wB97X-D	Thermochemical Recipes***	<table border="1"><tr><td>G3(MP2)</td></tr><tr><td>T1</td></tr></table>	G3(MP2)	T1
BP												
BLYP												
EDF1												
EDF2												
B3LYP												
MD6												
wB97X-D												
G3(MP2)												
T1												
Configuration Interaction	<table border="1"><tr><td>CIS</td></tr><tr><td>CIS(D)</td></tr></table>	CIS	CIS(D)									
CIS												
CIS(D)												

---

\* A variety of other functionals are available and may be entered directly into the **Options** box. A full listing is provided in **Appendix A**. Additionally, any combination of supported exchange and correlation functionals may be specified.

\*\* A number of other advanced correlated methods are available, but need to be entered directly into the **Options** box: CCSD(2), OD, OD(T), QCCD, QCCD(2), QCISD and QCISD (T).

\*\*\* G2 and G3 recipes are also available but need to be entered into the **Options** box.

The most important application of **molecular mechanics models** is conformational searching, in particular, for molecules with several degrees of conformational freedom. In these cases, even the simplest quantum chemical models may not be practical. The MMFF force field is known to properly assign equilibrium conformation in a variety of simple molecules for which experimental data are available. It also has been shown to provide a reasonable account of conformational energy differences in larger organic molecules as obtained from high-level quantum chemical calculations. Aside from conformational energy differences, molecular mechanics models are not suitable for thermochemical calculations.

**Semi-empirical models** are simplest methods based on quantum mechanics. They are applicable to molecules containing 100 - 200 atoms, but not to molecules containing thousands of atoms, for example, proteins. The MNDO, AM1, RM1, PM3 and PM6 models generally provide geometries in good accord with experimental structures. This makes them suitable for evaluation of such properties as polar surface area that depend solely on geometry. However, none of the present-generation semi-empirical models is suitable for the calculation of relative energies as might be required to conclude, whether a chemical reaction is weakly or strongly *exothermic*, thermoneutral or weakly or strongly *endothermic*, or whether one isomeric product of a reaction is likely to be more or less stable than another, or to identify the lowest-energy tautomer. Also, none of the methods reliably accounts for conformational energy differences, as would be required to establish whether a particular conformer has a good chance at actually being present or to calculate the average value of a property based on the relative (Boltzmann) populations of different conformers.

**Hartree-Fock models** follow from the **Schrödinger equation** by requiring that the electrons be independent particles. This is known as the **Hartree-Fock approximation**. Here, the motions of electrons in molecules (molecular orbitals) are approximated by a sum of the motions of electrons in atoms (atomic orbitals). A second approximation, termed the **LCAO** or **Linear Combinations of Atomic Orbital approximation**, distinguishes different Hartree-Fock models. Different models use different basis sets, that is, different numbers and kinds of atomic orbitals. For example, the model with the 3-21G split-valence basis set uses one atomic orbital to describe each non-valence electron and two atomic orbitals to describe each valence (core) electron, that is, it comprises one 1s atomic orbital and two sets of 2s and 2p atomic orbitals.

1s  
2s, 2p<sub>x</sub>, 2p<sub>y</sub>, 2p<sub>z</sub>  
2s', 2p'<sub>x</sub>, 2p'<sub>y</sub>, 2p'<sub>z</sub>

The model with the 6-31G\* polarization basis set increases the flexibility by including d type atomic orbitals that, while not occupied in the atom, are used in molecules.

1s  
2s, 2p<sub>x</sub>, 2p<sub>y</sub>, 2p<sub>z</sub>  
2s', 2p'<sub>x</sub>, 2p'<sub>y</sub>, 2p'<sub>z</sub>  
d<sub>xx</sub>, d<sub>xy</sub>, d<sub>xz</sub>, d<sub>yy</sub>, d<sub>yz</sub>, d<sub>zz</sub>

Larger basis sets involving even higher splitting of valence shells, including higher-order (f and g type) functions as well as diffuse functions can be and have been defined and are available in *Spartan*. Description of their makeup is provided in **Topics** available under the **Activities** menu and in *A Guide to Molecular Mechanics and Quantum Chemical Calculations* which is provided as a PDF under the **Help** menu.

3-21G and 6-31G\* and larger basis set Hartree-Fock models generally provide good accounts of the geometries of organic molecules. This extends to unusual or unstable systems for which experimental data may be limited (and hence may be poorly represented in parameterization of the semi-empirical models). The 6-31G\* model (and to a lesser extent, the 3-21G model) also properly accounts for thermochemistry, at least insofar as being able to say whether a particular reaction is weakly or strongly *exothermic* or weakly or strongly *endothermic*. It is also able to properly order the stabilities of isomeric products and to reliably identify the lowest-energy tautomer. Finally, 6-31G\* calculations provide a qualitatively correct account of the relative stabilities of different conformational arrangements. Hartree-Fock models with larger basis sets show similar behavior.

The most conspicuous and most important flaw of Hartree-Fock models is their inability to properly account for the energetics of bond dissociation. This is the case even with very large basis sets, and can be traced to the Hartree-Fock approximation which in effect replaces interactions between electrons by interactions between one electron and the average field created by all other electrons. As a consequence, electrons get in each other's way to a greater extent than they should. This leads to overestimation of the electron-electron repulsion energy, most significantly for electrons that are paired. Because bond dissociation reduces the number of electron-pairs (by one), the overestimation will

be greater for reactants than for products and the calculated dissociation energy will be too small.\*

Models that do away with the Hartree-Fock approximation, or at least lessen its effect, are termed **correlated models**. These divide into two broad categories, **density functional models** and **wavefunction-based models**. Density functional models explicitly introduce an approximate correlation term. They are not much more costly in terms of computation time than Hartree-Fock models, although it is not apparent how to improve on a particular choice of correlation correction. Wavefunction-based models, typified by **Møller-Plesset models** and **configuration interaction models** extend the flexibility of Hartree-Fock models by mixing ground-state and excited-state wavefunctions. They are significantly more costly than both Hartree-Fock and density functional models, but offer the advantage over the latter of a clear path to improvement. In the limit of complete mixing, so-called full configuration interaction and infinite order Møller-Plesset models, wavefunction-based models lead to the exact result, although in practice this limit cannot be reached.

Density functional and wavefunction-based correlated models make use of the same basis sets as Hartree-Fock models, except that minimal and split-valence basis sets do not yield satisfactory results. Among the most popular models are the B3LYP/6-31G\* density functional model and the MP2/6-31G\* (second-order) Møller-Plesset model. The primary use of these models (and other correlated models) is to obtain thermochemistry, including for reactions in which bonds are made or broken and including activation energies.

While geometry is well described with simple quantum chemical models, energy (or more precisely relative energy) has proven to be more difficult to calculate accurately. Hartree-Fock models properly describe the energies of *isodesmic* reactions, where formal bond count is maintained but are less successful for reactions in which reactants and products have markedly different bonding. Indeed, high-level correlated calculations with large underlying basis sets are needed to properly account for the energetics of bond breaking or to relate the energies of reactants to that of a transition state.

A number of thermochemical recipes or combinations of different quantum chemical methods have been formulated to reproduce experimental heats

---

\* Consider the bond dissociation energy of H<sub>2</sub>. The Hartree-Fock energy of the products (hydrogen atoms each with only one electron) is exact, while the energy of the reactant (H<sub>2</sub> with two electrons) is too large. Therefore, the bond dissociation energy will be too small.

of formation to within 4-8 kJ/mol. (Except for very simple molecules, experimental heats are seldom more accurate than this.) G3(MP2) is the simplest of the recipes that have been formulated. It reproduces experimental heats of formation for a wide variety of small molecules to within 6 kJ/mol (mean absolute error), with only a small number of systems showing errors twice this amount. The problem with G3(MP2) is its large computational cost, which limits the size of systems that can be dealt with. In practice, G3(MP2) may easily be applied to molecules comprising ten heavy (non-hydrogen) atoms and with great difficulty to molecules comprising fifteen or more heavy atoms.

Three different energies are required in the G3(MP2) recipe:

$E_A$  MP2/6-31G\* using the MP2/6-31G\* equilibrium geometry

$E_B$  MP2/6-311++G(2df, 2p) using the MP2/6-31G\* equilibrium geometry

$E_C$  QCISD(T)/6-31G\* using the MP2/6-31G\* equilibrium geometry

These are combined in order to take into account both the effects of increased basis set size and of electron correlation beyond the MP2 level. Note, that the two effects are treated independently in order to avoid having to perform a QCISD(T) calculation using a large basis set.

$$E = E_A + (E_B - E_A) + (E_C - E_A)$$

Two corrections are needed to bring the calculated energy in line with what is actually measured experimentally: a correction for zero-point energy;  $E_{\text{zero-point}}$ , and a correction to the enthalpy (energy) to account for finite temperature (298K),  $H(T)$ .

$$E_{\text{zero-point}} = \frac{1}{2} h \sum_i^{\text{vibrations}} \nu_i$$

$$\Delta H(T) = H_{\text{trans}}(T) + H_{\text{rot}}(T) + \Delta H_{\text{vib}}(T) + RT$$

$$H_{\text{trans}}(T) = \frac{3}{2} RT$$

$$H_{\text{rot}}(T) = \frac{3}{2} RT \quad (\text{RT for a linear molecule})$$

$$\Delta H_{\text{vib}}(T)^* = H_{\text{vib}}(T) - H_{\text{vib}}(T) - H_{\text{vib}}(0)$$

$$= Nh \sum_i^{\text{vibrations}} \frac{\nu_i}{(e^{h\nu_i/kT} - 1)}$$

$h$  is Planck's constant,  $k$  is the Boltzmann constant,  $R$  is the gas constant and  $N$  is Avogadro's number. Both zero-point and temperature corrections

require vibrational frequencies.\* These are obtained from Hartree-Fock/6-31G\* calculations and require the corresponding equilibrium geometry.

To summarize, the G3(MP2) recipe requires two equilibrium geometry calculations (at HF/6-31G\* and MP2/6-31G\*), a HF/6-31G\* frequency calculation, a MP2/6-311++G (2df, 2p) energy calculation and a QCISD(T)/6-31G\* energy calculation. The QCISD(T) calculation dominates with increasing molecular size. Next most costly is the MP2/6-31G\* geometry optimization, although both the large basis set MP2 energy and the HF/6-31G\* frequency may be significant.

T1 is a simplified recipe based on G3(MP2) and intended to closely reproduce G3(MP2) heats of formation. It eliminates the QCISD(T)/6-31G\* energy calculation ( $E_C$ ), the zero-point energy calculation ( $E_{\text{zero-point}}$ ) and the temperature correction to the enthalpy ( $\Delta H(T)$ ), and replaces the large basis set MP2 energy calculation ( $E_B$ ) with a dual-basis-set RI-MP2 calculation, where G3MP2large is the large basis set (the same basis set that is used in the G3(MP2) recipe) and 6-311G\* is the small basis set. To compensate, an empirical correction based on the Hartree-Fock and RI-MP2 energies and Mulliken bond orders is introduced. This follows from a linear regression analysis involving G3(MP2) heats of formation for  $\approx 1,050$  molecules. T1 uses HF/6-31G\* equilibrium geometries, instead of MP2/6-31G\* geometries.

An automated procedure has been developed to locate the best (lowest heat of formation) conformer, or at least a conformer that is very close to the best. This involves four steps: (i) examine up to 1,000 conformers using a restricted systematic search procedure and the MMFF molecular mechanics model; (ii) obtain equilibrium geometries for the four best conformers, making certain to exclude duplicate conformers and conformers that are nearly the same using the HF/6-31G\* model; (iii) perform (up to) two T1 calculations on the two best HF/6-31G\* structures; (iv) select the lower energy T1 conformer. The default numbers of structures examined in the second and third steps (four and two, respectively) may be altered by keywords typed into the **Options** box (**Appendix C**).

---

\* Note, that the expression for  $\Delta H_{\text{vib}}(T)$  goes to  $RT$  instead of  $1/2 RT$  as  $\nu_i$  goes to zero. This requires that contributions that would exceed  $1/2 RT$  be set to  $1/2 RT$ .

Except for semi-empirical methods and for thermochemical recipes, either **restricted** or **unrestricted scf methods** may be used. The former restricts paired electrons to the same orbital, whereas the latter allows the electrons to occupy different orbitals. The default (singlets with restricted methods and non-singlets with unrestricted methods) may be overridden using the **Scf** keyword from the **Options** box (**Appendix C**).

## Basis Sets

Except for molecular mechanics and semi-empirical methods and for thermochemical recipes, a basis set needs to be specified. Menus provide a number of choices.\* For Hartree-Fock methods:

```
STO-3G
3-21G
6-31G*
6-31G**
6-31+G*
6-311G*
6-311+G**
```

For all other methods:

```
6-31G*
6-31G**
6-31+G*
6-311G*
6-311+G**
6-311++G**
6-311++G(2df,2p)
cc-pVTZ
```

The 3-21G basis set for second-row and heavier main-group elements incorporates a set of d-type functions. While these are not occupied in the atomic ground state, they have been found necessary for proper description of bonding in molecules. 6-311++G (2df,2p) uses the 6-311++G (3d2f,2p) basis set for second-row elements.

The combination of a method, for example, B3LYP, and a basis set, for example, 6-31G\*, constitutes a **theoretical model** or more simply a **model**. The nomenclature is to separate method and basis set by a “/”, for example, B3LYP/6-31G\*. By convention, specification of Hartree-Fock

---

\* A number of other basis sets are available, and may be entered directly into the **Options** box: These include the Pople 6-31G or 6-311G basis sets followed by: (i)j(f) or (i)j(f, k)l(d) where *i* is the number of sets of d functions on non-hydrogen atoms, *j* is the number of sets of f functions on non-hydrogen atoms, *k* is the number of s functions on hydrogen atoms and *l* is the number of sets of p functions on hydrogen atoms. If 0 eliminate p, d, f; if 1 eliminate the number; 6-31G or 6-311G with “+” or “++” inserted before G; cc-pVDZ, cc-pVQZ and cc-pV5Z.



(HF) as a method is optional. That is, specification of basis set alone refers to Hartree-Fock models, for example, both HF/6-31G\* and 6-31G\* refer to combination of the Hartree-Fock method and the 6-31G\* basis set. Methods without basis sets, for example, MMFF (molecular mechanics) and PM3 (semi-empirical) are themselves models.

## Pseudopotentials

Pseudopotentials are employed only where all-electron basis sets are unavailable, and are used in conjunction with all-electron basis sets for elements where they are available. This is primarily for elements heavier than Kr, most importantly fourth-row main-group elements and second and third-row transition metals. For example, a B3LYP/6-31G\* calculation on ruthenium ethylene tetracarbonyl would employ a pseudopotential only for ruthenium and the all-electron 6-31G\* basis set for hydrogen, carbon and oxygen. The overall method is referred to as **B3LYP/6-31G\***. Pseudopotentials pertain to Hartree-Fock, density functional, MP2 and advanced correlated calculations only, but not to RI-MP2 calculations or calculations using the STO-3G and 3-21G basis sets.

## Solvent

A solvent is specified for Hartree-Fock and density functional calculations by replacing *vacuum* in the menu to the right of **in** by water or one of a selection of common organic solvents.\*

Vacuum
Water
Acetone
Diethyl Ether
DMF
DMSO
Ethanol
Methylene Chloride
THF
Toluene

The effect of the solvent on the energy and wavefunction (and hence on any properties derived from the wavefunction) is approximated by way of the SM8 model. This has been parameterized to reproduce experimental

\* Other solvents are available and can be entered via the **Solvent** keyword in the **Options** box.

energies of solvation and related quantities. It has not been parameterized to reproduce equilibrium geometries and properties other than the energy, and there is no solid calibration of its success or failure for anything other than energy. It should be used with caution.

The cost of SM8 in terms of computation time is roughly equal to that of an energy calculation. Therefore, solvated Hartree-Fock and density functional energy calculations are twice as costly as the corresponding calculation in vacuum while geometry calculations are roughly 30% more costly.

If *checked*, **Dual Basis** signifies that a Hartree-Fock, density functional, MP2 or RI-MP2 energy equilibrium and transition-state geometry calculations with the 6-311+G\*\*, 6-311++G\*\*, 6-311++G (2df, 2p), or cc-pVTZ basis set will be approximated by a calculation in which scf convergence is first achieved using the smaller 6-311G\* basis set and then corrected perturbatively for the effects of basis set extension. Other combinations of small and large basis sets may be specified under **Options** (see discussion following). Use of dual basis has only a small effect on calculated relative energies, but can reduce the cost of large-basis-set energy calculations by two to ten times (depending on the choice of the large basis set). A dual basis set is employed for calculations with the T1 thermochemical recipe.

### **Start from (Energy, Equilibrium Geometry)**

Energy and equilibrium geometry calculations (only) from one model may be preceded by equilibrium geometry calculations at another (simpler) model. Choices depend on the method selected:

**Molecular Mechanics**

not applicable

**Semi-Empirical**

Current  
MMFF  
MMFF Conformer  
MMFFaq Conformer

**Hartree-Fock**

Current  
MMFF  
AM1  
PM3  
MMFF Conformer  
MMFFaq Conformer

## Density Functional, Møller-Plesset, Advanced Correlated, Thermochemical Recipes

Current
MMFF
AM1
PM3
3-21G
6-31G*
MMFF Conformer
MMFFaq Conformer

Selection of **MMFF conformer** or **MMFFaq conformer** signifies that geometry optimization will follow a search for the lowest-energy conformer. These lead to an additional menu specifying the level of geometry optimization.

**Molecular Mechanics, Semi-Empirical**      not applicable

**Hartree-Fock**

MMFF
AM1
PM3
3-21G

**Density Functional, Møller-Plesset, Advanced  
Correlated, Thermochemical Recipes**

MMFF
AM1
PM3
3-21G
6-31G*

### Maximum (Equilibrium Conformer, Conformer Distribution)

**Maximum Conformers Examined** sets a maximum for the number of possible conformers that will be considered both for **Equilibrium Conformer** and **Conformer Distribution** tasks. Note, that if a systematic search has been requested (**SearchMethod** keyword; **Appendix C**) this invokes an algorithm whereby systematic conformers are randomly eliminated in order to enforce the conformer limit.

**Maximum Conformers Kept** sets a maximum for the number of conformers returned at the end of a **Conformer Distribution** task. It must be smaller or equal to the number of conformers examined.

### Subject to

*Spartan* allows calculations to be carried out in the presence of geometrical constraints or with atoms which have been frozen in place. In addition, *Spartan* attempts to make use of symmetry.

If *checked*, **Constraints** signifies use of any previously defined constraints on distances, angles and dihedral angles

into equilibrium and transition-state geometry optimization, conformation searching and generation of energy profiles. Does not apply to energy calculation, conformer library generation or similarity analysis. See **Constrain Distance**, **Constrain Angle** and **Constrain Dihedral** under the **Geometry** menu (**Chapter 15**) for information about constraining geometrical parameters.

NOE data are not treated as constraints, but rather as post-calculation filters in establishing conformer distributions.

If *checked*, **Frozen Atoms** signifies that the coordinates of any atoms that have previously been frozen will not be moved during equilibrium and transition-state geometry optimization, conformation searching and generation of energy profiles. Does not apply to energy calculation, conformer library generation or similarity analysis. See **Freeze Center** under the **Geometry** menu (**Chapter 15**) for information about freezing atoms. See also discussion in **Chapter 19** regarding the freezing of heavy atoms in structures resulting from a search of the Cambridge Structural Database.

If *checked*, **Symmetry** signifies that molecular symmetry is to be employed wherever possible to simplify the calculation. By default **Symmetry** is checked.

## Compute

Entries under this section request calculation of infrared, Raman, NMR and UV/visible spectra as well as a variety of descriptors used in QSAR analyses.

### IR

If *checked*, calculates vibrational frequencies and infrared intensities together with the corresponding vibrational modes. These are available in the output (**Output** under the **Display** menu; **Chapter 18**) along with selected thermodynamic

properties (entropies\* and Gibbs energies). Vibrational modes may be animated and an IR spectrum displayed from the **IR** dialog accessible from **Spectra** under the **Display** menu (**Chapter 18**). Second derivatives required for frequency calculation are evaluated analytically for molecular mechanics, semi-empirical, Hartree-Fock and density functional models, but involves numerical differentiation of analytical gradients for RI-MP2 and MP2 models, this is much more costly in terms of computation. Frequency calculations for higher-order correlated models are not available.

## Raman

If *checked*, calculates Raman intensities (**IR** must also be checked in order for vibrational frequencies to be available.) The Raman spectrum may be displayed from the **Raman** dialog accessible from **Spectra** under the **Display** menu (**Chapter 18**).

## NMR

If *checked*, calculates NMR chemical shifts. These are then available in the output (**Output** under the **Display** menu; **Chapter 18**) as well as from the **Atom Properties** dialog (accessible from **Properties** under the **Display** menu; **Chapter 18**) and as atom labels (**Configure...** under the **Model** menu; **Chapter 14**). <sup>13</sup>C (proton decoupled) and <sup>1</sup>H (shifts only) spectra from Hartree-Fock or density functional calculations may be displayed from the **NMR Spectra** dialog accessible from **Spectra** under the **Display** menu (**Chapter 18**)\*\*. <sup>13</sup>C spectra may also be displayed as DEPT spectra (carbons with one or three attached protons point up, carbons with two attached protons point down and quaternary carbons are

---

\* Note that the linear harmonic oscillator approximation incorrectly leads to a vibrational contribution to the entropy that goes to  $\infty$  as frequency goes to 0. Vibrational entropy contributions from frequencies below 300 cm<sup>-1</sup> should be treated with caution.

\*\* Chemical shifts for other nuclei are available in the **Output** dialog (**Output** under the **Display** menu) and may also be attached as labels (**Configure...** under the **Model** menu). All chemical shifts are given relative to standard compounds, for example, TMS in the case of both proton and carbon. Standards are listed in **Appendix M**.

absent). Finally, 2D COSY (proton vs. proton) and HSQC and HMBC (proton vs  $^{13}\text{C}$ ) spectra are available from the dialog.

$^{13}\text{C}$  chemical shifts calculated from the B3LYP/6-31G\* model will be empirically corrected for local topology (numbers and kinds of directly-bonded atoms), and  $^{13}\text{C}$  shifts calculated from the EDF2/6-31G\* model will be corrected using a more elaborate scheme involving both topology and bond orders to neighboring atoms. Uncorrected shifts are available in the printed output. Line intensities are assumed to be proportional to the number of equivalent carbons or hydrogens. HH and CH coupling constants are not calculated, but are estimated empirically from the molecular geometry.

## UV/vis

If *checked*, specifies that a series of excited-state energy calculations will be performed (following a ground-state calculation) and that state-to-state energy differences will be reported. A UV/visible spectrum may be displayed from the **UV/vis Spectra** dialog accessible from **Spectra** under the **Display** menu (**Chapter 18**). The default number of excited states is 6, but may be changed using the **UV States** keyword from the **Options** box (**Appendix C**). UV/visible spectra calculations are restricted to Hartree-Fock models (CIS for the excited state) and density functional models (time-dependent density functional theory for the excited state).


## QSAR

If *checked*, specifies calculation of a number of QSAR descriptors that would not otherwise be available. These include descriptors based on the electron density surface, the electrostatic potential map and the local ionization potential map, as well as polarizability and logP.\*


---

\* LogP is estimated according to the method of Ghose, Pritchett and Crippen (*J. Computational Chem.*, **9**, 80 (1988)). An alternative procedure due to Villar (*J. Computational Chem.*, **6**, 681 (1991); *Int. J. Quantum Chem.*, **44**, 203 (1992)) is available. See **Table 18-3**.

## Total Charge

Total charge. The default setting **Neutral** may be changed either by *clicking* on  to the right of the box, and selecting **Anion**, **Dianion**, **-3**, etc. from the menu which appears, or by replacing any menu entry by a numerical value. **Total Charge** is ignored for molecular mechanics calculations.

## Multiplicity

Spin multiplicity (number of unpaired electrons). The default setting **Singlet** may be changed either by *clicking* on  to the right of the box, and selecting instead **Doublet** or **Triplet** from the menu which appears, or by replacing any menu entry by a numerical value. Multiplicity is 1 for singlets (no unpaired electrons), 2 for doublets (one unpaired electron), 3 for triplets (two unpaired electrons), 4 for quartets (three unpaired electrons), etc. **Multiplicity** is ignored for molecular mechanics calculations.

## Print

If *checked*, writes the quantity to the output window. Text output may be seen by selecting **Output** from the **Display** menu (**Chapter 18**) and printed by first selecting the output dialog and then selecting **Print Output** from the **File** menu (**Chapter 12**). Additional printing may be requested from the **Options** box (see discussion following).

### Orbitals & Energies

If *checked*, writes the orbitals and energies to the output. Orbital energies are also available in the spreadsheet (**Chapter 18**). An orbital energy diagram is available from **Properties** under the **Display** menu (**Chapter 18**).

### Thermodynamics

If *checked*, writes thermodynamic quantities to the output. Requires that vibrational frequencies (**IR** under **Compute**) have been calculated. Selected thermodynamic quantities are also available in the spreadsheet (**Spreadsheet** under the

**Display** menu; **Chapter 18**).

### **Vibrational Modes**

If *checked*, writes vibrational frequencies and modes to the output. Requires that **IR** under **Compute** have been calculated.

### **Atomic Charges**

If *checked*, writes Mulliken, natural and electrostatic-fit charges to the output. Note that Mulliken and natural charges from semi-empirical models are identical. Atomic charges are also available from the **Atom Properties** dialog (accessible from **Properties** under the **Display** menu; **Chapter 18**), in the spreadsheet (**Spreadsheet** under the **Display** menu; **Chapter 18**) and as atom labels (**Configure...** under the **Model** menu; **Chapter 14**). Atomic charges are not available for calculations using advanced correlated models or for thermochemical recipes.

### **Options**

Program options may be specified using keywords entered into the **Options** box. Keywords may be either single words or expressions. Keyword=N indicates an integer argument, and keyword=F indicates a real argument. **Appendix C** contains a listing of commonly-used keywords.

### **Converge**

If *checked*, invokes procedures to assist in convergence of quantum chemical calculations. May lead to an increase in computation time.

### **Similarity**

*Spartan* provides a procedure to assess and quantify the similarity of one or more molecules in a document to one or more molecules in one or more different documents. Similarity may either be based on molecular structure (geometry) or on chemical function descriptors (CFD's). *Spartan* also assesses the similarity of one or more pharmacophores in a document to one or more molecules in one or more different documents.



Three different permutations of similarity analysis are available. Here, *template* refers to the selected document in the graphical user interface and *library* refers to one or more different documents specified in the **Calculations** dialog.

<i>template</i>	<i>library</i>
pharmacophores	molecular structures for which CFD's are assigned on-the-fly
molecular structures	molecular structures
CFD's from molecular structures	molecular structures for which CFD's are assigned on-the-fly

Two types of libraries are supported for molecular structures (and CFD's generated on-the-fly from molecular structures). In the first, each entry represents a single molecule. In the second, each entry represents a selection of different conformers of a single molecule and spanning all possible shapes that the molecule may take on. The latter are generated using **Conformer Library** under the **Calculate** menu (see discussion earlier in this chapter), and *may only be used as library entries in similarity analysis and not as templates*.

Of the three cases for similarity analysis, only the second makes direct use of molecular structures (selected atomic centers for the molecules in the template) as a similarity measure. The other two libraries make use of CFD's generated and (optionally) edited for the template and generated on-the-fly in the library. Selection of structure or CFD as the basis of similarity measure must be explicitly stated (see discussion following).

Seven different chemical function descriptors (CFD's) are recognized by the similarity analysis procedure, the first six of which are attributes of a molecule while the seventh derives from knowledge of a molecule bound to a protein or nucleic acid.

hydrophobe (sterically crowded region)  
 aromatic hydrophobe (aromatic ring)  
 hydrogen-bond donor  
 hydrogen-bond acceptor

positive ionizable site  
negative ionizable site  
excluded volume

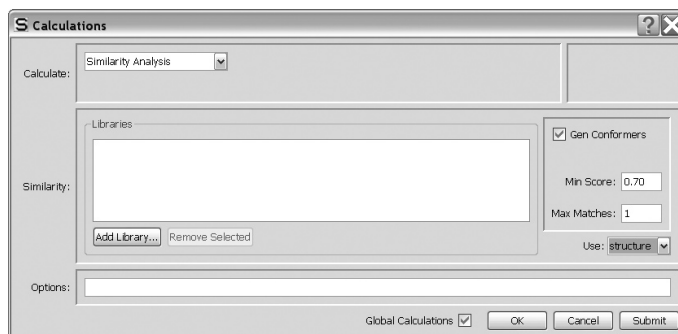
Discussion of CFD types has already been provided (see **CFD's** under the **Model** menu; **Chapter 14**).

The individual CFD's for molecules in the template (only) may either be used or ignored in the similarity procedure. This is specified from **Set Similarity Centers** under the **Geometry** menu (**Chapter 15**). The **CFD Properties** dialog (accessible from **Properties** under the **Display** menu; **Chapter 18**) allows for changing or extending the definition of a CFD. For example, the nitrogen in a dimethylamino group might initially be assigned as a hydrogen bond acceptor CFD, but could be changed to a positive ionizable site (as such a center could become positively-charged at biological pH) or could be extended to cover both possibilities.

*Spartan* assumes that any chiral molecules in a library are unresolved (they comprise an equal mixture of enantiomers). By default, similarity comparisons are carried out on both the library molecules as they are stored and (if the molecules are chiral) on their mirror images. The default procedure may be overridden (providing comparisons only to the stored molecules) by way of the **SingleEnantiomer** keyword entered in the **Options** box. See discussion later in this chapter and also **Appendix C**.

Note that only a single enantiomer of a chiral template molecule is considered. If comparisons involving both enantiomers are desired, the second enantiomer must be added to the query.

Selection of **Similarity** from the **Calculate** menu with either a molecule or pharmacophore (list of molecules/pharmacophores) leads to a dialog for setting up a similarity analysis.



One or more libraries need to be identified. This is accomplished by *clicking* on the **Add Library...** button, and gives rise to a file browser. Once selected, *clicking* on **Open** adds the library to the list of templates. A library may be deleted from the list by *clicking* on its name in the list and *clicking* on the **Remove Selected** button at the bottom of the dialog. **Shift** and **Ctrl** keys function in the usual manner in both library addition and removal.

Prior to starting a similarity analysis, you need to select either **structure** or **CFD** from the **Use** menu near the bottom of the **Calculations** dialog.

Similarity analysis is most valuable where the template molecule is conformationally flexible, and it is not apparent which (if any) of the different shapes it adopts might match a particular target in a library. It becomes indispensable where the target is also flexible and its conformation is not established. Here, the process must be repeated for every template conformer-target conformer pair, that is, thousands or tens of thousands of times.

This process can be shortened if the purpose is simply to identify which target molecules a particular template molecule might match, irrespective of conformer. In this case, the first good match concludes the analysis, and the user is left with the knowledge that template and target can be made to match (together with an example of a match).



The similarity algorithm incorporated in *Spartan* attempts to examine template-target pairs that are most likely to lead to matches before examining pairs that are less likely. It also allows early termination upon finding a single

good match, or a user-specified number of good matches\*. Taken together, these can lead to one or two orders of magnitude speedup in the analysis.

\* Criteria for early termination are specified using program options. See **Similarity Analysis Options** in **Appendix C**.

## Global Calculations

If *checked*, signifies that settings in the **Calculations** dialog are to be applied to all molecules in the document.

The **Calculations** dialog may be exited by clicking on **Submit**, **Cancel** or **OK** at the bottom right of the dialog, or on  at the top. (**Submit** and **OK** are not available if the job is already executing.) Clicking on **OK** or on **Submit** overwrites any previous information. Additionally, **Submit** enters the job in the execution queue (see discussion later this chapter). Clicking on **Cancel** or on  exits the **Calculations** dialog without saving any changes.

## Surfaces

*Spartan* allows graphical display of the HOMO and LUMO among other molecular orbitals, the electron density, the spin density for molecules with unpaired electrons, the electrostatic potential, the hydride potential and the local ionization potential.

The *electron density* is the number of electrons found at a point in space. It is the quantity measured in an X-ray diffraction experiment that is then used to locate atomic positions, that is, most electrons are closely associated with atoms. While the electron density is non-zero everywhere, it is possible to define surfaces of constant density. The most important of these contains most of a molecule's electrons and that roughly corresponds to a space-filling model, that is, a van der Waals surface. We will refer to this as the *electron density*. It is interesting because it reveals overall molecular size and shape and demarks the steric barrier seen by encroaching molecules. Another important surface, that we will refer to as the *bond density*, contains fewer electrons in total and demarks atomic connectivity.

The *spin density* is the difference in the number of electrons of  $\alpha$  and  $\beta$  spin at a point in space. It indicates the location of the unpaired electron in

a radical or unpaired electrons in a triplet or higher multiplicity state.

The **electrostatic potential** is the energy of interaction of a positive charge with a molecule. This assumes a fixed electron distribution for the molecule. It represents a balance between repulsive interactions involving the positively-charged nuclei and attractive interactions involving the negatively-charged electrons. Regions where the balance tips toward attraction are said to be electron rich (basic) and subject to attack by electrophiles, while regions where the balance tips toward repulsion are said to be electron poor (acidic) and subject to attack by nucleophiles. Electron-rich regions such as lone pairs are typically located outside the van der Waals surface. As such, they may be easily identified by constructing a surface of negative (attractive) electrostatic potential. While interesting electron-poor areas such as acidic hydrogens also lie outside the van der Waals surface, the electrostatic potential is also positive (repulsive) throughout the region inside this surface.

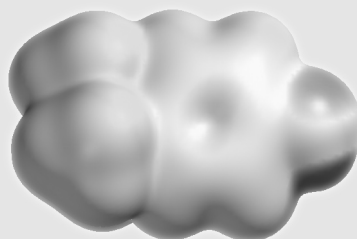
The **hydride potential** is the energy of interaction of hydride anion with a molecule. Like the electrostatic potential, this assumes a fixed electron distribution for the molecule and represents a balance between attractive and repulsive interactions. Regions where the balance tips toward attraction are said to be electron poor (acidic) and subject to attack by nucleophiles, while regions where the balance tips toward repulsion are said to be electron rich (basic) and subject to attack by electrophiles.

The **local ionization potential** indicates the ease or difficulty of electron removal (ionization). Like the negative regions of the electrostatic potential, regions of low local ionization potential are likely to be subject to attack by electrophiles.

Note that none of the last three quantities are experimental observables, although all relate to quantities that can be given clear chemical interpretation.


Additionally, any one of the quantities listed above (except the electron density) may be mapped onto any surface (except a molecular orbital surface). In practice, the only maps to have received widespread attention are those based on the electron density surface (depicting overall molecular size and shape). Most common are the electrostatic potential map, the hydride potential map, the local ionization potential map and the LUMO map.

The **electrostatic potential map** paints the value of the electrostatic potential onto an electron density surface. By convention, colors toward red depict negative potential, while colors toward blue depict positive potential, and colors in between (orange, yellow, green) depict intermediate values of the potential. Thus, an electrostatic potential map for *p*-*tert*-butylphenol will show oxygen to be red, its attached (acidic) hydrogen to be blue, the  $\pi$  faces of benzene to be orange or yellow and the *tert*-butyl group to be green.



The main advantages of this presentation relative to separate electron density and electrostatic potential surfaces are its clarity and its compactness. The main disadvantage is that it provides information only about the contact surface and does not reveal how far electron-rich and electron-poor areas extend beyond the surface.

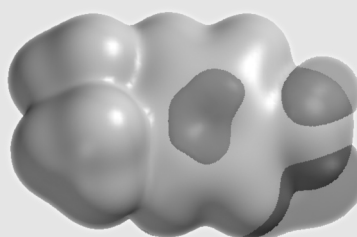An alternative to an electrostatic potential map, referred to as an **exposed electrostatic potential surface**, is a composite of three different surfaces: an electron density surface depicting overall molecular size and shape, a negative electrostatic potential surface identifying electron-rich regions and a positive electrostatic potential surface identifying electron-poor regions. These surfaces need to be generated and then displayed simultaneously. The electron density may either be displayed as an opaque solid or as a transparent solid (in order that the molecular skeleton may be seen inside). The two potential surfaces are best represented as transparent solids, the negative surface colored red and the positive surface colored blue. The exposed electrostatic potential surface for *p*-*tert*-butyl phenol is shown on the next page.




Note that the exposed electrostatic potential surface provides the same


information as the electrostatic potential map. Red areas in the map correspond to regions when the negative electrostatic potential surface is likely to protrude from the electron density while blue areas correspond to regions where the positive electrostatic potential surface is likely to stick out. There is broad similarity of this graphic to CFD representations discussed previously (see **CFD's** under the **Model** menu in **Chapter 14**). Hydrogen-bond acceptors and/or positive ionizable sites are associated with protruding negative potential surfaces, whereas hydrogen-bond donors and/or negative ionizable sites are associated with protruding positive potential surfaces.

The *hydride potential map* paints the value of the hydride potential onto an electron density surface. By convention, colors toward red depict negative potential, while colors toward blue depict positive potential, and colors in between (orange, yellow, green) depict intermediate values of the potential. Thus, a hydride potential map for

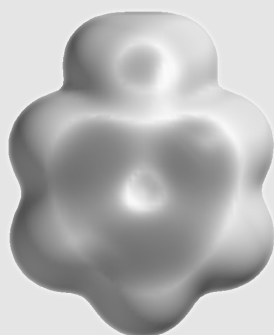The **local ionization potential map** paints the value of the local ionization potential onto an electron density surface. By convention, colors toward red indicate low ionization potential, while colors toward blue indicate high ionization potential. Thus, the local ionization potential map for aniline shows that the *ortho* and *para* ring positions have a lower ionization potential than the *meta* positions, consistent with the known directing ability of an amino group in electrophilic aromatic substitution.




The **LUMO map** paints the absolute value of the lowest-unoccupied molecular orbital (the LUMO) onto an electron density surface. By convention, colors near blue indicate high concentration of the LUMO, while colors near red indicate low concentration. Given that the LUMO designates space available for a pair of electrons, a LUMO map indicates where nucleophilic attack would likely occur. For example, a LUMO map for cyclohexenone shows concentration in two regions, one over the carbonyl carbon and the other over the  $\beta$  carbon, consistent with both carbonyl addition and Michael (conjugate) addition.

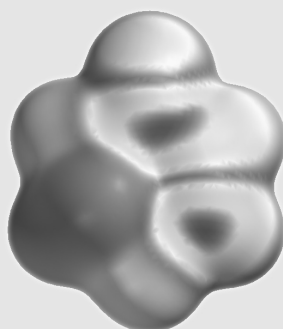

The **spin density map** paints the value of the spin density onto an electron density surface. By convention, colors near blue indicate high concentration



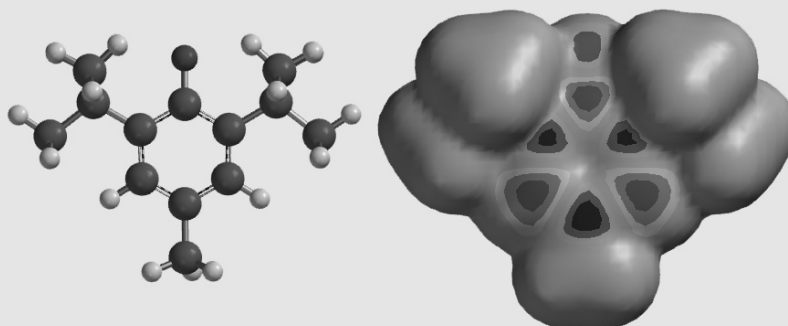of spin density, while colors near red indicate low concentration. For example, a spin density map for the radical resulting from loss of hydrogen from 3,5-di-*tert*-butylhydroxytoluene (BHT) shows that the spin has delocalized from oxygen onto the *ortho* and *para* ring positions.



This radical would be expected to be particularly stable, which explains why BHT acts as an antioxidant (scavenging less favorable localized radicals).

All of these maps use an electron density surface and to delineate overall molecular size and shape. Not all regions on this surface are accessible and therefore available for interaction with their environment (or with an incoming reagent). *Spartan* allows these regions to be identified.\*

Surfaces (including those underlying maps) connect points of equal value (they are isosurfaces), and may be displayed as an arrangement of dots, a mesh, or an opaque or translucent solid. Examples of graphical output in orthogonal projection are provided in **Figure 17-1**. Surfaces (and maps) may also be rendered in perspective (see **Chapter 14**) and in stereo (see **Chapter 2**).







Calculated quantities may also be displayed as two dimensional contour plots (slices). Unlike surfaces and maps, these can be translated, rotated and zoomed independently of the molecular skeleton. An example of a slice display is provided in **Figure 17-1**.

Several different surfaces, maps and slices may be simultaneously

---

\* A region on a density surface is designated as inaccessible if a sphere of radius 1.0 Å centered on a line normal to the surface and touching a point in the middle of the region, impinges on any other regions of the density surface. The sphere radius may be changed in the **Settings Preferences** dialog (**Preferences** under the **Options** menu; **Chapter 20**).

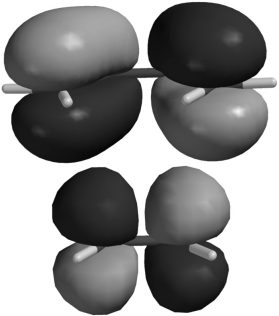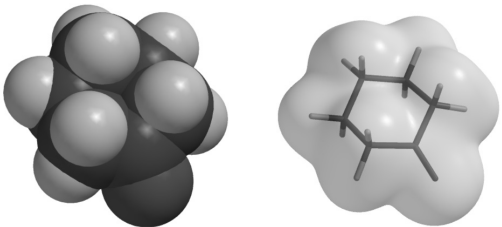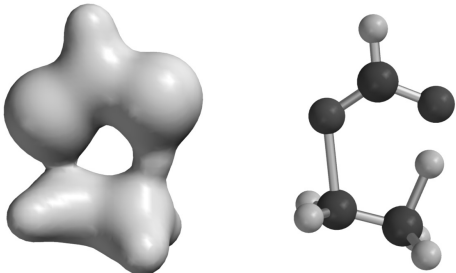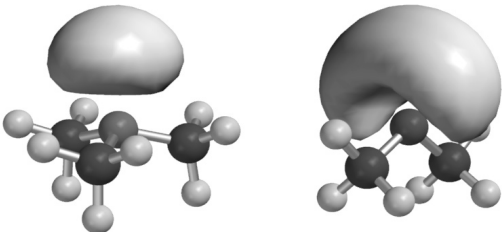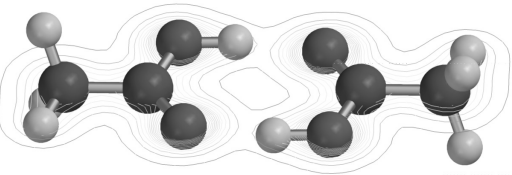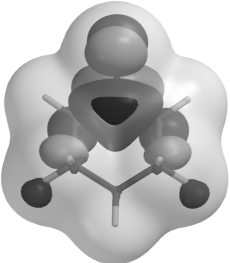
**Figure 17-1: Examples of Graphical Displays Available in *Spartan***

<p><b>Frontier orbitals for a symmetry-allowed Diels-Alder reaction,</b></p>  <p>showing interaction of the HOMO of 1,3-butadiene and the LUMO of ethylene.</p>	<p><b>Space-filling model and electron density surface of cyclohexanone,</b></p>  <p>showing overall molecular size and shape.</p>
<p><b>Electron density surface (0.08 electrons/au<sup>3</sup>) of transition structure for pyrolysis of ethyl formate,</b></p>  <p>showing bonding in the transition state.</p>	<p><b>Electrostatic potential surfaces (-40 kJ/mol) of trimethylamine (left) and dimethyl ether (right),</b></p>  <p>showing the lone pairs on nitrogen and oxygen, respectively.</p>
<p><b>Electron density slice for acetic acid dimer,</b></p>  <p>showing hydrogen bonding.</p>	<p><b>Simultaneous display of the LUMO and the electron density surfaces of cyclohexanone,</b></p>  <p>showing accessibility for nucleophilic attack.</p>

displayed. In addition, any of the usual structure models may be displayed along with the graphic. The total display can become very complex, and selective use of meshes and/or translucent solids (as opposed to opaque solids) may facilitate visualization.

Discussion of the utility of graphical models for describing molecular structure and chemical reactivity is provided in *Topics* available under the **Activities** menu and in *A Guide to Molecular Mechanics and Quantum Chemical Calculations*, provided as a PDF under the **Help** menu.

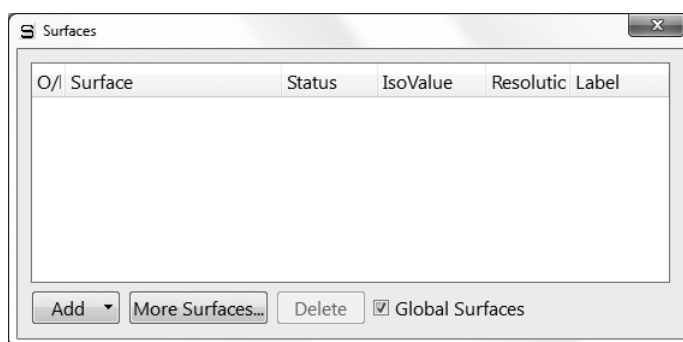
Selection of **Surfaces** from the **Setup** menu results in display of the **Surfaces** dialog.



This contains a box at the top for listing requested surfaces and property maps.

### Common Surfaces and Property Maps

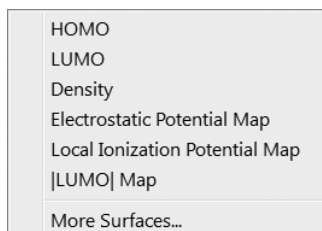**Add** at the bottom of the dialog is used to specify a number of commonly-used graphical surfaces and property maps\*. *Clicking* on it leads to a menu.




Selection of all but the last entry in the menu leads to a request for

---


\* The number of selections will increase if the molecule has unpaired electrons.

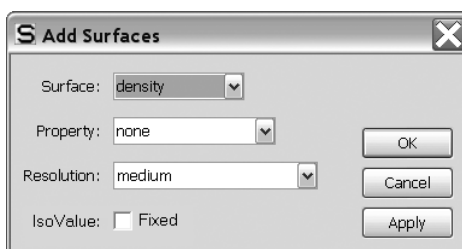the analogous surface or map. Calculated at medium resolution. Additional surfaces and maps or the same surfaces or maps at different resolution may be requested by selecting **More Surfaces...** from the menu (or by *clicking* on **More Surfaces...** at the bottom of the **Surfaces** dialog). This leads to the **Add Surfaces** dialog that contains three menus and a check box:



## Surface

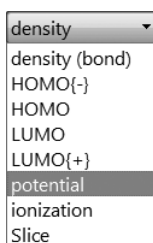Available surface types appear under the **Surface** menu.



**Density** portrays overall molecular size\*, **density (bond)** locates bonds, **HOMO{-}**, **HOMO**, **LUMO**, **LUMO{+}**, **SOMO\*\*** are molecular orbitals, **potential** is the electrostatic potential, **hydride** is the hydride potential, **ionization** is the local ionization potential and **spin\*\*** is the spin density.

Unlike the molecular orbitals, the electron density and the electrostatic potential, the local ionization potential does not go to zero with increasing distance. This makes its use as a surface problematic.

Selection of **HOMO{-}** and **LUMO{+}** results in display of a box to decrement the HOMO and increment the LUMO. This allows any molecular orbital to be specified.

\* In the case of molecular mechanics, the density will be based on van der Waals radii.


\*\* These menu entries appear only for molecules with unpaired electrons.



**Slice** designates that a plane will cut through the graphic defined by **Property** (see discussion following).

## Property

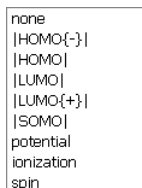Properties for maps appear in the **Property** menu.



Available properties are the molecular orbitals (**HOMO{-}**, **HOMO**, **LUMO**, **LUMO{+}**, **SOMO\***), the electrostatic potential (**potential**\*\*), the local ionization potential (**ionization**) and the spin density (**spin**\*). **none** indicates that no property is to be mapped onto the surface). As with **Surface** above, selection of **HOMO{-}** and **LUMO{+}** leads to a decrement (increment) box.

A **Spin** button will be displayed if **Multiplicity** (in the **Calculations** dialog) is set to a value other than **Singlet**, or if **Scf=unrestricted** has been specified in the **Options** box (**Calculations** dialog; see also **Appendix C**) and if **HOMO{-}**, **HOMO**, **LUMO** or **LUMO{+}** has been selected for **Surface** or for **Property**. *Clicking* on **Spin** toggles it between **Alpha** and **Beta**. **Alpha** designates that the molecular orbital either to be displayed as a surface or mapped as a property corresponds to  $\alpha$  spin; **Beta** designates that the molecular orbital corresponds to  $\beta$  spin.

---

\* These menu entries appear only for molecules with unpaired electrons.

\*\* In the case of molecular mechanics, the electrostatic potential will be based on charges.

## Resolution

Selection of surface resolution is from the **Resolution** menu.

low (8x Faster)
medium
intermediate (4x Slower)
high (8x Slower)

Medium resolution generally is intended for routine work, low resolution is used to get rough images very quickly, while intermediate and high resolution may be employed to obtain publication quality graphics. Both calculation time and disk storage increase significantly with increasing resolution.

## IsoValue

This allows a fixed isosurface to be generated in place of a three-dimensional grid of points (from which different isosurfaces, and/or slices can be constructed). Fixed isosurfaces are particularly attractive for property maps (that are generally based on a density surface representing molecular size and shape) as they require less computer time and disk storage. Note that a fixed density surface needs to be based on a specific isovalue (see below) and not on percentage enclosure of the total number of electrons.

*Checking* the box to the left of **Fixed** leads to display of the isovalue (0.002 electrons/au<sup>3</sup> in the case of density). This value may be changed prior to surface generation by editing the text box that appears to the right. **Fixed** may not be specified in conjunction with **Slice** as a **Surface**.

Following **Surface**, **Property**, **Resolution** and (optionally) spin selection, *clicking* on **OK** adds the requested graphic to the list and removes the (**Add Surfaces**) dialog. *Clicking* on **Apply** adds the requested graphic to the list but leaves the dialog on screen. *Clicking* on **Cancel** does not add a graphic to the list but removes the (**Add Surfaces**) dialog.

The process (*clicking* on **Add...**, surface and property selection, *clicking* on **OK** or **Apply**) may be repeated as required.

An existing graphic or a request for a graphic may be deleted from

the list by first highlighting (*clicking* on) it and then *clicking* on **Delete**.

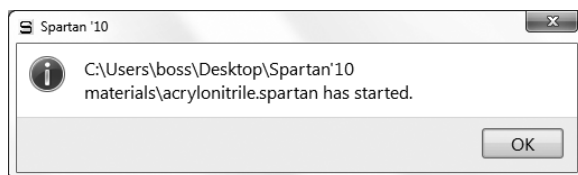
## Global Surfaces

If *checked*, **Global Surfaces** signifies that the requested surfaces will be calculated for all members of the list.

Only one copy of the **Surfaces** dialog may appear on screen, and any actions relate to the currently selected molecule. The dialog may be removed by *clicking* on .

## Submit

Following setup of a molecular mechanics or quantum chemical calculation, or a similarity analysis, as well as any requests for properties, spectra and/or graphical displays, the required calculations will begin when **Submit** is selected. If the job has not previously been named, selection of **Submit** triggers a request for a name. If the document being submitted comprises a single molecule and the molecule exists in the Spartan Molecular Database, the name in the database will be presented as a default name. Otherwise, the default name presented will be *spartan* for the first job and *spartanX* (X an integer starting with 1 for all successive jobs). After a name has been provided (or if a name already exists or if the default name is accepted) a dialog appears indicating that the job has actually been submitted.\*




*Click* on **OK** to remove it. After a job has started, and until it has completed, all associated files will be designated read only.

---

\* The job is submitted to a job queue and will be submitted for execution only when released from this queue. See **Monitor** under the **Options** menu (**Chapter 20**) for discussion.

In addition to the local (Windows, Macintosh or Linux) machine, *Spartan*'s computational codes may run on a remote (licensed) Linux server. The only exception is the similarity analysis code which must be run on the local machine. If remote submission has been set up (see **Hosts** from the **Preferences** dialog under the **Options** menu; **Chapter 20**), selection of **Submit** will also request a name for the destination machine.

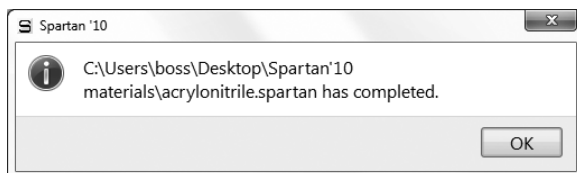Another dialog appears following completion of a calculation.



Click on **OK** to remove it.

Some types of jobs, including requests for a conformer distribution and an energy profile, lead to an additional document being created for each molecule in the original document. These new documents are named *document.task.identifier* where *document* is the name given to the original document, for example, *Conf* for a conformer distribution and *Prof* for energy profile, and *identifier* identifies the molecule inside the original document.

In the case where there is only one molecule in the original document, a query dialog is provided asking whether the resulting *Spartan* document is to be opened.






# Chapter 18

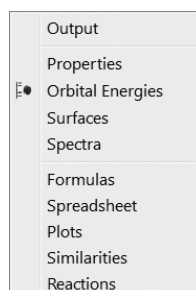
## The Display Menu

Functions available under the **Display** menu provide for text, dialog, spreadsheet and graphical displays following structure/property calculation, conformational analysis or similarity analysis. Functions are also available to query a variety of on-screen objects, display and compare IR, Raman, NMR and UV/visible spectra, animate vibrational motions, construct relationships among calculated properties or for a property between reactants and products of a chemical reaction, perform linear regression and prepare plots from spreadsheet data, display similarities between molecular structures and between molecular structures and pharmacophores and calculate reaction energies.



### Output

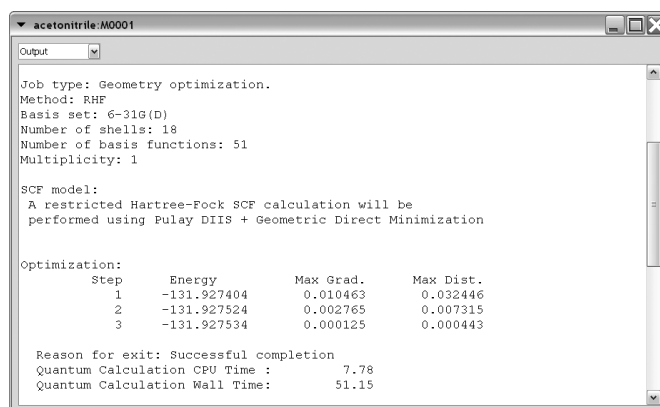Selection of **Output** opens a window.





acetonitrile-M0001

Output

```
Job type: Geometry optimization.  
Method: RHF  
Basis set: 6-31G(D)  
Number of shells: 18  
Number of basis functions: 51  
Multiplicity: 1  
  
SCF model:  
A restricted Hartree-Fock SCF calculation will be  
performed using Pulay DIIS + Geometric Direct Minimization  
  
Optimization:  
Step      Energy           Max Grad.        Max Dist.  
  1      -131.927404       0.010463        0.032446  
  2      -131.927524       0.002765        0.007315  
  3      -131.927534       0.000125        0.000443  
  
Reason for exit: Successful completion  
Quantum Calculation CPU Time :      7.78  
Quantum Calculation Wall Time:    51.15
```

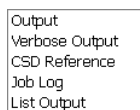A menu at the top left of the window selects the type of output.



**Output** provides normal output. **Verbose Output** contains more detailed output, but is eliminated upon normal completion unless **Keep Verbose** is checked in the **Setting Preferences** dialog (**Preferences** under the **Options** menu; **Chapter 20**). **CSD Reference** provides the literature reference for data retrieved from CSD. **Job Log** and **List Output** contain diagnostic information. The contents of the output window may be printed by selecting **Print Output...** from the **File** menu (this replaces **Print...** when an output window is selected), or copied to the clipboard by selecting **Copy** from the **Edit** menu (**Chapter 13**) when an output is selected. **Find...** and **Find Next** functions from the **Edit** menu are also available. An output window is associated with each document and pertains to the selected molecule. Several output windows may be open and may be closed by *clicking* on  at the top right.



Output is not accessible for jobs that are executing or are in the execution queue. Output for jobs that execute locally may be viewed using the **Monitor** under the **Options** menu (**Chapter 20**), and the current structure (in a geometry or transition-state optimization) can be examined.

## Properties

*Spartan* provides specialized dialogs for reporting and (in some cases) changing the properties of molecules, atoms, bonds, surfaces and constraints, and for changing default plot styles and limits, the fitting functions, and for editing CFD's and substituent libraries. Only one **Properties** dialog may be open, and this refers either to the selected molecule (**Molecule Properties**), or to the selected component (atom, bond, etc.) or attribute (spectra, graphical surface, substituent, etc.) of the selected molecule (**Atom Properties**, **Bond Properties**, **Surface Properties**, etc.), or to a plot (**Plot Properties**, **Curve Properties**, etc.), or fitting function (**Regression Properties**). Selection of a

different molecule leads to the **Molecule Properties** dialog for that molecule. Dialogs that refer to components/attributes of the (newly selected) molecule follow by *clicking* on the component/attribute.

Dialog selection operates in toggle mode, with the **Molecule Properties** dialog being in the neutral position. With the **Molecule Properties** dialog on screen, *clicking* on a component/attribute brings up the appropriate **Properties** dialog. For example, *clicking* on an atom brings up the **Atom Properties** dialog. *Clicking* on a different component/attribute brings up the appropriate **Properties** dialog. *Clicking* a second time on the same component reverts back to the **Molecule Properties** dialog.\*

Most **Properties** dialogs have an associated **Utilities** or **Style** dialog. For example, associated with the **Molecule Properties** dialog is a **Molecule Utilities** dialog. These access additional information about the molecule and its components/attributes, or provide style and color controls. This is useful for highlighting (or dehighlighting) a particular molecule, component or attribute. **Utilities/Style** dialogs are reached by *clicking* on  at the bottom right of the appropriate **Properties** dialog. Return to the **Properties** dialog follows from *clicking* on  at the bottom right of the associated **Utilities/Style** dialog.

The **Properties** (or **Utilities/Style**) dialog may be removed from the screen by *clicking* on .

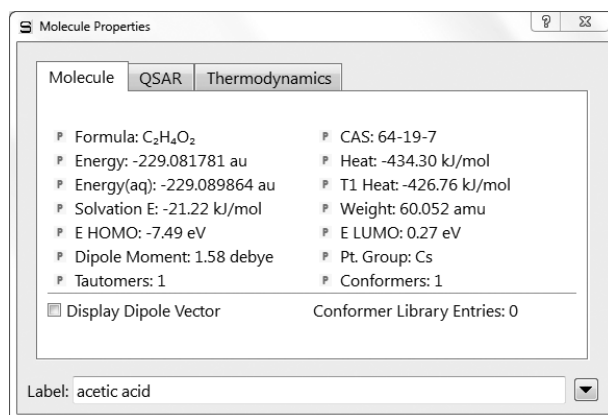
### **Molecule Properties, QSAR Descriptors and Thermodynamics**

The **Molecular Properties** dialog is divided into three displays **Molecule**, **QSAR** and **Thermodynamics**, controlled by tabs at the top left. Entries under the **Molecule** tab relate to common molecular properties that (if applicable to the selected level of calculation) are automatically calculated.

---

\* The exception to this involves *clicking* on a graphical surface or property map, for example, *clicking* on a property map to obtain the value of the property at a particular surface location. *Clicking* a second time on the surface or map will report a new value of the property. *Clicking* on the background leads to the **Molecule Properties** dialog.

## Molecule




These include the energy (in au), the heat of formation (in kJ/mol), or the sum of the strain energy (in kJ/mol) as appropriate, the energy (heat of formation, strain energy) corrected for an aqueous environment according to the SM5.4 model (in the same units as the non-corrected quantity), the HOMO and LUMO energies (in eV), dipole moment (in debyes), molecular weight (in amu), the point group, the predicted number of tautomers and the predicted number of conformers, and (if available) the experimental heat of formation (in kJ/mol), the heat of formation from the T1 thermochemical recipe (in kJ/mol) and the CAS number. These may be posted to the spreadsheet using the **P** buttons, or *dragged* into the spreadsheet, into the formula editor (see **Formulas** in this chapter), or onto any of the dialogs associated with the Spartan Molecular Database or Spartan Spectra and Properties Database (see **Open Database** under the **Search** menu; **Chapter 19**). **Label** identifies the molecule in a document and is the name given to the molecule in the Spartan Molecular Database or the Spartan Spectra and Properties Database.\* It appears in the first column of the spreadsheet (see **Spreadsheet** later in this chapter). The dialog also displays the number of entries in a conformer library (see discussion under **Calculations** from the **Setup** menu; **Chapter 17**). Finally, the dialog permits display of the dipole moment vector by *checking* the box to the right of **Display Dipole Vector**.

\* This assumes that the molecule exists in one of the databases. If it does not, the name will be **M000X** (X=1,2,...) unless altered by the user.

## QSAR

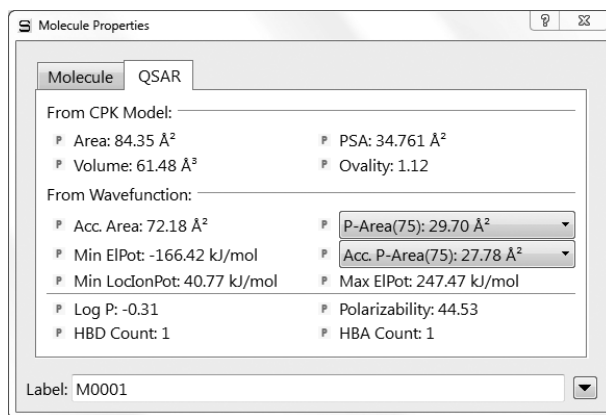Entries under the **QSAR** tab provide additional properties, some of which may be particularly valuable in QSAR type analyses.




These include: the area, volume, polar surface area (PSA)\* and ovality obtained from a space-filling model, and other structure-dependent indicators: LogP, polarizability and the number of hydrogen-bond donor (HBD) and acceptor sites (HBA).\*\* All of these are provided automatically. Additional quantities based on the electron density surfaces as well as on electrostatic potential maps are also available: the accessible area, the polar area corresponding to absolute values of the electrostatic potential greater than 75, 100 and 125 kJ/mol, the accessible polar area corresponding to each of these, the minimum and maximum values of the electrostatic potential (as mapped onto an electron density surface) and the minimum value of the local ionization potential (as mapped onto an electron density surface). These quantities need to be requested by checking **QSAR** inside the **Calculations** dialog (**Calculations** under the **Setup** menu; **Chapter 17**).



\* By default, polar surface area is defined as the area due to nitrogen and oxygen and any hydrogens attached to nitrogen and oxygen. Polar surface areas corresponding to arbitrary alternative definitions are available for posting into the spreadsheet using the **PAREA** function. See **Table 18-3**.


\*\* Counts of hydrogen-bond donors and acceptors and positive and negative ionizable sites derive from CFD assignments. These in turn follow from atomic connectivity together with special case designations for common organic functional groups. CFD's assignments may be modified on a case-by-case basis using the **CFD Properties** dialog (**Properties** under the **Display** menu; this chapter).

## Thermodynamics

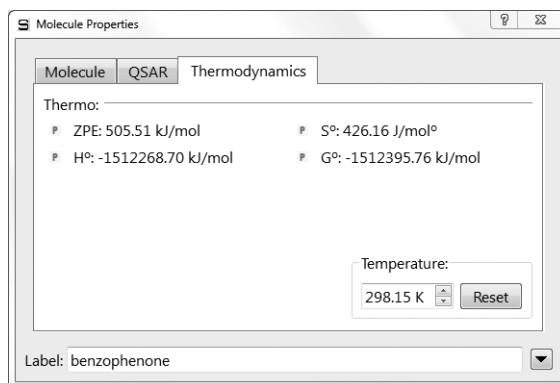Entries under the **Thermodynamics** tab provide the zero-point energy, the enthalpy of temperature T, the entropy and the Gibbs energy at temperature T. All require vibrational frequencies.



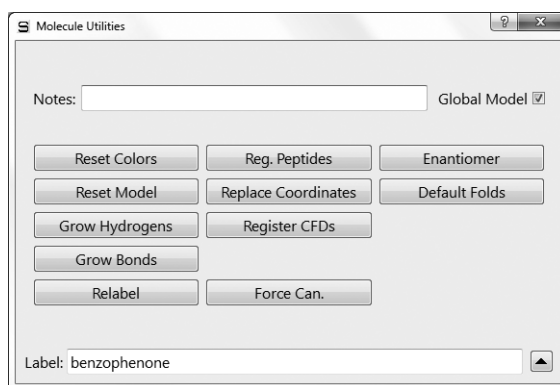
Clicking on  at the bottom right of the **Molecule Properties** dialog brings up the **Molecule Utilities** dialog (clicking on  returns to the **Molecule Properties** dialog).





**Notes** is reproduced in the output. Controls reset model color and style, add missing hydrogens and bonds, provide information about amino acids in polypeptides and about CFD's, replace coordinates by standards based only on atomic connectivities, change enantiomers and reset default conformer selections.


## Thermodynamics

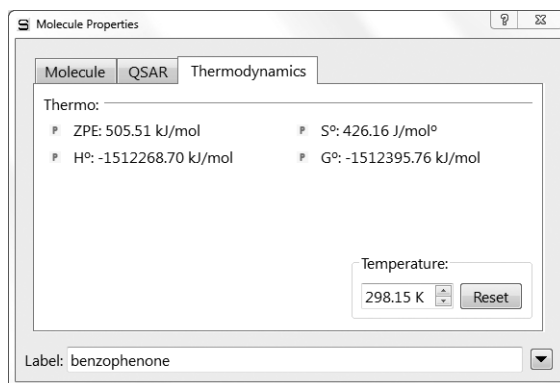Entries under the **Thermodynamics** tab provide the zero-point energy, the enthalpy of temperature T, the entropy and the Gibbs energy at temperature T. All require vibrational frequencies.



## Molecule Utilities

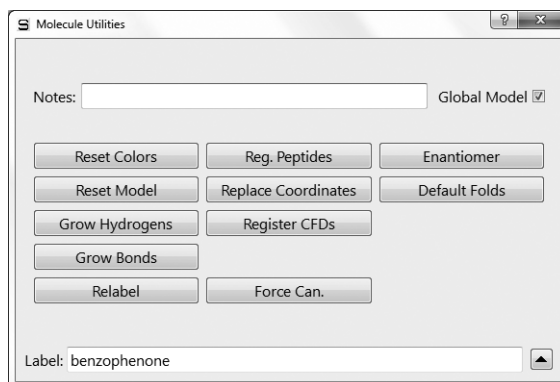Clicking on  at the bottom right of the **Molecule Properties** dialog brings up the **Molecule Utilities** dialog (clicking on  returns to the **Molecule Properties** dialog).




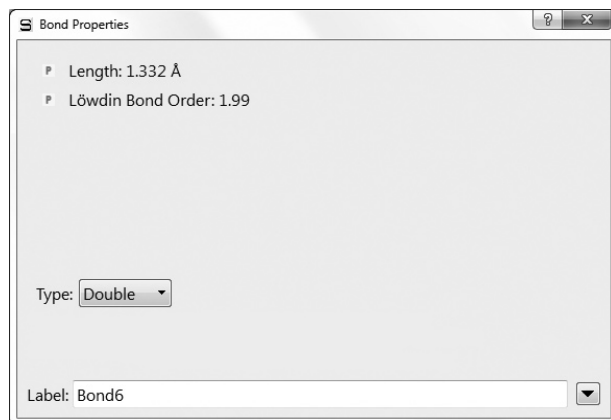**Notes** is reproduced in the output. Controls reset model color and style, add missing hydrogens and bonds, provide information about amino acids in polypeptides and about CFD's, replace coordinates by standards based only on atomic connectivities, change enantiomers and reset default conformer selections.


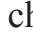


## Bond Properties

Selection of a bond with a **Properties** dialog on screen, or selection of **Properties** following selection of a bond leads to the **Bond Properties** dialog.

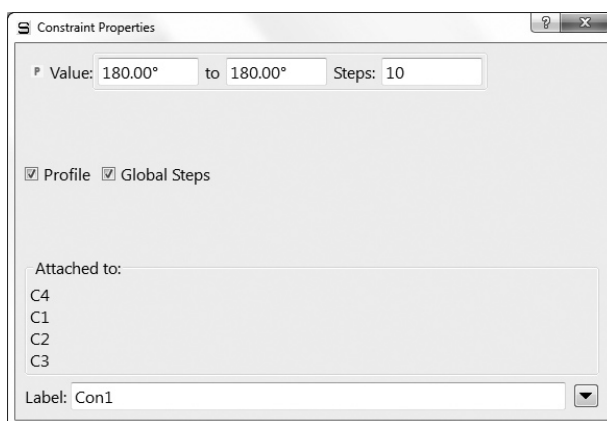


This displays the bond length, bond type (and allows changing the bond type) and Löwdin bond order (the bond order used in NMR shift corrections).



*Clicking* on  at the bottom right of the **Bond Properties** dialog brings up the **Bond Style** dialog (not shown). This contains controls for changing bond color. *Clicking* on  at the bottom right of this dialog returns to the **Bond Properties** dialog.

## Constraint Properties



Selection of a constraint marker with a **Properties** dialog on screen, or selection of **Properties** following selection of a constraint marker, leads to the **Constraint Properties** dialog (not shown). This allows setting the value of a constraint, posting it to the spreadsheet and changing the default constraint label. *Checking Profile* leads to an expanded dialog.



This allows specifying a sequence of constraints for an energy profile (see **Calculations...** under the **Setup** menu; **Chapter 15**). The value of the starting constraint is given in the lefthand box to the right of **Value**, and the value of the ending constraint is given in the righthand box to the right of **Value**. The number of steps in the profile is given in the box to the right of **Steps**. Initially, the numbers in both boxes to the right of **Value** will be the same, and **Steps** will be set to **1**. These may be altered by typing the desired numbers into the appropriate boxes and then *pressing* the **Enter** key (**return** key on Mac). This functionality may also be accessed from **Constrain Distance (Angle, Dihedral)** under the **Geometry** menu (**Chapter 15**). An energy profile may involve constraints on more than one geometrical variable (changed in lock-step with one another). In this case, the constraint ranges for the individual variables need to be selected. If **Global Steps** is *checked*, the independent variables are moved in concert, meaning that the number of steps must and will be the same for each variable. If *unchecked*, the variables are moved independently and the number of steps may differ from one to another. In this case, the total number of steps in the profile is the product of the number of steps for each variable.


*Clicking* on  at the bottom right of the **Constraint Properties** dialog brings up the **Constraint Style** dialog (not shown). This contains controls for selecting the color of the constraint marker. *Clicking* on  returns to the **Constraint Properties** dialog.

## Point and Plane Properties

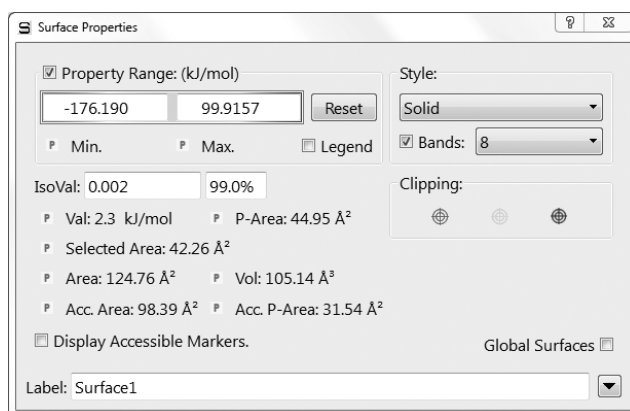Selection of a user-defined point or plane with a **Properties** dialog on screen, or selection of **Properties** from the **Display** menu following selection of a point or plane, leads to the **Point Properties** or **Plane Properties** dialog (not shown). These allow changing point or plane labels. *Clicking* on  at the bottom right of the **Point (Plane) Properties** dialog leads to the **Point (Plane) Style** dialog (not shown), with controls for adjusting point (plane) colors. *Clicking* on  returns to the **Point (Plane) Properties** dialog.

## Surface Properties

Selection of a graphical display with a **Properties** dialog on screen, or selection of **Properties** following selection of a graphical display, leads to the **Surface Properties** dialog.



This allows changing display style, isovalue (and for electron density surfaces), percentage of the electrons contained inside the surface, turning on mapped properties, selecting between continuous and banded displays, setting the property range of maps, displaying accessible area of surfaces and changing default surface labels. The dialog reports (and optionally posts) the area, volume and accessible area\* of the surface, the polar area of an

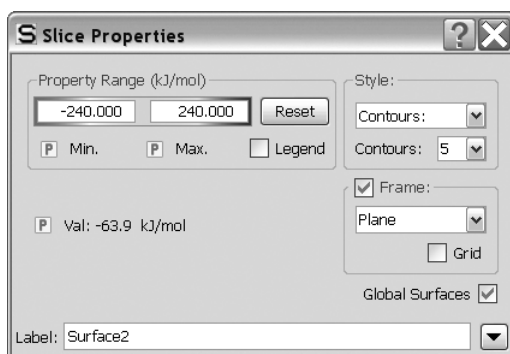
---

\* A region on a density surface is designated as inaccessible if a sphere of radius 1.0 Å centered on a line normal to the surface and touching a point in the middle of the region, impinges on any other regions of the surface. The default radius (**Accessible Area Radius**) may be changed in the **Settings Preferences** dialog (**Preferences** under the **Options** menu; **Chapter 20**).



electrostatic potential map<sup>\*</sup>, maximum and minimum value of the mapped property and its value at the cursor position<sup>\*\*</sup> and the area enclosed in a selected band. If *checked*, **Legend** displays a color scale. If *checked*, **Global Surfaces** designates that the settings apply to all molecules in the document.

The dialog also provides for up to three independent clipping planes which operate in toggle mode. Following selection, a plane may be rotated, translated or zoomed with the usual mouse operations.

Where the graphic is a slice, the **Slice Properties** dialog appears.



This contains controls to account for the number of contours displayed and to the type of slice (plane, cylinder or sphere). Check boxes allow for a frame around the slice and for a grid.

*Clicking* on  at the bottom right of the **Surface Properties** dialog brings up the **Surface Style** dialog (not shown). This contains controls for selecting surface color. *Clicking* on  at the bottom right of this dialog returns to the **Surface Properties** dialog.

## Curve and Plot Properties

Plots created either using data from the spreadsheet or from spectra may be edited using a set of closely-related dialogs. These refer

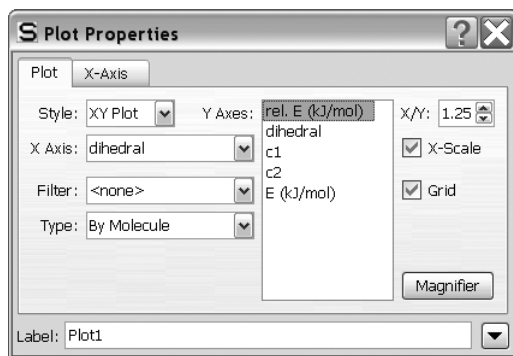
---

\* This is defined as that part of the surface area for which the absolute value of the electrostatic potential is  $> 100$  kJ/mol. The cutoff (**Polar Area Range**) may be changed in the **Settings Preferences** dialog (**Preferences** under the **Options** menu; **Chapter 20**).

\*\* To determine property value at another position *click* on it. To bring up the **Molecule Properties** dialog, *click* on the background.


to curves (a relationship between variables), and plots (a set of relationships that share variables).

Selection of a plot frame with a **Properties** dialog on screen, or selection of **Properties** following selection of a plot frame, leads to the **Plot Properties** dialog (XY plot case shown).





The **X Axis** menu and the **Y Axes** list duplicate functionality under the **Plots** dialog (see **Plots...** later in this chapter) and allow specification of a molecular property to serve as the X axis of the plot and one or more properties to serve as the Y axes. The contents of the **Y-Axes** list depends on selection from the **Type** menu: **By Molecule** and **Free Form** reference explicit entries in the spreadsheet (see discussion later in this chapter) while **By Atom** references *Spartan's* file of atom-based properties (the *Property Archive*). Initial selections may be altered as desired. In addition, the aspect ratio of the plot may be changed, the numerical scale of the X (horizontal) axis turned on or off, a grid added to the plot, and the default label changed. **Magnifier** brings up a *magnifying glass* that may be moved independently. The magnifier may be removed from the plot using **Delete** (**Build** menu) or the **Delete** key.

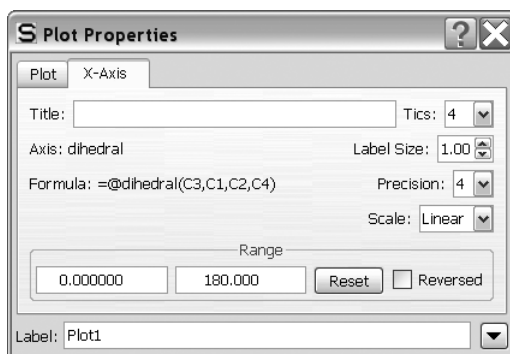
The scale along the X (horizontal) axis may be altered from default settings by *clicking* on the **X-Axis.** tab at the left of the dialog. This brings up the **Plot X-Axis Properties** dialog (*clicking* on the **Plot** tab returns to the **Plot Properties** dialog).



To alter the scale, edit the values inside the **Range** boxes, and *press* the **Enter** key (**return** key on Mac). Controls are also available to alter the number of tic marks on the X axis, the label size and precision and change the scale from linear to logarithmic.

*Clicking* on  at the bottom right of either the **Plot Properties** or **Plot X-Axis Properties** dialog brings up the **Plot Style** dialog (not shown). This contains controls for selecting the color of the plot frame, as well as the label associated with the X axis. *Clicking* on  at the bottom right of this dialog returns to either the **Plot Properties** or **Plot X-Axis Properties** dialog.

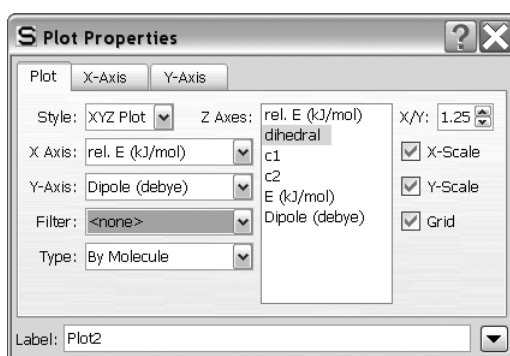
The **Style** menu in the **Plot Properties** dialog allows shifting between 2D (XY) and 3D (XYZ) plots. Selection of **XYZ Plot** leads to an extended **XYZ Plot Properties** dialog.



This is identical to the **XY Plot Properties** dialog, except that it includes an additional menu (**Y Axis**) from which one property needs to be selected. In addition, the **Y Axes** list in the **XY Plot Properties** dialog has been renamed to **Z Axes**.

Scales along X and Y axes may be adjusted by *clicking* on **X-Axis** and **Y-Axis** tabs to the left of the dialog. This leads to the **Plot X-Axis (Plot Y-Axis) Properties** dialog (*clicking* on the **Plot** tab at the left of this dialog returns to the **Plot Properties** dialog). These are analogous to the (**Plot X-Axis Properties**) dialog discussed for XY plots.

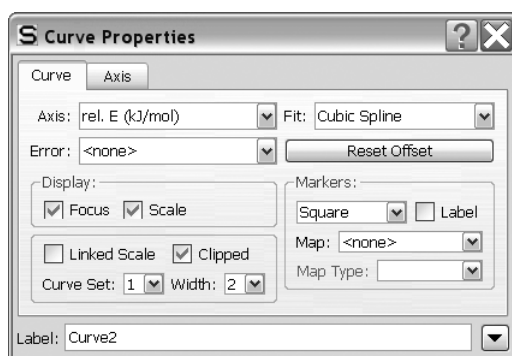Selection of an individual curve in a plot with a **Properties** dialog on screen, or selection of **Properties** from the **Display** menu following selection of an individual curve, leads to the **Curve Properties** dialog.




This allows changing the quantity that is plotted, providing a scale and error bars (down from the column in the spreadsheet designated in the **Error** menu) and turning clipping on. Controls are available to designate the kind of markers for the data points, request marker labels, and designate the kind of fit to the data points and the line width of the fitting curve\*. If *checked*, **Focus** marks the selected molecule on the curve.



The scale along the Y axis of an XY plot or Z axis of an XYZ plot may be altered from the default (minimum to maximum value of the Y or Z) by *clicking* on the **Axis** tab. This brings up the **Curve Axis Properties** dialog.

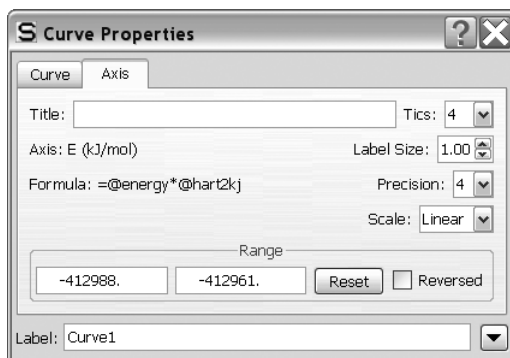---

\* Fitting functions include linear, quadratic and cubic least square, a Fourier series, a Lorentzian or a Gaussian. The data points may also be connected by straight lines (point-to-point) or by a smooth curve (cubic spline) or be presented as a skyline. Finally, the fitting curve can be removed.




To alter the scale, type in new values inside the boxes underneath **Range**. You need to *press* the **Enter** key (**return** key on Mac) following each entry. Controls are also available to alter characteristics of the Y (XY plot) or Z axis (XYZ plot) labels.

*Clicking* on  at the bottom right of either the **Curve Properties** or **Curve Axis Properties** dialog brings up the **Curve Style** dialog (not shown). This contains controls for selecting the color of the curve (and the label associated with the Y (or Z) axis). *Clicking* on  at the bottom right returns to the original dialog.

## Substituent Properties

Selection of a substituent icon with a **Properties** dialog on screen, or selection of **Properties** following selection of a substituent icon, leads to the **Substituent Properties** dialog.

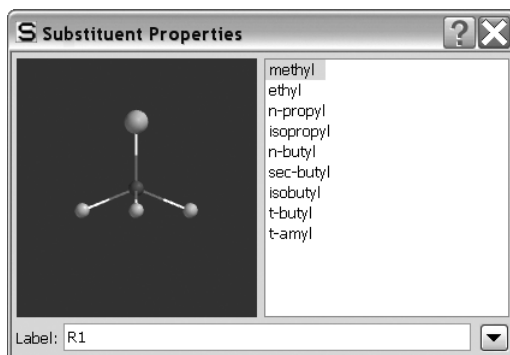



This allows substituent libraries already attached to a molecule to be edited, and the label on the substituent icon to be changed. It does not allow for attachment points to be reassigned. The dialog




comprises a scroll box (listing the substituent names), a viewing window (showing a model of the substituent) and a label field.

To delete a substituent from the library, *click* on its name in the scroll box, then *right click* and select **Delete** from the menu that appears. To add a new substituent that has been copied to the clipboard to the library, *right click* inside the scroll box and select **Paste** from the menu that appears. Note, that if the clipboard contains several molecules, all will be added to the substituent library. A substituent (or list of substituents) in a *Spartan* document can be added by *right clicking* inside the scroll box, selecting **Append** from the menu that appears and choosing a document from the browser, or *dragging* the document that contains the substituent(s) from a window into the scroll box.

*Clicking* on  at the bottom right of the **Substituent Properties** dialog brings up the **Substituent Style** dialog (not shown). This contains controls for selecting substituent (icon) color. *Clicking* on  returns to the **Substituent Properties** dialog.



## CFD Properties

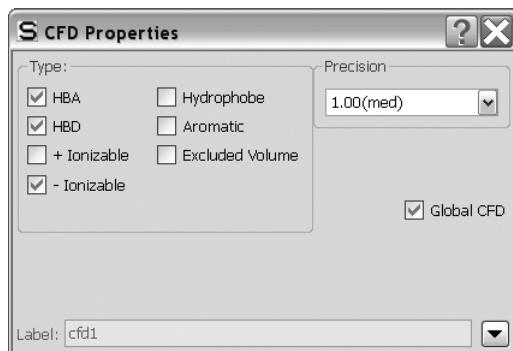Selection of a chemical function descriptor (CFD) or of an element of a pharmacophore with a **Properties** dialog on screen, or selection of **Properties (Display menu)** following selection of a CFD or pharmacophore element, leads to the **CFD Properties** dialog.




This allows changing the definition of the CFD, its label and the precision of the fit; the bigger the radius of the CFD, the less precise the fit). If *checked*, **Global CFD** specifies that CFD definitions

will apply to all molecules in the document.\*

Clicking on  in the dialog brings up the **CFD Style** dialog (not shown). This contains controls for selecting CFD color. Clicking on  in this dialog returns to the **CFD Properties** dialog.

## Regression Properties

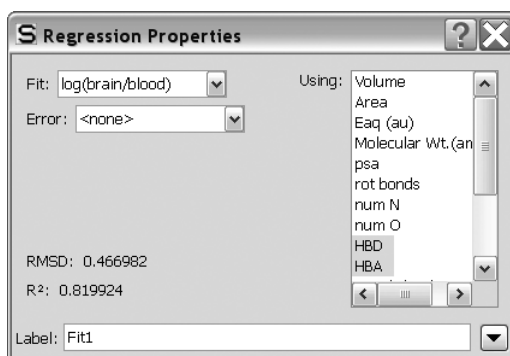Following a linear regression analysis, a new row, labeled **Fit1\*\***, appears near the bottom of the spreadsheet. This contains information about the fit. Clicking on this line with a **Properties** dialog on screen, or selecting **Properties (Display menu)** after clicking on the line, leads to the **Regression Properties** dialog.



This reports RMSD and  $R^2$ , as well as allows for changing what is to be fit (**Fit**) and what it is to be fit to (**Using**). The error statistics will immediately update.

## “Selected” Style

The “**Selected**” **Style** dialog is accessed either by defining a selection box (see **Mouse and Keyboard Operations in Chapter 2**) with a **Properties** dialog on screen, or by selecting **Properties** from the **Display** menu after defining a selection box. This contains controls for changing the color and model type of whatever is included in the selection box.

---

\* This is generally applicable only where the molecules are conformers, or the CFD labels have been adjusted to associate positions on different molecules. Label reassignment is accomplished using the **Atom Properties** dialog (see discussion earlier in this chapter).

\*\* More precisely, a row will be written for each fit, and labelled **Fit1**, **Fit2**, . . .

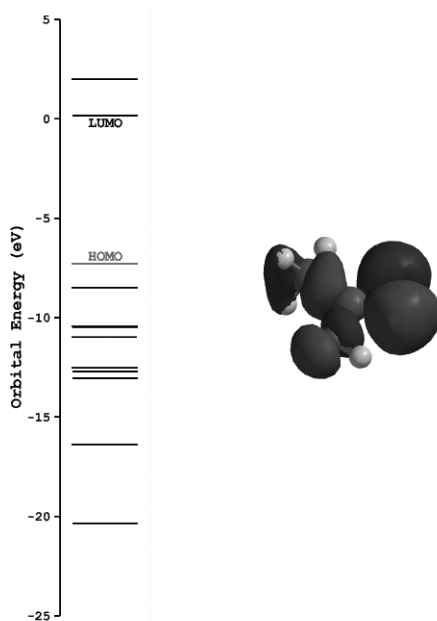
## Orbital Energies

Given that a wavefunction is available, *Spartan* allows display of an orbital energy diagram.

This comprises up to twelve occupied molecular orbitals and two unoccupied molecular orbitals, the highest-occupied (HOMO) and lowest-unoccupied (LUMO) being explicitly designated.

*Clicking* on an energy level in the diagram leads to display of the corresponding molecular orbital. This may be manipulated in the usual way and the display style altered (from the menu at the bottom right of the screen).

Only one orbital at a time may be displayed and this is not saved.



## Surfaces

This accesses the same dialog already described in **Chapter 17**.

## Spectra

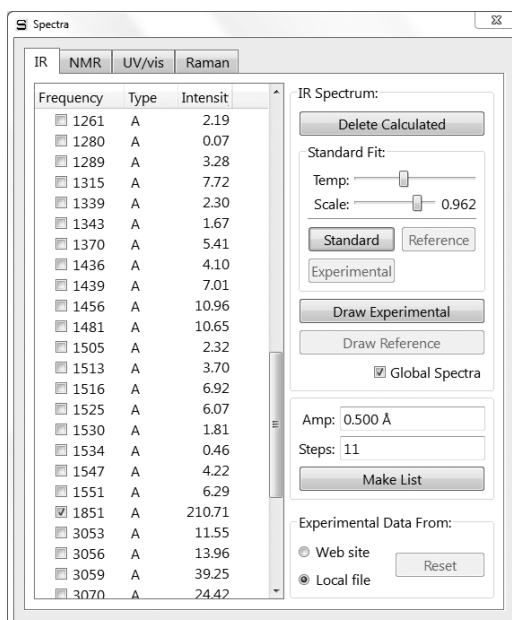
*Spartan* displays calculated IR, Raman, NMR and UV/visible spectra. (Spectra need to have been previously requested from the **Calculations** dialog under the **Setup** menu; **Chapter 17**.) In addition, it allows on-line access and display of experimental IR, NMR and

UV/visible spectra from publicly available collections, facilitating comparison with calculated spectra. Calculated IR and Raman spectra may be fit to measured spectra, using a linear scaling parameter and a peak width parameter.

Access to spectra display is via a series of tabbed dialogs (**IR**, **Raman**, **NMR** and **UV/vis**).

**IR, Raman** (these two dialogs are nearly identical and will be discussed together)

*Clicking* on the **IR (Raman)** tab leads to the **IR Spectra (Raman Spectra)** dialog. This displays the calculated and (if available) experimental infrared (Raman) spectra, and allows for animated display of the vibrational modes as well as for generation of a sequence of structures along a vibrational coordinate.



At the left-hand side of the dialog is an ascending list of frequencies (in  $\text{cm}^{-1}$ )\* together with (infrared) intensities and symmetry labels. A frequency is selected for display (animation) by *clicking* on it. *Clicking* again deselects it (stops the animation).

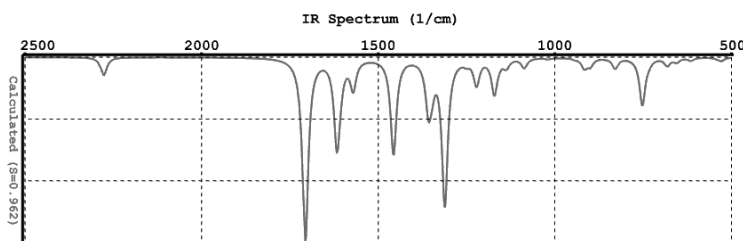
\* Imaginary frequencies, for example, corresponding to the reaction coordinate for a transition state, will appear first in the list, and will be designated by the letter **i** in front of the frequency value.

The amplitude of vibrational motion (the maximum displacement away from equilibrium of any pair of atoms) may be changed from the default amplitude (0.5Å which is much larger than appropriate for vibrational motion at room temperature) by altering the contents of the box to the right of **Amp** at the lower right of the dialog. The default number of steps that make up the display (11) may be changed by altering the contents of the box to the right of **Steps**\* at the lower right of the dialog. The greater the number of frames, the slower will be the vibration.

*Clicking* on **Make List** creates an unnamed *Spartan* document containing a sequence of structures corresponding to motion along the selected vibrational mode. Sequences (corresponding to different vibrations) can be created.

Note that the information required to produce vibrational motions (animations) and to construct lists of structures along different vibrational coordinates *is not available* as part of an entry in the Spartan Molecular Database. It *is available* as part of the entry in the Spartan Spectra and Properties Database (see discussion in **Chapter 19**). Only frequencies and intensities are provided in the Spartan Molecular Database, but vibrational motions are provided in the Spartan Spectra and Properties Database.


The calculated spectrum may be drawn by *clicking* on **Draw Calculated** at the top right of the dialog (IR spectrum shown).



The experimental spectrum (if available) may be drawn with or without the calculated spectrum by *clicking* on **Draw Experimental** at the center right of the dialog (IR spectrum shown).

---

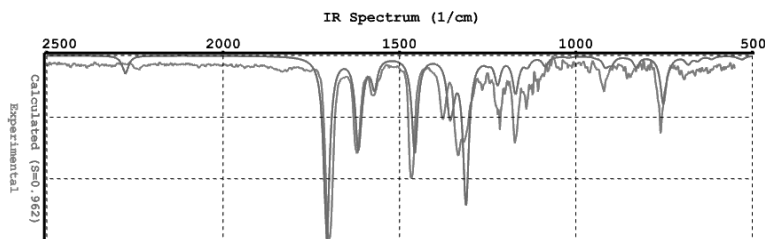\* It is recommended that the number of steps be odd ensuring that the center point corresponds to the actual equilibrium or transition-state structure.



Finally, a reference spectrum (a spectrum shared by all molecules in the document; see discussion of the Spartan Infrared Database in **Chapter 19**) if available, may be displayed along with the calculated and/or experimental spectra by *clicking* on **Draw Reference**.

Draw operations for both calculated and experimental IR and Raman spectra, apply to all molecules in the document (not just the selected molecule). Note, however, that the reference spectrum is common to all molecules in the document.

Experimental IR spectra are accessed from the NIST (National Institute of Standards and Technology) website (<http://webbook.nist.gov/chemistry>). This comprises approximately 7,000 IR spectra, primarily for organic molecules. *Spartan* keeps an index to facilitate retrieval but the actual data are on the website. Alternatively, the user may supply a spectra file (see **Appendix I**). Selection (website vs. file) is made at the bottom right of the dialog.

The default range for calculated, experimental and (for infrared only) reference IR and Raman spectra is from  $4000\text{ cm}^{-1}$  to  $500\text{ cm}^{-1}$ . This corresponds to the range commonly measured and reported. This range may be changed by selecting **Properties** from the **Display** menu, *clicking* on the horizontal axis to bring up the **Plot Properties** dialog, *clicking* on the **X Axis**. tab at the left of the dialog and editing the values under **Range**. The spectrum is treated like any other graphical object, and may be translated and scaled using the appropriate mouse operations. It may not be rotated.

The calculated spectrum may be removed by *clicking* on **Delete Calculated** (that has replaced **Draw Calculated**) in the dialog.

Similarly, the experimental spectrum (infrared only) may be removed by *clicking* on **Delete Experimental** (that has replaced **Draw Experimental**) and the reference spectrum may be removed by *clicking* on **Delete Reference** (that has replaced **Draw Reference**). The last spectrum (calculated, experimental or reference) removed will also cause removal of the plot axes. These operations apply to all molecules in the document (not just the selected molecule). An individual (calculated or experimental) spectrum or entire plot may also be removed by selecting **Delete** from the **Build** menu and then *clicking* on a spectrum or on the plot axes, respectively. Holding down the **Delete** key while *clicking* on the spectrum or plot axes has the same effect.

Calculated IR and Raman frequencies exhibit systematic errors. Frequencies obtained from (limiting) Hartree-Fock models are typically 10-12% larger than measured frequencies, while frequencies obtained from density functional and MP2 models are typically 3-4% larger than measured frequencies. Systematic errors may be revealed (and quantified) by first simultaneously displaying calculated and experimental spectra and then scaling one of the spectra to provide a best visual fit. To scale the calculated spectrum, use the slider bar to the right of **Scale** near the top right of the dialog. This is limited in range (from 0.8 to 1.0). If it is necessary to scale outside this range, *click* on the spectra (not on an axis) to select and move the mouse up and down while simultaneously depressing the **Shift** and **Alt** keys in addition to the right button. The scale factor (multiplying the frequency) will be displayed on the vertical axis. Scaling applies uniformly to all molecules in the list (not just the selected molecule).


A scaling factor may be applied to frequencies, and to all properties that make use of frequencies, prior to reporting using the **FreqScale** keyword typed into the **Options** box in the **Calculations** dialog (**Calculations...** under the **Setup** menu, **Chapter 17**).

A second slider bar marked **Temp** controls peak width. This is loosely connected with the temperature at which the experimental measurement is carried out. Low “temperature” (slider to the

left) will produce sharp spikes, whereas high “temperatures” (slider to the right) will produce broad bands. The default setting is intended to give calculated spectra that broadly resemble experimental infrared spectra. Both **Scale** and **Temp** slider bars actually control the fit of the calculated frequencies and intensities to a Lorentzian function. The default settings can be restored by *clicking* on the **Standard** button. A best fit to an experimental spectrum superimposed onto the calculated spectrum can be obtained by *clicking* on the **Experimental** button. A best fit to an unknown reference spectrum (associated with all the molecules in the document rather than a single molecule) can be obtained by *clicking* on the **Reference** button.

## NMR

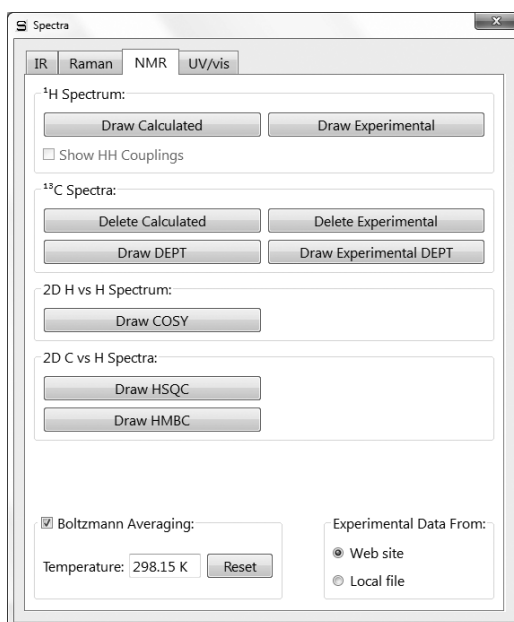Spectra are divided into four categories\*:  $^1\text{H}$  spectra,  $^{13}\text{C}$  and DEPT spectra, COSY spectra and HSQC and HMBC spectra. Access to experimental 1D spectra is provided either on-line or from a local file. *Clicking* on the **NMR** tab leads to the **NMR Spectra** dialog.




---


\* Chemical shifts for other nuclei may be displayed as labels (see **Configure...** under the **Model** menu; **Chapter 14**).



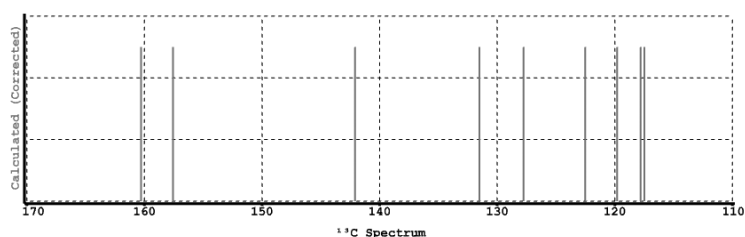Calculated proton,  $^{13}\text{C}$  and DEPT spectra may be drawn by *clicking* on their respective **Draw Calculated** buttons. Carbon shifts (only) for both  $^{13}\text{C}$  and DEPT spectra from the B3LYP/6-31G\* model are corrected based on topology and from the EDF2/6-31G\* model are corrected based on topology and bond orders. Carbon shifts from other models and proton shifts are not corrected.




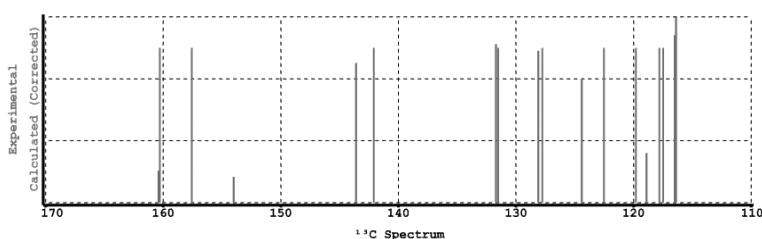Experimental  $^{13}\text{C}$  and DEPT spectra (if available) may be drawn together with or independently of the calculated spectra by *clicking* on **Draw Experimental**. Draw operations apply to all molecules in the list.



Experimental proton and  $^{13}\text{C}$  spectra are accessed from <http://nmrshiftdb.ice.mpg.de>. This comprises approximately 15,000 experimental spectra, all for organic molecules. Alternatively, the user may supply a spectra file (**Appendix I**). Selection (website vs. file) is made at the bottom right of the dialog.

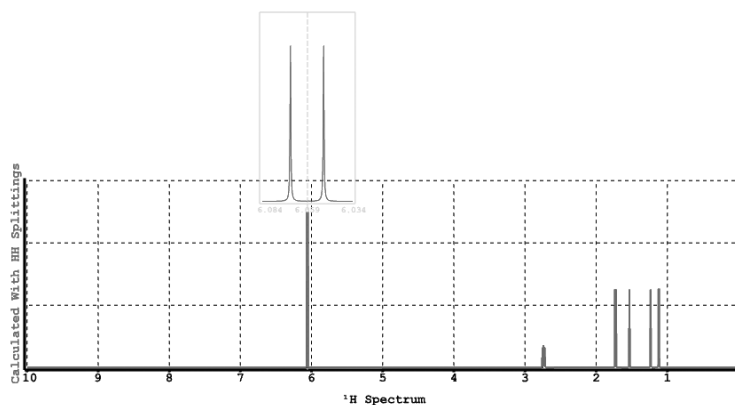
Proton spectra may be drawn to reflect spin-spin coupling due to neighboring nuclei. Coupling constants are not calculated, but are estimated empirically based on three-dimensional geometry. A magnifier is available under the **Plot Properties** dialog (see discussion earlier in this chapter).



Default ranges for proton spectra (0 to 10 ppm) and  $^{13}\text{C}$  spectra (0 to 150 ppm) may be changed by selecting **Properties** from the **Display** menu, *clicking* on the horizontal axis to bring up the **Plot Properties** dialog, *clicking* on the **X Axis** tab at the left of the dialog and editing the values under **Range**. Spectral plots may be translated and scaled.

A calculated spectrum may be removed by *clicking* on **Delete Calculated** (that has replaced **Draw Calculated**) and an experimental spectrum removed by *clicking* on **Delete Experimental**. The last spectrum removed will also cause removal of the plot axes. These operations apply to all molecules in the document. An individual (calculated or experimental) spectrum or entire plot may also be removed by selecting **Delete** from the **Build** menu (or holding down the **Delete** key) and then *clicking* on the spectrum or on the plot axes, respectively.

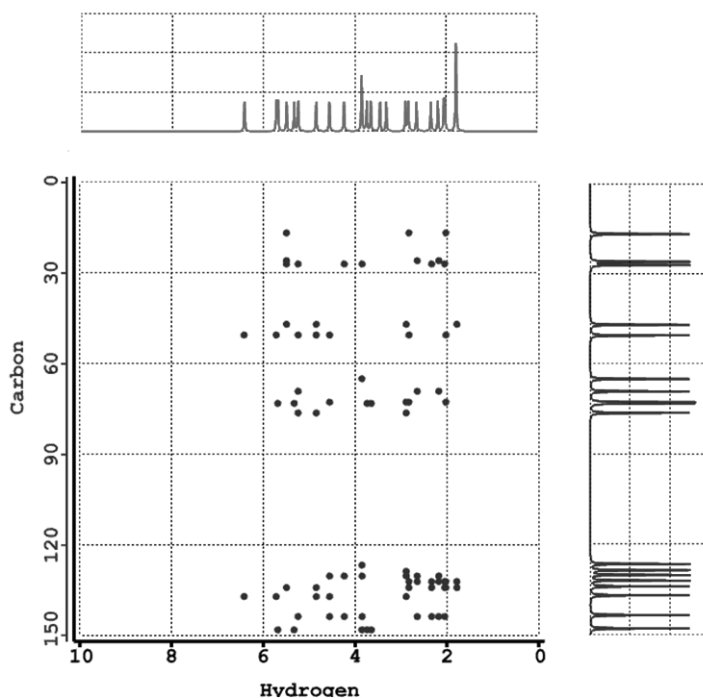
2D spectra are available from **Draw Calculated** buttons. COSY spectra identify protons three bonds removed, HSQC spectra identify protons directly attached to carbons and HMBC spectra identify protons that are two or three bonds removed. The corresponding proton and carbon spectra axes may be displayed along with 2D spectra by selecting them in addition to the 2D spectra. 2D spectra (and the 1D axes spectra) may be removed using the **Delete Calculated** buttons (that have replaced the **Draw Calculated** buttons). Experimental 2D spectra displays are not yet available. An HMBC spectrum is shown below.




Uncorrected spectra (all proton spectra and  $^{13}\text{C}$  and DEPT spectra from models other than B3LYP/6-31G\* and EDF2/6-31G\*) are likely to show significant systematic errors in chemical shifts. To scale either the calculated or experimental spectra in order to provide a best visual fit, *click* on the spectra (not on an axis) to select and move the mouse up and down while simultaneously depressing the **Shift** and **Alt** (**option** on Mac) keys in addition to the right button. The scale factor (multiplying the chemical shifts) will be displayed on the vertical axis. Scaling applies uniformly to all molecules in the document.

## UV/vis

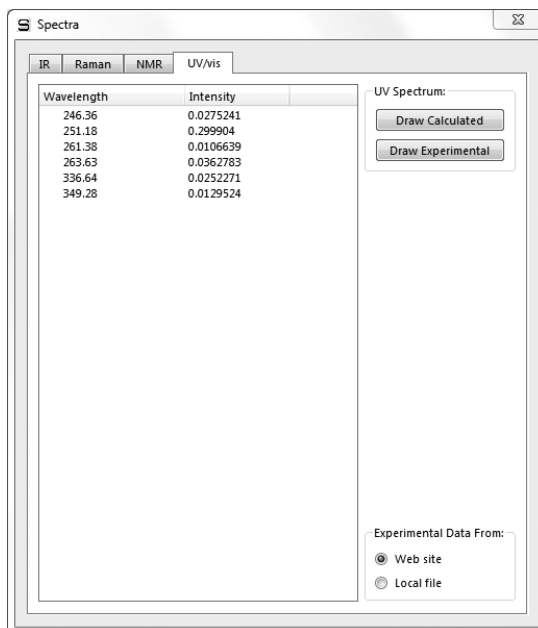
Clicking on the **UV/vis** tab leads to the **UV/vis Spectra** dialog that provides for display of UV/visible spectra.



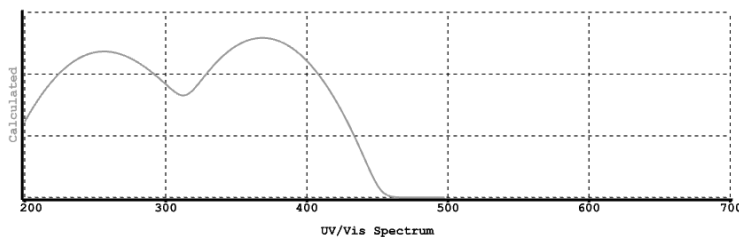This lists wavelengths and intensities. Click on **Draw Calculated** to draw the calculated spectrum.



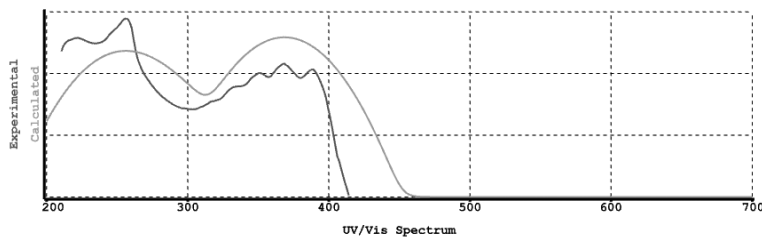Click on **Draw Experimental** to draw the corresponding experimental spectrum. Draw operations for both calculated and experimental UV/visible spectra and apply to all molecules in the list.




Experimental UV/visible spectra are accessed from the NIST (National Institute of Standards and Technology) website (<http://webbook.nist.gov/chemistry>). This comprises approximately 1,500 UV/visible spectra, primarily for organic molecules. *Spartan* keeps an index to facilitate retrieval but the actual data are on the website. Alternatively, the user may supply a spectra file (see **Appendix I**). Selection (website vs. file) is made at the bottom right of the dialog.

The default range for both calculated and experimental UV/visible spectra is from 200 nm to 700 nm. This corresponds to the range commonly measured and reported. This range may be changed by selecting **Properties** from the **Display** menu, *clicking* on the horizontal axis to bring up the **Plot Properties** dialog, *clicking* on the **X Axis** tab at the left of the dialog and editing the values under **Range**. The spectrum is treated like any other graphical object, and may be translated and scaled using the appropriate mouse operations.

The calculated spectrum may be removed by *clicking* on **Delete Calculated** (that has replaced **Draw Calculated**) in the dialog. Similarly, the experimental spectrum may be removed by *clicking* on **Delete Experimental** (that has replaced **Draw Experimental**). The last spectrum removed will also cause removal of the plot axes. These operations apply to all molecules in the document (not just the selected molecule). An individual (calculated or experimental) spectrum or entire plot may also be removed by selecting **Delete** from the **Build** menu (or holding down the **Delete** key) and then *clicking* on the spectrum or on the plot axes, respectively.


To help reveal similarities and differences between calculated and experimental UV/visible spectra, one or the other may be scaled. To scale, *click* on the spectra (not on an axis) to select and move the mouse up and down while simultaneously depressing the **Shift** and **Alt** (**option** on Mac) keys in addition to the right button. The scale factor will be displayed on the vertical axis. Scaling applies to all molecules in the document.

Only a single copy of one of the dialogs under **Spectra** may appear on screen; scaling will relate to the currently selected molecule. The dialog may be removed by *clicking* on .

## Formulas

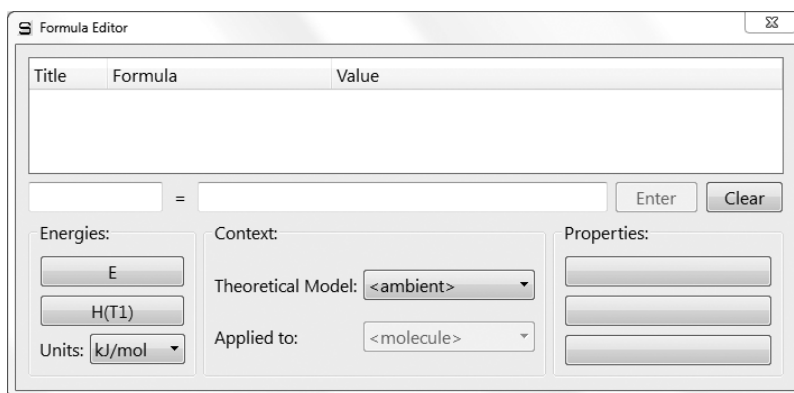The molecular mechanics and quantum chemical methods available in *Spartan* are able to provide a wide variety of molecular and atomic properties and QSAR descriptors. Quantities related to structure are available through functions under the **Geometry** menu (**Chapter 15**), while most other properties may be accessed from **Properties** dialogs or from the spreadsheet. In some cases, however, it is not the properties or descriptors themselves that are of interest, but rather *relationships* among them. For example, while the energies of the HOMO and LUMO of a molecule may not be of interest, the difference in their energies (the HOMO-LUMO gap) may be. Another important kind of relationship is that involving a change in property between reactants and products of a chemical reaction, in particular, the change in energy.

*Spartan* provides a tool to help define relationships among properties of a single molecule or of a single property between molecules that formally correspond to the reactants and products of a chemical reaction. The former will be discussed here and the latter will be discussed in the next chapter in the context of *data mining*. The **Formula Editor** dialog is accessed by selecting **Formulas** from the **Display** menu.



The bottom half of the dialog contains controls to build formulas, while the top half is a box to store finished formulas. As the most important property is arguably the energy, it is provided its own control (button). So too is the heat of formation obtained from the

T1 model. Both are found under **Energies** at the lower left of the **Formula Editor** dialog along with a menu to specify units. Other properties need to be accessed by *dragging* the property from its on-screen display (for example, a bond distance from its display at the bottom of the screen following selection of **Measure Distance** from the **Geometry** menu) or from one of the **Properties** dialogs (for example, the dipole moment in the **Molecule Properties** dialog) onto one of the three buttons under **Properties** at the lower right of the **Formula Editor** dialog. Following this, the button takes the name of the property. The property obtained from *Spartan* is now available for use in a formula.

The only other setting required is the **Theoretical Model** (under **Context** in the lower center of the **Formula Editor** dialog). The choices reflect the theory levels stored in databases available to *Spartan* (see discussion in **Chapter 19**). **Ambient** refers to the theoretical model actually used for the calculation of the selected molecule, and should not be selected if a calculation has not actually been run. The **Applied to** menu under **Context** will be grayed out if the molecule that *Spartan* is pointing to does not contain reaction information. Use of the formula editor to relate information for the reactants and products of chemical reactions is discussed in **Chapter 19**.

Pressing either the **E** or **H(T1)** buttons or a **Property** button onto which a property has been *dragged*, builds a formula and places it inside the second box at the center of the dialog. The default name (**F<sub>x</sub>**, x=0, 1, 2...) appears in the first box. Either the formula or its name may be edited. Pressing **Enter** (**return** on Mac) stores the formula with the default or edited name. Pressing **Clear** removes the formula. *Double clicking* on a formula in storage brings it back to the boxes at the center of the formula editor (but does not remove it from storage). The formula can then be edited and returned to storage (by *clicking* on **Enter**). If the name of the formula has not been changed, it will replace the original in storage. If it has been changed, it will be treated as a new formula. Once in storage, a formula can be *dragged* to wherever it is needed, in particular, into the spreadsheet, the **Plots** dialog or any of the dialogs associated with the *Spartan* Molecular

Database or Spartan Spectra and Properties Database. To remove a stored formula, *right click* on it and select **Delete** from the menu.

Formulas use the usual set of arithmetic operations and parenthesis to connect different properties. In addition, Boolean operators and standard mathematical functions designated in **Table 18-1** may be used. These need to be typed into the box at the center of the dialog. For example, to construct a formula for the HOMO-LUMO gap, first *drag* the LUMO energy from the **Molecule Properties** dialog onto one of the properties buttons in the **Formula Editor** dialog. *Pressing* the button leads to the text string **@LUMO\*@hart2eV** at the center of the dialog. Insert (type) “-” at the end of the string. Next, drag the HOMO energy from the **Molecule Properties** dialog onto a properties button in the editor and *press* the button. The text string is now **@LUMO\*@hart2eV-@LUMO\*@hart2eV**. Note that while orbital energies are evaluated in atomic units (hartrees), the formula includes a conversion factor **@hart2eV** to bring them into electron volts (eV). You can change units by manually editing the text string using the expressions provided in **Table 18-2** (preceded by **@**). Change the title to the left of the string from **Fx** to **HOMO/LUMO gap**, and *press* the **Enter (return)** key to the right of the string. The full formula is moved to the storage area at the top of the dialog. Note that the formula will be evaluated and numerical results presented if the appropriate data for the selected molecule with the selected theoretical model are available in the database (in this case, HOMO and LUMO energies).

A number of the quantities that are calculated by *Spartan* are not available from the graphical interface, but must be extracted from the archive of properties maintained by *Spartan*. Examples include orbital energies other than those corresponding to the HOMO and LUMO, and the polar surface area defined to include sulfur and any attached atoms in addition to nitrogen, oxygen and any attached hydrogens (the default). However, these and other quantities may be used in formulas by manually editing the text string. The expressions provided in **Table 18-3** need to be preceded by a **@**. For example, **@HOMO(-1)** accesses the energy of the molecular orbital immediately



**Table 18-1: Arithmetic and Boolean Operations and Mathematical Functions**

<b>arithmetic operations</b>		<b>boolean operations</b>	
+	addition	>	greater than
-	subtraction	>=	greater than or equal to
*	multiplication	<	less than
/	division	<=	less than or equal to
^	raise to a power	==	equal to
		!=	not equal to
			or
		&	and
<b>mathematical functions</b>			
ABS(x)	absolute value	LN(x)	natural logarithm
ACOS(x)	inverse cosine	LOG(x)	log (base 10)
ASIN(x)	inverse sine	SIN(x)	sine
ATAN(x)	inverse tangent	SQRT(x)	square root
COS(x)	cosine	TAN(x)	tangent
EXP(x)	exponential		

**Table 18-2: Conversion Factors and Constants**

ANGS2AU	Ångstroms to atomic units
AU2ANGS	atomic units to Ångstroms
EV2HART	eV to atomic units (hartrees)
EV2KCAL	eV to kcal/mol
EV2KJ	eV to kJ/mol
HART2KCAL	atomic units (hartrees) to kcal/mol
HART2EV	atomic units (hartrees) to eV
HART2KJ	atomic units (hartrees) to kJ/mol
KCAL2EV	kcal/mol to eV
KCAL2HART	kcal/mol to atomic units (hartrees)
KCAL2KJ	kcal/mol to kJ/mol
KJ2EV	kJ/mol to eV
KJ2HART	kJ/mol to atomic units (hartrees)
KJ2KCAL	kJ/mol to kcal/mol
PI	$\pi$

**Table 18-3: Calculated Quantities**

AREA	area of a user-defined plane ( $\text{\AA}^2$ )
ELECTRONEGATIVITY	QSAR descriptor $-(\text{HOMO}+\text{LUMO})/2$ (eV)
HARDNESS	QSAR descriptor $-(\text{HOMO}-\text{LUMO})/2$ (eV)
HOMO(-n)	energy of $n^{\text{th}}$ orbital below the HOMO (eV)
HOMOBETA(-n)	energy of the $n^{\text{th}}$ orbital below the $\beta$ HOMO (eV)
HYPERPOLARIZABILITY	beta hyperpolarizability parameter
INERTIA(i)	principle movements of inertia from largest (i=1) to smallest (i=3)
ISOTOPE(i)	mass number of atom i
LENGTH (i)	length of bond i ( $\text{\AA}$ )
LOGPC	LogP from Crippen model
LOGPV	LogP from Villar model
LUMO(+n)	energy of $n^{\text{th}}$ orbital above the LUMO (eV)
LUMOBETA(+n)	energy of $n^{\text{th}}$ orbital above the $\beta$ LUMO (eV)
PAREA (i, j, . . .)	surface area of space-filling model due to atoms of type i, j, ..., (atomic numbers or elemental symbols). i:, j:, ... signifies that attached hydrogens are also included. The surface area of nitrogen, oxygen and bonded hydrogens is the polar surface area. ( $\text{\AA}^2$ )
POLARIZABILITY	alpha polarizability parameter

**Table 18-4: Specialty Functions**

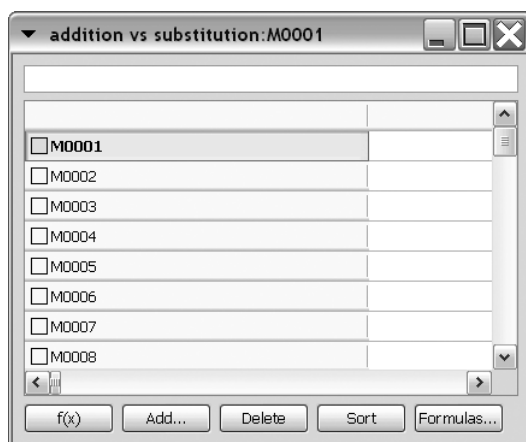
AVG (column name)	average of values in column
FITVAL (fit name)	column of fit values from regression analysis
MIN (column name)	minimum of values in column
MAX (column name)	maximum of values in column
NUM (column name)	number of defined entries in column
ROW	the number of the row in the spreadsheet
ROW(molecule name)	the number of the row of molecule
REF(i, x)	the value of the x referenced to row i
STDEV (column name)	standard deviation of values in column
SUM (column name)	sum of values in column

below the HOMO, while *@PAREA(7,8,16,7:,8:,16:)* accesses the polar area extended to sulfur and any attached hydrogens.




Formulas are used in conjunction with spreadsheet manipulations and plotting (both discussed later in this chapter) and with data mining (discussed in **Chapter 19**).

## Spreadsheet

Associated with each *Spartan* document (including documents comprising a single molecule) is a spreadsheet. This may be displayed by selecting **Spreadsheet**.



The spreadsheet comprises a series of rows (corresponding to different molecules in the document) and columns (corresponding to different molecular properties). This gives rise to cells, the number of which is the product of the number of rows (molecules) and the number of columns (molecular properties). The spreadsheet may be expanded or contracted by positioning the cursor at one of the corners, *pressing* the left mouse button and *dragging* the mouse.

Only one molecule from one document may be selected (although several molecules may be simultaneously displayed). Molecule selection follows either by *clicking* on the spreadsheet cell containing the molecule label or identifier (leftmost column), or by using the  and  buttons or the scroll bar at the bottom left of the screen. Molecules may be *animated* (stepped through in succession) using the  button at the bottom left of the screen. Animation speed may


be adjusted from the **Settings Preferences** dialog (**Preferences** under the **Options** menu; **Chapter 20**). Selection of a new molecule in the document results in deselection of the previously selected molecule. A molecule may be designated for permanent display by *checking* the box to the left of its identifier in the spreadsheet. The molecules in a document may either be translated and rotated in concert or manipulated independently. This is controlled by **Coupled** (which operates in toggle mode) under the **Model** menu (**Chapter 14**). By default (**Coupled** checked) molecules move in concert. *Uncheck Coupled* to move them independently.

Upon initial entry, all columns of the spreadsheet except the leftmost column, are blank. The leftmost column contains a label that may be changed either by directly typing a new label into the spreadsheet or into the **Label** box in the **Molecule Properties** dialog (see discussion earlier in this chapter). Additionally, default identifiers (**M0001**, ...) may be replaced by chemical names if the molecule exists in the Spartan Molecular Database (SMD) by *clicking* on **Label** (the leftmost header cell), then *right clicking* and selecting **Rename Selected using SMD** from the menu that results. To replace all identifiers, *right click* inside the header cell of the leftmost column and again select **Rename Selected using SMD**. Renaming also applies to molecules in the Spartan Spectra and Properties Database (SSPD) in which case the menu entry is **Rename Selected Using SSPD**.


Information may be added to the spreadsheet in several different ways:

### **From the Add Dialog**

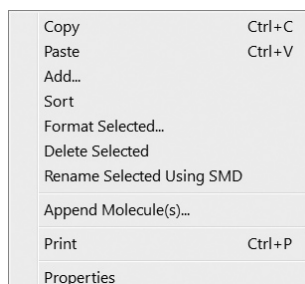A selection of molecular properties may be entered into the spreadsheet by first *clicking* on the header cell of an empty column, and then *clicking* on **Add...** at the bottom of the spreadsheet. Alternatively, *right click* inside the header cell and then to select **Add...** from the menu that results.



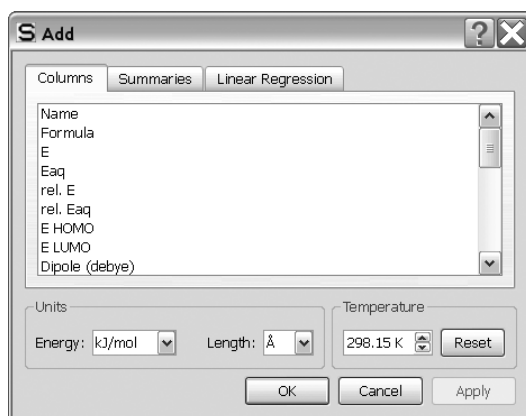This leads to the **Add Columns** dialog.




Name	molecule name as it appears in SMD or SSPD
Formula	molecular formula
E	energy (heat of formation, strain energy)
E <sub>aq</sub>	aqueous energy (heat of formation, strain energy) based on the SM5.4 model
rel. E	energy (heat of formation, strain energy) relative to selected molecule
rel. E <sub>aq</sub>	aqueous energy (heat of formation, strain energy) relative to selected molecule based on the SM5.4 model
E HOMO	energy of highest-occupied molecular orbital
E LUMO	energy of lowest-occupied molecular orbital
Dipole	dipole moment (in debyes)
Molecular Weight	molecular weight (in amu)
Boltzmann Distribution	Boltzmann distribution based on energy

Boltzmann Distribution (aq)	Boltzmann distribution based on aqueous energy as estimated from SM5.4 model
Alignment Scores	$1 - R^2/N$ , where $R^2$ is the root mean square distance and N is the number of alignment centers. 1 is a perfect score

One or more properties may be added to the spreadsheet by *clicking* on their entries, then specifying appropriate units from the **Energy** and **Length** menus, varying temperature from the **Temp** menu (as applied to the Boltzmann distributions), and finally *clicking* on **OK** or on **Apply**. In the former case, the dialog is dismissed and in the latter it is left on screen. *Clicking* on **Cancel** or  removes the dialog without further affecting the contents of the spreadsheet. (Any previous additions to the spreadsheet as a result of *clicking* on **Apply** are maintained.)

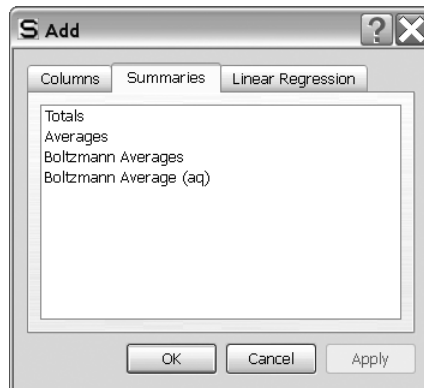
Column totals, averages and Boltzmann-weighted averages may be accessed by *clicking* on the **Summaries** tab. This leads to the **Add Summaries** dialog.



Totals	sum of column values
Averages	average of column values
Boltzmann Averages	Boltzmann weighted average of column values
Boltzmann Averages (aq)	Boltzmann weighted average of column values using energies corrected for aqueous environment based on the SM5.4 model

*Clicking* on one or more of these followed by *clicking* on **OK** or **Apply**, leads to the requested summaries as rows at the bottom of the spreadsheet, identified as **Totals**, **Averages**, etc. (one row for each summary). Note, that this dialog can be brought up directly by *clicking* on an empty row header (instead of an empty column header) inside the spreadsheet prior to *clicking* on **Add...**

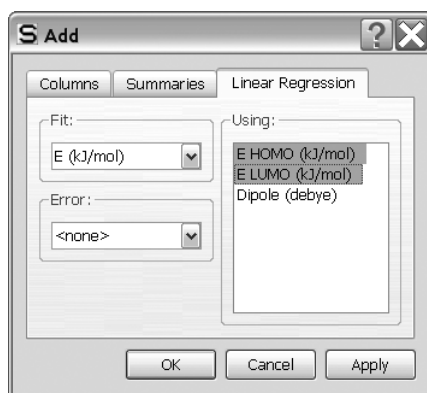Finally, a linear regression analysis may be performed on the data in the spreadsheet. First *click* on the **Linear Regression** tab to bring up the **Add Linear Regression** dialog.



Select one entry from the **Fit** menu and one or more entries from the list under **Using**. *Clicking* on **OK** or **Apply** performs the linear regression analysis and places the results in a row at the bottom of the spreadsheet identified by **Fit**. As many regression analyses as desired may be performed on the data in the spreadsheet. The individual results will be entered as separate rows in the spreadsheet, with names **Fit1**, **Fit2**, etc. Additional information about the regression analyses is available from the **Regression Properties** dialog (see discussion earlier in this chapter).

### **From Post (P) Buttons**

Post buttons (P) found in a number of properties dialogs provide an alternative method to the **Add** dialog for entering calculated properties into the spreadsheet. Note that some properties may require user specification. These include individual bond distances, angles and dihedral angles (available from **Measure**

**Distance**, **Measure Angle** and **Measure Dihedral** under the **Geometry** menu; **Chapter 15**), bond distance, angle and dihedral angle constraints (available from **Constrain Distance**, **Constrain Angle** and **Constrain Dihedral** under the **Geometry** menu; see discussion in **Chapter 15**), atomic charges, chemical shifts (available from the **Atom Properties** dialog; this chapter), the accessible area of an electron density surface, the polar area and accessible polar area of an electrostatic potential map, the area of a selected region (band) of a banded property map, minimum and maximum property values on a map and the value of the property at a specific location on a property map (available from the **Surfaces Properties** dialog; this chapter). With the exception of the property value on a map and the area of a selected band, post generates an entire column. Where atom labels are involved, for example, in defining a specific distance, post can be expected to yield consistent results for all molecules in a document only where the molecules are closely related, for example, molecules resulting from a conformational search, or where labels have been explicitly reassigned\*. The property value and the area of a selected band on a map is posted only for the selected molecule. Post buttons are also available for CAS numbers and experimental heats of formation contained in *Spartan's* internal database, and for T1 heats of formation contained in both SMD and SSPD databases.

### From the Clipboard

Properties of one or more molecules in a document may be copied into the clipboard and then pasted into individual spreadsheet cells. These include (but are not restricted to) bond distances, angles and dihedral angles (**Measure Distance**, **Measure Angle** and **Measure Dihedral** under the **Geometry** menu), bond distance, angle and dihedral angle constraints (**Constrain Distance**, **Constrain Angle** and **Constrain Dihedral** under the **Geometry** menu), atomic charges and chemical shifts (**Atom Properties** dialog), and the value of a property on a property

---

\* Label reassignment is accomplished using the **Atom Properties** dialog (see discussion earlier in this chapter).



map and the area of a selected band (**Surface Properties** dialog). To copy the spreadsheet, first highlight the numerical value of the property in the appropriate screen location (distances, etc.) or dialog (charges, etc.), then select **Copy** from the **Edit** menu, then *click* on the appropriate (destination) cell in the spreadsheet, and finally select **Paste** from the **Edit** menu.

## User-Defined Expressions

An expression may be entered either into a header cell (in which case it refers to all molecules) or into an individual cell (in which case it refers only to a single molecule). Expressions in the column header take the form *name=formula*, where *formula* is made up of arithmetic operations (**Table 18-1**), specialty functions (**Table 18-4**), calculated quantities (**Table 18-3**), conversion factors and constants (**Table 18-2**) in addition to numerical values. References to specialty functions, calculated quantities and conversion factors and constants must be preceded by @. For example, *mu=@DIPOLE* typed into a header cell gives the dipole moment. Some functions have arguments, for example, *c1* and *c2* in the expression *c12=@DISTANCE(c1,c2)* refer to atoms *c1* and *c2*, while *3* in the expression *orbitale=@HOMO(-3)* designates the energy of the molecular orbital three orbitals below the HOMO. It is necessary to *press* the **Enter** key (**return** key on Mac) following entry of the expression into a cell. The leading *name=* is optional for individual entries.

Formula construction inside the spreadsheet uses the same rules and references the same set of functions described in **Formulas** earlier in this chapter. It is more convenient to construct expressions in the formula editor and *drag* them into the spreadsheet. Note that spreadsheet expressions refer to the molecules in the spreadsheet and not to their reaction products.

## Numerical Data

Numerical data may be entered by typing directly into the spreadsheet. A column header first needs to be specified. *Double*

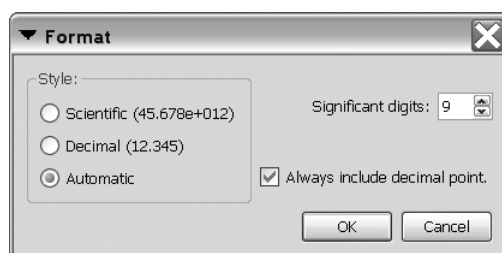
*click* on an empty column header cell, type in a name and *press* the **Enter** key (**return** key on Mac). Then, type the data into individual cells of the new column (*press* the **Enter** key following each data entry). Alternatively, use the ↓ key.

Each row in a spreadsheet corresponds to a molecule in a document, and new rows are automatically added in response to adding new molecules to the document. New molecules may be added by building (**New Molecule** under the **File** menu; **Chapter 12**), by appending one or more existing documents each containing one or more molecules using either **Append Molecule(s)...** under the **File** menu (**Chapter 12**), or by *right clicking* inside the header cell of the first available row and selecting **Append** from the menu that appears, by pasting from the clipboard, or by *dragging* from the file system. To copy a molecule into the clipboard, first select (*click* on) it, and then select **Copy** from the **Edit** menu, or *click* on its identifier (left most column) in its spreadsheet, and then select **Copy** from the **Edit** menu. Alternatively *right click* either on the molecule or on its identifier in the spreadsheet and select **Copy** from the menu that appears. Use of the clipboard permits several molecules to be selected (and copied) at once using the **Shift** and **Ctrl** keys in the usual manner. To copy the contents of the clipboard to its destination, *click* on an empty row header in the spreadsheet (for the destination document), and then select **Paste** from the **Edit** menu. An alternative to the two-step **Copy-Paste** procedure is to *drag* the molecule or set of molecules from one spreadsheet to another.

A row (molecule) may be deleted from a spreadsheet, either by first selecting the molecule and then selecting **Delete Molecule** from the **File** menu (**Chapter 12**), or by first *clicking* on its identifier in the spreadsheet (leftmost column) and then either *clicking* on the **Delete** button at the bottom of the spreadsheet, or by *right clicking* on its identifier in the spreadsheet and then selecting **Delete Selected** from the menu that appears. In all cases, a warning is provided prior to deletion. An entire column in the spreadsheet may be deleted by first *clicking* inside its header cell and then *clicking* on the **Delete** button (or **Delete Selected** from the menu).

Rows in the spreadsheet may be sorted according to the numerical values in any column either by first *clicking* inside the header cell and then *clicking* on the **Sort** button at the bottom of the spreadsheet or by *right clicking* inside the header cell and selecting **Sort** from the menu that appears. The rows are placed in ascending order, the smallest (least positive) value of the selected property at the top, largest (most positive) value at the bottom. To sort in descending order, hold down the **Shift** key before *clicking* on the **Sort** button or selecting **Sort** from the menu.


Information in one or more columns of the spreadsheet may be formatted by *right clicking* inside the header cell(s) and selecting **Format Selected** from the menu that appears.



Format as desired and *click* on **OK** to remove the dialog. The full contents of the spreadsheet may be formatted by *right clicking* inside the header cell for the left most column and then selecting **Format Selected** from the menu.

A button at the bottom right of the spreadsheet toggles between numerical representation of data, **f(x)**, and formula presentation, **=?**.

The spreadsheet may be printed by selecting **Print Spreadsheet...** from the **File** menu (that replaces **Print...** when a spreadsheet is selected).


As many spreadsheets as desired (corresponding to the same or to different documents) may be open on screen. A spreadsheet is removed when the associated document is closed and may also be removed by *clicking* on .

The contents of the spreadsheet may be brought into Excel™ using the clipboard. Select whatever cells are to be copied, select **Copy** from the **Edit** menu. Alternatively, *right click* with the proper cells selected and select **Copy** from the menu that appears. **Paste** into Excel.


The contents of an Excel spreadsheet may be brought into *Spartan*. Copy whatever information is to be transferred to the clipboard, move into *Spartan*, then *click* on the appropriate cell and select **Paste** from the **Edit** menu (or *right click* on the appropriate cell and select **Paste** from the menu that appears). Note, that information on the clipboard that goes beyond the number of rows in *Spartan*'s spreadsheet will be ignored.

## Plots...

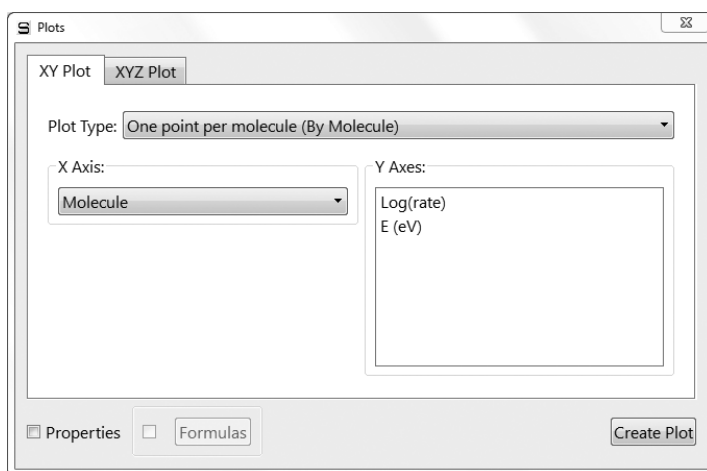Plots may be constructed from data in a spreadsheet and a variety of simple curves fit to these data. Selection of **Plots...** from the **Display** menu leads to the **Plots** dialog.



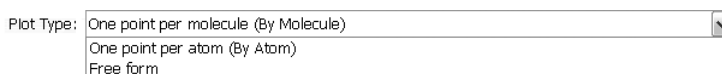A tab at the top of the dialog selects between **XY Plot** (two dimensions) or **XYZ Plot** (three dimensions). Common to both is the **Plot Type** menu.



This distinguishes between plots based on the entire set of molecules in a document (**By Molecule**), to those based on the atoms in a single molecule (**By Atom**). The latter is the mode used to plot calculated


IR, Raman and NMR spectra. The third choice (**Free Form**) is used to construct single molecule plots, for example, calculated, UV/vis spectra as well as plots of experimental IR, Raman, NMR and UV/visible spectra.

By atom and free form plots make use of vector capabilities in *Spartan's* spreadsheet. See the FAQ *Vector Spreadsheet Operations* available under **Help** from the **Help** menu (**Chapter 21**).

Also common to both is the **Properties** checkbox at the bottom left of the dialog. Unchecked, this provides access for the **X-Axis** (X-Axis and Y-Axis for an XYZ plot) menu and **Y-Axes** (**Z-Axes** for an XYZ plot) only to quantities that have been entered in the spreadsheet. Checked, it provides access to all quantities in the *Spartan's* Property Archive.

## XY Plots

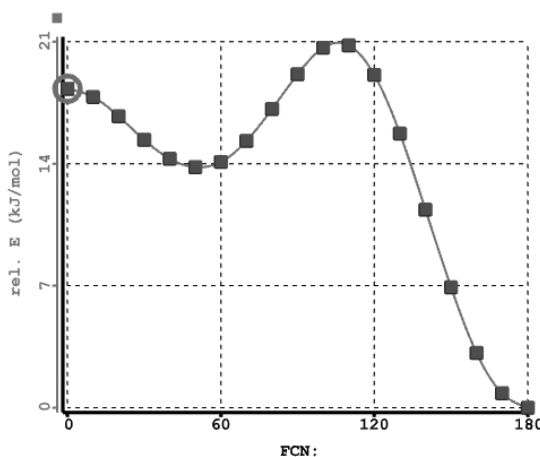
The XY Plot incorporates an **X Axis** menu designating the molecular property to be displayed among the X (horizontal) axis, and a list of properties to be displayed along the Y (vertical) axis. These properties correspond one to one to the columns in the associated spreadsheet. To construct a 2D plot, select an item from the **X Axis** menu, then *click* on one or more items from the **Y Axes** list and finally *click* on **OK**. (Repeated *clicking* on a property in the **Y Axes** turns the list on and off.) The dialog is removed from the screen and a plot appears.




The plot may be moved about the screen. Select (*click* on) either the frame or on one of the curves (it will turn gold in response), then hold down the right mouse button and *drag* the mouse. The plot may also be scaled (expanded or shrunk) by first selecting either the frame or one of the curves and *dragging* the mouse while holding down both the right button and the **Shift** key. The plot may not be rotated.

## XYZ Plots

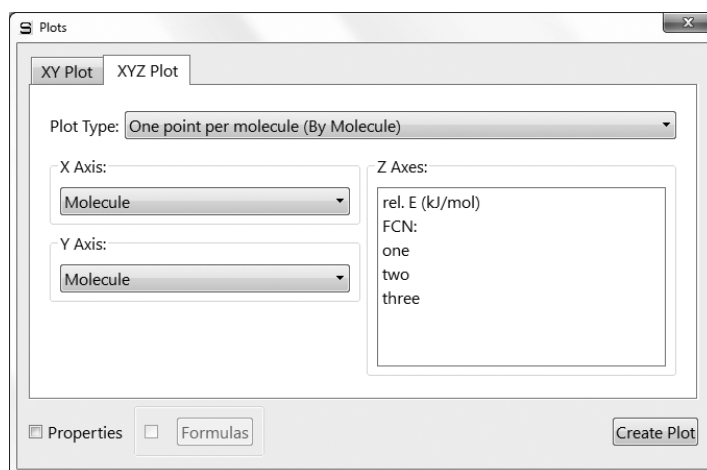The XYZ Plot is reached by clicking on the **XYZ Plot** tab at the top of the dialog. This leads to a new dialog.




This is very similar to the XY Plot, except that it incorporates **X Axis** and **Y Axes** menus to designate quantities to be displayed along the X and Y axes, respectively, as well as a list of properties to be displayed along the Z axis. To construct a 3D plot, select one item from each of the **X Axis** and **Y Axis** menus, then *click* on one or more items from the **Z Axes** list and finally *click* on **OK**. A plot appears.

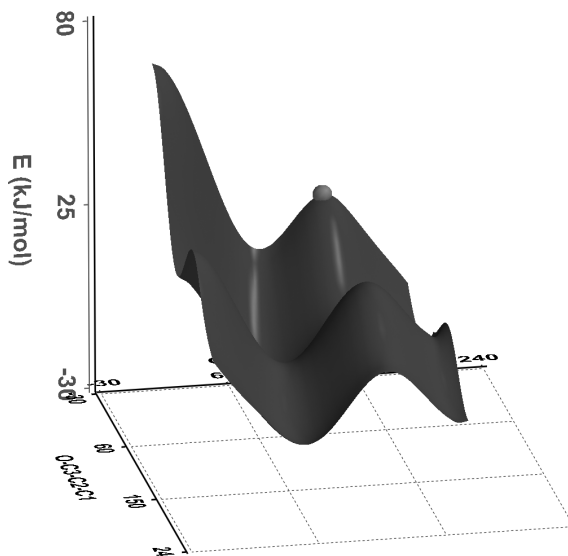

The plot can be moved about the screen. Select (*click* on) either the frame or on one of the curves (it will turn gold in response), then hold down the right mouse button and *drag* the mouse. The plot may also be scaled (expanded or shrunk) by first selecting either the frame or one of the curves and *dragging* the mouse while holding down both the right button and the **Shift** key. Finally, 3D plots (unlike 2D plots discussed earlier) may be rotated, by first selecting either the frame or one of the curves and *dragging* the mouse while holding down the left button.

Plots initially presented connect the data points with cubic splines. A variety of curve fits (linear, least squares, Fourier, etc.) as well as a variety of different presentation formats are available under the **Plot Properties** and associated dialogs (see **Properties** earlier in this chapter).

A plot may be deleted by first selecting (*clicking* on) **Delete** from the **Build** menu, or *clicking* on , or holding down the **Delete** key, and then *clicking* on either the plot frame or on an individual curve. *Clicking* on the frame removes it and all associated curves, while *clicking* on a curve removes only this curve. However, the last curve deleted also deletes the frame. IR, Raman, NMR and UV/vis spectral plots may also be removed from the corresponding **Spectra** dialog.

## Similarities...

The results of a similarity analysis based either on selected atoms in a structure or on a set of CFD's as a template may be viewed using the **Similarities** dialog.

Template Molecule	Library Molecule	Score	Collisions	Conformer	Enantiomer
M001	terfenadine	0.79	0	176	Inverted
M001	terfenadine	0.78	0	194	Inverted


The box at the right of the dialog identifies each hit in terms of the identity of the **Template Molecule**, the **Library Molecule**, and in the case where a library molecule has been expanded in terms of a set of diverse conformers, the identity of the conformer. It also provides a score\* and, in the case where excluded volume CFD's are present in the template, the number of collisions between these elements and the library molecule. Finally, it indicates whether the enantiomer for the molecule included in the library has been inverted. Entries in the dialog may be sorted by score or number of collisions (as well as by the names of the template and library molecules) by *clicking* on the appropriate column header at the top of the box.

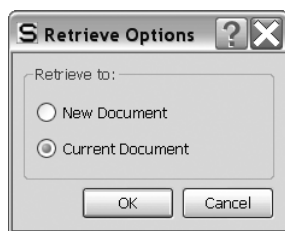
A minimum score can be set from the **Minimum Score** menu at the bottom right of the dialog. The number of matches that satisfy this minimum is reported.

\* The score is defined as  $[(1-R^2)/N]$ -penalty, where  $R^2$  is the root mean square distance between template and library molecule centers and  $N$  is the number of similarity centers. Whereas molecule alignment as a prelude to similarity scoring is based on atomic, CFD and pharmacophore centers **only**, penalties are given for alignments that lead to unfavorable steric interactions between template and library entries or (in the case of molecule alignment by CFD's) to incorrect orientation of hydrogen bond donor or acceptor CFD's. These are subtracted from the score based on  $R^2$ . 1 is a perfect score.



As hits are selected by *clicking* on them in the box, the library entry superimposed onto the template will be displayed in the window at the left of the dialog. This graphic may be manipulated (rotated, translated, scaled) with the usual mouse/keyboard commands (the cursor needs to be positioned inside the window). The display style cannot be altered, nor can any measurements be made.

The library entry associated with a particular hit may be retrieved (brought into *Spartan's* main window) by first selecting it and then *clicking* on the **Retrieve** button at the bottom of the **Similarities** dialog. The template cannot be retrieved. Retrieval can either be to a New Document or to the Current Document depending on the setting in the **Retrieve Options** dialog accessed by *clicking* on  to the right of the **Retrieve** button.



Retrieval to the Current Document (that with the template) allows both the template and one or more hits from the library to be displayed simultaneously, with a full range of model styles and display options.


The **Similarities** dialog is closed by *clicking* on the **Done** button.

## Reactions...

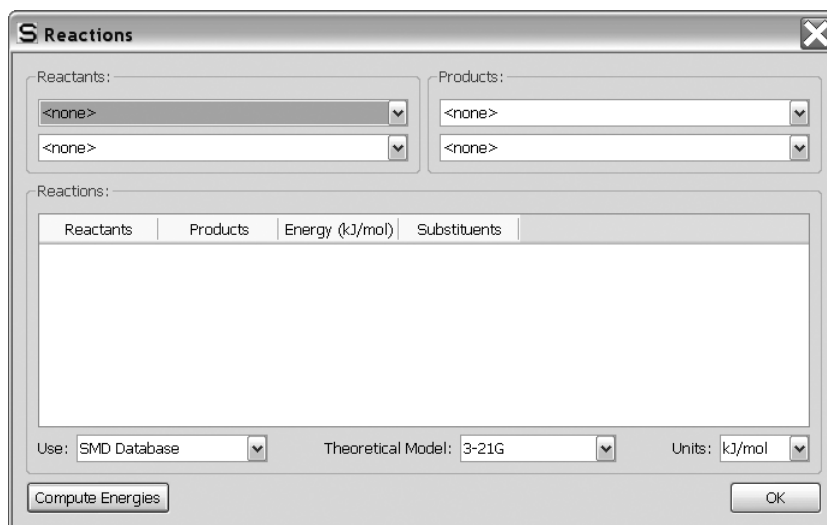Data entered in a *Spartan* document may be used to calculate reaction energies including activation energies.


$$\Delta E = E_{\text{product1}} + E_{\text{product2}} - E_{\text{reactant1}} + E_{\text{reactant2}}$$

Selection of **Reactions...** from the **Display** menu leads to the **Reactions** dialog.




Two menus under **Reactants:** and two menus under **Products:** identify the reactants and products of reaction. They correspond to the labels (identifiers) of the molecules in the document, plus a null entry **<none>**. Note that the overall reaction is mass balanced.

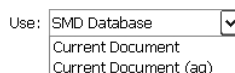The **Use** menu identifies the source of the energies to be used in the reaction energy calculation.




**SMD Database** or **SSPD Database** refers to energies drawn from entries in the Spartan Molecular Database (SMD), or Spartan Spectra and Properties Database (SSPD), **Current Document** refers to energies in document and **Current Document (aq)** refers to energies in the document that have been empirically corrected for aqueous solvent using the SM5.4 model.

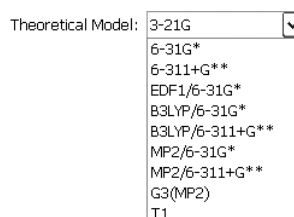Selection of **SMD Database** or **SSPD Database** (from the **Use** menu) requires specification of theoretical models. This is made from the **Theoretical Model** menu.




Note that while T1 is not included in either SMD or SSPD, T1 heats of formation have been made available.

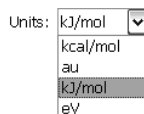The **Units** menu allows for selection of units.



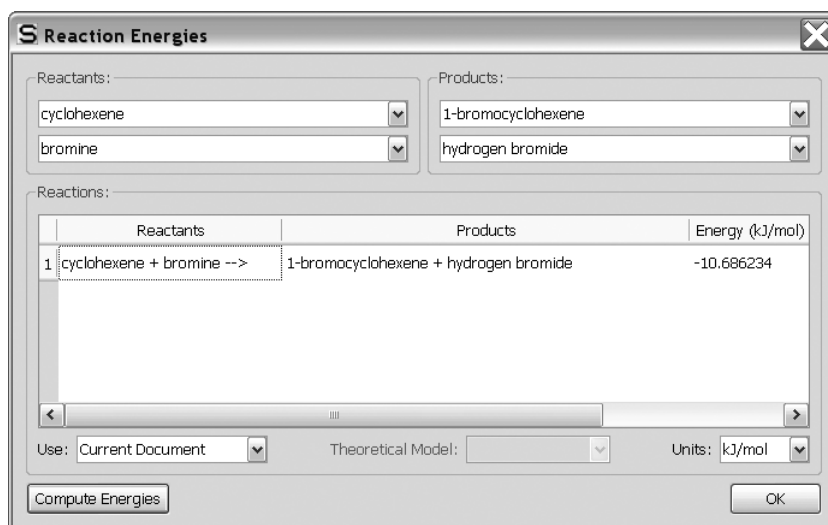A reaction energy is computed by *clicking* on **Compute Energies** at the bottom left of the dialog.



Substituted molecules can be used in place of real molecules as reactants and products. By accessing SMD or SSPD, this means that the energies for a series of related (by substitution) reactions can be computed at once.\*

The results of a reaction energy calculation may be printed by *right clicking* inside the display area of the **Reactions** dialog and selecting **Print** from the menu that results.

The **Reactions** dialog is closed by either *clicking* on **OK** or on .

---


\* An example is provided in the tutorial *Positional Selectivity of Substituted Naphthalenes* in **Chapter 5**.

# Chapter 19
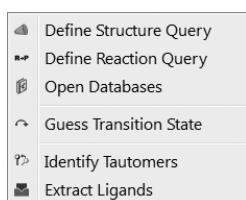
## The Search Menu

---

Functions available under the **Search** menu provide for defining structure and reaction queries to be used in searches of databases attached to **Spartan**. These include the Cambridge Structural Database (CSD) of experimental X-ray structures, the Spartan Molecular Database (SMD) of calculated structures, energies, and properties, the Spartan Spectra and Properties Database (SSPD) of calculated spectra, properties and QSAR descriptors, the Spartan Infrared Database (SIRD) of calculated infrared spectra and the Spartan Reaction Database (SRD) of calculated transition states. CSD and SRD are searchable by substructure, whereas SMD and SSPD are searchable by substructure, name, formula, molecular weight and isomer and SIRD is searched based on matching to an unknown infrared spectrum. The **Search** menu also accesses procedures for automatically guessing transition states for reactions based on their similarity to transition states in SRD, for identifying tautomers, and for extracting ligands from PDB files.




### Define Structure Query ( )

This designates one or more free valences on a structure as attachment points. *Clicking* on the free valence marks it with an orange cone and designates it as an attachment point. This means that anything directly bonded to the valence (including hydrogen) is an acceptable hit for a search on a database. For example, *clicking* on the *para* free valence in toluene designates that any *para*-substituted toluene is an acceptable hit. *Clicking* on both free valences of hydrogen peroxide means that all peroxides including hydrogen peroxide and

all hydroperoxides) are acceptable hits. Free valences not designated as attachment points are assumed to be hydrogens. *Clicking* on a free valence previously designated as an attachment point removes the cone and its designation as an attachment point.

### Define Reaction Query ( )

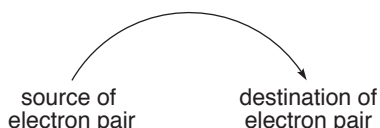This serves either to identify the product or products of a chemical reaction from the reactants to be used for mining the Spartan Molecular Database or Spartan Spectra and Properties Database, or to provide an entry into the Spartan Reaction Database in order to access a previously calculated transition state, or to identify a similar transition state. All three functions make use of a set of curly arrows familiar to generations of organic chemists and formally corresponding to the movement of a pair of electrons.



There are two possible sources of an electron pair and three possible destinations, leading to six combinations:

lone pair → lone pair	move lone pair
lone pair → bond	use lone pair to increase bond order
lone pair → space between atoms	use lone pair to create new (single) bond
bond → lone pair	decrease bond order to make lone pair
bond → bond	decrease order of one bond to increase order of another bond
bond → space between atoms	decrease order of one bond to make a new (single) bond

The first of these is a null operation, and its inclusion has no effect.

The reaction query procedure presently does not distinguish between a lone pair and a radical center. Examples are provided at the end of this section.



Selecting **Define Reaction Query** results in a message at the bottom left of the screen.

Select atom or bond as tail.

The tail of the arrow corresponds to the source of the electron pair. If the source is a lone pair, then select the atom that holds the lone pair. If the source is a bond, then select the bond. *Clicking* on an atom or bond highlights (colors gold) the atom or bond and leads to a new message at the bottom left of the screen.

Select atom, bond, or two atoms as head. If two atoms hold SHIFT key.

The head of the arrow corresponds to the destination of the electron pair. If the destination is an atom (leading to a lone pair), then select the atom that will hold the lone pair. If the destination is an existing bond (leading to an increase in bond order from single  $\rightarrow$  double or double  $\rightarrow$  triple), then select the bond. If no bond presently exists, *press* the **Shift** key and select the two atoms that will become bonded upon reaction. These operations result in a curved arrow being drawn on the reactant structure. This extends from an atom, or the center of a bond to an atom, or the center of a bond, or the center of a dotted line that has been drawn between atoms that are to be bonded. The original message returns to the bottom left of the screen.


The process (tail selection followed by head selection) is repeated as necessary to fully define the reaction. Incorrect reaction arrows may be removed by selecting **Delete** from the **Build** menu () followed by *clicking* on the arrow to be deleted. You need to again select **Define Reaction Query** () in order to continue arrow selection. More simply, *click* on the arrow(s) to be deleted while holding down the **Delete** key.

Once defined, reaction queries can be used to relate the properties of the products of chemical reactions to those of the products for the purpose of mining the Spartan Molecular Database or Spartan Spectra and Properties Database. Full discussion is provided in **Open Databases** later in this chapter. Reaction queries can also be used to search the Spartan Reaction Database for transition-state structures associated with the defined reaction (also discussed in **Open Databases**). In

both of these cases, structure queries (discussed in **Define Structure Query** earlier in this chapter) will also normally be employed. Finally, reaction queries may be employed to automatically provide a guess at a transition state based on similarity to an entry in the Spartan Reaction Database (see discussion in **Guess Transition State**) later in this chapter. Here, no structure queries would be defined.


### Examples of Reaction Queries
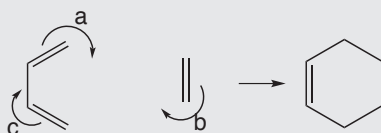
#### *Diels-Alder reaction of 1,3-butadiene and ethylene*




- a, b. double bond to empty space leading to single bond and single bond
- c. double bond to single bond leading to single bond and double bond

#### *S<sub>N</sub>2 reaction of chloride and methyl iodide*

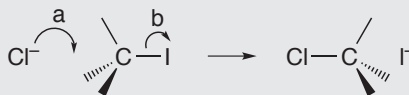


- a. lone pair on Cl<sup>-</sup> to empty space leading to ClC bond
- c. CI bond leading to lone pair on I

#### *Ene reaction of 1-pentene*

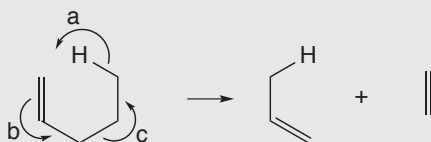

- a. CH bond to empty space leading to CH bond
- b. double bond to single bond leading to single bond and double bond
- c. single bond to single bond leading to empty space and double bond

### Ring closure of 1-hexenyl radical to methylcyclopentyl radical

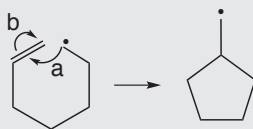

- radical center to empty space (between C<sub>2</sub> and C<sub>6</sub>) leading to single bond
- double bond to terminal carbon leading to single bond and radical center

## Open Databases (📖)

This provides access to the Cambridge Structural Database (CSD) of experimental X-ray crystal structures, the Spartan Molecular Database (SMD) of calculated structures, energies, spectra and molecular properties, the Spartan Spectra and Properties Database (SSPD) of spectra, properties and QSAR descriptors, the Spartan Reaction Database (SRD) of calculated transition states and the Spartan Infrared Database (SIRD) of calculated infrared spectra.\* CSD, SMD and SRD may be searched based on substructure. SMD and SSPD may also be searched by name, formula, molecular weight and isomer. If a selected molecule on screen exists in SMD and/or in SSPD, it may be replaced by a database entry and, even if it is not replaced, the molecule name in SMD may be used. SIRD is searched by comparison with an input infrared spectrum.

## Cambridge Structural Database (CSD)

The Cambridge Structural Database (CSD) is a collection of nearly 500,000 experimental X-ray crystal structures for organic and organometallic molecules.\*\* It is maintained by the Cambridge Crystallographic Data Centre (CCDC) and grows approximately 10% per year. CSD not only contains information about the molecular geometry, but also about the manner in which molecules pack in

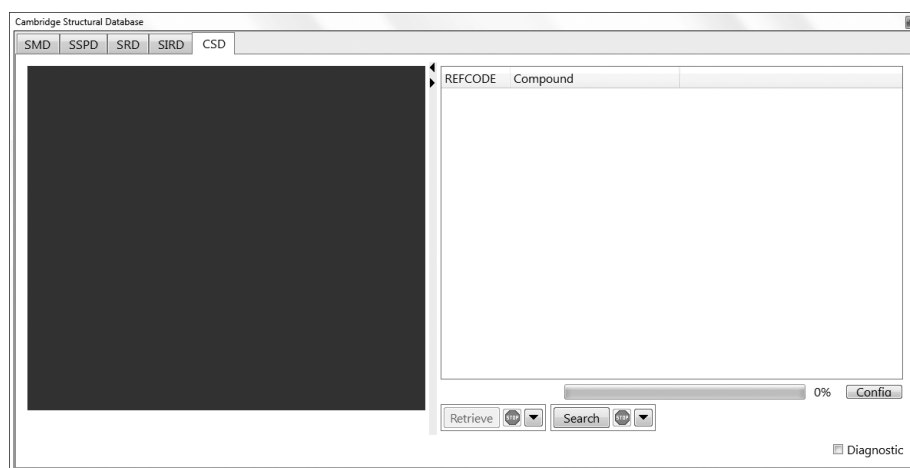
\* Databases may be located on the local machine or on a remote server. Setup details are provided in **Appendix F** for CSD and in **Preferences** under the **Options** menu in **Chapter 20** for SMD, SSPD, SRD and SIRD.

\*\* CSD is not included with *Spartan* but is available by subscription from CCDC (<http://www.ccdc.cam.ac.uk/>). For installation instructions, see **Appendix F**.



the solid state\*. In short, the CSD represents a gold mine of detailed experimental structural information, and also serves to identify molecules that can be (and have been) synthesized and purified through crystallization.

Access to the data in CSD is via a substructure searching procedure. A molecule is first constructed and one or more attachment points are identified. *Clicking* on the **CSD** tab at the top of the **Database** dialog leads to the **CSD** dialog.



This contains a window at the left for previewing structures, a box at the right for listing hits on CSD and buttons at the bottom to control aspects of the search and subsequent transfer of CSD data into the file system. Both the preview window and the text box may be scaled independently. A substructure search is set up by specifying one or more of the following\*\*:

(i) **Attachment Points**

A free valence is designated as an attachment point by first selecting **Define Structure Query** from the **Search** menu and then *clicking* on it. It will be marked with an orange cone. *Clicking* again removes the designation and the cone.

\* *Spartan*'s interface to CSD accesses only molecular information. Additional information related to packing in the crystal may be accessed using CCDC's program ConQuest, that is supplied as part of the CSD subscription.

\*\* If nothing is specified, a search for an exact match is carried out assuming that all atoms and bonds are exactly as given in the query and that all free valences are assumed to be hydrogens.


Anything may be grown off an attachment point (including hydrogen), although substituents are not permitted. Free valences not designated as attachment points are assumed to be hydrogens.

(ii) **Wild-Card Atoms**


An atom may be designated as a wild-card (meaning that element type is unimportant) by first selecting **Define Structure Query** from the **Search** menu and then *clicking* on it. An orange ball will surround the atom. *Clicking* again removes the wild-card designation and the ball. Use of wild-card atoms will result in structures that incorporate variants of the original substructure with different atoms at designated positions. For example, wild-card designation of an unsubstituted position of a substituted benzene will allow substituted pyridines to be located.

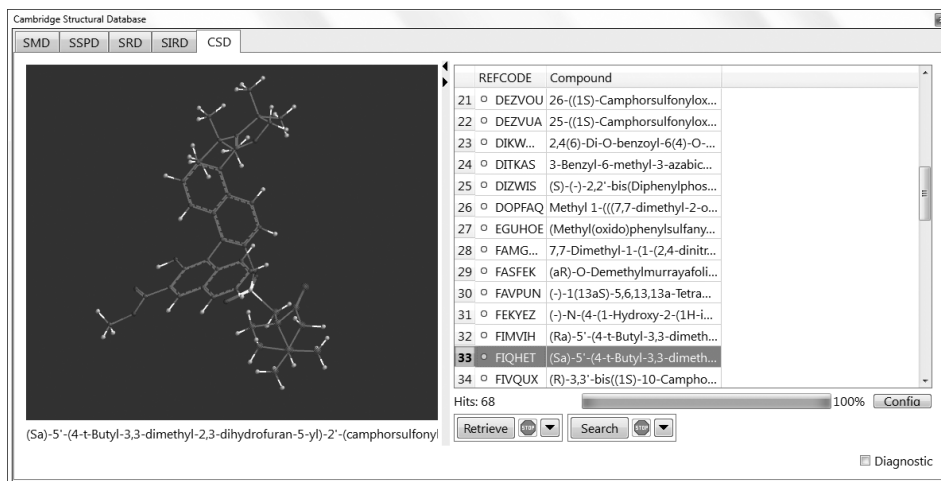
(iii) **Wild-Card Bonds**

A bond may be designated as a wild-card (meaning that bond type is unimportant) by first selecting **Define Structure Query** from the **Search** menu and then *clicking* on it. An orange cylinder will be drawn around it. *Clicking* on this cylinder removes the wild-card designation and the cylinder. This is important in certain heterocycles where bond typing may be ambiguous.

To limit a CSD search either to organic molecules (about half the total collection) or to molecules containing transition metals, lanthanides or actinides, bring up the **Search Options** dialog by *clicking* on  to the right of the **Search** button at the bottom right of the dialog, and select **Organics** or **Inorganics** respectively.




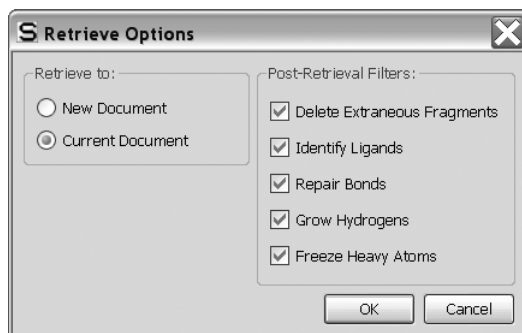
After a substructure has been specified, a search of CSD is initiated by *clicking* on the **Search** button at the bottom right of the dialog. If there are hits to be found, these will appear in the scroll box at the right of the dialog. The search can be terminated at any time by *clicking* on  to the right of the **Search** button at the bottom of the dialog. Once completed (or stopped), the CSD reference codes (**REFCODE**) and the names for the hits are displayed in the box and the number of hits is indicated immediately below. Structural data are available only for entries preceded by a filled yellow circle, although experimental references are available for all entries (see **Output** under the **Display** menu; **Chapter 18**). A ball-and-wire model for any entry for which a structure is available may be displayed in the window at the left of the dialog by *clicking* on its name in the box.



The model can be manipulated (rotated, translated, scaled) inside the window with the usual mouse/keyboard commands (you need to position the cursor inside the window).

The structure is *exactly as it appears in CSD*. In particular, it may contain two or more different molecules, two or more different conformations of the same molecule and/or extraneous molecules (most commonly solvent molecules and counterions). Hydrogens will be present only where they have been assigned in the experimental structure. The structure may also contain errors in bonding due to uncertainties in the original data. Modifications to the CSD entry to

address these issues may be made prior to transfer to the file system using the **Retrieve Options** dialog. This is accessed by *clicking* on  to the right of the **Retrieve** button at the bottom of the **CSD** dialog.



(i) **Delete Extraneous Fragments**

If *checked*, deletes any detached molecules that do not contain the original substructure.

(ii) **Identify Ligands**

Identifies ligands in structures that are returned from CSD (primarily for organometallic compounds) and replace whatever bonding is supplied by a representation in which a dashed line is drawn between the geometric center of the ligand and the atom to which it connects.

(iii) **Repair Bonds**

If *checked*, attempts to fix bonding errors based on the actual geometry in CSD and on normal valence rules.

(iv) **Grow Hydrogens**

If *checked*, adds hydrogens wherever they appear to be missing, based on the actual geometry and on normal valence rules. This function is also available under **Grow Hydrogens** in the **Molecule Utilities** dialog (**Chapter 18**).


(v) **Freeze Heavy Atoms**


If *checked*, freezes all atoms except for hydrogens. This allows experimental X-ray hydrogen positions that are often poorly determined, to be refined using molecular mechanics or quantum chemical methods, while maintaining the heavy

atom skeleton.\*

(vi) **Retrieve to:**

Selection of **New Document** means that the search results are to be stored in a new document and leads to a request for a document name, whereas selection of **Current Document** results in the retrieved molecules being appended to the end of the current document.

The **Retrieve Options** dialog is dismissed with all changes to existing settings maintained by *clicking* on **OK**. *Clicking* on **Cancel** or on  also dismissed the dialog but any changes to settings are lost.

Retrieval is accomplished by first selecting (*clicking* on) one or more hits in the scroll box and then pressing the **Retrieve** button at the bottom of the **CSD** dialog. The **Shift** and **Ctrl** keys may be used in the usual way to specify retrieval of multiple hits. Retrieval can be stopped at any time by *clicking* on  to the right of the **Retrieve** button at the bottom of the **CSD** dialog.

### **Spartan Molecular Database (SMD) and Spartan Spectra and Properties Database (SSPD)**

*Spartan* connects to two different collections of molecules for which quantum chemical calculations have already been performed and the results of these calculations are available. The first is the Spartan Molecular Database (SMD). This comprises more than 500,000 records of calculations performed using one or more of the following quantum chemical models: Hartree-Fock with 3-21G, 6-31G\*, 6-311+G\*\* basis sets, B3LYP and EDF1 density functional and MP2 with 6-31G\* and 6-311+G\*\* basis sets and G3(MP2). The LACVP pseudopotential is used in place of the 6-31G\* basis set for main-group elements heavier than krypton and for second and third-row transition metals. Thus, a B3LYP/6-31G\* calculation on benzene tungsten tricarbonyl ( $C_6H_6WCO_3$ ) uses the LACVP

---

\* Hydrogen positions may be refined using molecular mechanics using **Minimize** under the **Build** menu; **Chapter 16**. Refinement using quantum chemical models is specified in the **Calculations** dialog (**Calculations...** under the **Setup** menu; **Chapter 17**) by way of **Frozen Atoms** (**Freeze Center** under the **Geometry** menu; **Chapter 15**).

pseudopotential for tungsten and the 6-31G\* basis set for hydrogen, carbon and oxygen. Except for the G3(MP2) entries, the geometry of each molecule has been optimized, assuming the best (lowest-energy) conformer obtained from the MMFF molecular mechanics model. G2(MP2) entries use MP2/6-31G\* equilibrium geometries. Included in the SMD entry is the geometry, the energy, the energy corrected empirically for aqueous solvation, the dipole moment and Mulliken, natural and electrostatic atomic charges. Properties related to geometry (area, volume and polar surface area) are calculated on-the-fly from the geometry.

The following recipe has been used to create the individual entries in SMD:

1. Identify the lowest-energy conformer using the MMFF molecular mechanics model. Molecules with fewer than 200-300 conformers are searched systematically (guaranteeing that the global minimum will be found), while Monte-Carlo sampling is used for molecules with more conformers.
2. Starting from the best conformer, perform a geometry optimization using the appropriate quantum chemical model. This implies that all quantum chemical models for a particular molecule will refer to the same conformer.

The Spartan Spectra and Properties Database (SSPD) presently comprises approximately 60,000 molecules. This number is expected to double by the end of 2011. SSPD differs from SMD in several important ways. The first difference is that all calculations have been obtained from a single theoretical model, specifically the EDF2/6-31G\* density functional model. This represents a credible attempt to get around the commonly-asked and difficult to answer question “which model is best to use?”. The EDF2/6-31G\* density functional model consistently yields reliable structures, reaction energies, spectra and properties, yet is simple enough to be routinely applied to reasonable-size molecules, that is, up to 300-500 amu.

The second difference is that, *wherever possible*, SSPD entries use the best conformation assigned from the T1 thermochemical recipe rather than from MMFF molecular mechanics. T1 has been shown

to reliably account for conformational energy differences in small molecules (where experimental data are available), which in turn suggests that it should be able to reliably identify the lowest-energy conformer. Note, however, that T1 is presently defined only for uncharged, closed-shell molecules incorporating H, C, N, O, F, Si, P, S, Cl and Br. The conformations of molecules outside of these limits (perhaps 2% of the entries in SSPD) fall back to use of the MMFF method for assigning the lowest-energy conformer.

The following recipe has been used to create the individual entries in SSPD:

1. Identify a set of up to 100 diverse conformers using the PM3 semi-empirical model. This is done with the so-called restricted systematic search algorithm described under **Calculations** in the **Setup** menu (**Chapter 17**).
2. Perform T1 calculations on all conformers and select the one with the lowest (T1) heat of formation.
3. Starting from this conformer, perform a geometry optimization using the EDF2/6-31G\* model.

A simpler and automated procedure for providing T1 heats of formation for conformationally flexible molecules has previously been described in **Chapter 17**.

The third difference is that the vast majority of SSPD entries include both the infrared and (proton and  $^{13}\text{C}$ ) NMR spectra. (The future plan is for some entries to also include the Raman spectrum.) The vibrational motions associated with selected lines in the calculated infrared spectrum can be animated. In addition, COSY, HSQC and HMBC NMR spectra (2D spectra) can be generated on-the-fly from the information in SSPD.

The fourth difference is that, in addition to the properties available from SMD, each SSPD entry includes the T1 heat of formation and information useful for QSAR type analyses, specifically, information based on electrostatic potential and local ionization potential maps.

The fifth difference is that each SSPD entry includes the wavefunction, allowing on-the-fly generation of any graphical surface and property

maps available to *Spartan*.



Because of inclusion of the wavefunction (and to a lesser extent inclusion of the vibrational motions), SSPD has a significantly larger disk footprint than SMD (3Gb vs. 1.5 Gb). However, the in-memory footprint of SSPD is actually smaller than that for SMD.

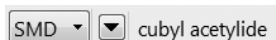
SMD and SSPD are configured separately (**Paths** from **Preferences** under the **Options** menu; **Chapter 20**), and can either be used independently or in concert. SMD facilitates comparisons among different models, whereas SSPD provides simple access to a wide selection of data from a single high quality model.


### Accessing and Processing Data From SMD and SSPD

There are two ways to access data from SMD and SSPD. The first way is to replace a structure that has been built (or read in) with a matching database entry. (Note that it is possible to replace all structures in list of molecules for which there are matches in SMD or SSPD in a single operation.) The second way is to search SMD or SSPD for entries that share particular structural features.


### Molecule Replacement

For the purpose of molecule replacement, SMD and SSPD are treated as independent databases. This gives the user the choice between the diversity of theoretical models represented in SMD and the simplicity brought on by the single model in SSPD. Selection is made from the menu to the left of the  at the bottom of *Spartan*'s main screen. Once selected, the existence of the selected molecule in SMD *or* in SSPD is evidenced by its name appearing to the right of .




SMD  cubyl acetylide

Both collections may be placed under either the SMD or SSPD heading. See **Paths** from **Preferences** under the **Options** menu; **Chapter 20**.

The structure of any SMD or SSPD entry can be seen by *clicking* on . This brings up a dialog with a viewing screen on the left and a



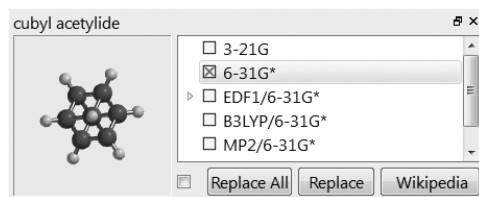
list of quantum chemical models for which database entries exist on the right. The ▲ changes to a ▼.



By default, the SMD entry of a 6-31G\* calculation (if available) is checked and a ball-and spoke model of its structure appears in the viewing screen. A different SMD entry (corresponding to a different quantum chemical model) may be selected by *checking* the box to the left of the name of the model in the dialog. Note, that there is only one entry corresponding to the EDF2/6-31G\* model for SSPD. The model in the viewing window may be manipulated using the usual mouse/keyboard commands (you need to position the cursor inside the viewing area). Model style may not be changed.

The selected model may be replaced by the selected database entry by *clicking* on **Replace** at the bottom of the dialog. (Replacement can be undone by selecting **Undo** from the **Edit** menu; **Chapter 13**.) Replacement destroys all information except the molecule identifier (leftmost column of the spreadsheet (see **Spreadsheet** under the **Display** menu; **Chapter 18**). If desired, the molecule identifier can be changed to the molecule name (from the database) by first bringing up the spreadsheet, *clicking* on **Label** (the header cell corresponding to the leftmost column), *right clicking* and selecting **Rename Selected from SMD** or **Rename Selected from SSPD** from the menu that results.

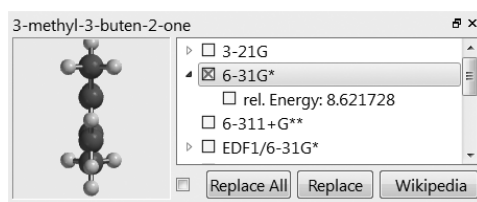Some SMD and SSPD entries also include structures for higher-energy conformers (in addition to that for the lowest-energy conformer). These entries are designated by a ▶ to the left of the check box in the preview dialog, *clicking* on which leads to an expanded entry. This accesses a higher-energy conformer and provides their energies (in kJ/mol) relative to the best conformer.



Entries are selected as before, after which they may be examined in the viewing screen and used to replace the on-screen molecule.

In the event that the selected (on-screen) molecule is part of a multi-molecule document, it is possible to replace all molecules in the document for which database entries for the specified theoretical model are available.\* This is accomplished by *clicking* on **Replace All** instead of **Replace**. A warning message is provided prior to the replacement. Note that the lowest-energy conformer for all other (non-selected) molecules in the document is used, even though a higher-energy conformer may be indicated for the selected molecule. (Replacement of document entries with non-lowest energy conformers needs to be done one-by-one using **Replace** instead of **Replace All**.) Molecules in the document for which there are no database entries at the specified theoretical model are not affected. As with single molecule replacement, all information except the molecule identifier is replaced by that in the database. Molecule names from the selected database for all molecules in a document may be substituted for the original molecule identifiers by *clicking* on **Label** (the header cell of the leftmost column in the spreadsheet) *right clicking* and selecting **Rename Selected from SMD** or **Rename Selected from SSPD** from the menu that results.

The **Wikipedia** button at the bottom right of the database preview dialog brings up a Wikipedia page based on the name of the selected molecule.

### Substructure and Other Modes of Searching SMD and SSPD

Access to data from SMD and SSPD follows selecting **Open Databases**, and then *clicking* on either the **SMD** or **SSPD** tab (at the

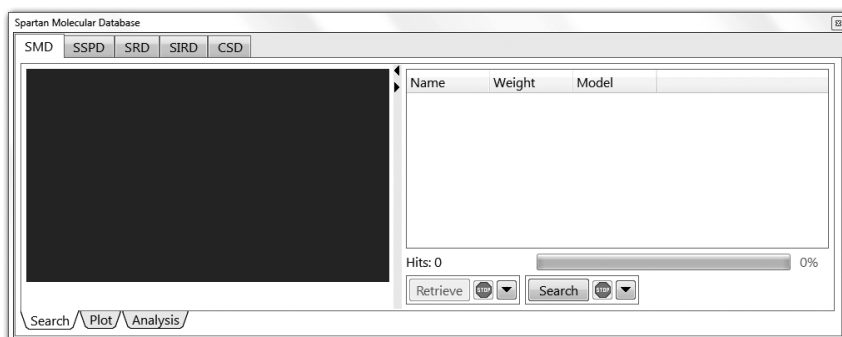
---

\* For 32 bit implementation, the practical replacement limit is 1,500-2,000 molecules for SMD and 500-800 for SSPD.

top left of the resulting window). (Aside from the title, the resulting displays are identical.) Each comprises three different panels, **Search**, **Plot** and **Analysis**, selection of which is controlled by tabs at the bottom left of the window. A search of SMD or SSPD is carried out from the **Search** panel, plotting of the search results from the **Plot** panel and regression analysis of search results from the **Analysis** panel.

Substructure searching of SMD and SSPD is similar to that already described for substructure searching of the Cambridge Structural Database (CSD). The only required input is a three-dimensional structure. While one or more growth points (see **Define Structure Query** earlier in this chapter) will normally be specified, they are optional and in their absence an exact search will be carried out. Wild-card atoms and bonds used in searches of CSD are not supported. One or more substituents may be attached to the structure (see **Substituent Model Kit** under the **Build** menu; **Chapter 16**), in effect allowing multiple starting structures. Name searching requires a name (or partial name), formula searching requires an exact formula, weight searching requires a molecular weight and isomer searching requires a three-dimensional structure without any growth points.


*Clicking* on the **Search** tab leads to the **Search** dialog.



This contains a box on the right for listing the search hits and providing values for any requested properties (at the outset only the molecule name is provided), and a window on the left for viewing and manipulating the individual hits. Controls along the bottom right access the **Search Options** dialog which specifies the databases to be searched, the kind of search, and (for SMD only) the theoretical model(s) to be searched. They are also used to initiate and terminate

a search and to designate the location for any results that are retrieved from the search.

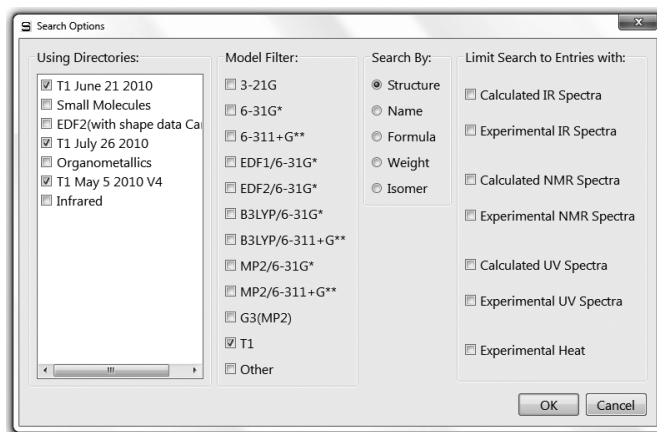With the SMD tab selected, *clicking* on  to the right of the **Search** button leads to the **Search Options** dialog.



**Using Directories** lists the collections that make up SMD. **SMD (sample)** is the sample of ~5000 molecules supplied with all copies of *Spartan* while **SMD** is the full collection. If **SMD** is installed, **SMD (sample)** should be removed (**Paths** from **Preferences** under the **Options** menu; **Chapter 20**). **Method Filters** specifies the theoretical models to be included in the search. **Other** accesses data with models that are not explicitly listed. **Search By** designates the type of search: **Structure**, **Name**, **Formula**, **Weight** or **Isomer**. **Limit Search to Entries with** designates information that must be in the hit.

(i) **Name**

Selection followed by *clicking* on **OK** in the **Search Options** dialog leads to a box **By Name** at the bottom right of the dialog. The search will return all entries that contain the string that is entered into the box. For example, *typing* toluene will not only result in toluene, but also in molecules such as *para*-toluenesulfonic acid and 4-chloro-2-fluorotoluene.

(ii) **Formula**

Selection followed by *clicking* on **OK** in the **Search Options** dialog leads to a box **By Formula** at the bottom right of


the dialog. A formula of the form  $C_xH_yHet_z$ , where “Het” corresponds to elements other than carbon and hydrogen (each element needs to be specified explicitly).  $x$ ,  $y$ ,  $z$  are integers and are required. For example, the formula for nitrobenzene needs to be written as **C6H5N1O2**. The order of atoms is irrelevant. The search will return all entries with this formula. For example, the formula **C6H5N1O2** will return not only nitrobenzene but also molecules such as phenyl nitrite and pyridine-4-carboxylic acid.

(iii) **Weight**

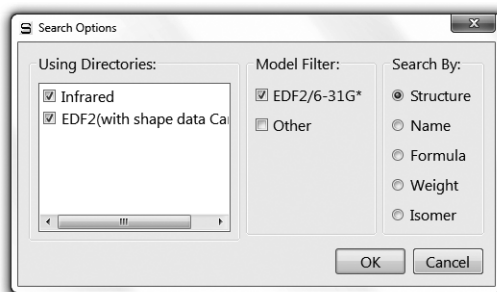
Selection followed by *clicking* on **OK** in the **Search Options** dialog leads to a box **By Weight** at the bottom right of the dialog. The search will contain all entries with the specified molecular weight  $\pm 0.05$  amu.


(iv) **Isomer**

Selection followed by *clicking* on **OK** in the **Search Options** dialog specifies that all isomers of the selected molecule will be returned. The result is the same as specification of a molecular formula.

*Clicking* on **OK** removes the **Search Options** dialog with all settings maintained. *Clicking* on **Cancel** or  removes the dialog without maintaining the settings.

The **Search Options** dialog associated with SSPD is simpler in that the **Methods Filters** and **Limit Search** fields are missing. Only a single theoretical model (the EDF2/6-31G\* density functional model) is available and most entries contain calculated IR and NMR spectra (and none contain UV spectra). **Using Directories** lists the collections that make up SSPD. **SSPD (sample)** is the sample of ~5000 molecules supplied with all copies of Spartan while **SSPD** is the complete collection. If **SSPD** is installed, then **SSPD (sample)** should be removed (**Paths** from **Preferences** under the **Options** menu; **Chapter 20**). The same search methods available for SMD are available here.




A search is carried out by *clicking* on the **Search** button at the bottom right of the search dialog, following which hits will appear at the right. The search may be terminated at any time by *clicking* on  to the right of the **Search** button. Any hits that have already been identified will be kept. A ball-and spoke model of a hit may be examined at any time by *clicking* on its name in the box. It may be manipulated in usual way, although the model style cannot be changed.


Properties may be added to the hit list by *dragging* them (names not numerical values) from one of the properties dialogs (see **Properties** under the **Display** menu; **Chapter 18**), the formula editor (**Formulas** under the **Display** menu; **Chapter 18**) or from the main *Spartan* screen, primarily entries associated with the **Geometry** menu; **Chapter 15**. The hit list can be sorted using the property, by *right clicking* on the column header cell and selecting **Sort** from the menu.


The results of a search on SMD or SSPD as well as any subsequent plotting and analyses are saved in the *Spartan* document corresponding to the original structure query or to the blank document that was used for a name or spectra query. This behavior is different from previous versions of *Spartan*.


### **Plotting Data Resulting from an SMD or SSPD Search**

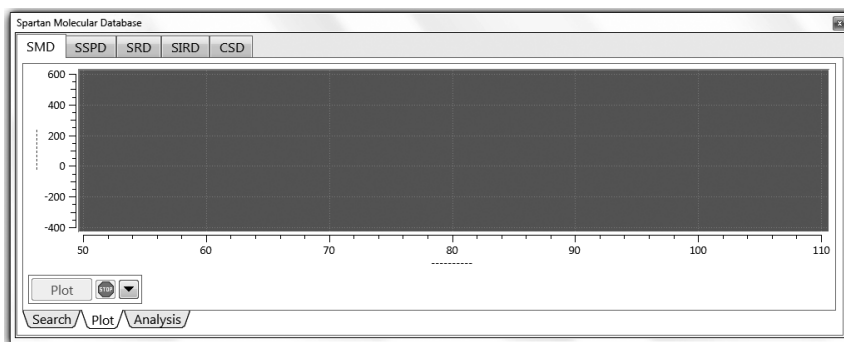
Data resulting from a search on SMD or SSPD may be displayed either as a histogram (property value vs. number of occurrences) or as a 2D scatter plot (property value vs. property value). Access to plotting is via the **Plot** tab at the bottom left of the dialog.



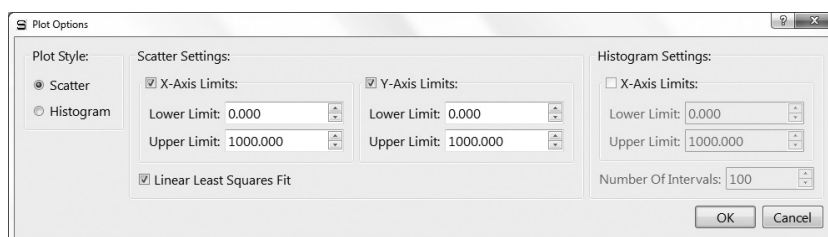Selection of plot type and variable ranges is from the **Plot Options** dialog accessible from  to the right of the **Plot** button.





Default plot ranges are from minimum to maximum property values. *Checking **Linear Least Squares Fit (Scatter Settings)*** fits a straight line to the data. ***Number of Intervals (Histogram Settings)*** sets the number of bars in the histogram. *Clicking on **OK*** removes the dialog with all settings maintained. *Clicking on **Cancel*** or  removes the dialog without maintaining the settings.

A plot is drawn by *clicking* on the **Plot** button at the bottom left of the dialog. The points will not appear instantaneously as the data need to be extracted from disk. Plotting may be terminated at any time by *clicking* on  to the right of the **Plot** button. If a least-squares line has been requested, it will be based on the accumulated data. At anytime, the name of the molecule associated with a point on a scatter plot can be displayed by positioning the cursor over that point. The point can be deleted by *right clicking* and selecting **Delete** from the resulting menu. The plot will rescale as appropriate.

## Regression Analysis of Data Resulting from an SMD or SSPD Search

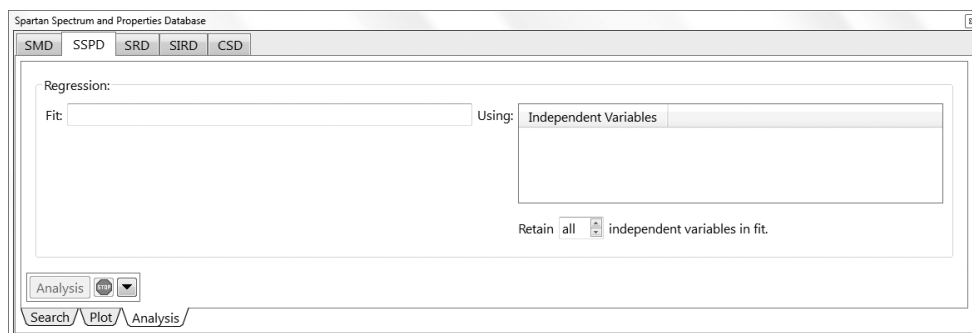A regression analysis of data resulting from a search may be carried out by selecting the **Statistics** tab at the bottom left of the dialog.




The name (not the numerical value) of the property to be fit is *dragged* from one of the properties dialogs, the formula editor or from the main *Spartan* screen, into the box to the right of **Fit**. The names of one or more properties to be used in the fit are *dragged* into the box to the right of **Using**. The fitting procedure (singular value decomposition) will use only variables that have been selected (*clicked* on) from among the list, up to a maximum number of independent variables (the most important in terms of their contribution to the fit from among this list) specified in the menu at the bottom right of the dialog. Note that the number of variables actually used in the fit may be smaller than the maximum number selected (from the menu at the lower right of the dialog), as variables that are linearly dependent and/or do not contribute significantly to the fit will be discarded. The analysis is started by *clicking* on the **Analysis** button at the bottom left of the dialog. As with plotting, this does not occur instantaneously as the underlying data needs to be extracted from the disk. Upon completion, the fitting function will be displayed at the lower left of the dialog, and can be copied and then pasted in the **Plot** dialog. Analysis may be terminated at any time by *clicking* on  to the right of the **Analysis** button, following which results from a fit to a portion of the data will be displayed.



## Relationships Involving Chemical Reactions

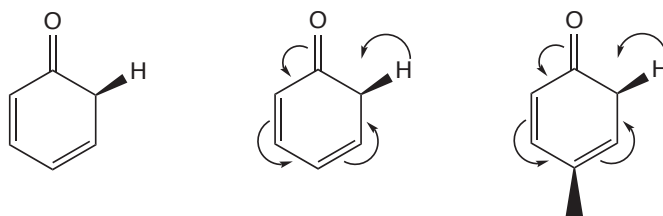We have already discussed how the formula editor may be used to construct relationships between among properties for a single molecule (see **Formulas** under the **Display** menu; **Chapter 18**). The formula editor can also be employed to construct relationships between the reactants and the products of a chemical reaction of the reaction, because the products of a reaction are completely and unambiguously defined given the structure of the reactant (or reactants) and a set of curly arrows (see **Define Reaction Query** earlier in this chapter). Note that the same information can also be used to guess a transition state for the reaction (see **Guess Transition State** later in this chapter).

For example, consider the difference in energies of 2,4-cyclohexadieneone and phenol. Whereas the normal preference is for keto forms be favored over enols, here the reverse is certainly true as phenol benefits from aromatic stabilization. Three steps are needed to use SMD or SSPD to ask to what extent substituents alter this preference. First, a structure query (the reactant) must be built, arrows placed to designate the product of reaction and one or more growth points added to define the limits of the search. For example, use 2,4-cyclohexadieneone and limit substituents to the 4 position.




Next, a formula for the energy of reaction needs to be constructed. Inclusion of arrows extends the formula from the query itself (a 4-substituted 2,4-cyclohexadieneone), to the product of reaction (a *para*-substituted phenol) and to the difference between product and reactant. Selection of which is to be employed in the formula follows from the **Applied to** menu under **Context** in the formula editor (highlighted only if arrows are present). In this case, the choices are (R) 2,4-cyclohexadieneone, (P1) phenol or (P-R) reaction. The formula for the reaction energy, `@energy*@hart2kj`, appears in the

box at the center of the dialog following selection of **(P-R) reaction** from the **Applied to** menu and *clicking* on **E** under **Energies**. The default name may be changed by editing the text inside the box to the left of the formula. *Clicking* on **Enter** enters the (renamed) formula **@rea(@energy\*@hart2kj)** into the list at the top of the formula editor.

A search may now be carried out. After it has completed, dragging the formula for the reaction energy onto an empty header cell of the hit list will result in reaction energies for all entries to appear.


As a second example, consider the energy for addition of hydrogen fluoride to propene, leading to 1-fluoropropane. Here there are two reactants (hydrogen fluoride and propene), but only the latter will be provided growth points. Thus, the query is for information about the energy of HF addition as a function of substitution on the alkene.

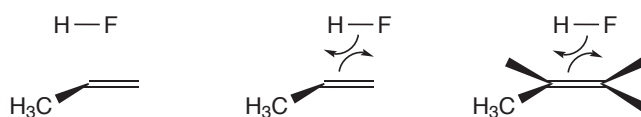The query starts from an intramolecular complex of propene and hydrogen fluoride, oriented such that the C=C and H-F bonds parallel with fluorine on the unsubstituted side of the alkene. Two arrows, one from the CC bond to the space between the unsubstituted carbon and fluorine and the other from the HF bond to the space between the substituted carbon and hydrogen, define the product (1-fluoropropane). Growth points attached to all three open valences extend the search to all available substituted propenes (for which analogous substituted 1-fluoropropanes are also available).




A formula for the energy of reaction may be constructed using the same sequence as in the previous example, selecting **(P-R) reaction** from the **Applied to** menu under **Context** in the formula editor, *clicking* on **E**, changing the default name and *clicking* on **Enter**.

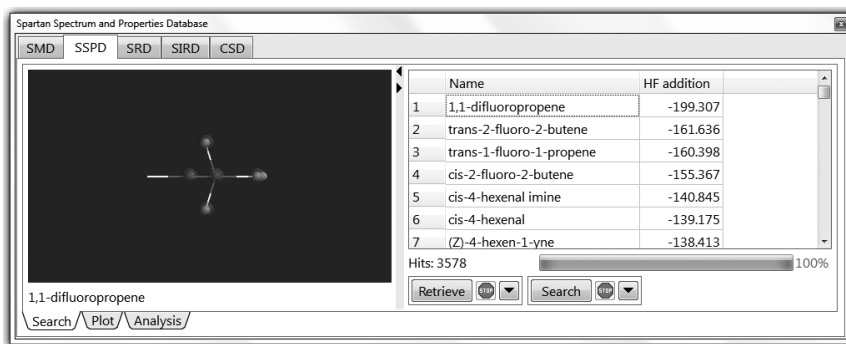
Following a search using the query, dragging the formula for the reaction energy onto an empty header cell of the hit list adds a column of reaction energies to the list.




It is possible to apply different queries to the same set of search results, for example, to ask about the energy of addition of HCl (instead of HF) to a substituted propene. Whereas the query changes (the reaction of HCl and propene rather than of HF and propene), the underlying search (on propene) does not. The result is two sets of energies. (Results designated *No DB Entry* indicate that the required information is not available.)

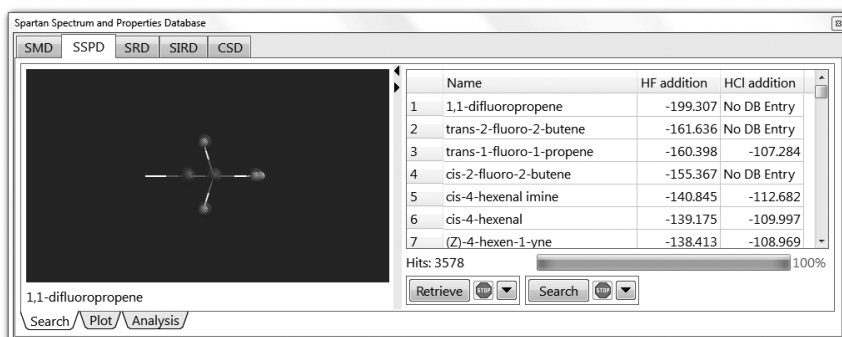


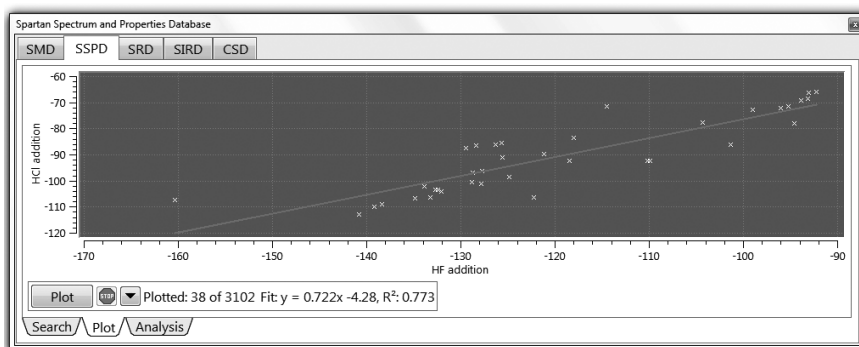
A plot of HCl vs. HF reaction energies shows the extent to which they are actually correlated.




## Retrieving Data from an SMD or SSPD Search

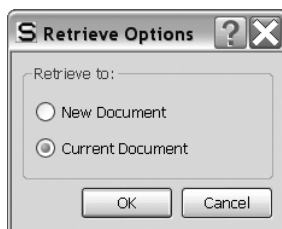Prior to retrieval, it is necessary to specify a destination from the **Retrieve Options** dialog. *Click* on  to the right of the **Retrieve** button.



Selection of **New Document** results in request for a file name; selection of **Current Document** results in the retrieved molecules being appended to the end of the current document. The **Retrieve Options** dialog is dismissed by *clicking* on **OK**. *Clicking* on **Cancel** or on  also dismisses the dialog but any changes to settings are lost.

To retrieve, select (*click* on) one or more hits in the box and *press* the **Retrieve** button. The **Shift** and **Ctrl** keys may be used to specify retrieval of multiple hits. Retrieval may be stopped at any time by *clicking* on  to the right of the **Retrieve** button.


## Spartan Reaction Databases (SRD)

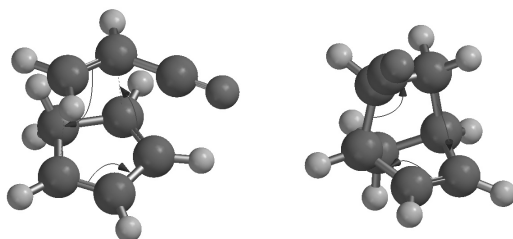The Spartan Reaction Database (SRD) is a collection of transition states for approximately 1,500 organic and organometallic reactions obtained from one or more of the following theoretical models: AM1, PM3, HF/3-21G and HF/6-31G\*. Each entry contains the infrared spectrum, which should have a single imaginary frequency corresponding to the reaction coordinate.


An entry in SRD is in terms of a three-dimensional structure (corresponding to the reactants or the product), together with an appropriate set of curly arrows (connecting the reactants to the product or vice versa) and designation of theoretical model.

Instructions for drawing reaction arrows has been provided earlier in this chapter (**Define Reaction Query**).

For example, an entry for the Diels-Alder reaction of cyclopentadiene and acrylonitrile is given on the left, while that for the corresponding retro Diels-Alder is given on the right.



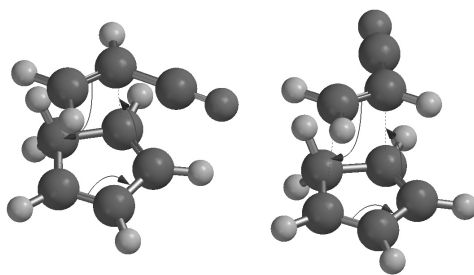
Note, that stereoisomers are distinguished by the way the reactants are positioned. For example, reactions leading to *endo* and *exo* Diels-Alder adducts are distinct.



SRD is open-ended and users may construct their own data collections with additional reactions and/or results from additional theoretical models on existing reactions (**Appendix G**).


## Substructure Searching

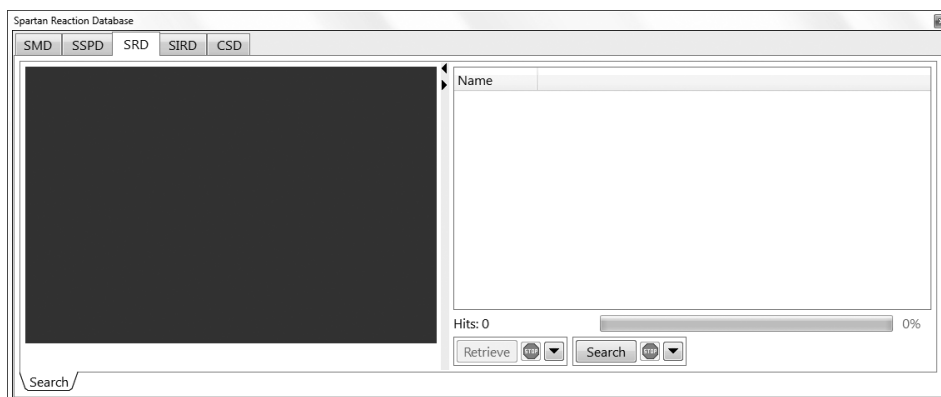Substructure searching from SRD is very similar to that from SMD. Access to SRD is provided by *clicking* on the **SRD** tab at the top of the **Database** dialog. This gives rise to the **SRD** dialog.




The **SRD** dialog contains a window on the left for previewing structures and a box on the right for listing hits on SRD, as well as buttons at the bottom to control various aspects of the search and subsequent transfer of data into *Spartan*'s file system.

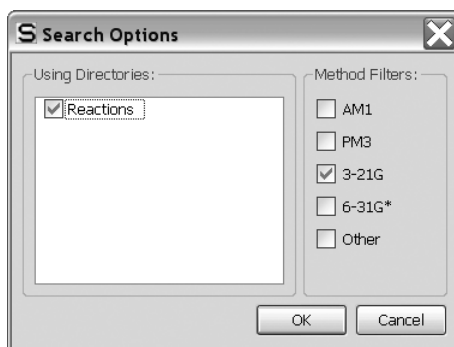A substructure search is set up by specifying one or more attachment points on the reactant structure (see **Define Structure Query** earlier in this chapter) that has been augmented with reaction arrows (see **Define Reaction Query** earlier in this chapter). Prior to starting the actual search, it is necessary to specify which theoretical model or set of models are to be included. This is done from the **Search Options** dialog, and is reached by *clicking* on  to the right of the **Search** button at the bottom right of the dialog.




**Using Directories** indicates the libraries to be searched. One is supplied with *Spartan* and the users are free to replace or supplement this with their own collections. *Checking* one or more boxes under the **Method Filters** specifies the theoretical model(s) to be searched. *Clicking* on **OK** removes the dialog. (*Clicking* on **Cancel** or  also removes the dialog, but selections are not made.)

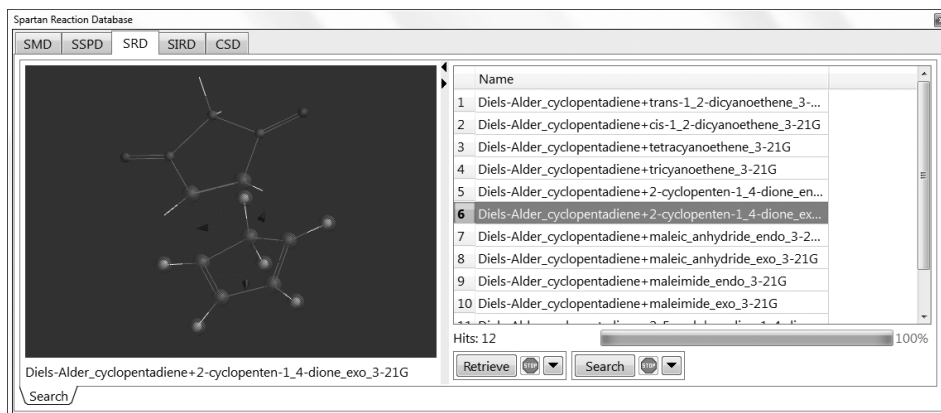
Following specification of structure/reaction arrows, attachment points and designation of theory level(s), a search on SRD is carried out by *clicking* on the **Search** button at the bottom right of the dialog. If there are hits, they will begin to appear in the box at the right of the dialog after a few seconds. The search may be terminated at any time by *clicking* on  to the right of the **Search** button at the bottom right of the dialog.





Database hits are displayed in the box at the right of the dialog. They may be sorted by name, molecular weight and theoretical model by *clicking* on the appropriate tab. A 3D structure for a particular hit may be previewed by *clicking* on its name in the box. A ball-and-wire model of the transition state with bonding consistent with the reactant and associated reaction arrows is displayed in the window at the left of the dialog, and its name appears immediately below the window. The model may be manipulated inside the window in the usual way, but the model style cannot be changed.

## Retrieving Data from an SRD Search

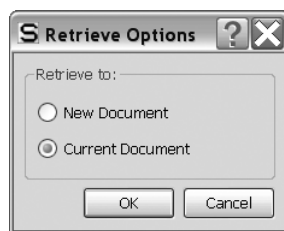Prior to retrieval of one or more database hits, it is necessary to specify a destination. This is accomplished with the **Retrieve Options** dialog, brought up by *clicking* on  to the right of the **Retrieve** button.



Selection of **New Document** results in request for a file name; selection of **Current Document** results in the retrieved molecules being appended to the end of the current document. The **Retrieve Options** dialog is dismissed by *clicking* on **OK**. *Clicking* on **Cancel** or on  also dismisses the dialog but any changes to settings are lost.

Retrieval is accomplished by first selecting (*clicking* on) one or more hits in the box and then pressing the **Retrieve** button at the bottom of the **SMD** dialog. The **Shift** and **Ctrl** keys may be used to specify retrieval of multiple hits. Retrieval may be stopped at any time by *clicking* on  to the right of the **Retrieve** button.

## Spartan Infrared Database (SIRD)

The **Spartan Infrared Database (SIRD)** is actually a second access point into the **Spartan Spectra and Properties Database (SSPD)** discussed earlier in this chapter. It comprises infrared spectra for ~60,000 molecules\* calculated from the EDF2/6-31G\* model and (where the T1 model is defined) based on the lowest-energy conformer from T1. Unlike the other databases of calculated information attached to *Spartan*, SIRD is searched not on substructure but rather on a measured (or calculated) infrared spectrum. This is input in JCAMP format (.dx).\*\* The search is carried out in two steps: a fast search on the full database is performed first to identify a small set of possible matches which are then examined more thoroughly.

### Spectra Searching

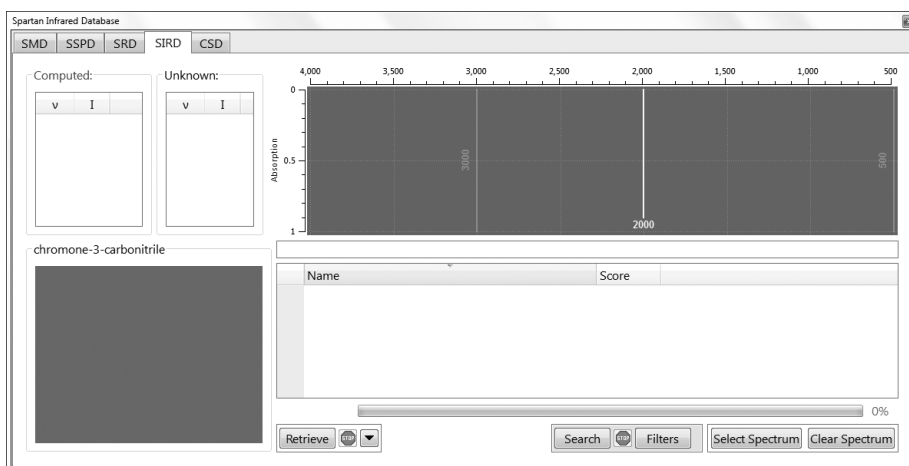
Access to SIRD is provided by *clicking* on the **SIRD** tab at the top of the **Database** dialog. A *Spartan* document must be open even if it is an empty document (**New** from the **File** menu). This gives rise to the **SIRD** dialog.

---

\* As with SMD, a small subset of ~5000 molecules is provided with *Spartan*. The full SSPD database may be separately licensed.


\*\* A sample .dx file is provided in **Appendix I**.





While this dialog is quite unlike those associated with other databases accessible from *Spartan*, it does, share two aspects with these dialogs, a window at the lower left for examining the 3D structure of the selected hit and a box at the lower right for listing all hits. It also contains a window at the top right for examining the match between the unknown (measured) infrared spectrum and that of a hit in the database, as well as boxes at the top left listing the frequencies and intensities for both unknown and database spectra. Finally, buttons at the bottom of the dialog access the unknown spectrum, provide (optional) filters for its search against spectra in the database and control transfer of any search results into *Spartan*'s file system.

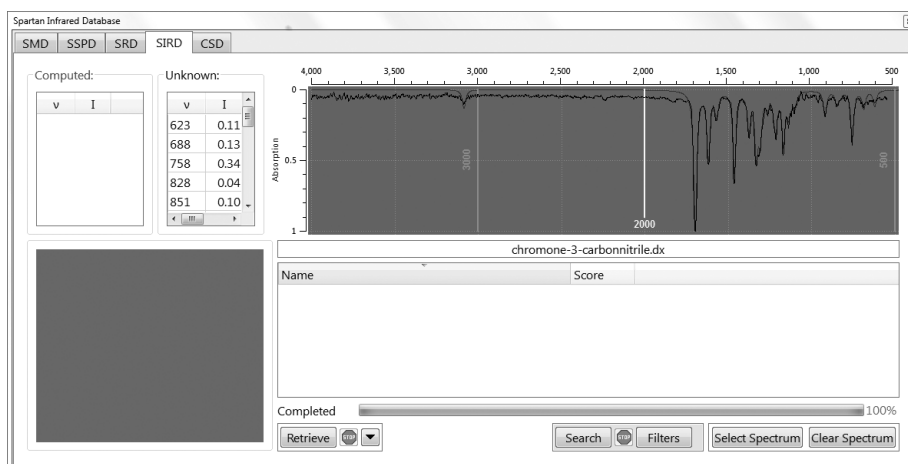
A search is initiated by *clicking* on **Select Spectrum** at the bottom right of the dialog. This leads to a file browser that has been set up to recognize only .dx files. Selection of a file in the usual manner followed by *clicking* on **Open** (or *double clicking* on the name of the .dx file) leads to the display of the unknown spectrum in the window at the top right of the dialog.



The name of the (.dx) file is provided underneath the window. The unknown spectrum has been fit on-the-fly to a Lorentzian function allowing extraction of a set of vibrational frequencies. These are displayed in the box marked **Unknown:** immediately to the left of the spectrum, along with the corresponding intensities. Features in the unknown spectrum with very low intensity cannot be fit.


## Filters

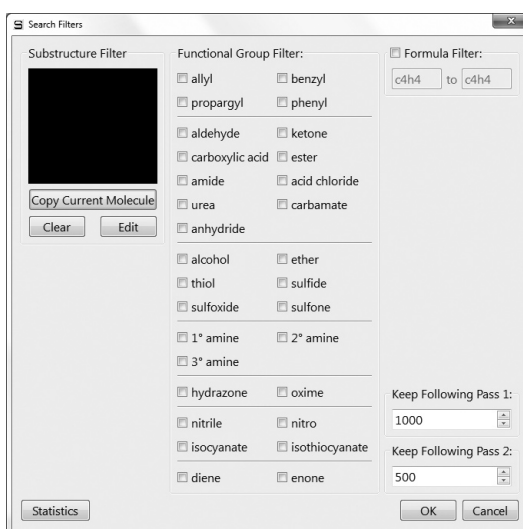Conditions under which the search is carried out may be specified by *clicking* on the **Filters** button at the bottom of the dialog. (This step is optional.)



The **Filters** dialog contains a window for constructing a substructure filter on the left, a listing of common functional groups in the middle, and controls to limit the molecular formula of hits on the database on the top right. Controls at the bottom right govern the number of preliminary hits that are moved from the first pass of the search to the (more costly) second pass and the number of hits that are finally kept.

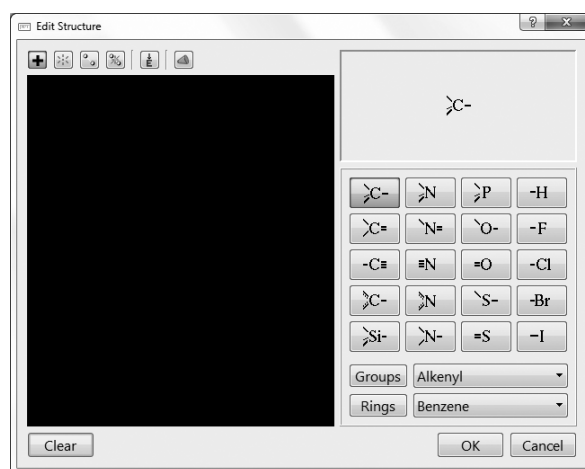This allows a substructure to be constructed (or edited from the currently selected molecule) and used as a substructure filter. To build a substructure, *click* on **Edit** below the (empty) **Substructure Filter** window at the top left of the screen. Alternatively, you can start with the currently selected molecule by *clicking* on **Copy Current Molecule** prior to *clicking* on **Edit**. A builder panel appears.



This contains most of the functionality in the organic model kit (see **Chapter 16**) and in addition, the ability to define one or more attachment points (see **Define Structure Query** earlier in this chapter). Whatever substructure is built needs to be in the hit. Note that if no attachment points are designated, the molecule that is built is the only acceptable hit. *Click* on **OK** to exit the build panel.

**Functional Group Filters** allow the search to be restricted to molecules that incorporate one or more of the functional groups. *Check* to select and *uncheck* to deselect.

**Formula Filter** allows selection of the minimum and maximum molecular formula (or both) to be considered in the search. The


format of the input is  $C_cN_nO_o \dots$ , where C, N, O, etc. are atomic symbols and c, n, o are the number of each of the designated atoms. Neither hydrogens nor any elements that are not explicitly designated are considered.

**Keep Following Pass 1** controls the number of entries that are kept following the initial pass of *Spartan*'s spectra matching algorithm. The larger the number, the more computer time will be needed, but the more likely it will be to produce a good hit (assuming that the unknown or a closely-related molecule is actually in the database).

**Keep Following Pass 2** controls the number of hits that are returned to the user. Does not affect the time required for the search.

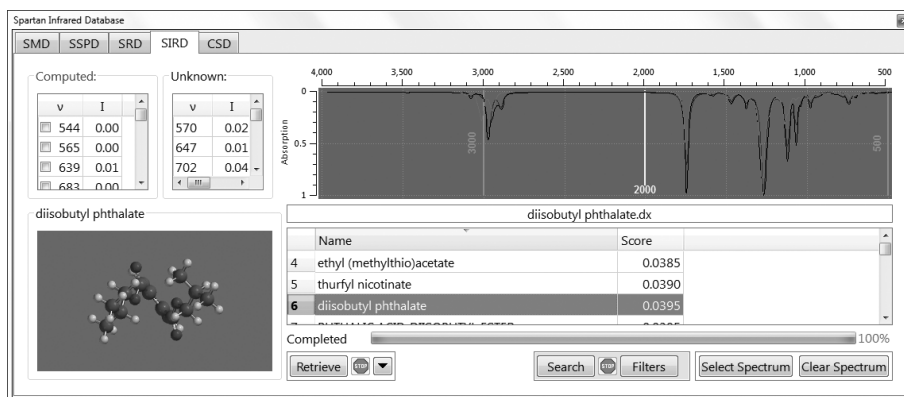
*Clicking* on **OK** exits the dialog. *Clicking* on **Cancel** or  also dismisses the dialog but any changes to the settings are lost.

In addition, the search may be restricted to a specific frequency range. This is controlled by the two grey slider bars displayed inside the dialog containing the spectrum of the unknown. Frequencies below 500-600  $\text{cm}^{-1}$  (which typically correspond to torsional motions and therefore may depend on conformation) require special instrumentation to measure and are seldom reported. While the experimental spectrum is usually recorded and reported to 4500  $\text{cm}^{-1}$ , the region beyond  $\sim 2800 \text{ cm}^{-1}$  is dominated by CH stretching vibrations and may be too crowded to be of value in spectral searching. A reasonable range is from 700-2200  $\text{cm}^{-1}$ , but this will depend on the details of the unknown spectrum.


A search on the reference spectrum restricted by whatever conditions are specified in the **Search Options** dialog is carried out by *clicking* on the **Search** button at the bottom of the **SIRD** dialog. The search is carried out in two passes. The first pass may require several tens of seconds depending on the size of the database and the filters that have been set. The second pass may require several minutes (depending on the number of hits kept following the initial pass as specified in the **Search Options** dialog). Searching may be stopped at any time by *clicking* on  to the right of the **Search** button at the bottom of the **SIRD** dialog. However no results will


be returned if the search is stopped prematurely.

When the search has completed a listing of hits (the number of which is specified by **Keep Following Pass 2** in the **Search Options** dialog) ordered by their least squares errors (lowest is best) appears in the box at the right of the **SIRD** dialog. Hits may be sorted by least squares error as well as name and molecular weight by *clicking* on the appropriate tab. As each hit is selected (*clicked* on the list), its spectrum will be superimposed onto the spectrum of the unknown.





## Retrieving Data from an SIRD Search

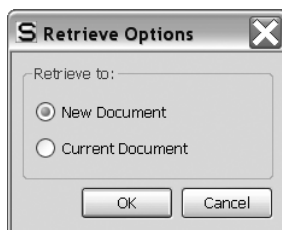
Prior to retrieval on one or more database hits, it is necessary to specify a destination. This is accomplished with the **Retrieve Options** dialog, brought up by *clicking* on  to the right of the **Retrieve** button.



Selection of destination is accomplished by *toggling* between **New Document** and **Current Document**. Selection of the former creates a new document which will need to be named later, while selection of the latter results in the retrieved molecules

being appended to the end of the current document. The dialog is dismissed by *clicking* on **OK**. *Clicking* on **Cancel** or  also dismisses the dialog but any changes are lost.

Retrieval is accomplished by first selecting (*clicking* on) one or more hits in the box and then *pressing* the **Retrieve** button at the bottom of the **SIRD** dialog. The **Shift** and **Ctrl** keys may be used to specify retrieval of multiple hits. Retrieval may be stopped at any time by *clicking* on () to the right of the **Retrieve** button.

Whichever database dialog is open may be closed by *clicking* on , or alternatively by selecting **Close Databases** (that has replaced **Open Databases**) in the **Search** menu.

### Guess Transition State ()

*Spartan* provides a facility for automatically guessing transition state geometries based on the similarity with one or more entries in *Spartan*'s reaction database\*. An exact hit replaces the guess with transition state for the same reaction with the particular theoretical model entered in the database\*\*. More commonly, *Spartan* will attempt to provide as close a match as possible with a database entry. This will generally involve a less substituted system or one in which substituents differ. Here, those parts of the structure of the transition state in the database that are common will be used, and the remaining parts will be optimized using molecular mechanics. No conformational searching is performed. Note that it may be essential for the reactants to be properly oriented to reflect the desired stereochemical outcome of the reaction.

Where a reaction is completely unknown to the database, a fallback technique similar to the linear synchronous transit method is automatically invoked.


Input to *Spartan*'s transition-state guessing procedure will be familiar to organic chemists, in that it is based on curly arrows, where each

---

\* The Spartan Reaction Database (SRD) may also be directly searched by substructure. See discussion earlier in this chapter.

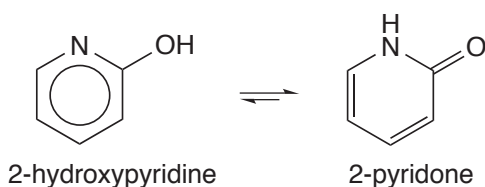
\*\* At the present time, the majority of database entries derive from very simple (semi-empirical and Hartree-Fock) quantum chemical models.

arrow identifies the movement of one electron pair. The difference is that arrows are superimposed onto a three-dimensional structure rather than a two-dimensional drawing. The procedure has already been described and illustrations provided in **Define Reaction Query** earlier in this chapter.

After all reaction arrows have been properly designated, *clicking* on  at the bottom right of the screen to replace the reactant with a guess at the transition state. In the event that the guess is unreasonable, this operation may be undone (select **Undo** from the **Edit** menu). This allows you to review your assignment of arrows and make changes as needed.


### Identify Tautomers ()

*Spartan* provides a facility to automatically identify tautomers, that is, isomers arising from labile transfer of hydrogens among heteroatoms. Tautomer identification is especially important in dealing with heterocyclic compounds where two or more different tautomers may exist in equilibrium and the identity of the dominant tautomer may not be obvious. For example, while 2-hydroxypyridine is an aromatic molecule and would be expected to be unusually stable, the non-aromatic 2-pyridone molecule is actually thermodynamically favored.





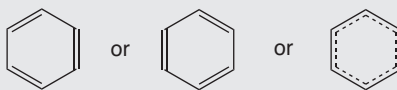This, of course, has consequences in trying to locate a particular molecule in a database where it may be represented by a different tautomer.


In contrast, *Spartan* almost always recognizes different ways to draw the same structure. The most important example is benzene, which can be represented as one of two structures with alternating double and single bonds or a single delocalized structure in which all bonds are the same and intermediate in length between single and double bonds.




All of these translate into the same molecule.

The procedure incorporated into *Spartan* is limited to tautomers involving nitrogen, oxygen, phosphorous and sulfur. Tautomers involving carbon have intentionally been excluded. Also excluded are zwitterion tautomers, for example,  $\text{H}_3\text{N}^+\text{CH}_2\text{CO}_2^-$  as a tautomer of glycine. Within these limits, all possible tautomers will be identified.

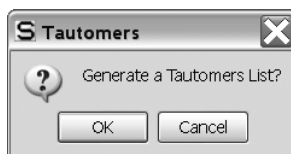The existence of tautomers is signaled by the word **Tautomer** displayed at the bottom right of the screen. If tautomers exist, individual structures may be examined or a complete list of structures generated by selecting **Tautomers** from the **Search** menu (or *clicking* on the  icon on the top of the screen). Following this, *clicking* on the step buttons at the bottom right of the screen moves through the list of tautomers. Any one of these may be used in place of the original. The full list of tautomers may be generated by *clicking* on  (at the bottom right of the screen). This leads to a dialog.



*Clicking* on **OK** dismisses the dialog and leads to a new (unnamed) document containing the full list of tautomers. (The original document is unaffected.) *Clicking* on either **Cancel** or  dismisses the dialog without generating a tautomer list.

As tautomers are generated without consideration of steric crowding, it is important to minimize each in turn (**Minimize** under the **Build** menu; **Chapter 16**) prior to using them for quantum chemical calculations.



## Extract Ligands ( )

*Spartan* provides a facility for extracting ligands from PDB files\* either as molecular structures or as footprints (pharmacophores) that these structures leave behind. The latter is separated into three parts, the first two of which follow from the molecular structure of the ligand and the third which follows from the immediate (protein or nucleotide) environment:

- (i) hydrogen-bond donor/acceptor sites and positive/negative ionizable sites
- (ii) hydrophobes
- (iii) excluded volumes

A pharmacophore can be constructed by requiring (i) to account for hydrogen-bonded and electrostatic interactions and (optionally) *either* (ii) or (iii) to account for steric interactions. Selection of (ii) leads to a pharmacophore that is visually simpler, but selection of (iii) may provide a more realistic picture.

Selection of **Extract Ligands** leads to a message at the bottom of the screen.

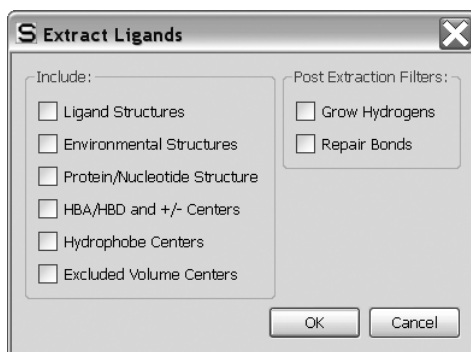
Select Ligands.

Ligands in PDB files are represented as a transparent space-filling model. *Clicking* on a ligand selects it and displays the PDB **HET name\*\*** at the bottom of the screen (if a HET name is available). Multiple ligands may be selected if desired by holding down on the **Ctrl** key. All ligands may be selected by *clicking* on the **Select All Ligands** button at the bottom right of the screen. After one or more ligands have been selected, *clicking* on the **Extract Ligands** button at the bottom right of the screen brings up the **Extract Ligands** dialog.

---

\* This capability is not restricted to files in PDB, but extends to files written in PDB format.

\*\* This is the code given to small molecules that are associated with proteins/nucleotides in PDB files.



This allows selection of whether **Ligand Structures** or one or more pharmacophore elements (**HBA/HBD and +/- Centers**, **Hydrophobe Centers** including aromatic centers, **Excluded Volume Centers**) or *both* structure and one or more pharmacophore elements are to be abstracted from the PDB file, and whether amino acid and/or nucleotide residues in close hydrogen-bonding or charge-charge contact with the ligand (**Environmental Structures**) are to be identified. It also allows the **Protein/Nucleotide Structure** (minus all ligands) to be extracted.

In the case of extraction of **Ligand Structures**, the **Extract Ligands** dialog also provides a set of utilities:

(i) **Grow Hydrogens**

If *checked*, adds hydrogens wherever they appear to be missing, based on the actual geometry and on normal valence rules. This function is also available under **Grow Hydrogens** in the **Molecule Utilities** dialog (**Properties** under the **Display** menu; **Chapter 18**).

(ii) **Repair Bonds**

If *checked*, attempts to fix bonding errors based on the actual geometry in the PDB entry and on normal valence rules.

*Clicking* on **OK** at the bottom of the **Extract Ligands** dialog extracts the ligand(s) or pharmacophore(s), placing them in a single new (unnamed) document. *Clicking* on **Cancel** dismisses the dialog without ligand (pharmacophore) extraction.




# Chapter 20
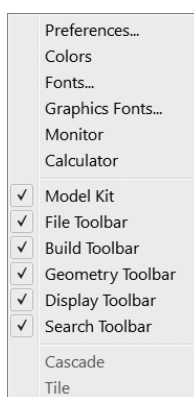
## The Options Menu

---

*Functions under the **Options** menu\* set default colors, fonts, user preferences and van der Waals radii, locate databases, identify remote servers, engage parallel computing, set program queue's, monitor executing jobs, and sets icon displays.*




### Preferences...

This sets up preferences relating to the graphical user interface (**Settings**), and to molecule displays (**Molecules**). It permits changes to default van der Waals radii used for space-filling models as well as for calculating molecular surface areas and volumes (**VDW Radii**). It also specifies the locations of databases (**Paths**), identifies remote servers (**Hosts**), sets job queues (**Jobs**), (optionally) specifies parallel operations (**Parallel**), specifies platform dependent features (**X-Platform**) and specifies which icons are to be displayed (**Toolbars**). Selection results in one of nine dialogs, depending on which tab has been selected in the previous entry. *Clicking* on a tab brings up the associated dialog. To exit a **Preferences** dialog *click* on **OK**. *Clicking* on **Cancel** or  exits the dialog without instituting any changes.

---

\* **Preferences** are also under the *Spartan'10* menu on the Macintosh version.
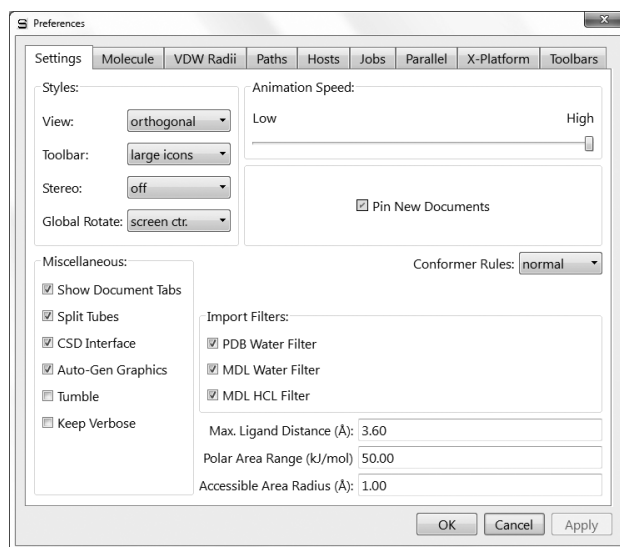
## Settings



- (i) **View: Orthogonal/Perspective**  
Controls the view of structural models and graphics.
- (ii) **Toolbar Style: text/small icons/medium icons/large icons**  
Controls presentation style (text vs. icons) and size of the icons.
- (iii) **Stereo: off/red-cyan**  
Turns **red-cyan** stereographic display on and off.
- (iv) **Global Rotate: screen ctr./molecule ctr.**  
**Screen ctr.** rotates all molecules about a common center; while **molecule ctr.** rotates each molecule about its own center.
- (v) **Animation Speed**  
A slider bar controls the maximum speed for animations. This has become an important control as the performance of graphics cards has increased.
- (vi) **Pin New Documents**  
If *checked*, defaults to display of any new documents (from building or brought in from the **File** menu) irrespective of whether or not they are selected. Does not affect the

status of existing documents. Has no meaning unless **Show Document Tabs** (see discussion following) is *checked*.

(vii) **Show Document Tabs**

If *checked*, this displays a tab at the bottom of the screen for each open document. Only if a tab is *checked* will the molecule display if it is not the selected document. If **Show Document Tabs** is not *checked*, the display corresponds to previous versions of *Spartan*.

(viii) **Split-Tubes**

If *checked*, the tubes in tube and ball-and-spoke models are split to designate multiple bonds.

(ix) **CSD Interface**

If *checked*, displays the **CSD** tab in the **Database** dialog (**Database** under the **Search** menu; **Chapter 19**). *Uncheck* if CSD is not installed.

(x) **Auto-Gen Graphics**

If *checked*, performs graphics calculations in the interface. Does not pertain to orbital displays initiated from **Orbital Energies** under the **Display** menu. Note that graphics calculations on documents containing more than 25 molecules will not be performed in the interface irrespective of this setting.

(xi) **Tumble**


If *checked*, allows automatic tumbling of molecule. To tumble a molecule, select it, *press* the left mouse button, move the mouse and release the button.

(xii) **Keep Verbose**

If *checked*, keeps extended (verbose) output. Normally discarded upon successful completion, this may be useful for identifying the source of problems for calculations that have not successfully completed or have led to suspicious results. Note that keeping verbose output significantly increases the size of the *Spartan* document.

- (xiii) **Conformer Rules**  
Sets rules for conformational searching. **Normal** should be used for **Equilibrium Conformer** and **Conformer Distribution** where Monte Carlo methods are involved. **Trimmed** eliminates less important degrees of conformational freedom and should be used for **Conformer Library** where systematic searches are carried out.
- (xiv) **PDB Water Filter**  
If *checked*, removes water molecules from imported PDB files.
- (xv) **MDL Water Filter**  
If *checked*, removes water molecules from SD files.
- (xvi) **MDL HCl Filter**  
If *checked*, removes HCl molecules from SD files.
- (xvii) **Max-Ligand Distance**  
Set maximum distance (in Å) for an amino acid in a protein to be recognized as bonded to a ligand.
- (xviii) **Polar Area Range**  
Set potential (in kJ/mol) for calculating polar area from electrostatic potential map. Values  $> | \text{range} |$  contribute to polar area. Default is 100 kJ/mol.
- (xix) **Accessible Area Radius**  
Set sphere radius (in Å) for determining accessible area. Default is 1.0.

## Molecule

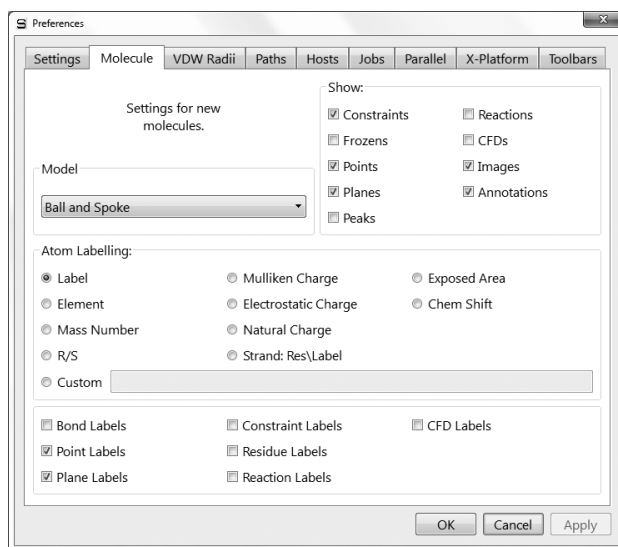

This specifies default settings for model appearance. These settings may be overridden for a specific molecule (or list of molecules) using entries under the **Model** menu and for specific portions of a molecule using **Utilities/Style** dialogs associated with **Properties** dialogs (**Properties** under the **Display** menu; **Chapter 18**).

- (i) **Model: Wire/Ball and Wire/Tube/Ball/Spoke/Space Filling/Line**  
Controls default model style.
- (ii) **Show: Constraints/Frozens/Points/Planes/Peaks/Reactions/CFD's/Images/Annotations**  
If *checked*, constraints and frozen markers, points and planes, reaction arrows and CFD's will always be shown as part of the model. Otherwise, they will be shown only in the appropriate mode.
- (iii) **Atom Labeling: Label/Element/Mass Number/Mulliken Charge/Electrostatic Charge/Natural Charge/Strand: Residue/Label/R/S/Exposed Area/Chem Shift/Custom**  
Controls default label type. **Custom** allows the user to enter an expression designating a label. These are of the same form as used in the spreadsheet. See the section on **User-**



**Defined Expressions** under **Spreadsheet** in the **Display** menu (**Chapter 18**).

(iv) **Bond Labels**

If *checked*, bond labels will be shown.

(v) **Point Labels**

If *checked*, point labels will be shown.

(vi) **Plane Labels**

If *checked*, plane labels will be shown.

(vii) **Constraint Labels**

If *checked*, constraint labels will be shown.

(viii) **Residue Labels**

If *checked*, residue labels will be shown.

(ix) **Reaction Labels**


If *checked*, reaction arrow labels will be shown.

(x) **CFD Labels**

If *checked*, CFD labels will be shown.

## VDW Radii

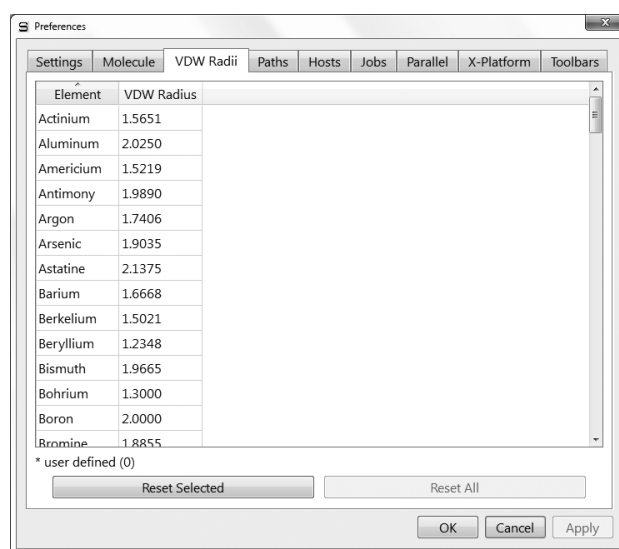
This provides a list of van der Waals radii



To order the list by element name *click* on **Element**, and by atomic radius *click* on **VDW Radius**. Individual entries may be changed from default values by first *clicking* on the entry and then entering a new value. The currently selected entry may be returned to its default radius by *clicking* on **Reset Selected** at the bottom of the dialog, and the full set of radii may be returned to their default values by *clicking* on **Reset All** at the bottom of the dialog.

## Paths

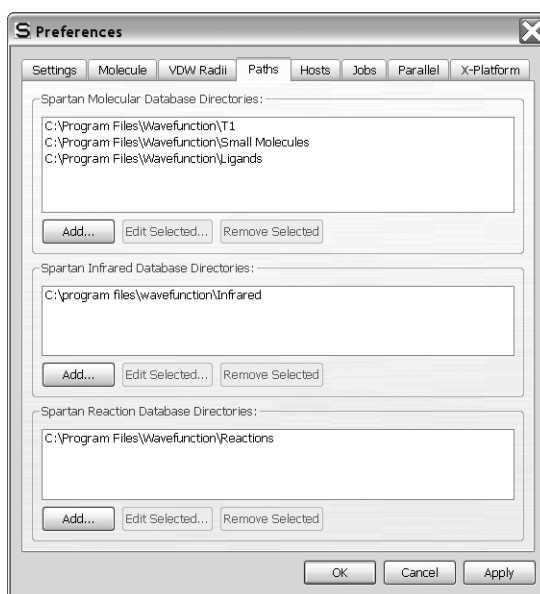This allows setting up of paths for SMD, SSPD/SIRD and SRD database directories.



**Molecular Database Directories** defines the content of SMD, SSPD/SIRD and SRD in terms of one or more directories. Multiple directories can be accessed at once. Specification of which directory (or directories) is to be searched in the **Search Options** dialogs associated with SMD, SSPD, SIRD and SRD (**Databases** under the **Search** menu; **Chapter 19**). Note that SIRD is the same database as SSPD, and only the method of access (spectral matching vs. substructure matching) are different.


All copies of *Spartan* are shipped with ~5,000 molecule subsets

of both SMD and SSPD databases (referred to as **SMD (sample)** and **SSPD (sample)**, respectively. The full databases (**SMD** and **SSPD**) are available from Wavefunction.

Other collections may be added, for example, collections of diverse conformations for use in similarity analyses, or user-supplied collections. **Appendix G** discusses creation of user-defined databases.

## Hosts

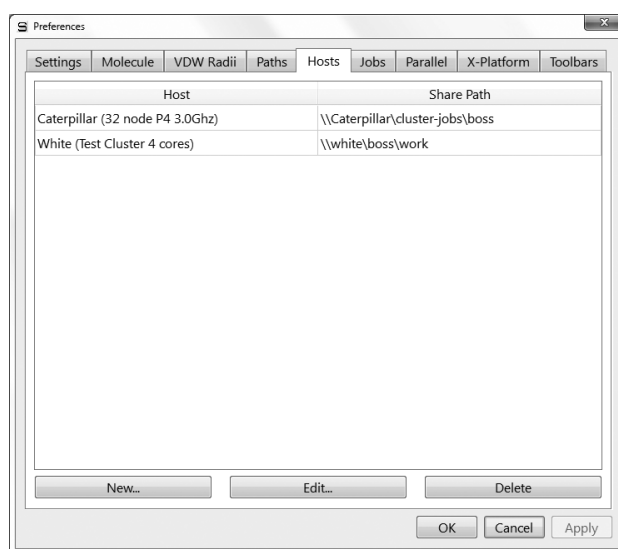
This allows specification of one or more computers that may act as servers.



To add a new host, *click* on the **New...** button and follow the directions provided in the **Host** dialog that appears.

## Jobs

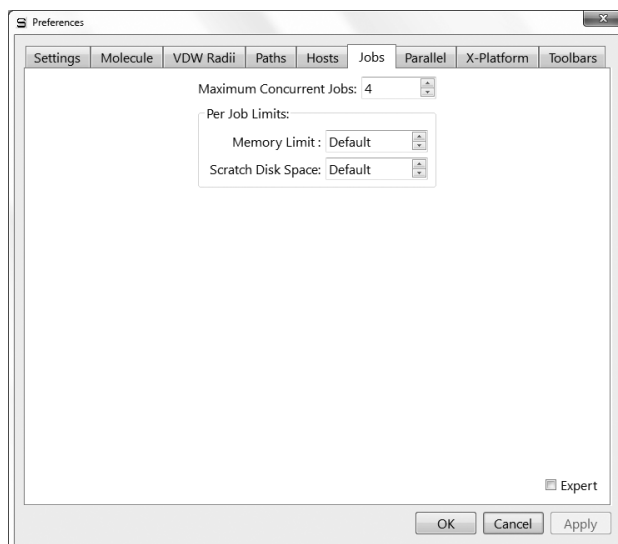This allows setting of job limits.




**Maximum Concurrent Jobs** designates the maximum number of jobs that can be run at one time and **Memory** and **Scratch Disk Space** maximum amounts of memory and disk for a job.

You need to use caution that the maximum number of jobs does not significantly exceed the number of processors (cores) on your computer and also that you have sufficient memory to support this maximum number.

## Parallel

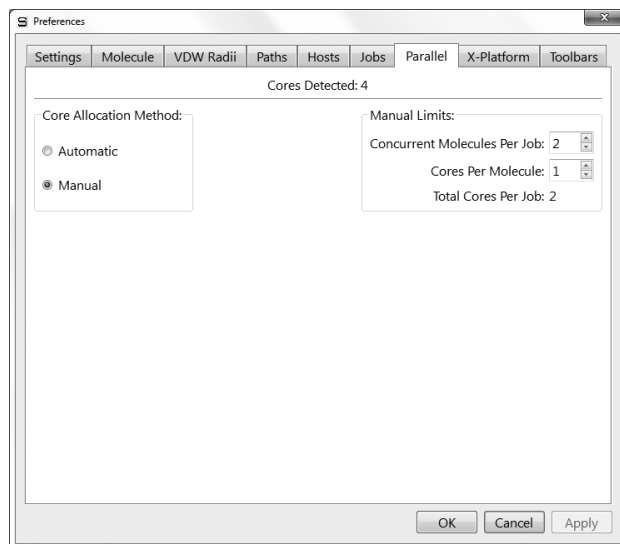This controls the use of a multi-core machine. The number of available cores is indicated at the top of the dialog.



### (i) Core Allocation Method

Toggles between **Automatic** and **Manual**. **Automatic** will assign all available cores to a single molecule (single molecule inside a *Spartan* document) if there is only one molecule in the document and will divide the cores among the molecules in the document if there is more than one. **Manual** will follow settings under **Manual Limits**.

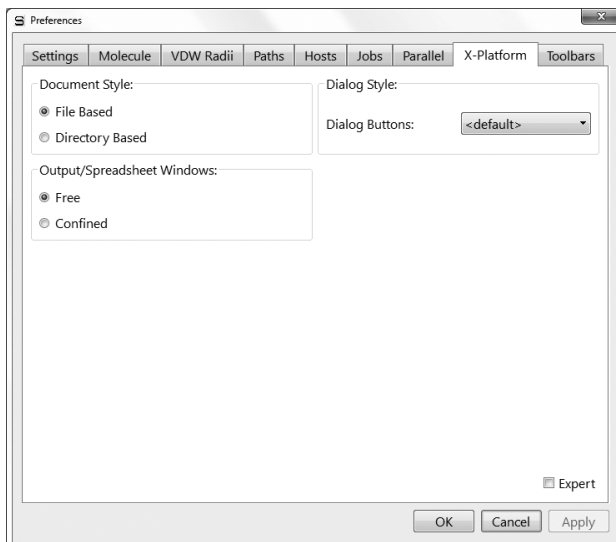
### (ii) Manual Limits

**Concurrent Molecules Per Job** sets the number of jobs (molecules inside a *Spartan* document) that will run simultaneously and **Cores Per Molecule** sets the number of cores used on a single molecule.

If you are simultaneously working with two or more documents, you need to use caution in specifying more than one core/job. It is very easy (and, depending on available memory, potentially very inefficient) to demand more cores than are actually available. This is not an issue if you are working only with a single document even if it contains more than one molecule.

## X-Platform

Items here refer to cross-platform applications and should not be of concern to users of a single platform.



(i) **Document Style**

Toggles between **File Based** and **Directory Based**. The latter is for compatibility between Windows and previous Macintosh and Unix/Linux versions of *Spartan*.

(ii) **Dialog Style**

Menu selects between several dialog formats: **Windows**, **Macintosh**, and **KDE** and **Gnome** Linux. **Default** uses the style associated with the machine on which *Spartan* is installed.

(iii) **Gzip Compression**


(iv) **Archive Compression**

## Toolbars

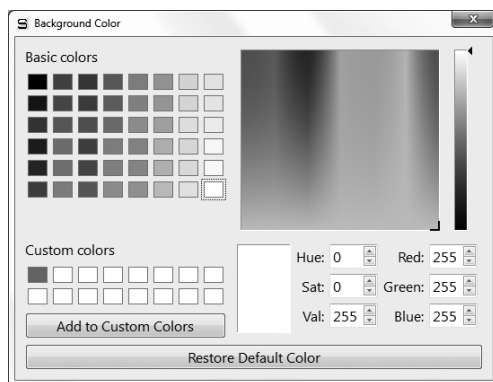A listing of all icons available for entries in the **File**, **Geometry**, **Build**, **Display** and **Search** menus is provided. If *checked*, the icon will appear at the top of the *Spartan* screen.

## Colors


This alters default colors. Selection leads to the **Colors** dialog.




After selecting an object, its color may be set by choosing from the palette, moving the cursor inside the window of colors, or by selecting either a set of hue, saturation and values, or red, green and blue settings. The default color may be reset by *clicking* on **Restore Default Color**. Color selection applies to all objects of the same type, for example, all carbon atoms, and not just to the selected carbon. Further control of colors is available from **Utilities/Style** dialogs associated with **Properties** dialogs (**Properties** under the **Display** menu; **Chapter 18**). *Clicking* on  removes the dialog.

## Fonts/Graphic Fonts

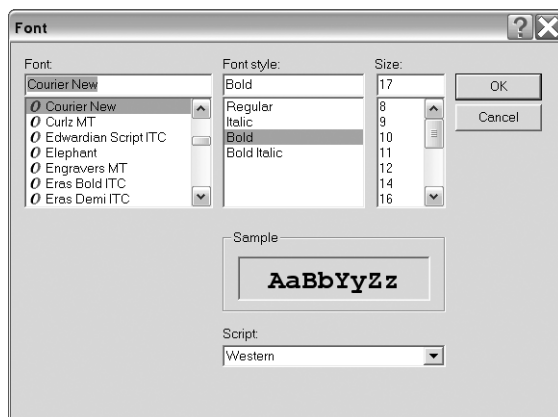This selects fonts, style and size of labels attached to molecules (**Labels** and **Configure...** under the **Model** menu; **Chapter 14**), and plots (**Plots...** under the **Display** menu; **Chapter 18**). Selection leads to the **Fonts** dialog.



Selections are made from the **Font**, **Font Style** and **Size** menus. *Clicking* on **OK** dismisses the dialog with selections kept. *Clicking* on **Cancel** or on  dismisses the dialog but selections are lost.

## Monitor


Job Name	Status	Wall Time	CPU Time
4-methyl-1,2,3-thiadiazole-5-carbonyl chloride	Queued		
4-methylthiazole-2-thiol	Queued		
5-amino-3-methoxy-1,2,4-thiadiazole	Queued		
5-chloro-3-methylsulfonyl-1,2,4-thiadiazole:M00...	SCF Calculation	22:17	20:03
5-chloro-3-methylsulfonyl-1,2,4-thiadiazole:M00...	SCF Calculation	22:17	19:57
5-fluoro-2-methoxy-pyridine	Queued		
6-amino-3-hyoxypyridine	Queued		
bis(4-n,n-dimethyl aminophenyl)-(2-methyl indol...	\\White\Spartan10		
wako910a:2-(Trifluoromethyl)cinnamic acid	SCF Calculation	18:03	09:32
wako910at:Butyl 2,2,2-Trichloroacetimidate	Gradient Calculation	10:10	07:52

Output	Energy	Force	Displacement
6	-826.240399	0.003579	0.128518
7	-826.240425	0.008281	0.073099
8	-826.240875	0.002948	0.046622
9	-826.241003	0.002603	0.185009
10	-826.241234	0.002063	0.189794
11	-826.241545	0.003858	0.188532
12	-826.241957	0.007231	0.155193
13	-826.242655	0.007728	0.180917
14	-826.243261	0.011207	0.190586

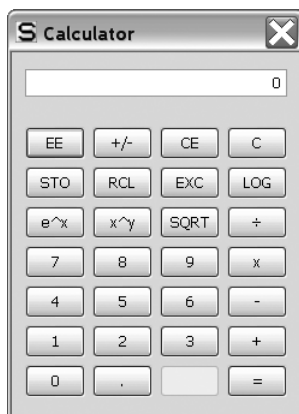
This provides a listing of all executing/queued jobs and their status. To see accumulated output for an executing job, *click* on its name. To kill a job, *click* on its name, and then select **Terminate** from the **Actions** menu at the top of the dialog. To start a queued job (irrespective of the imposed queue limits; see previous discussions under **Jobs** and **Parallel**), *click* on its name and select **Start** from the **Actions** menu.




The **Monitor** may be removed either by selecting **Exit** from the **File** menu or by *clicking* on  at the top of the dialog.

## Calculator

Selection brings up a **Calculator**.



This functions the same way as a normal pocket calculator. The **Calculator** is removed by *clicking* on .

## Model Kit

If *checked*, this signifies that a model kit (organic, inorganic, peptide, nucleotide or substituent) is to remain on screen. Does not apply to ChemDraw.

## File Toolbar, Geometry Toolbar, Build Toolbar, Display Toolbar, Search Toolbar

If *checked*, this provides display of toolbars that access functions contained in the **File**, **Build**, **Geometry**, **Display** and **Search** menus, respectively.

## Cascade, Tile

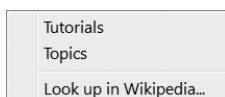
Arranges open windows in a cascade or as tiles on top of *Spartan*'s main window (**NAM**).

# Chapter 21

## The Activities Menu


---

*The **Activities** menu permits on-screen display of the full set of **Spartan** tutorials and a series of topics of practical relevance to molecular modeling. It also allows a Wikipedia page to be brought up (external to **Spartan**).*




### Tutorials/Topics

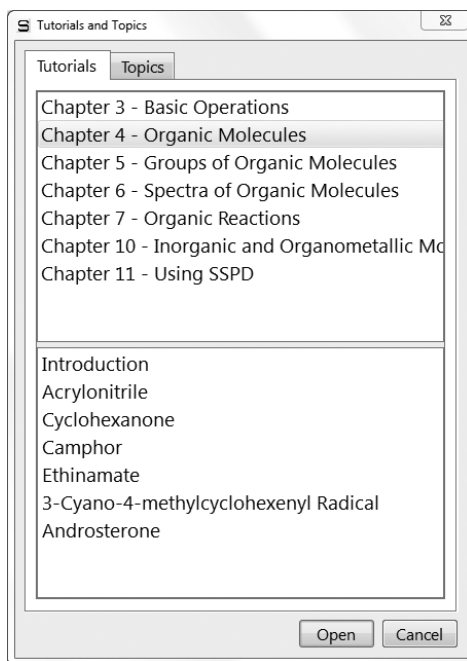Selection of **Tutorials** or **Topics** brings up a dialog with either the **Tutorials** or **Topics** tab selected (**Tutorials** tab shown).



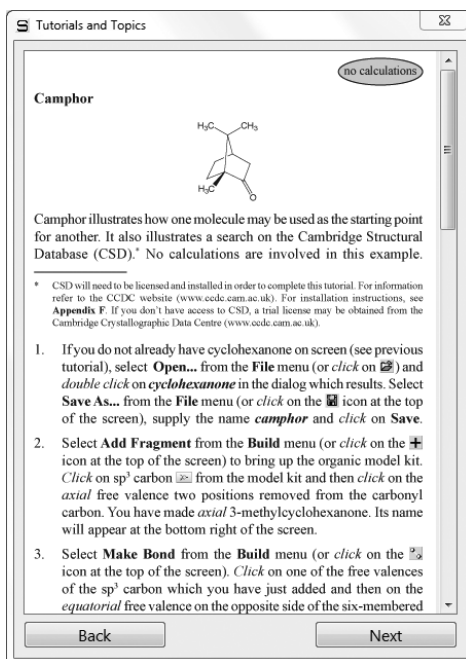
Selection of a *Chapter* (in this case, corresponding to a chapter in this User's Guide) from the upper half of the dialog, followed by *clicking* on **Open** leads to a scrollable window containing the appropriate text. This may be moved about the screen but may not be sized.



Topics are specified by name only (in the upper half of the dialog).

## Look Up in Wikipedia...

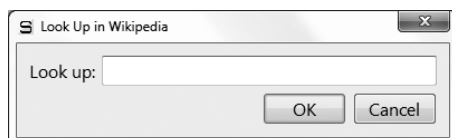Selection results in a dialog.




Entering a query followed by *clicking* on **OK** leads to a Wikipedia page. This occupies a window that is external to *Spartan*.

# Chapter 22
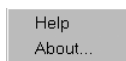
## The Help Menu

---

*This chapter describes help.*




### Help

This provides information relating to application of computational methods available in *Spartan*, as well as technical details regarding the program's operation.

A number of topics are dealt with under **Help**; including *Mouse/Keyboard Functions*, *General Operating Features*, *Units*, *Differences between Spartan'08 and '10*, *Selecting a Theoretical Model*, *Relative Computational Times*, *Memory Usage*, *Calculation Options*, and a number of FAQ documents (*Energy Questions*, *Convergence Questions*, *Questions about Basis Sets*, *Problems with Memory*, *Questions about the Property Module*, *Questions about Conformations and Energy Profiles*, *MMFF94 Compliance*, *Units*) related to computational questions. **Help** also provides a link to Wavefunction's website. Help files are HTML documents and require that Internet Explorer is installed to access them.

The *Spartan Tutorial and User's Guide* (this document) and *A Guide to Molecular Mechanics and Quantum Chemical Calculations* are also available (as PDF files) under **Help**.

Finally, note that several dialogs, in particular, the **Calculations** dialog incorporate imbedded help messages. *Clicking* on  at the upper right, followed by *clicking* on a menu, button, etc. in the dialog gives rise to a brief informative message about the object queried.

## About...\*

Provides information about the user's release of *Spartan*.



**Spartan<sup>10</sup>**  
*for Windows, Macintosh, and Linux*

**WAVEFUNCTION**  
WAVEFUNCTION, INC.  
18401 Von Karman Ave.,  
Suite 370  
Irvine, CA 92612  
sales@wavefun.com • support@wavefun.com  
www.wavefun.com

**Q-CHEM**  
The Design Center, Suite 690  
5001 Baum Blvd., Pittsburgh, PA 15213

The image features a large, stylized molecular model of a complex organic molecule, possibly a protein or a large pharmaceutical compound, rendered in a 3D ball-and-stick style with semi-transparent surfaces. The model is positioned on the right side of the advertisement, partially overlapping the text area.

---

\* About is located under the **Spartan 10** menu on the Macintosh version.

# Appendix A

## Capabilities and Limitations

---

### **Molecular Mechanics Models<sup>1</sup>**

The molecular mechanics module calculates the energy (a combination of strain energy and intramolecular interaction energy), equilibrium geometry and vibrational frequencies, as well as provides for conformational searching\*. SYBYL (Tripos, Inc.) and MMFF (Merck Pharmaceuticals) force fields are available. Energies from the latter may be corrected for the effects of aqueous solvent (MMFFaq). Structures are unaffected. There are no atom limits.

### **Semi-Empirical Models<sup>1</sup>**

The semi-empirical module calculates the heat of formation, wavefunction, equilibrium and transition-state geometries and vibrational frequencies, as well as provides for conformational searching. MNDO, AM1, RM1, PM3 and PM6 models are supported. MNDO/d replaces MNDO for second-row (and heavier) main-group elements, and PM3 has been parameterized for transition metals.\*\* MNDO has been parameterized for H, He, Li-F, Al-Cl, Ca, Zn, Ge, Br, Cd, Sn, I, Hg and Pb; AM1 for H, B-F, Al-Cl, Zn, Ge, Br, Sn and I; RM1 for H, C-F, P, S, Cl, Br and I; PM3 for H-Ne, Mg-Ar, Ca, Ti-Br, Zr, Mo-Pd, Cd-I, Hf-Pt and Hg-Bi and Gd; PM6 for H-Bi except for Ce-Yb. All semi-empirical models involve a minimal valence basis set. Preset limits for semi-empirical calculations are 300 atoms.

### **Hartree-Fock Models<sup>1-3</sup>**

The Hartree-Fock module calculates the energy and wavefunction, equilibrium and transition-state geometries and vibrational frequencies,

---

\* MMFF is the only method available for generation of conformer libraries.

\*\* Parameters for all available transition metals except Zn, Cd and Hg include d-type orbitals in the basis set.

as well as provides for conformational searching. Both closed-shell and open-shell (either ROHF or UHF) calculations are supported. Available basis sets are enumerated later in this appendix. Hartree-Fock models have been implemented for use on multi-core processors. Preset limits for Hartree-Fock calculations are 200 atoms and 2,000 basis functions.

### **Density Functional Models<sup>1,2</sup>**

The density functional module calculates the energy and wavefunction, equilibrium and transition-state geometries and vibrational frequencies, as well as provides for conformational searching. The module supports local density calculations and a wide variety of gradient corrected and hybrid functionals. Only a few functionals are provided in the menu of functionals, but other selections can be entered into the **Options** line of the **Calculations** dialog (**Setup** menu; **Chapter 17**). The **Options** line may also be used to specify combinations of exchange and correlation components to be used in building a total functional. For example, the line **EXCHANGE = SLATER: 0.9, BECKE: 0.1** indicates a combination of 9 parts Slater exchange and 1 part Becke exchange. Similarly, the line **CORRELATION = VWN: 0.19, LYP: 0.81** indicates a combination of 19 parts VWN and 81 parts LYP correlation.

A full listing of available functionals and instructions for constructing custom functionals from available exchange and correlation components is provided in **Table A-1**. Available basis sets are enumerated later in this appendix. Density functional models have been implemented for use on multi-core processors.

Preset limits for density functional calculations are 200 atoms and 2,000 basis functions.

### **Møller-Plesset Models<sup>1,2</sup>**

The Møller-Plesset module calculates the energy and wavefunction and equilibrium and transition-state geometries as well as provides for conformational searching using either the MP2 or RI-MP2 models. Vibrational frequencies are also available from both models, but

**Table A-1: Available Exchange, Correlation and Combination Functionals**

<b>exchange functionals</b>	
HF	Fock exchange
Slater	Slater (Dirac 1930)
Becke	Becke (1988)
Gill96	Gill (1996)
GG99	Gilbert, Gill (1999)
B (EDF1)	Becke (uses EDF1 parameters)
PW91	Perdew, Burke, Ernzerhof (1996)
<b>correlation functionals</b>	
VWN	Vasko, Wilk, Nusair #5
LYP	Lee, Yang, Parr
PW91	GGA91, Perdew
LYP	(EDF1)
P86	Perdew (1986)
PZ81	Perdew, Zungner
PBE	Perdue, Burke, Ernzerhof (1996)
<b>combination functionals</b>	
B3PW91	
B3LYP	
B3LYP5 (Gaussian's version of B3LYP)	
EDF1	
EDF2	
BMK	
M05	
M05-2X	
M06	
M06-2X	
M06-HF	
M06-L	
$\omega$ B97	
$\omega$ B97X	
$\omega$ B97X-D	



require numerical differentiation of analytical gradients. Available basis sets are enumerated later in this appendix. RI-MP2 models have been implemented for use on multi-core processors.

Depending on basis set, practical considerations limit application of these techniques to molecules comprising 20-30 heavy (non-hydrogen) atoms, although program limits are set much larger.

MP3 and MP4 calculations are supported under **Advanced Correlated Models**.

### **Advanced Correlated Models<sup>1,2</sup>**

The advanced correlated module calculates only the energy using the MP3, MP4, CCSD, CCSD(T), CCSD(2), OD, OD(T), QCISD, QCISD(T), QCCD and QCCD(2) correlated models. Available basis sets are enumerated later in this appendix.

Depending on basis set, practical considerations limit application of these techniques to molecules comprising 10-20 heavy (non-hydrogen) atoms, although program limits are set much larger.

### **Thermochemical Recipes**

The G2, G3 and G3(MP2) recipes are available for highly-accurate heat of formation calculations. Practical considerations generally limit application of the first two recipes to molecules comprising less than 10 heavy atoms and the G3(MP2) recipe to molecules comprising less than 15 heavy atoms.

The T1 recipe<sup>4</sup> has been developed to closely reproduce G3(MP2) heats of formation. It makes use of an HF/6-31G\* geometry, instead of the MP2/6-31G\* geometry, replaces the large basis set MP2 calculation in G3(MP2) by a dual-basis set RI-MP2 calculation, eliminates both the QCISD(T) calculation and the vibrational frequency calculation (needed to obtain zero-point energy and to correct the energy for finite temperature) and introduces an empirical correction based on the atom counts and Mulliken bond orders. The result is that T1 requires 2-3 orders of magnitude less computation time than G3(MP2). Calculations on molecules with 30-50 heavy atoms are practical.

T1 is limited to closed-shell, uncharged molecules comprising H, C, N, O, F, Si, P, S, Cl and Br only, and that can be properly described in terms of a conventional valence structure. Development of a “T1 like” recipe for transition states is underway.

### **Excited-State Models**

Two classes of models are available for calculations of excited states: configuration interaction (CI) models and time-dependent density functional (TDDFT) models. The former comprises three models: CIS (configuration interaction singles), CISD (configuration interaction singles with doubles correction) and RI-CIS(D), while the latter supports the same functionals or combinations of functionals available for ground state density functional calculations. At the present time, only the CIS and RI-CIS(D) models are available for geometry optimization and frequency calculations, all other models being restricted to energy calculations. Available basis sets are enumerated later in this appendix.

Preset limits for excited-state calculations are 200 atoms and 2,000 basis functions. The practical limit is much smaller than this.

### **Basis Sets and Pseudopotentials**

Available basis sets for Hartree-Fock, density functional, MP2 (RI-MP2) and advanced correlated models for ground states and CIS, CISD (RI-CIS(D)) and time-dependent density functional (TDDFT) models for excited states include the full range of Pople basis sets: minimal STO-3G, split-valence 3-21G, 6-31G, 6-311G and polarization 6-31G\*, 6-31G\*\*, 6-311G\* and 6-311G\*\*. These may be supplemented with additional polarization functions and/or with diffuse functions. Also available are the cc-pVDZ, cc-pVTZ, cc-pVQZ, cc-pV5Z and aug-cc-pCVQZ basis sets and the G3Large basis set.

STO-3G is available for H–I, 3-21G for H–Cs, 6-31G and basis sets built on it for H–Kr, 6-311G and basis sets built on it for H–Kr and I, cc-pVDZ, cc-pVTZ, cc-pVQZ, and G3 Large for H–Ca and Ge–Kr and aug-cc-pCVQZ for H–Ar.

Use of minimal and split-valence basis sets for any but Hartree-Fock calculations generally leads to poor results and is not recommended.

Dual-basis set techniques are supported for ground-state energy, equilibrium and transition-state geometry calculations for Hartree-Fock, density function and MP2 models.

Pseudopotentials are available for calculations on molecules incorporating heavy atoms (for which all-electron basis sets are not available). Their use is automatic.

## Solvent Models

*Spartan'10* supports two different empirical solvation models. The SM5.4 model<sup>5</sup> (available for H, C–F, S–Cl, Br and I) estimates the aqueous solvation energy. This in turn may be added to the gas-phase energy obtained from any molecular mechanics or quantum chemical model. Because the computational cost of SM5.4 is very low (on the order of that for an energy using a semi-empirical model), this solvent correction is provided automatically for all quantum chemical calculations.

The SM8 model<sup>6</sup> (available for H, C–F, S–Cl, Br and I) is applicable only to Hartree-Fock and density functional models. It leads to a change in the wavefunction, meaning that, in addition to energy, equilibrium or transition-state geometry, properties, infrared spectra and graphical displays may be obtained in the presence of solvent. SM8 has been parameterized for water as well as a variety of organic solvents. A few are available in a menu (acetone, diethyl ether, DMF, DMSO, ethanol, methylene chloride, THF and toluene) while others may be specified using the **Solvent** keyword (**Appendix C**). The computational cost of SM8 is roughly that of an energy with the corresponding Hartree-Fock or density functional model.

## Properties and Spectra

The properties module is automatically invoked following completion of a molecular mechanics or quantum chemical calculation. By default, it provides for text output printing, population analyses and

atomic charge calculation (Mulliken, natural bond orbital and based on fits to electrostatic potentials), evaluation of thermodynamic quantities (enthalpy, entropy and free energy) if vibrational frequencies are available, and calculation of the dipole moment.

Several properties based on “graphical displays” have been implemented for the first time in *Spartan*’10. These include the area, exposed area and volume of an electron density surface, the polar area of an electrostatic potential map for cutoff values of 75, 100 and 125 kJ/mol, the maximum and minimum values of the electrostatic potential on an electrostatic potential map and the minimum value of the local ionization potential on a local ionization potential map. These require calculation of the respective surfaces and property maps, but only a small set of numerical values and not the surfaces and maps themselves are actually stored. These properties need to be explicitly requested by *checking* the **QSAR** box inside the **Calculations** dialog (**Calculations** under the **Setup** menu; **Chapter 17**). They are either available from the **QSAR** tab of the **Molecule Properties** dialog (**Properties** under the **Display** menu; **Chapter 18**) or from the **Formula Editor** dialog (**Formulas** under the **Display** menu; **Chapter 18**).

The properties module is also responsible for calculating quantities related to infrared/Raman spectra (vibrational frequencies and infrared/Raman intensities), NMR spectra (chemical shifts) and UV/visible spectra (excited-state transition energies). IR spectra calculations may be carried out with molecular mechanics models, semi-empirical models, Hartree-Fock models, density functional models and MP2 and RI-MP2 models. NMR spectra calculations are limited to Hartree-Fock and density functional models and UV/visible spectra calculations are limited to Hartree-Fock models (CIS for excited states) and density functional models (TDDFT for excited states).

## Graphical Models

The graphics module provides for data preparation associated with the display as surfaces, property maps and slices of molecular orbitals, electron densities, spin densities, electrostatic potentials, hydride

potentials and local ionization potentials. The sizes of electron density surfaces (and property maps based on electron density surfaces) may be chosen either using a specific value of the density or a value that encloses a specific percentage of the total number of electrons. Accessible and inaccessible regions may be distinguished for electron density surfaces and all property maps based on electron density surfaces.

**Spartan'10** displays an orbital energy diagram (12 highest-occupied molecular orbitals and 2 lowest-unoccupied molecular orbitals). Molecular orbital graphics may be generated on-the-fly.

## Similarity Analysis

The similarity analysis module evaluates and quantifies the extent to which molecules are similar based either on selected atomic centers or on selected chemical function descriptors (CFD's). In addition, it quantifies the degree of fit of a molecule into a chemical environment (a pharmacophore). Similarity analyses are carried out between one or more molecules or pharmacophores in a document (the template) and one or more molecules or pharmacophores in one or more documents (the library). Molecules in the library (only) may either correspond to single conformers or to collections of diverse conformers.

- 
1. For a general discussion and assessment of the techniques and methods available in **Spartan**, see: W.J. Hehre, *A Guide to Molecular Mechanics and Quantum Chemical Calculations*, Wavefunction, Inc., Irvine, CA 2003. This is available as a PDF under the **Help** menu.
  2. For a review of the quantum chemical methods utilized in **Spartan** (except for semi-empirical methods) with emphasis on recent developments, see: Y. Shao, L.F. Molnar, Y. Jung, J. Kussmann, C. Ochsenfeld, S.T. Brown, A.T.B. Gilbert, L.V. Slipchenko, S.V. Levchenko, D.P. O'Neill, R.A. DiStasio Jr., R.C. Lochan, T. Wang, G.J.O. Beran, N.A. Besley, J.M. Herbert, C.Y. Lin, T. Van Voorhis, S.H. Chien, A. Sodt, R.P. Steele, V.A. Rassolov, P.E. Maslen, P.P. Korambath, R.D. Adamson, B. Austin, J. Baker, E.F.C. Byrd, H. Dachsel, R.J. Doerksen, A. Dreuw, B.D. Dunietz, A.D. Dutoi, T.R. Furlani, S.R. Gwaltney, A. Heyden, S. Hirata, C-P. Hsu, G. Kedziora, R.Z. Khalliulin, P. Klunzinger, A.M. Lee, M.S. Lee, W.Z. Liang, I. Lotan, N. Nair, B. Peters, E.I. Proynov, P.A. Pieniazek, Y.M. Rhee, J. Ritchie, E. Rosta, C.D. Sherrill, A.C. Simmonett, J.E. Subotnik, H.L. Woodcock III, W. Zhang, A.T. Bell, A.K. Chakraborty, D.M. Chipman, F.J. Keil, A. Warshel, W.J. Hehre, H.F. Schaefer, J. Kong, A.I. Krylov, P.M.W. Gill and M. Head-Gordon, *Phys. Chem. Chem. Phys.*, **8**, 3172 (2006).
  3. For an older account see: W.J. Hehre, L. Radom, P.v.R. Schleyer and J.A. Pople, *Ab Initio Molecular Orbital Theory*, Wiley, New York, 1986.
  4. W.S. Ohlinger, P.E. Klunzinger, B.J. Deppmeier, W.J. Hehre; *J. Phys. Chem. A.*, **113**, 10, 2165 (2009).
  5. C.C. Chambers, G.D. Hawkins, C.J. Cramer and D.G. Truhlar, *J. Phys. Chem.*, **100**, 16385 (1996).
  6. A.V. Marenich, R.M. Olsen, C.P. Kelly, C.J. Cramer and D.G. Truhlar, *J. Chem. Theory Comput.*, **3**, 2011 (2007).

# Appendix B

## Menus

---

### *Spartan Screen*

#### **File**

---

<u>N</u> ew	Brings up a model kit for molecule building or substitution
<u>O</u> pen...	Opens (imports) a molecule or set of molecules
<u>C</u> lose	Closes a molecule or set of molecules
<u>S</u> ave	Saves (exports) a molecule
Save <u>A</u> s...	Saves a molecule as a <i>Spartan</i> document under a user-specified name; creates files for use in the Spartan Molecular Database, Spartan Reaction Database, Spartan Infrared Database and Spartan Spectra and Properties Database
<u>N</u> ew Molecule	Adds a molecule to an existing list; brings up a model kit for molecule building or substitution
<u>D</u> elete Molecule	Deletes a molecule (or molecules) from a list
Append <u>M</u> olecule(s)...	Appends molecules to an existing list
Access PDB Online...	Accesses the online Protein Data Bank (PDB)
<u>P</u> rint...	Prints on-screen display; also prints contents of output window and the <b>Spreadsheet</b>
Embedded Data	Associates external files with a <i>Spartan</i> document, either at molecule or document level (or both)
Start/Stop QuickTime Recording	Starts and stops QuickTime recording of contents of main <i>Spartan</i> screen (only)
<u>E</u> xit	Exits <i>Spartan</i>

## **Edit**

---

<u>U</u> ndo	Undoes previous operations
<u>C</u> ut	Moves the current molecule or contents of the selection box to the clipboard
<u>C</u> opy	Copies the current molecule or contents of the selection box to the clipboard
<u>P</u> aste	Pastes contents of the clipboard to the screen
<u>F</u> ind...	Locates a text string in the output dialog or an on-screen molecular fragment
Find <u>N</u> ext	Locates next occurrence of a text string or molecular fragment
<u>C</u> enter	Centers the molecule on screen; applies to all molecules in a document
<u>C</u> lear	Clears the selected molecule

## **Model**

---

<u>W</u> ire	Displays structure as wire-frame model
<u>B</u> all and Wire	Displays structure as ball-and-wire model
<u>T</u> ube	Displays structure as tube model
<u>B</u> all and Spoke	Displays structure as ball-and-spoke model
<u>S</u> pace Filling	Displays structure as space-filling model
<u>L</u> ine	Displays structure as a line drawing
<u>H</u> ide	Hides structure model from view
<u>G</u> lobal Model	Applies model type and labels of current molecule to all molecules in the document
<u>C</u> oupled	Couples motions of all molecules in the document
<u>H</u> ydrogens	Toggles hydrogens “on” and “off”
<u>L</u> abels	Toggles labels “on” and “off”
<u>R</u> ibbons	Toggles ribbons “on” and “off”
<u>R</u> amachandran Plot	Display a Ramachandran plot for a protein that has been brought in from the Protein Data Bank (PDB)

Hydrogen <u>B</u> onds CFD's	Toggles hydrogen bonds “on” and “off” Toggles chemical function descriptors (hydrophobes, hydrogen-bond donors/acceptors, positive/negative centers) “on” and “off”
<u>C</u> onfigure...	Labels atoms, bonds, etc., provides information about polypeptides/polynucleotides residues and designates ribbon displays

## **Geometry**

---

Measure <u>D</u> istance	Displays and/or sets bond distance
Measure <u>A</u> ngle	Displays and/or sets bond angle
Measure <u>D</u> ihedral	Displays and/or sets dihedral angle
<u>F</u> reeze Center	Freezes selected atomic positions
Set <u>T</u> orsions	Selects bond rings for conformational searching
Set Similarity Centers	Identifies atoms /CFD's for similarity analyses
<u>C</u> onstrain Distance	Constrains bond distance
<u>C</u> onstrain Angle	Constrains bond angle
<u>C</u> onstrain Dihedral	Constrains dihedral angle
Define <u>P</u> oint	Defines a point as a geometric mean of a set of atoms; defines ligand point as a position that is perpendicular to the centroid of a plane made by three or more atoms
Define <u>P</u> lane	Defines a plane made by three or more atoms
Define CFD	Defines the location of a CFD
<u>A</u> lign	Aligns molecules in a document according to selected atoms, structure, or CFD's

## **Build**

---

<u>V</u> iew	Removes the model kit (unless designated for permanent display)
Add <u>F</u> ragment	Brings up a model kit that contains libraries of atomic fragments, functional groups, rings, ligands and substituents; accesses



	user-defined libraries; provides access to ChemDraw ( <b>NAL</b> , <b>NAM</b> )
<u>D</u> elete	Deletes atoms, bonds, points, planes, etc.
<u>M</u> ake Bond	Makes bonds between free valences or atoms
<u>B</u> reak Bond	Breaks a bond
<u>M</u> inimize	Performs energy minimization using MMFF molecular mechanics

## **Setup**

---

<u>C</u> alculations...	Sets up molecular mechanics and quantum chemical calculations and similarity analyses; specifies calculation of IR, NMR and UV/vis spectra
<u>S</u> urfaces	Sets up generation of and displays graphical surfaces
<u>S</u> ubmit	Submits job to the execution queue

## **Display**

---

<u>O</u> utput	Displays text output
<u>P</u> roperties	Displays molecule, bond and atom properties as well as information about CFD's, geometrical constraints, graphical surfaces, plots (including spectral plots) and statistical analyses
<u>O</u> rbital Energies	Displays an orbital energy diagram and allows on-the-fly generation and display of molecular orbitals
<u>S</u> urfaces	Sets up generation of and displays graphical surfaces (same as entry in <b>Setup</b> menu)
<u>S</u> pectra	Displays IR, Raman, NMR and/or UV/visible spectra, animates vibrational modes (IR and Raman), and accesses on-line experimental spectral databases; fits calculated spectra to experimental spectra and (for IR and Raman) to reference spectrum

Formulas	Constructs expressions for use in spreadsheet and database searching/mining operations
<u>S</u> preadsheet	Displays spreadsheet
<u>P</u> lots...	Creates plots from the data in the spreadsheet
<u>S</u> imilarities	Displays results of similarity analyses
<u>R</u> eactions	Calculates reaction (activation) energies using data either from current document or from the Spartan Molecular Database

## **Search**

---

Define Structure Query	Add “growth points” to a structure defining an underlying substructure
Define Reaction Query	Add “curly arrows” to structure query (“reactant”) in order to relate it to the product of a chemical reaction or to the transition state leading to that product
<u>O</u> pen/ <u>C</u> lose <u>D</u> atabases	Performs substructure search on Cambridge Structural Database (CSD), Spartan Molecular Database (SMD) and Spartan Reaction Database (SRD); performs spectral search of the Spartan Infrared Database (SIRD)
<u>G</u> uess Transition <u>S</u> tate	Provides transition-state guess based on reaction database or, lacking a database entry, based on linear synchronous transit
<u>I</u> dentify <u>T</u> automers	Identifies tautomers
<u>E</u> xtract Ligands	Extracts ligands and pharmacophores based on protein-ligand binding from PDB files

## **Options**

---

<u>P</u> references...	Sets various run-time and labeling preferences and specifies execution preferences and run-time queue; invokes parallel calculation; sets up icon displays
<u>C</u> olors	Sets screen and model colors

<u>F</u> onts...	Sets fonts for labels and plot displays
<u>G</u> raphics <u>F</u> onts...	Sets fonts for graphical displays
<u>M</u> onitor	Monitors and allows for killing executing jobs and bypassing queue to start jobs
<u>C</u> alculator	Pocket calculator
<u>M</u> odel <u>K</u> it	Keeps a model kit on screen
<u>F</u> ile, <u>B</u> uild, <u>G</u> eometry, <u>S</u> earch <u>T</u> oolbar	Toggles toolbars “on” and “off”
<u>C</u> ascade	Arranges open windows in a cascade
<u>T</u> ile	Arranges open windows as tiles ( <b>NAL</b> , <b>NAM</b> )

## **Activities**

---

<u>T</u> utorials	Brings up <i>Spartan’10</i> tutorials
<u>T</u> opics	Brings up selection of topics relevant to calculations performed in <i>Spartan</i> .
<u>P</u> roblems	Brings up a selection of problems suitable for <i>Spartan’10</i>
<u>L</u> ook up in Wikipedia	Brings up an html page pointing to a Wikipedia page

## **Help**

---

<u>H</u> elp	Provides information about the performance and timing of computational methods in <i>Spartan</i> ; provides information about using graphical models in <i>Spartan</i> ; bulletin board for FAQs about <i>Spartan</i>
<u>A</u> bout...	Provides program version information for citation and support

## *Contextual*

### **Main Screen**

---

Copy	Copies selected molecule to the clipboard
Paste	Pastes the contents of the clipboard into the selected document
Delete Selected	Deletes selected molecule from document
Properties	Brings up the <b>Molecular Properties</b> dialog

### **Spreadsheet**

---

Copy	Copies text of selected cell or cells to the clipboard. If leftmost cell (or cells) selected, copies molecule(s) to the clipboard
Paste	Pastes the contents of the clipboard into selected cells. If leftmost cell (or cells) selected, either pastes text or molecule(s) depending on menu choice
Add	Brings up the <b>Add</b> dialog (spreadsheet) for adding calculated quantities into the spreadsheet
Sort	Sorts the column from low to high. <i>Pressing the <b>Shift</b> key prior to menu selection sorts from high to low</i>
Format Selected	Formats selected cell(s), selected column(s) if selection is in a header cell, or entire spreadsheet if selection is header cell of leftmost column
Delete Selected	Deletes selected molecule(s) from document
Append	Appends the contents of <i>Spartan</i> document(s) to the spreadsheet (corresponding to the selected document)
Rename Selected Using SMD/SSPD	Rename selected molecule(s) with names in the Spartan Molecular Database or Spartan Spectra and Properties Database; appears only when leftmost cell(s) selected

Properties                      Brings up the **Molecular Properties** dialog

### **Substituent Model Kit**

---

Copy	Copies selected molecule(s) from a custom library onto the clipboard
Paste	Pastes molecule(s) from the clipboard into a custom library
Delete	Deletes selected molecule(s) from the custom library
Append	Appends the contents of <i>Spartan</i> document(s) to the custom library; brings up a file browser

### **Reactions**

---

Copy	Copies selected text to the clipboard
Print	Prints selected text

### **Output Window**

---

Copy	Copies selected text to the clipboard
Print	Prints selected text

# Appendix C

## Commonly-Used Program Options

---

This appendix describes a number of commonly-used program options. These are specified using keywords input into the **Options** box in the **Calculations** dialog (**Setup** menu). In addition to the above, any method, functional or combination of functionals and/or basis set supported in *Spartan* (and not included explicitly in the menus) may be typed into the **Options** box. For example, typing **EDF2 6-311G(2df,2pd)** will result in density functional calculation using the EDF2 functional and the 6-311G(2df,2pd) basis set. (See **Chapter 17** for further details.)

Keywords are case insensitive. An equals sign (=) separates a keyword from its integer, real or character string value. Keywords may either be single words or expressions. **keyword=N** indicates an integer argument, **keyword=C** indicates a character argument, and **keyword=F** indicates a floating point argument. Real values can optionally include an **eNN** or **e-NN** (floating point power of ten). Default values for keywords are indicated in **bold**, and alternative values in *italics*.

### Conformational Search

keyword	value	description
<b>SearchMethod</b>	<b>C</b>	conformational searching method employed; overrides default choice which depends on complexity of system may be overridden with one of the following: <i>Systematic</i> use systematic method <i>MonteCarlo</i> use Monte-Carlo method <i>Sparse</i> use systematic method but randomly eliminate conformers

to stay within set conformer limits. Used in conformer library generation.

<b>MaxEnergy</b>	<b>F</b>	sets the maximum energy (in kJ/mol) at which a conformer will be kept to (minimum energy + F), where minimum energy is the energy of the lowest energy conformer encountered at that point in the search; default is <b>40</b> kJ/mol
<b>StartTemperature</b>	<b>F</b>	sets initial temperature (in K) for a Monte Carlo search (only) to F; high temperatures will result in the molecule further exploring global space, while low temperatures will result in the molecule further exploring local space; default is <b>5000K</b>

### Similarity Analysis

<b>keyword</b>	<b>value</b>	<b>description</b>
<b>SingleEnantiomer</b>	<b>-</b>	Limits similarity analysis to the enantiomer actually in specified library documents. The default is to consider both enantiomers
<b>QuitSimilarity</b>	<b>N</b>	Quit similarity analysis involving a conformer list (in the graphical user interface) and/or one or more conformer libraries when N good matches have been identified; default is <b>1</b>
<b>QuitCriterion</b>	<b>F</b>	R <sup>2</sup> criterion for a good match; default is <b>0.9</b>

### Geometry Optimization

<b>keyword</b>	<b>value</b>	<b>description</b>
<b>GeometryCycle</b>	<b>N</b>	set maximum number of equilibrium and transition-state geometry optimization cycles to N; default is <b>number of independent variables +20</b> for equilibrium geometry optimization and <b>3x number of independent variables</b> for transition-state geometry optimization
<b>GradientTolerance</b>	<b>F</b>	set convergence criterion in equilibrium and transition-state) geometry optimization for the maximum gradient component (in hartrees/bohr) to F; default is <b>4.5 x 10<sup>-4</sup></b>

<b>DistanceTolerance</b>	<b>F</b>	set convergence criterion in equilibrium and transition-state geometry optimization for the maximum change in a bond length (in Å) to F; default is <b>1.8 x 10<sup>-3</sup></b>
<b>Hess</b>	<b>C</b>	Choose Hessian (matrix of second derivatives) for start of geometry (transition-state geometry) optimization (in the absence of a Hessian from a previous calculation, the default is <b>MMFF</b> molecular mechanics for geometry optimization and <b>PM3</b> semi-empirical for transition-state geometry optimization; may be overridden with one of the following:
		<i>unit</i> unit matrix
		<i>MMFF</i> MMFF molecular mechanics
		<i>AM1</i> AM1 semi-empirical
		<i>STO-3G</i> STO-3G Hartree-Fock
		<i>3-21G</i> 3-21G Hartree-Fock
		<i>6-31G*</i> 6-31G* Hartree-Fock

## Molecular Mechanics Calculations

<b>keyword</b>	<b>description</b>
<b>MMFF94s</b>	Use MMFF94s method instead of MMFF94. The only significant difference between the two is that MMFF94 forces nitrogens attached to aromatic rings to be planar (as in aniline). While this is widely perceived to be the case, in fact the nitrogen in aniline and related molecules are pyramidal not planar.

## Quantum Chemical Calculations

<b>keyword</b>	<b>value</b>	<b>description</b>
<b>Diis</b>	<b>N</b>	switch on diis all the time in SCF procedure. N is the size of the iterative subspace; it should be an integer between 2 and 10; default is <b>5</b>
<b>NoDiis</b>		turn off diis SCF convergence accelerator
<b>ScfCycle</b>	<b>N</b>	set maximum number of SCF iterations to N; default is <b>50</b>
<b>Guess</b>	<b>C</b>	choose initial wavefunction guess; in the absence of a guess from a previous calculation, the



default is either **sad** or **PM3** depending on the calculation type, and may be overridden with one of the following:

*core* diagonalize the core Hamiltonian  
*sad* superposition of atomic densities  
*PM3* PM3 semi-empirical

**Scf** C SCF procedure; default is **restricted** for closed-shell systems and **unrestricted** for open-shell systems, either of which may be specified for both closed and open-shell systems

**SmallBasis** C specifies the small basis set to be used in dual basis set calculations. The large basis set is that indicated in the menu.

**Core** C use of frozen core approximation in Møller-Plesset and advanced correlated calculations; default is **frozen** which may be changed to *thawed*

**ScfTolerance** C control of all tolerances in SCF procedure and to affect the precision of the algorithm; default is **normal** which may be changed to *high* or *veryhigh*

F set SCF energy convergence criterion (in hartrees) to F; default is **1.0 x 10<sup>-6</sup>**

**BigGrid** – uses very large grid in density functional calculations

**Mix** – specifies that  $\alpha$  and  $\beta$  HOMO's in the guess wavefunction should be constructed according to:

$$\text{HOMO}_\alpha = \frac{\text{HOMO} + \text{LUMO}}{\sqrt{2}}$$

$$\text{HOMO}_\beta = \frac{\text{HOMO} - \text{LUMO}}{\sqrt{2}}$$

useful for generating a guess wavefunction for singlet diradical

**Exchange** E specifies custom exchange functional or combination of exchange functionals for density functional calculation. See **Appendix A** for details.

<b>Correlation</b>	<b>C</b>	specifies custom correlation functional or combination of correlation functionals for density functional calculation. See <b>Appendix A</b> for details.
<b>Energy Only</b>		skips calculation of the density matrix for an MP2 energy calculation, reducing the time by 50%

## Semi-Empirical Calculations

<b>keyword</b>	<b>value</b>	<b>description</b>
<b>NoAmideFix</b>	–	turns “off” molecular mechanics amide correction in PM3 calculations
<b>hh</b>	<b>C</b>	turns “on” and “off” hydrogen-hydrogen repulsion term in PM3 calculations; default is <b>on</b> only where transition metals are involved and <b>off</b> otherwise
<b>Polarizability</b>	–	calculates and prints polarizabilities and hyperpolarizabilities

## Property and Spectra Calculations

<b>keyword</b>	<b>value</b>	<b>description</b>
<b>NBO</b>	<b>C</b>	specify details of NBO population analysis; default is <b>normal</b> which may be overridden for <i>ionic</i> and <i>3C</i> (three center) systems
<b>NoElCharge</b>	–	turns “off” calculation of electrostatic charges
<b>ElCharge</b>	<b>N</b>	adjusts size of grid used to calculate charge from electrostatic potential; default is <b>1</b> point/atomic unit
<b>Temperature</b>	<b>F</b>	temperature (in K) used in calculation of thermodynamic properties; default is <b>298°K</b>
<b>TEMPRANGE</b>	<b>Fi,Fj,Fk</b>	temperature range applied to thermodynamic properties when IR frequencies are calculated, where i=starting temperature, j=ending temperature, and k=increment, all values in K
<b>Pressure</b>	<b>F</b>	pressure (in atm) used in calculation of thermodynamic properties; default is <b>1</b> atm
<b>Moments</b>	–	calculates and prints principal moments of inertia

<b>NoSolvent</b>	–	turns “off” calculation of aqueous solvation energy using SM5.4 model
<b>Solvent</b>	<b>C</b>	specifies solvent for SM8 model. See <b>FAQ (Properties)</b> available under <b>Help (Help menu; Chapter 19)</b> for a list of available solvents.
<b>NoAreaVolume</b>	–	turns “off” calculation of surface area and volume
<b>MaxPropWeight</b>	<b>N</b>	sets limit of molecular weight for electrostatic charge, solvent and area/volume calculations; default is <b>500</b> amu
<b>Polarizability</b>	<b>F</b>	calculates and prints either static (F not specified) or frequency dependent (frequency F in $\text{cm}^{-1}$ ) polarizabilities and hyperpolarizabilities
<b>FreqScale</b>	<b>F</b>	scales calculated vibrational frequencies by F; default is <b>1</b>
<b>NumericalFreq</b>	–	uses numerical differentiation (of analytical gradients) for frequency calculation
<b>UVstates</b>	<b>N</b>	specifies number of states to be examined in UV/visible spectra calculation; default is <b>5</b>
<b>IRCSteps</b>	<b>N</b>	maximum number of steps allowed in constructing IRC; default is <b>40</b>
<b>IncludeTriplets</b>		include excited triplet states in a UV spectra calculation where the ground state is a singlet. Ignore where the ground state is a triplet.
<b>IRCStepSize</b>	<b>N</b>	maximum step size (in mass-weighted atomic units x 100) for IRC; default is <b>150</b>
<b>DropVibs</b>	<b>N</b>	ignore all vibrational modes with frequencies below N $\text{cm}^{-1}$ in calculating thermodynamic products
<b>ClampThermo</b>		do not allow individual vibrational contributions to the temperature correction to the enthalpy to exceed 1/2 RT, to the entropy to exceed 1/2 R, and for the heat capacity to exceed 1/2 R.
<b>AvgMass</b>		use average mass instead of mass of dominant isotope in calculating thermodynamic properties

---

## Miscellaneous and Printing

keyword	value	description
<b>3DPlot</b>	–	treat dynamic constraints independently (rather than in concert). This allows data for 3D (XYZ) plot to be collected. Two dynamic constraints only.
<b>Execute</b>	<b>C</b>	either skip execution ( <b>skip</b> ) or force execution regardless of history ( <b>force</b> )
<b>PrintVibCoords</b>	–	print coordinates corresponding to normal mode vibrational frequencies
<b>PrintCoords</b>	–	prints Cartesian coordinates
<b>PrintLev</b>	<b>C</b>	print level; default is <b>normal</b> which may be changed to <i>verbose</i>
<b>Print OrbE</b>	–	print energies of molecular orbitals
<b>PruneVirtual</b>	<b>N,C</b>	keep N virtual orbitals; default is <b>10</b> ; <b>None</b> keeps all virtual orbitals



# Appendix D

## Units

---

### **Geometries**

Cartesian coordinates are given in Ångstroms (Å), and in atomic units (au).

Bond distances are given in Å and in au. Bond angles and dihedral angles are given in degrees (°).

Surface areas, accessible surface areas and polar surface areas are available in Å<sup>2</sup> and volumes in Å<sup>3</sup>, and in au<sup>2</sup> (au<sup>3</sup>).

$$1 \text{ \AA} = 0.1 \text{ nm} = 1.889762 \text{ au}$$

### **Energies, Heats of Formation and Strain Energies, Zero-Point Energies, Enthalpies and Gibbs Energies and Entropies**

Total energies from Hartree-Fock calculations are available in au, kcal/mol, kJ/mol and electron volts (eV).

Experimental heats of formation as well as those from semi-empirical calculations and from thermochemical recipes are available in kJ/mol, au, kcal/mol and eV.

Strain energies from molecular mechanics calculations are available in kJ/mol, au, kcal/mol and eV.

Energies, heats of formation and strain energies corrected empirically for the effects of aqueous media are given in the same units as the corresponding gas-phase quantities.

Zero-point energies, enthalpies and Gibbs energies available in kJ/mol, kcal/mol and au/mol. Entropies are available in kJ/mol·degree, kcal/mol·degree and au/mol·degree.

### **Orbital Energies**

Orbital energies are available in eV, kcal/mol, kJ/mol and au.

## Energy Conversions

	au	kcal/mol	kJ/mol	eV
1 au	-	627.5	2625	27.21
1 kcal/mol	1.593 (-3)	-	4.184	4.337 (-2)
1 kJ/mol	3.809 (-4)	2.390 (-1)	-	1.036 (-2)
1 eV	3.675 (-2)	23.06	96.49	-

a) exponent follows in parenthesis, e.g., 1.593 (-3) = 1.593 x 10<sup>-3</sup>

### Electron Densities, Spin Densities, Dipole Moments, Charges, Electrostatic Potentials and Local Ionization Potentials

Electron densities and spin densities are given in electrons/au<sup>3</sup>.

Dipole moments are given in debyes.

Atomic charges are given in electrons.

Electrostatic and hydride potentials are given in kJ/mol.

Local ionization potentials are given in eV.

### Vibrational Frequencies, Infrared and Raman Intensities

Vibrational frequencies are given in wavenumbers (cm<sup>-1</sup>). Intensities are given relative to the most intense peak.

### Chemical Shifts, Coupling Constants

Chemical shifts are given in parts-per-million (ppm) relative to appropriate standards: hydrogen, tetramethylsilane; carbon, tetramethylsilane; nitrogen, nitromethane; fluorine, fluorotrichloromethane; silicon, tetramethylsilane, phosphorous, phosphoric acid. These are available for Hartree-Fock, BP, BLYP, B3LYP, EDF1 and EDF2 models with the 6-31G\*, 6-31+G\*, 6-31G\*\*, 6-31+G\*\*, 6-311G, 6-311+G\*, 6-311+G\*\* and 6-311++G\*\* basis sets. Coupling constants are in ppm.

### UV/visible Spectra

$\lambda_{\max}$  is given in wavenumbers (cm<sup>-1</sup>).

# Appendix E

## Citation

---

The proper citation for *Spartan'10* is as follows:

*Spartan'10*  
Wavefunction, Inc.  
Irvine, CA

Except for molecular mechanics and semi-empirical models, the calculation methods used in *Spartan* have been documented in: Y. Shao, L.F. Molnar, Y. Jung, J. Kussmann, C. Ochsenfeld, S.T. Brown, A.T.B. Gilbert, L.V. Slipchenko, S.V. Levchenko, D.P. O'Neill, R.A. DiStasio Jr., R.C. Lochan, T. Wang, G.J.O. Beran, N.A. Besley, J.M. Herbert, C.Y. Lin, T. Van Voorhis, S.H. Chien, A. Sodt, R.P. Steele, V.A. Rassolov, P.E. Maslen, P.P. Korambath, R.D. Adamson, B. Austin, J. Baker, E.F.C. Byrd, H. Dachsel, R.J. Doerksen, A. Dreuw, B.D. Dunietz, A.D. Dutoi, T.R. Furlani, S.R. Gwaltney, A. Heyden, S. Hirata, C-P. Hsu, G. Kedziora, R.Z. Khalliulin, P. Klunzinger, A.M. Lee, M.S. Lee, W.Z. Liang, I. Lotan, N. Nair, B. Peters, E.I. Proynov, P.A. Pieniazek, Y.M. Rhee, J. Ritchie, E. Rosta, C.D. Sherrill, A.C. Simmonett, J.E. Subotnik, H.L. Woodcock III, W. Zhang, A.T. Bell, A.K. Chakraborty, D.M. Chipman, F.J. Keil, A. Warshel, W.J. Hehre, H.F. Schaefer, J. Kong, A.I. Krylov, P.M.W. Gill and M. Head-Gordon, *Phys. Chem. Chem. Phys.*, **8**, 3172 (2006).

A discussion and assessment of commonly-used calculation methods is found in: W.J. Hehre, *A Guide to Molecular Mechanics and Quantum Chemical Calculations*, Wavefunction, Inc., Irvine, CA, 2003. This is available as a PDF under the **Help** menu.





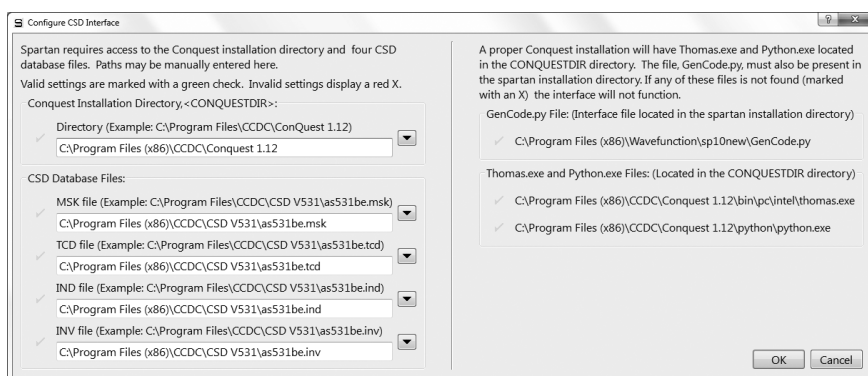
# Appendix F

## Installing the Cambridge Structural Database

Access to the Cambridge Structural Database (CSD) from *Spartan* should be automatic, following successful installation of the ConQuest program (supplied with CSD) according to the detailed instructions provided. If access is denied, the following steps need to be taken:

### Windows:

- (i) *Click* on **Config** near the bottom right of the dialog. This brings up the **Configure CSD Interface** dialog,



This will list any files which are missing. *Click* on **Setup Paths** at the bottom of the dialog to bring up the **CSD File Paths** dialog,

This dialog may show one or more red lights.

- (ii) Make certain that the path to ConQuest is correctly specified. The most common error will be that the version number is not current. (At the time of writing, the current version is 1.12.) You can get the current version number from your installation of ConQuest.

- (iii) Given the proper version number for ConQuest, update the file names in the boxes below **MSK file**, **TCD file**, **IND file** and **INV file**. You can get the current file names from your installation of ConQuest. As you enter the current file names, the red lights turn to green lights. When all lights are green, *click* on **OK** to dismiss the **CSD File Paths** dialog and then *click* on **OK** to dismiss the **Configure CSD** dialog.

**Linux:**

- (i) **CSD HOME** must be defined in order for *Spartan* to locate the CSD program.
- (ii) Please refer to the “Cambridge Structural Database System Release and Installation Notes” for instructions on setting up **CSD HOME**.

# Appendix G

## Constructing Custom Databases

---

Custom collections may be added to SMD, SIRD and SRD. These may either include the geometry, energy (heat of formation, strain energy), HOMO and LUMO energies, dipole moment, atomic charges, chemical shifts and vibrational frequencies obtained from quantum chemical calculations, or (for SMD and SIRD only) conformer libraries generated from molecular mechanics calculations. Entries for SRD (transition states) need to derive from *Spartan's* transition state guessing facility (**Guess Transition States** under the **Search** menu; **Chapter 19**) in order that they include the original set of reaction arrows.

Molecules (or lists of molecules) that are to be placed in a user database are first saved using **Save As** under the **File** menu with **Save as type** set to **Spartan Database**. This results in two files, a **.spentry** and a **.spindex** file for each molecule (or list). These should then be put in a single directory, the path to which must be entered into the **Paths Preferences** dialog under **Preferences...** in the **Options** menu (**Chapter 20**).

Note that it is important to specify a label for each molecule. Label specification is done either in the **Molecule Properties** dialog (see **Properties** under the **Display** menu in **Chapter 18**) or in the spreadsheet (see **Spreadsheet** under the **Display** menu in **Chapter 18**).



# Appendix H

## Pharmacophore Input

---

The following file format is used for input of a pharmacophore (illustrated for a pharmacophore with four hydrophobes and one positive ionizable center) :

```
Blank line
Blank line
Blank line
ENDCART
BEGINCFD
HPHOBE|RING_AROM  8.050 -4.180 -0.400  0  2.000  1.000  1
HPHOBE|RING_AROM  9.630 -0.640 -3.590  1  2.000  1.000  1
HPHOBE|RING_AROM 10.280 -0.380  3.200  2  2.000  1.000  1
HPHOBE|RING_AROM  8.130 -2.300  4.330  3  2.000  1.000  1
POS                14.370 -1.700  0.280  4  3.000  1.000  1
ENDCFD
```

The lines between BEGINCFD and ENDCFD define the individual pharmacophore elements (one line/element):

type1|type2... (| separates types for multiple definitions), Cartesian (X,Y,Z) coordinates (in Angstroms), unique number, radius of element (in Angstroms), weight of element in fit, use the element in fit (1) or ignore the element (5)

The following element types are available:

HPHOBE	hydrophobe
RING_AROM	aromatic ring
POS	positive ionizable center
NEG	negative ionizable center
HBA_CENTER	hydrogen bond acceptor
HBD_CENTER	hydrogen bond donor
EXCL_VOL	excluded volume



# Appendix I

## Input of Experimental Spectra

---

The databases of infrared and UV/visible spectra accessible online from *Spartan* are in JCAMP format. Specification of both and examples can be found at the NIST website (NIST.gov).

The database of  $^{13}\text{C}$  spectra is in CML (Chemical Markup Language) format. A description may be found at <http://www.ch.ic.ac.uk/omf/cml>. A sample file is provided below.

```
<?xml version="1.0" encoding="ISO-8859-1"?>
<cml>
  <spectrum type="NMR">
    <metadataList>
      <metadata name="nmr:OBSERVENUCLEUS" content="13C"/>
    </metadataList>
    <peakList>
      <peak xValue="125.1399993896484" peakHeight="0.458" />
      <peak xValue="128.4499969482422" peakHeight="0.558" />
      <peak xValue="129.0099945068359" peakHeight="0.658" />
      <peak xValue="21.64999961853027" peakHeight="0.758" />
      <peak xValue="137.8399963378906" peakHeight="0.858" />
    </peakList>
  </spectrum>
</cml>Blank line
```

This contains five shifts (xValue) and five associated intensities (peakHeight).





# Appendix K

## <sup>13</sup>C Chemical Shift Corrections

---

Two different empirical schemes have been developed to correct calculated <sup>13</sup>C chemical shifts. The first, implemented in *Spartan*'08 for the B3LYP/6-31G\* density functional model, is based on the number and kinds of directly-bonded atoms, specifically, -H, -C(sp<sup>3</sup>), -C(sp<sup>2</sup>), -C(sp), -N(sp<sup>3</sup>), -N(sp<sup>2</sup>), -N(sp), -F, -S(sp<sup>3</sup>), -S(sp<sup>2</sup>) and -Cl. A total of 29 parameters was involved in a linear regression to an experimental dataset involving ~800 sp<sup>3</sup> carbons, ~1500 sp<sup>2</sup> and aromatic carbons and ~100 sp carbons. Because this scheme depends only on topography and not on three-dimensional structure, *changes* in <sup>13</sup>C shifts between stereoisomers are unaffected.

Corrected <sup>13</sup>C chemical shifts (<sup>13</sup>C) for sp<sup>3</sup>, sp<sup>2</sup> and sp carbons are described by the following linear functions, where <sup>13</sup>C is the chemical shift from the B3LYP/6-31G\* model, *nH*, *nC*, *nN*, *nO*, *nS*, *nF* and *nCl* are the number of single bonds to hydrogen, carbon, nitrogen, oxygen, sulfur, fluorine and chlorine, respectively, *n=C*, *n=N*, *n=O* and *n=S* are the number of double bonds to carbon, nitrogen, oxygen and sulfur.

$$\text{sp}^3: {}^{13}\text{C} = -2.5752 + 1.0497 \cdot {}^{13}\text{C} + 0.6081 \cdot n\text{H} - 0.3631 \cdot n\text{C} + 0.2715 \cdot n\text{N} \\ + 0.0613 \cdot n\text{O} - 3.9978 \cdot n\text{S} - 1.6804 \cdot n\text{F} - 9.7005 \cdot n\text{Cl}$$

$$\text{sp}^2: {}^{13}\text{C} = 4.7483 + 1.0113 \cdot {}^{13}\text{C} + 0.3977 \cdot n\text{H} + 0.0783 \cdot n\text{C} + 1.3293 \cdot n\text{N} \\ + 0.9851 \cdot n\text{O} - 4.1232 \cdot n\text{S} - 1.5531 \cdot n\text{F} - 8.8051 \cdot n\text{Cl} - 0.4127 \cdot n=\text{C} + \\ 1.7721 \cdot n=\text{N} + 4.5154 \cdot n=\text{O} + 4.4296 \cdot n=\text{S}$$

$$\text{sp}: {}^{13}\text{C} = 18.3286 + 0.8740 \cdot {}^{13}\text{C} - 2.4100 \cdot n\text{H} + 1.1251 \cdot n\text{C} - \\ 1.1723 \cdot n=\text{C} + 6.9164 \cdot n=\text{N}$$

This leads to an RMS errors of 2.0 ppm (reduced from 4.0 ppm), 2.4 ppm (reduced from 4.0 ppm), 2.6 ppm (reduced from 4.0 ppm) for sp<sup>3</sup>, sp<sup>2</sup> and sp for <sup>13</sup>C shifts, respectively.

The second scheme, implemented in *Spartan*'10 for the EDF2/6-31G\* density functional model, is based both on bond counts and on

bond orders. It individually parameterizes 19 important  $sp^2$  and  $sp$  carbon centers. This scheme involves 43 independent parameters for the general scheme (up from 29 due primarily to the inclusion of bonds to silicon, phosphorus and bromine) and 78 parameters (in 19 independent sets) for the special cases (15  $sp^2$  cases and 4  $sp$  cases). The dataset used for the linear regression is much larger, comprising ~8000  $sp^3$  carbons, ~6200  $sp^2$  carbons (including aromatic carbons) and ~450  $sp$  carbons.

The general case formula

$$^{13}\mathbf{C}_i = ^{13}\mathbf{C}_i \cdot \text{Scale} + \sum_k [X^{(0)}_k + X^{(1)}_k (Q_{i,k} - 1) + X^{(2)}_k (Q_{i,k} - 1)^2]$$

$^{13}\mathbf{C}$  is the corrected chemical shift and  $^{13}C$  is the uncorrected chemical shift and *Scale* is a scaling factor. Summation is carried out over all bonds (1 or 2 for  $sp$  carbons, 3 for  $sp^2$  carbons and 4 for  $sp^3$  carbons).  $X^{(0)}_k$ ,  $X^{(1)}_k$  and  $X^{(2)}_k$  are parameters that depend on the atom bonded to carbon (H, C, N, O, F, Si, P, S, Cl or Br) and  $Q_{i,k}$  are natural bond orders. Different expressions apply to  $sp$ ,  $sp^2$  and  $sp^3$  carbons. The equation for  $sp$  carbon involves four  $X^{(0)}$  parameters, four  $X^{(1)}$  parameters and one  $X^{(2)}$  parameter (11 parameters in total), that for  $sp^2$  carbon involves all ten  $X^{(0)}$  parameters, six  $X^{(1)}$  parameters and four  $X^{(2)}$  parameters (21 parameters in total), and that for  $sp^3$  carbon involves all ten  $X^{(0)}$  parameters but no  $X^{(0)}$  or  $X^{(1)}$  parameters (11 parameters in total). As there are no special cases for  $sp^3$  carbons, this means that the new correction scheme for  $sp^3$  carbons, like the previous scheme, depends only on topography and not on three-dimensional structure. Changes in chemical shifts between stereoisomers involving  $sp^3$  centers obtained from the EDF2/6-31G\* model are unaffected.

The special case formula is:

$$^{13}\mathbf{C}_i = ^{13}\mathbf{C}_i \cdot \text{Scale}_j + X^{(0)}_j + X^{(1)}_j \cdot \sum_k (Q_{i,k} - 1) + X^{(2)}_j \cdot \sum_k (Q_{i,k} - 1)^2]$$

The subscript *j* designates the particular special case.  $^{13}\mathbf{C}$  is the corrected chemical shift and  $^{13}C$  is the uncorrected chemical shift. *Scale* is a scaling factor unique to the particular special case and  $X^{(0)}_j$ ,  $X^{(1)}_j$  and  $X^{(2)}_j$  are parameters that depend only on the particular special case. Summations are carried out over all bonds (1 or 2 for  $sp$  carbons and 3 for  $sp^2$  carbons) and  $Q_{i,k}$  are natural bond orders.  $X^{(1)}_j$  or

both  $X^{(1)}_j$  and  $X^{(2)}_j$  are absent in some of the special case formulas.

RMS errors from the general case and special case formulas are 1.8 and 1.7 ppm, respectively.



# Appendix L

## Installing a Network HASP

---

*Spartan'10* requires a hardware license key (HASP) which is normally installed on each machine. However, a network hardware key (NetHASP) can be obtained from Wavefunction and installed on a server. This allows *Spartan'10* to be run on a specified number of client machines that are on the network. A one-time installation is required.\*

### On the Server:

1. The License Manager software must be installed. To do so, please run the **lmsetup.exe** file.

a. (Recommended). The latest release of this file can be obtained by visiting the following link at the Aladdin Web Site:

<http://www.aladdin.com/support/hasp/enduser.asp#lm>

b. The **lmsetup.exe** file may also be obtained directly from a subdirectory in the *Spartan* folder if you install the program directly onto your server . . . if you opt to install the program on your server, run through the installation procedure guided via the installation wizard. Upon completion, **lmsetup.exe** can be found in the following subdirectory:

SPARTAN'10\*\* : C:\...\Wavefunction\Spartan10Vxxx\Network\SERVERS\

2. It is suggested that you choose to configure the program to run as a 'service' during the installation procedure. When the program is finished installing, restart the server and verify that the NetHASP

---

\* The network HASP license management option requires a Windows machine to act as the license manager, client machines may be either Windows or Macintosh (NAL).

\*\* *Spartan* has multiple Windows release versions, the pathname depends on which version you are running (Spartan'10 v.1.0.0 program directory = Spartan10V100) – Each contains the Network and SERVERS folders. Your server must be running Windows 2000 or higher.

License Manager has loaded the necessary protocols for your network.

3. If you are running a Firewall, open ports 475 and 1947. (Used by default, but also configurable).

For additional server options, you can also configure an NHSRV.ini file. Such options include specifying which client machines can access Spartan by IP range.

### **On Each Client Computer:**

1. Install software locally on the network clients that will access the NetHASP through the network. Follow the installation procedure as guided by the installation wizard.

A network configuration file (nethasp.ini) must be added to the directories described below. Several variations of the file can be downloaded from the Aladdin website or you can use the template that is installed with Spartan:

Aladdin: <http://www.aladdin.com/support/hasp/enduser.asp#ini>

**Spartan:** The template nethasp.ini can be found in:

C:\...\Wavefunction\Spartan10Vxxx\Network\

2. Locate the NETHASP.ini file and open with a text editor (ex. Notepad). The line that requires modification is the “;;NHSEVER\_ADDR = <Addr1>, <Addr2>”. Delete the two semicolons (;) and delete the <Addr1>, <Addr2>. Insert the IP address of your server on the right of the =. If you have more than one NETHASP installed, subsequent server IP addresses may be entered by placing a comma after the first address, followed by a space. After these modifications are completed, save the file.

Next, right Click on this nethasp.ini file and select ‘Copy’. Move up one level to the **Spartan** folder. Right click on the background and select ‘Paste’. Open both the P2 and P4 folder and also right click and ‘Paste’ the NETHASP.ini file in these folders as well.

Need additional help? Visit the Windows or Macintosh FAQs, [www.wavefun.com](http://www.wavefun.com), support section.